



# **MICROSCOPY CATALOGUE**

# **OPTIKA** Milestones

V.		OPTIKA OPTIM
	1996	First <b>LED brightfield</b> microscope in the world
	2002	First <b>2-in-1 gemological</b> microscope in the world
	2002	First <b>LED darkfield</b> illuminator for stereomicroscopes in the world
	2003	First LED fluorescence microscope in the world
	2004	Infinifix (infinity optical system with 160mm objectives)
	2005	First <b>LED metallurgical</b> microscope
11	2006	X-LED illumination system - the original
	2007	First <b>LED darkfield</b> illuminator <b>for blood analisys</b> in the world
	2008	First <b>Digital Imaging Unit with integrated microscope camera</b> in the world
	2009	First full <b>Automatic Light Control</b> (ALC) in the world
	2010	Full motorized microscope controlled by mouse in the world
	2011	Full motorized microscope controlled by touch-screen monitor in the world
	2013	First digital microscope with integrated Digital Imaging Unit in the world
	2015	First <b>rechargeable microscope with Li-Ion battery</b> in the world
	2017	First educational microscope with oil/water 100x objective in the world
	2023	X-OEF (Optika EasyFocus System) Electronically controlled triple focus adjustment

# **INDEX**

1 EDUCATIONAL Microscopes	page	7
2 LABORATORY Microscopes	page	79
3 INSPECTION & INDUSTRIAL Microscopes	page	279
4 CAMERAS & DIGITAL Solutions	page	355
5 POLARIMETRY AND REFRACTOMETRY	page	437











# **Complete** Icon List

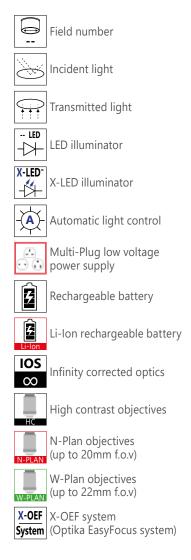

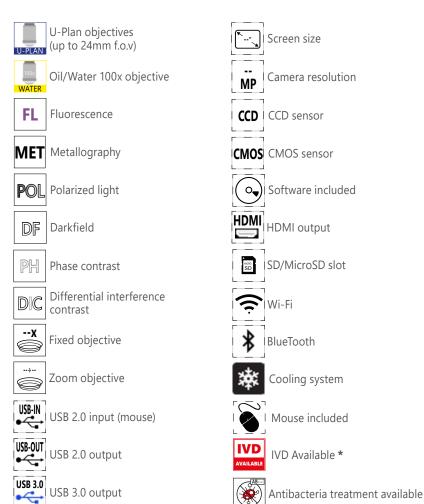

<sup>\*</sup> The IVD code must be requested at order

# **OPTIKA X-OEF** System



Optika is proud to introduce the first routine microscope with fully electronically controlled triple focus adjustment: Macro, Semi-Micro and Micro.

The new **X-OFF** system (Optika EasyFocus system) replaces the classic mechanical focusing system with an innovative solution that allows more functions.

EasyFocus allows the selection of three different resolutions (shift for each step):

#### Macro (250 μm); Semi-Micro (8 μm); Micro (1 μm)

The single knob allows the focus movement and the selection of the resolution is easily controlled by pressing it. Each focus resolution (Macro, Semi-Micro and Micro) is clearly displayed to the user by three dedicated blue LEDs.

The presence of a Semi-Micro adjustment with 8  $\mu$ m resolution makes finding the focus point with medium-low magnification objectives (10x,20x) exceptionally convenient and quick. The **X-OFF** system includes a stage retraction function that facilitates and speeds up immersion oil laying and slide changing operations.


B-820 is the result of OPTIKA Microscopes' long experience in the field of optical microscopy, offering an extremely valuable product for routine and research applications in brightfield and phase contrast.

Personal requirements are met in terms of application as well as ergonomics thanks to the microscope's extensive modularity.

The innovative B-820 model is a particularly high-performance and robust solution, offering 22mm field of view, the unique **X-LED**<sup>3</sup> illumination source (3.6 W) and a sturdy dye-cast structure for high stability, combined with a wide choice of heads, objectives and condensers.



Two different stands and mechanical stages are available for both right and left-handed users.

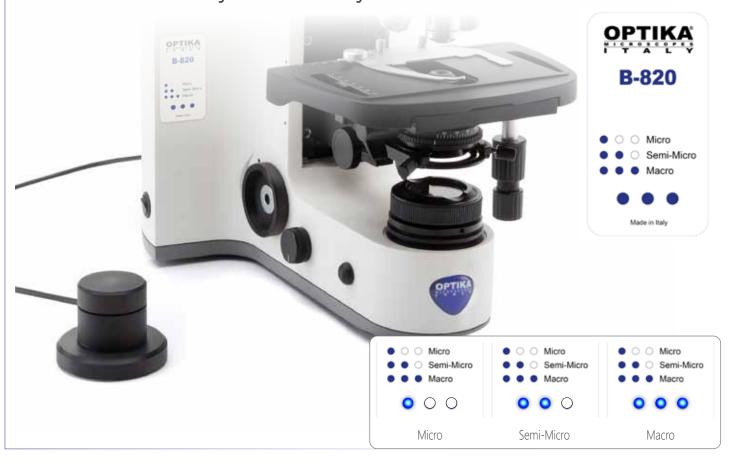

# **New IM-300 Series**



OPTIKA is pleased to present the brand new range of *IM-300* inverted microscopes.

These models, completely redesigned, offer exceptional working comfort and routine use of the highest level.

Inverted microscopes are useful for observing living cells or organisms at the bottom of a large container (for example, a tissue culture flask) under more natural conditions than on a glass slide, as is the case with a conventional microscope.

The *IM-300 series* includes a version for the simultaneous brightfield and phase contrast method, designed and engineered to be the ideal solution for fast and reliable routine inspections.

The round glass stage surface allows for optimal visual access to the objective turret.

A particularly simple and ingenious optical design enables stable alignments and smooth, precise movements.

The *IM-300 range* is completed with digital versions equipped with a large 4k touch screen monitor and a very high sensitivity camera.



# **OPTIKA AB-Shield**



Optika **AB-Shield** is a coating that allows you to treat different surfaces and guarantee the elimination and prevention of bacterial growth in hygiene-sensitive areas.

This treatment guarantees you a strong antibacterial effect due to the highly effective micro-silver additive and Destruction of 99.9% of surface bacteria.

The antibacterial effect remains for the lifetime of the coating. Also eliminates bad smells caused by bacteria. Proven effective against resistant bacteria such as MRSA, E-Coli, EHEC etc.



Optika **AB-Shield** is available as an add-on to your new microscope, thus making its use healthier for yourself and when it needs to be shared among several people.



The effectiveness was successfully tested using the international standard:

#### ISO 22196

"Measurement of antibacterial activity on plastics and other non-porous surfaces" against these pathogens:

Staphylococcus aureus DSM 346/ATCC 6538P Escherichia coli DSM 1576/ATCC 8739.

The method has been tested for bacteriostatic properties (growth inhibitor) and bactericidal (killing bacteria).






# **EDUCATIONAL Microscopes**

# **EDUCATIONAL** Microscopes

# Biological Microscopes ECOVISION SERIES - Entry-Level Biological Microscopes For Beginners page 9 B-60 SERIES - Entry-Level Biological Microscopes For Students page 17 B-150 SERIES - Middle-Level Biological Microscopes For Students page 25 B-190 SERIES - Advanced Biological Microscopes For Students And Teachers page 49 Stereomicroscopes MS/SFX SERIES - Entry-Level Monoscopes & Stereomicroscopes For Students page 61 SLX SERIES - Stereomicroscopes For Students And Teachers page 71

# Icons

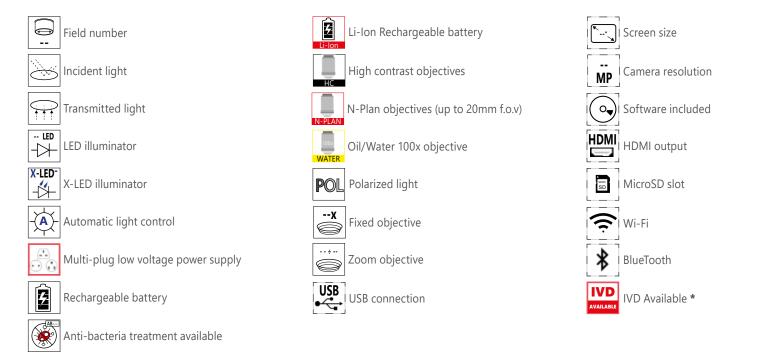

<sup>\*</sup> The IVD code must be requested at order




# ECOVISION Series



**Entry-Level Biological Microscopes For Beginners** 

# A Range Of Quality Microscopes For Beginners

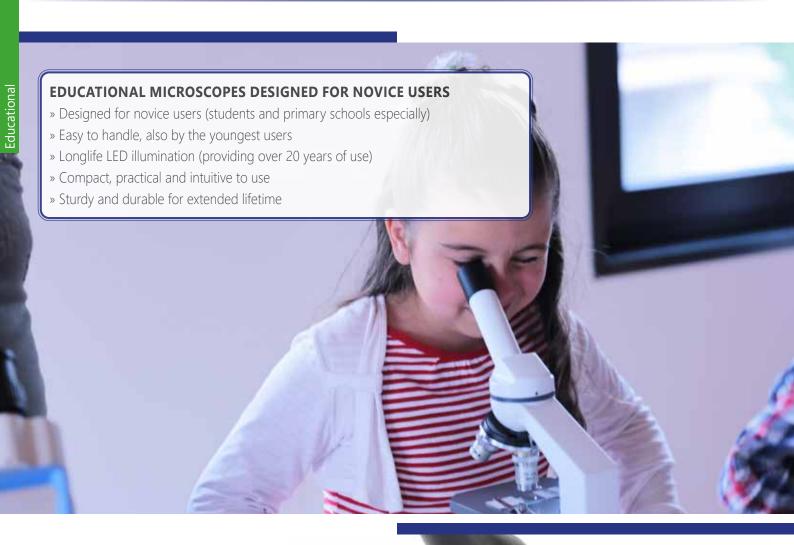

#### **COMFORTABLE, INTUITIVE & RELIABLE SOLUTIONS**

- » 18 mm field number for a wide observation area
- » Achromatic optics ensuring good contrast and quality images
- » Pre-aligned illumination and condenser to simplify operations
- » Cordless use, totally independent from mains/batteries connection
- » External power supply for enhanced safety and convenient servicing





# **ECOVISION** Series

A range of mainly cordless monocular microscopes ideal for students and mainly primary schools with achromatic lenses, FN 18 eyepiece, finite optical system, coaxial focusing, fixed or mechanical stage and 1 W LED illumination.

Slim and easy to carry, all the models are equipped with long lasting LED illumination to provide over 20 years of use.



#### Easy to handle, also by the youngest users

Extreme compactness and portability to ensure easy transportation in the classroom and outdoor, with slim body and useful handle for a facilitated and enjoyable teaching activity.

#### Cordless use, totally independent from mains connection

Most of the models work with or without the batteries in place. For outdoor use (4-hour autonomy, at medium intensity), three NiMH rechargeable batteries must be used.



# External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit.

#### Longlife LED illumination (providing over 20 years of use)

Money & energy saving thanks to LED long lifetime (65.000 hours, 22 years in case of 8 hours/day) which is more than 20 times compared to a standard halogen bulb.

#### (1)

# **ECOVISION** Series - B-20 Models

## **B-20R**















Cordless monocular microscope ideal for students and mainly primary schools, with achromatic lenses (400x), FN 18 eyepiece, finite optical system, round fixed stage and 0.3 W LED illumination with rechargeable batteries (not provided).

Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

Eyepiece: WF10x/18 mm.

Nosepiece: Triple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** X-Y moving and 360° rotating, 90 mm diameter, with sample clips.

**Focusing:** Separate coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Diffusing filter with rotating diaphragm wheel.

**Illumination:** 0.3 W LED, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-20CR**















Cordless monocular microscope ideal for students and mainly primary schools, with achromatic lenses (400x), FN 18 eyepiece, finite optical system, coaxial focusing, mechanical stage and 0.5 W LED illumination with rechargeable batteries (not provided).

Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

**Nosepiece:** Triple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** Double layer, 105x95 mm, moving range 50x15 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** N.A. 0.65 with iris diaphragm.

**Illumination:** 0.5 W LED, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **ECOVISION** Series - M-100 Models

## M-100FLed















Cordless monocular microscope ideal for students and mainly primary schools, with achromatic lenses (400x), FN 18 eyepiece, finite optical system, coaxial focusing, fixed stage and 0.5 W LED illumination with rechargeable batteries (not provided).

Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

Nosepiece: Triple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** Fixed stage, 120x110 mm, with sample clips.


**Focusing:** Separate coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

Condenser: N.A. 0.65 with iris diaphragm.


**Illumination:** 0.5 W LED, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **ECOVISION** Series - Comparison chart


Model	Head	Eyepieces	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
B-20R	Monocular, 45° inclined, 360° rotating	WF 10x/18	Triple, reversed	Achromatic 4x, 10x, 40x	X-Y moving, 360° rotating, 90 mm diameter, with sample clips	Separate coarse and fine	Diffusing filter with rotating diaphragm wheel	0.3 W LED, with brightness control, rechargeable batteries (not provided)
B-20CR	Monocular, 45° inclined, 360° rotating	WF 10x/18	Triple, reversed	Achromatic 4x, 10x, 40x	Double layer, 105x95 mm, moving range 50x15 mm	Coaxial coarse and fine	N.A. 0.65, with iris diaphragm	0.5 W LED, with brightness control, rechargeable batteries (not provided)
M-100FLed	Monocular, 45° inclined, 360° rotating	WF 10x/18	Triple	Achromatic 4x, 10x, 40x	Fixed, 120x110 mm, with sample clips	Separate coarse and fine	N.A. 0.65, with iris diaphragm	0.5 W LED, with brightness control, rechargeable batteries (not provided)

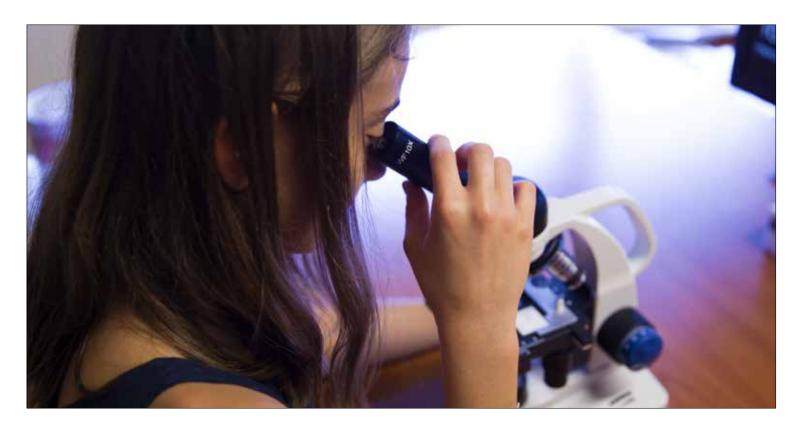


# **ECOVISION** Series - **Zoom comparison**

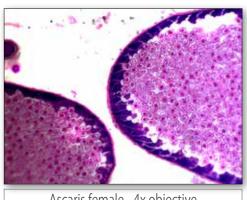

Ascaris female - 4x objective

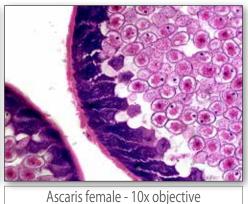

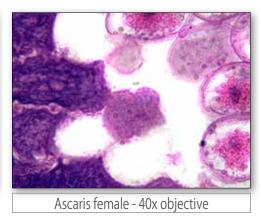

# **ECOVISION** Series - Accessories

#### **ACCESSORIES FOR B-20R / B-20CR**

#### **Eyecups & Eyepieces**

M-002.2 WF10x/18 eyepiece (23mm Ø)

M-003.2 WF15x/12 eyepiece (23mm Ø)

M-004.2 WF10x/18 micrometric eyepiece (23mm Ø)

M-008.2 WF10x/18 eyepiece, pointer (23mm Ø)

M-162 WF20x/10 eyepiece (23mm Ø)

#### **Additional Lenses**

M-114 0.35x C-Mount projection lens

M-115 0.5x C-Mount projection lens

M-118 0.75x C-Mount projection lens

#### Miscellaneous

15104 Cleaning kit

M-005 Micrometric slide, 26x76mm, 2 scales (1mm/100 & 10mm/100)

M-069 Solar charger

DC-001 Plastic dust cover, small, 340(l)x400(h) mm

AB-010 Antibacterial surface treatment, only for newly purchased microscope

#### **ACCESSORIES FOR M-100FLed**

#### **Eyecups & Eyepieces**

M-001 Huygens 5x eyepiece (23mm Ø)

M-002.2 WF10x/18 eyepiece (23mm Ø)

M-004.2 WF10x/18 micrometric eyepiece (23mm Ø)

M-008.2 WF10x/18 eyepiece, pointer (23mm Ø)

M-003.2 WF15x/12 eyepiece (23mm Ø)

M-162 WF20x/10 eyepiece (23mm Ø)

#### **Objectives**

M-131 Achromatic objective 4x/0.10

M-132 Achromatic objective 10x/0.25

M-133 Achromatic objective 20x/0.40

M-134 Achromatic objective 40x/0.65

M-135 Achromatic objective 60x/0.85

M-136 Achromatic objective 100x/1.25 (oil)

#### **Camera Adapters**

M-115 0.35x C-Mount projection lens

M-114 0.5x C-Mount projection lens

M-118 0.75x C-Mount projection lens

#### Stages

M-040 Attachable mechanical stage

#### **Condensers & Filters**

M-099 Polarising set (filters and rotating stage)

#### Miscellaneous

15008 Immersion oil, 10ml

15009 Immersion oil, 100ml

15104 Cleaning kit

DC-001 Plastic dust cover, small, 340(l)x400(h) mm

M-005 Micrometric slide, 26x76mm, 2 scales (1mm/100 & 10mm/100)

M-069 Solar charge

AB-010 Antibacterial surface treatment, only for newly purchased microscope



#### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. Output voltage: 5 Vdc. -

Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>).

Charging models: with solar panel (12h), with external USB power supply (2.5h).



#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.



 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# B-60 Series



**Entry-Level Biological Microscopes For Students** 

# Cordless Educational Microscopes, Ideal To Start Exploring

#### PERFECT FOR STUDENT'S FIRST EXPERIENCES

- » Designed for novice users (students and primary schools especially)
- » Easy to handle, also by the youngest users
- » Longlife LED illumination (providing over 20 years of use)
- » Compact, practical and intuitive to use
- » Sturdy and durable for extended lifetime

#### **COMFORTABLE, INTUITIVE & RELIABLE SOLUTIONS**

- » 18 mm field number for a wide extended observation area
- » StagErase™ eraseable stage to reduce scratches
- » Arm/wrist rest support to reduce the fatigue during use
- » Cordless use, totally independent from mains/batteries connection
- » External power supply for enhanced safety and convenient servicing





(1)



# **B-60** Series

A wide range of cordless, modern microscopes ideal for students and mainly primary schools with achromatic lenses, FN 18 eyepieces, finite optical system, coaxial focusing, StagErase™ eraseable fixed or mechanical stage and 1 W LED illumination with rechargeable batteries (not provided). Slim and easy to carry, all the models are equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use



#### Arm/wrist rest support to reduce the fatigue during use

Students get relaxed and stay relaxed when using the microscope! Effective in preventing fatigue during operation, increasing the ergonomy and the performance as a result

#### StagErase™ eraseable stage to remove scratches

Here's something you've never seen before! This new, revolutionary stage is coated with a special painting to reduce accidental scratches to the minimum and facilitate their removal



#### Cordless use, totally independent from mains/batteries connection

All models work with or without the batteries in place and are for outdoor use (4-hour autonomy, at medium intensity)

# External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage multiplug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit

#### Longlife LED illumination (providing over 20 years of use)

Money & energy saving thanks to LED long lifetime (65.000 hours, 22 years in case of 8 hours/day) which is more than 20 times compared to a standard halogen bulb

### U

# **B-60** Series - Models

## **B-61**















**B-62** 

Cordless, modern monocular microscope ideal for students and mainly primary schools with achromatic lenses (400x), FN 18 eyepiece, finite optical system, coaxial focusing, StagErase™ fixed stage and 1 W LED illumination with rechargeable batteries (not supplied).

Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase<sup>™</sup> eraseable fixed stage, 120x110 mm, with sample clips.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

Condenser: N.A. 0.65 with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.















Cordless, modern monocular microscope ideal for students and mainly primary schools with achromatic lenses (400x), FN 18 eyepiece, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

years of use.

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase™ eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-60** Series - Models

## **B-63**















**B-65** 













Cordless, modern monocular microscope ideal for students and mainly primary schools with achromatic lenses (600x), FN 18 eyepiece, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.
- Achromatic 60x/0.85, W.D. 0.13 mm.

All with anti-fungus treatment

**Specimen stage:** StagErase™ eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

Cordless, modern monocular microscope ideal for students and mainly primary schools with achromatic lenses (1000x), FN 18 eyepiece, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 45° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.
- Achromatic 100x/1.25 (Oil), W.D. 0.13 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase™ eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-60** Series - Models

## **B-66**



Cordless, modern binocular microscope ideal for students and mainly primary schools with achromatic lenses (400x), FN 18 eyepieces, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

**Head:** Binocular, 30° inclined; 360° rotating. Dioptric adjustment: Left eyepiece.

**Eyepiece:** WF10x/18 mm, secured by screw

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase<sup>™</sup> eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

**B-67** 



Cordless, modern binocular microscope ideal for students and mainly primary schools with achromatic lenses (600x), FN 18 eyepieces, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

**Head:** Binocular, 30° inclined; 360° rotating. Dioptric adjustment: Left eyepiece.

**Eyepiece:** WF10x/18 mm, secured by screw

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.
- Achromatic 60x/0.85, W.D. 0.13 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase<sup>™</sup> eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-60** Series - Models

## **B-69**



Cordless, modern binocular microscope ideal for students and mainly primary schools with achromatic lenses (1000x), FN 18 eyepieces, finite optical system, coaxial focusing, StagErase™ mechanical stage, Abbe condenser and 1 W LED illumination with rechargeable batteries (not supplied). Slim and easy to carry, it is equipped with arm/wrist rest support to reduce the fatigue during use and long lasting LED illumination to provide over 20 years of use.

**Head:** Binocular, 30° inclined; 360° rotating. Dioptric adjustment: Left eyepiece.

**Eyepiece:** WF10x/18 mm, secured by screw

Nosepiece: Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- Achromatic 4x/0.10, W.D. 10.6 mm.
- Achromatic 10x/0.25, W.D. 7.0 mm.
- Achromatic 40x/0.65, W.D. 0.5 mm.
- Achromatic 100x/1.25 (Oil), W.D. 0.13 mm.

All with anti-fungus treatment.

**Specimen stage:** StagErase™ eraseable mechanical stage, 125x125 mm, 62x24 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** 1 W LED, with brightness control, rechargeable batteries. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **Digital Video Bundles**



Two models of the B-60 series can be equipped with a camera and a 7" LCD screen, both in high definition. In these bundles, the normal head (supplied with the instrument) can be replaced with the digital system in a few minutes. This solution provides a system suitable for viewing specimens by several students at the same time, without removing the possibility of using the microscope in the classical way through the eyepiece.

LCD screen: High definition 7" LCD.

**Camera:** 1920x1080 pixels, 30fps (video). Up to 1844x1080 pixels (photo).

Storing capacity: On Micro Sd card


Video recording: Yes


**Measuring funcion:** Yes, simple line measurement

#### Models available:

**-B-61V:** Same fatures than standard B-61, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.** 

**-B-62V:** Same fatures than standard B-62, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.** 

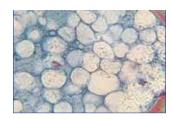

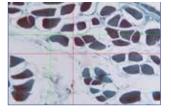

# **B-60** Series - Comparison chart

Model	Head	Eyepiece(s)	Nosepiece	Objectives	Stage	Focusing	Condenser	Illuminator
B-61	Monocular, 360° rotating, 45° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x	Fixed, 120x110 mm	Coaxial coarse and fine focusing	N.A. 0.65 with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-62	Monocular, 360° rotating, 45° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-63	Monocular, 360° rotating, 45° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x, 60x	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-65	Monocular, 360° rotating, 45° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x, 100x (oil)	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-66	Binocular, 360° rotating, 30° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-67	Binocular, 360° rotating, 30° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x, 60x	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries
B-69	Binocular, 360° rotating, 30° inclined	Wide Field 10x/18mm	Quadruple, reversed	Achromatic 4x, 10x, 40x, 100x (oil)	Mechanical, 125x125 mm with 62x24 mm X-Y moving range	Coaxial coarse and fine focusing	N.A. 1.25 Abbe type with iris diaphragm	1 W LED, brightness control, rechargeable batteries

## **B-60** Series - Accessories

#### **Eyecups & Eyepieces**

M-001 Huygens 5x eyepiece (23mm Ø) M-002.2 WF10x/18 eyepiece (23mm Ø)

M-004.2 WF10x/18 micrometric eyepiece (23mm Ø)

M-008.2 WF10x/18 eyepiece, pointer (23mm Ø)

M-003.2 WF15x/12 eyepiece (23mm Ø)

#### **Condensers & Filters**

M-155.2 Polarising set (filters only)

#### **Camera Adapters**

M-115 0.35x C-Mount projection lens M-114 0.5x C-Mount projection lens

M-118 0.75x C-Mount projection lens

#### Miscellaneous

<u>15008</u> Immersion oil, 10ml 15009 Immersion oil, 100ml <u>15104</u> Cleaning kit DC-001 Plastic dust cover, small, 340(l)x400(h) mm M-005 Micrometric slide, 26x76mm, 2 scales (1mm/100 & 10mm/100) M-069 Solar charger M-970 Plane-concave mirror, with base (only for B-61) Antibacterial surface treatment, only for newly purchased microscope AB-010

#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.



#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. Output voltage: 5 Vdc. -

Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>).

Charging models: with solar panel (12h), with external USB power supply (2.5h)



#### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA® S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

OPTIKA<sup>®</sup> North America OPTIKA<sup>®</sup> Central America OPTIKA<sup>®</sup> Africa namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# B-150 Series



**Middle-Level Biological Microscopes For Students** 

1

# The Most Comprehensive Series Dedicated To Students

#### A VARIETY OF CONFIGURATIONS TO MEET EVERY NEEDS

- » Designed to fulfill primary/secondary schools and educational labs
- » 18 mm field number for a wide observation area
- » Cordless use, totally independent from mains/battery connection (R-PL Line)
- » Sturdy and durable for extended lifetime; compact and intuitive
- » External power supply for enhanced safety and convenient servicing






# Inventors Of A New Way To Teach Microscopy

#### 100X OIL/WATER OBJECTIVE - ONLY AVAILABLE AT OPTIKA

- » Same objective for oil and water use
- » Oil represents the best media for high numerical aperture
- » Water combines results with convenience for educational purposes
- » Save time forget about tedious cleaning and maintenance
- » Save money no additional expenses due to inappropriate cleaning

100xOil

100x Water

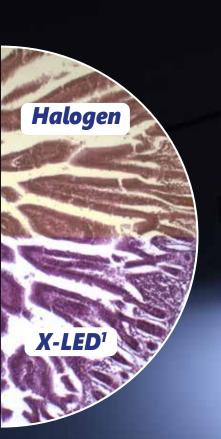

#### X-LED¹ FOR 65,000 HOURS OF OPERATION - ONLY AVAILABLE AT OPTIKA

- » State-of-the-art illumination system for incomparable light intensity
- » Exclusive lens & collector design, unmatched uniformity & brightness
- » Excellent color fidelity, constant pure-white color temperature
- » Money & energy saving, cutting electricity bills by 90%
- » Simple polarization versions with polarizer and analyzer (P Line)



# 1

# Plenty Of Smart & Innovative Light-Related Technologies

#### **AUTOMATIC LIGHT CONTROL - ONLY AVAILABLE AT OPTIKA (ALC LINE)**

- » Choose the light intensity according to your preference
- » Press the ALC button and the light will be automatically re-adjusted
- » When another objective is used
- » When the diphragm aperture changes
- » When processing samples with different opacity











#### STEP 3

#### Forget about the illumination!

The microscope will automatically adjust the brightness for you, in case of:

- Another objective is used
- The diaphragm aperture is changed
- Another specimen with different opacity is processed



Regulation of diaphragm aperture

# Plenty Of Smart & Innovative Light-Related Technologies

#### LI-ION BATTERIES PROS (on B-150R models):

- » **Reliable:** Significantly lower self-discharge rate than NiMH
- » Faster recharge: Li-lons can be charged in about 6 hours
- » Number of charges: approx. 2,000 times (+100% than NiMH batteries)
- » No "memory effect": can be charged at any time, without effects
- » Temperature tolerance to low temperature (more than NiMH batteries)



1

# **B-150** Series

A very comprehensive range of modern microscopes ideal for students and primary/secondary schools, available in brightfield or polarized light. Provided with achromatic or PLAN achromatic lenses, FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, fixed or mechanical stage and powerful, uniform, white color temperature 1 W X-LED1 illumination. Slim and easy to carry, all the models are equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

#### A variety of configurations to meet every needs

Configurations for every taste, including regular brightfield and the one ready for polarization analysis (P Line), automatic light control (ALC Line), with built-in cameras for image acquisition (D Line) and cordless versions with advanced features (R-PL Models)

#### High eyepoint eyepieces for glasses wearers

These eyepieces are designed in such a way that the exit pupil is further away from the eye lens than standard eyepieces, being are well suited for eyeglasses wearers

#### X-LED¹ - State-of-the-art illumination system for incomparable light intensity

Provided with an exclusive lens & collector design, OPTIKA X-LED technology ensures unmatched uniformity & brightness (more than a 20 W halogen lamp) for excellent color fidelity with constant pure-white color temperature



1



#### 100x oil/water objective: same objective for dual use

This new, revolutionary objective is something you've never seen before! Oil ensures the best performance achievable; water represents the most convenient solution as eliminates tedious cleaning

#### Incomparable comfort with the exclusive Automatic Light Control (ALC Line)

Light intensity is automatically adjusted by the microscope: no matter if the aperture of the diaphragm changes, if another objective is used, and if the opacity of the sample is different...the microscope will set the light!



#### N-PLAN objectives combined with exclusive Li-ion battery (R-PL line)

Laboratory grade optics meets the latest technology in terms of battery, for unparalleled lifetime (2000 charges), extended autonomy (15 hours/charge) and incredibly fast recharging time (6 hours)



#### External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit

#### GET THE MOST OUT OF OUR ACCESSORIES

#### M-974 - Blue filter

Increase the colour temperature of light (toward the blue)

#### M-976 - Green filter

Optimize the resolution of phase contrast

#### M-978 - Yellow filter

Decrease the colour temperature of light (toward the red)

#### M-988 - Frosted glass filter

Increase the uniformity of illumination, even further

#### U

# **B-150** Series - Standard Models

## **B-151**















Monocular microscope ideal for students and primary schools, with three achromatic lenses (400x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, fixed stage and powerful, uniform, white color temperature settable 1 W *X-LED*<sup>1</sup> illumination.

Slim and easy to carry, the LED illumination will provide over 20 years of use.

Head: Monocular, 30° inclined; 360° rotating.

Eyepiece: WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.

All with anti-fungus treatment.

Specimen stage: Fixed stage, 130x120 mm. With sample clips.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** N.A. 0.65, pre-centered, fixed with iris diaphragm.

**Illumination:** X-LED1 with white 1 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

## **B-151ALC**

















Monocular microscope ideal for students and primary schools, with three achromatic lenses (400x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, fixed stage and powerful, uniform, white color temperature settable 1 W **X-LED**<sup>†</sup> illumination. Slim and easy to carry, the LED illumination will provide over 20 years of use of use. The exclusive ALC will automatically adjust the brightness according to your preferences

**Head:** Monocular, 30° inclined; 360° rotating (when ALC cable is

**Eyepiece:** WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.

All with anti-fungus treatment.

**Specimen stage:** Fixed stage, 130x120 mm. With sample clips.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** N.A. 0.65, pre-centered, fixed with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. With **ALC** for automatic light control. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-150** Series - Standard Models

## **B-151R-PL**

















Cordless monocular microscope ideal for students and primary schools, with three PLAN achromatic lenses (400x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, fixed stage and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, the LED illumination will provide over 20 years of use of use. The Li-Ion battery ensures unparalleled duration and fast recharge.

Head: Monocular, 30° inclined; 360° rotating.

Eyepiece: WF10x/18 mm, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN 4x/0.10, W.D. 15.2 mm.
- N-PLAN 10x/0.25, W.D. 5.5 mm.
- N-PLAN 40x/0.65, W.D. 0.45 mm.


All with anti-fungus treatment.

Specimen stage: Fixed stage, 130x120 mm. With sample clips.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

Condenser: N.A. 0.65, pre-centered, fixed with iris diaphragm.


**Illumination:** X-LED¹ with white 1 W LED and light intensity control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

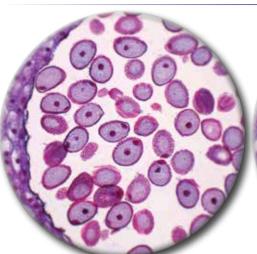

Lily Anther, Mature Pollen Grains, c.s Conventional Achromatic Objective

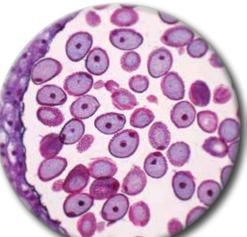

Lily Anther, Mature Pollen Grains, c.s **OPTIKA HC** Objective

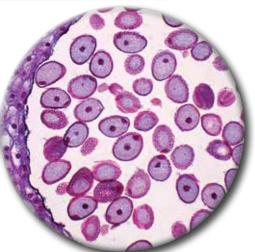

Lily Anther, Mature Pollen Grains, c.s **OPTIKA N-PLAN** Objective



N-PLAN 160 HC

**N-PLAN** 

**OPTIKA HC**: This series of objectives ensures a versatile and reasonably priced entry-level solution for brightfield and simple polarization applications. They are specifically designed to achieve optimal contrast and thus maximize yield on an instrument intended for education on F.N. 18.

**OPTIKA N-PLAN**: In addition to the advantages of the HC objectives, the total flatness of the field and an even greater contrast are achieved with the N-PLAN series.

#### $\bigcirc$

# **B-150** Series - Standard Models

## B-152 / B-153















Monocular microscope ideal for students and primary schools, with four achromatic lenses (600x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 30° inclined: 360° rotating.

Eyepiece: WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.
- Achromatic HC type 60x/0.85, W.D. 0.45 mm (only in B-153 model). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-152ALC / B-153ALC**

















Monocular microscope ideal for students and primary schools, with four achromatic lenses (600x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color 1 W *X-LED*<sup>7</sup>. Slim and easy to carry, equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED. The exclusive ALC will automatically adjust the brightness according to your preferences.

**Head:** Monocular, 30° inclined; 360° rotating (when ALC cable is unplugged).

**Eyepiece:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.
- Achromatic HC type 60x/0.85 (only in B-153ALC model) All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. With **ALC** for automatic light control. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-150** Series - Standard Models

# B-152R-PL / B-153R-PL

















Cordless monocular microscope ideal for students and primary schools, with three or four PLAN achromatic lenses (400x on B-152R-PL or 600x on B-153R-PL), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color 1 W **X-LED**<sup>1</sup>. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination. The Li-lon battery ensures unparalleled duration and fast recharge.

Head: Monocular, 30° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10, W.D. 15.2 mm.
- N-PLAN plan achromatic 10x/0.25, W.D. 5.5 mm.
- N-PLAN plan achromatic 40x/0.65, W.D. 0.45 mm.
- N-PLAN plan achromatic 60x/0.85, W.D. 0.45 mm. (only for B-153R-PL). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.


**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.




Lily Anther, Mature Pollen Grains, c.s Conventional Achromatic Objective

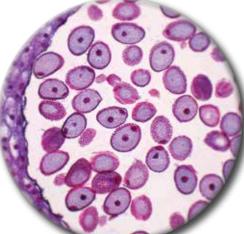

Lily Anther, Mature Pollen Grains, c.s **OPTIKA HC** Objective

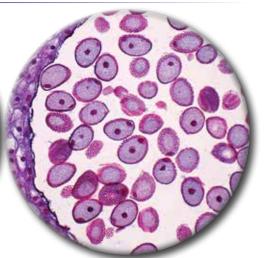

Lily Anther, Mature Pollen Grains, c.s **OPTIKA N-PLAN** Objective

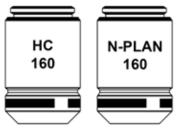

HC

**N-PLAN** 

**OPTIKA HC**: This series of objectives ensures a versatile and reasonably priced entry-level solution for brightfield and simple polarization applications. They are specifically designed to achieve optimal contrast and thus maximize yield on an instrument intended for education on F.N. 18.

**OPTIKA N-PLAN**: In addition to the advantages of the HC objectives, the total flatness of the field and an even greater contrast are achieved with the N-PLAN series.

#### $\bigcirc$

# **B-150** Series - Standard Models

## **B-155**

















**B-155ALC** 



Monocular microscope ideal for students and primary schools, with four achromatic lenses (1000x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use.

Head: Monocular, 30° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10
- Achromatic HC type 10x/0.25
- Achromatic HC type 40x/0.65
- Achromatic HC type 100x/1.25 (oil/water).

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

Monocular microscope ideal for students and primary schools, with four achromatic lenses (1000x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup>. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope with long lasting LED illumination. The ALC will automatically adjust the brightness according to your preferences.

**Head:** Monocular, 30° inclined; 360° rotating (when ALC cable is unplugged).

Eyepiece: WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10
- Achromatic HC type 10x/0.25
- Achromatic HC type 40x/0.65
- Achromatic HC type 100x/1.25 (oil/water).

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. With **ALC** for automatic light control. Multi-plug 100-240Vac/5Vdc external power supply.

## **B-155R-PL**



Cordless monocular microscope ideal for students and primary schools, with four PLAN achromatic lenses (1000x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color 1 W **X-LED**<sup>1</sup>. Slim and easy to carry, equipped with all the main controls to start learning how to use an advanced microscope with long lasting LED illumination. The Li-lon battery ensures unparalleled duration and fast recharge.

Head: Monocular, 30° inclined; 360° rotating.

**Eyepiece:** WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10, W.D. 18 mm.
- N-PLAN plan achromatic 10x/0.25, W.D. 7 mm.
- N-PLAN plan achromatic 40x/0.65, W.D. 0.53 mm.
- N-PLAN plan achromatic 100x/1.25, W.D. 0.13 mm (oil/water).


All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.


**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

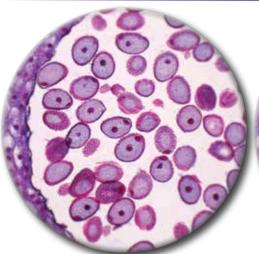

Lily Anther, Mature Pollen Grains, c.s Conventional Achromatic Objective

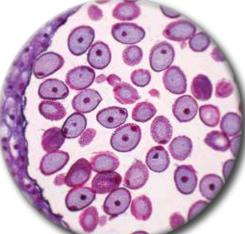

Lily Anther, Mature Pollen Grains, c.s **OPTIKA HC** Objective

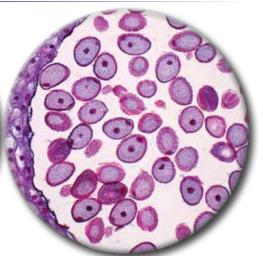

Lily Anther, Mature Pollen Grains, c.s **OPTIKA N-PLAN** Objective



**OPTIKA HC**: This series of objectives ensures a versatile and reasonably priced entry-level solution for brightfield and simple polarization applications. They are specifically designed to achieve optimal contrast and thus maximize yield on an instrument intended for education on F.N. 18.

**OPTIKA N-PLAN**: In addition to the advantages of the HC objectives, the total flatness of the field and an even greater contrast are achieved with the N-PLAN series.

## B-156 / B-157



Binocular microscope ideal for students and primary schools, with four achromatic lenses (600x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use.

Head: Binocular, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.
- Achromatic HC type 60x/0.85, W.D. 0.45 mm (only in B-157 model). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# **B-156ALC / B-157ALC**



Binocular microscope ideal for students and primary schools, with four achromatic lenses (600x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W X-LED¹ illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The exclusive **ALC** will automatically adjust the brightness according to your preferences.

Head: Binocular, 30° inclined; 360° rotating.

**Dioptric adjustment:** Left eyepiece.

Eyepieces: WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm
- Achromatic HC type 60x/0.85, W.D. 0.45 mm (only in B-157ALC model). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

Condenser: Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. With **ALC** for automatic light control. Multi-pluq 100-240Vac/5Vdc external power supply.

## B-156R-PL / B-157R-PL



Cordless binocular microscope ideal for students and primary schools, with four PLAN achromatic lenses (600x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The Li-Ion battery ensures unparalleled duration and fast recharge.

Head: Binocular, 30° inclined; 360° rotating.

**Dioptric adjustment:** Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**


- N-PLAN plan achromatic 4x/0.10, W.D. 15.2 mm.
- N-PLAN plan achromatic 10x/0.25, W.D. 5.5 mm.
- N-PLAN plan achromatic 40x/0.65, W.D. 0.45 mm.
- N-PLAN plan achromatic 60x/0.85, W.D. 0.45 mm (only in B-157R-PL model). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.


**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

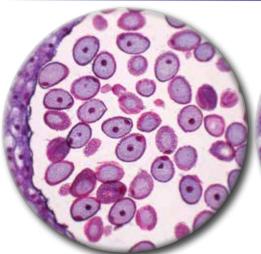


Lily Anther, Mature Pollen Grains, c.s Conventional Achromatic Objective

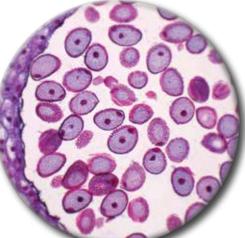

Lily Anther, Mature Pollen Grains, c.s **OPTIKA HC** Objective

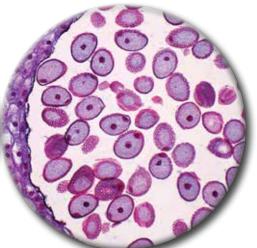

Lily Anther, Mature Pollen Grains, c.s **OPTIKA N-PLAN** Objective

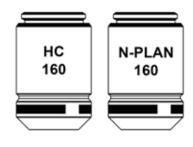

HC

**N-PLAN** 

**OPTIKA HC**: This series of objectives ensures a versatile and reasonably priced entry-level solution for brightfield and simple polarization applications. They are specifically designed to achieve optimal contrast and thus maximize yield on an instrument intended for education on F.N. 18.

**OPTIKA N-PLAN**: In addition to the advantages of the HC objectives, the total flatness of the field and an even greater contrast are achieved with the N-PLAN series.

## 1)

## **B-150** Series - Standard Models

## **B-159**



Binocular microscope ideal for students and primary schools, with four achromatic lenses (1000x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W *X-LED*<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use.

Head: Binocular, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.
- Achromatic HC type 100x/1.25, W.D. 0.13 mm (oil/water).

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

## **B-159ALC**



Binocular microscope ideal for students and primary schools, with four achromatic lenses (1000x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W X-LED¹ illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The exclusive **ALC** will automatically adjust the brightness according to your preferences.

**Head:** Binocular, 30° inclined; 360° rotating (when ALC cable is unplugged).

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- Achromatic HC type 4x/0.10, W.D. 18 mm.
- Achromatic HC type 10x/0.25, W.D. 7 mm.
- Achromatic HC type 40x/0.65, W.D. 0.53 mm.
- Achromatic HC type 100x/1.25, W.D. 0.13 mm (oil/water).

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

Condenser: Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. With **ALC** for automatic light control. Multi-plug 100-240Vac/5Vdc external power supply.

## **B-159R-PL**



Cordless binocular microscope ideal for students and primary schools, with four PLAN achromatic lenses (1000x), FN 18 high eyepointeyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W *X-LED*<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The Li-lon battery ensures unparalleled duration and fast recharge.

Head: Binocular, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw. **Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10, W.D. 15.2 mm.
- N-PLAN plan achromatic 10x/0.25, W.D. 5.5 mm.
- N-PLAN plan achromatic 40x/0.65, W.D. 0.45 mm.
- N-PLAN plan achromatic 100x/1.25, W.D. 0.13 mm (oil/water).


All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED¹ with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

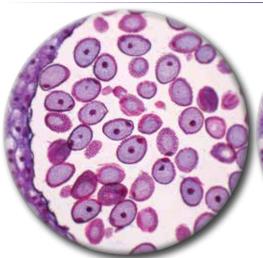


Lily Anther, Mature Pollen Grains, c.s Conventional Achromatic Objective

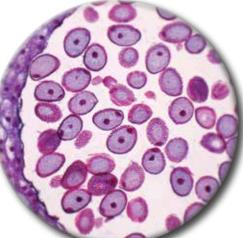

Lily Anther, Mature Pollen Grains, c.s **OPTIKA HC** Objective



Lily Anther, Mature Pollen Grains, c.s **OPTIKA N-PLAN** Objective

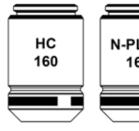

N-PLAN 160 HC

**N-PLAN** 

**OPTIKA HC**: This series of objectives ensures a versatile and reasonably priced entry-level solution for brightfield and simple polarization applications. They are specifically designed to achieve optimal contrast and thus maximize yield on an instrument intended for education on F.N. 18.

**OPTIKA N-PLAN**: In addition to the advantages of the HC objectives, the total flatness of the field and an even greater contrast are achieved with the N-PLAN series.

## (1)

# **B-150** Series - Polarizing Models

## B-150P-MRPL



















## B-150P-BRPL



Cordless monocular polarized light microscope ideal for students and primary schools, with three PLAN achromatic lenses (400x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, rotating stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination. Rotating swing-out polarizer and sliding-out fixed analyzer included. The Li-Ion battery ensures unparalleled duration and fast recharge.

**Head:** Monocular, 30° inclined; 360° rotating.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10
- N-PLAN plan achromatic 10x/0.25
- N-PLAN plan achromatic 40x/0.65.

All with anti-fungus treatment.

**Specimen stage:** Rotatable round stage, 120 mm diameter, with sample clips.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** N.A. 1.25, pre-centered, fixed, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

**Polarizing filters:** Rotating polarizer (swing-out) and fixed analyzer (sliding-out).

Cordless binocular polarized light microscope ideal for students and primary schools, with three PLAN achromatic lenses (400x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, rotating stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. Rotating swing-out polarizer and sliding-out fixed analyzer included. The exclusive Li-lon battery ensures unparalleled duration and fast recharge

Head: Binocular, 30° inclined; 360° rotating.

**Dioptric adjustment:** Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10
- N-PLAN plan achromatic 10x/0.25
- N-PLAN plan achromatic 40x/0.65. All with anti-fungus treatment.

**Specimen stage:** Rotatable round stage, 120 mm diameter, with sample clips.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

Condenser: N.A. 1.25, pre-centered, fixed, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

**Polarizing filters:** Rotating polarizer (swing-out) and fixed analyzer (sliding-out).

# **B-150** Series - Digital Models

## B-150D-MRPL



Cordless digital monocular microscope ideal for students and primary schools, with three PLAN achromatic lenses (400x), FN 18 high eyepoint eyepiece, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W **X-LED**<sup>1</sup> illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The exclusive Li-lon battery ensures unparalleled duration and fast recharge

**Head:** Monocular with integrated 1.3 MP camera, 30° inclined; 360° rotating.

Eyepieces: WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10
- N-PLAN plan achromatic 10x/0.25
- N-PLAN plan achromatic 40x/0.65


**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.

## B-150D-BRPL

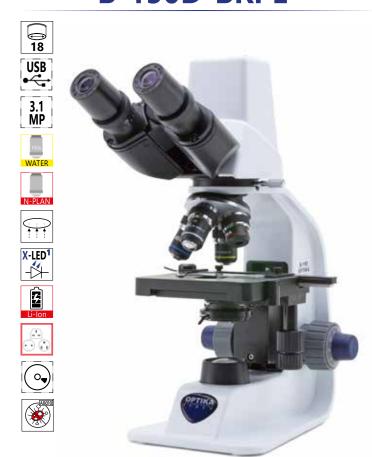

Cordless digital binocular microscope ideal for students and primary schools, with four PLAN achromatic lenses (1000x), FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 1 W X-LED¹ illumination. Slim and easy to carry, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use. The exclusive Li-lon battery ensures unparalleled duration and fast recharge

**Head:** Binocular with integrated 3.2 MP camera, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece.

#### **Objectives:**

- N-PLAN plan achromatic 4x/0.10
- N-PLAN plan achromatic 10x/0.25
- N-PLAN plan achromatic 40x/0.65
- N-PLAN plan achromatic 100x/1.25 (oil/water). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x116 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>1</sup> with white 1 W LED and brightness control. Color temperature: 6,300 K. Li-lon battery for long lasting operation. Multi-plug 100-240Vac/5Vdc external power supply.


# B-150 Series - B-150D Camera specifications

	B-150D-MRPL / B-150D-BRPL		
Camera resolution (n° of pixels: W x H)	2048x1536		
Color / Monochrome	Color		
Sensor Size	1/2.5"		
Sensor Technology	CMOS		
Sensor Type	APTINA CMOS		
Image Format	4/3		
Pixel Size (mm)	2.2x2.2		
Acquisition technology	Rolling shutter		
Sensitivity	0.53 V/lux-second 40.5		
Signal/Noise Ratio (DB)			
Dynamic range (DB)	66.5		
ADC conversion	12 Bit		
Color Depth	8 Bit		
Exposure time	0.1 - 1.5 sec		
Binning	1x1		
IR filter	Yes		
IR filter range (nm)	380-650 (IR CUT)		
Camera resolution (MP)	3.1		
USB standard	USB2.0		
Frame rate full resolution (fps)	4@2048x1536		
Frame rate other resolution (fps)	8@1280x1024; 30@640x480		

# **B-150** Series - Optical performance


Eyepiece			10x (M-002.1) 18		16x	(M-003)
Field number (mm)						12
Objective	N.A.	W.D. (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
4x	0.1	18	40x	4.5	64x	3
10x	0.25	7	100x	1.8	160x	1.2
20x	0.4	2	200x	0.9	320x	0.6
40x	0.65	0.53	400x	0.45	640x	0.3
60x	0.8	0.45	600x	0.3	960x	0.2
100x	1.25 (oil/water)	0.13	1000x	0.18	1600x	0.12

# **B-150** Series - Zoom comparison

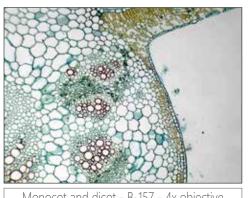

Monocot and dicot - B-157 - 4x objective



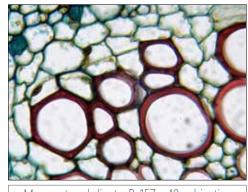

# **B-150** Series - Digital Video Bundles



Five models of the B-150 series can be equipped with a camera and a 7" LCD screen, both in high definition. In these bundles, the normal head (supplied with the instrument) can be replaced with the digital system in a few minutes. This solution provides a system suitable for viewing specimens by several students at the same time, without removing the possibility of using the microscope in the classical way through the eyepiece.

**LCD screen:** High definition 7" LCD.

Camera: 1920x1080 pixels, 30fps (video). Up to 1844x1080 pixels (photo).


**Storing capacity:** On Micro Sd card

Video recording: Yes


Measuring funcion: Yes, simple line measurement

#### Models available:

- **-B-151V:** Same fatures than standard B-151, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.**
- **-B-153V:** Same fatures than standard B-153, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.**
- **-B-151R-PLV:** Same fatures than standard B-151R-PL, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.**
- **-B-152R-PLV:** Same fatures than standard B-152R-PL, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.**
- **-B-159R-PLV:** Same fatures than standard B-159R-PL, but delivered as bundle together with the 7" screen with built in camera. **Optical head with eyepiece included.**

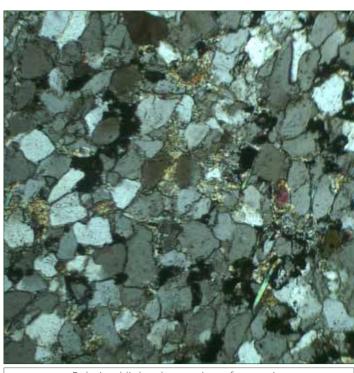

Polarized light observation of quartzite with B-150P-MRPL and 10x objective.

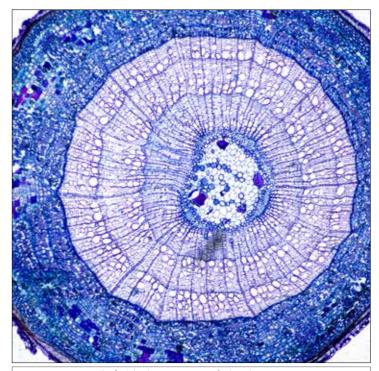

Brightfield observation of tilia three-year stem with B-159 and 20x objective.

# B-150 Series - Comparison charts

-150 - Star	ndard Mode	els, with H	IC Objectiv	res				
odel	Head	Eyepiece(s)	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
151	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x	Fixed, 130x120 mm, with sample clips	Coaxial coarse and fine, limit stop	N.A. 0.65, iris diaphragm, fixed	1 W X-LED <sup>1</sup> , manual brightness contro
·152 ·153*	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 60x*	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro
155	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro
156 157*	Binocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 60x*	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro
159	Binocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro
150 - ALC	Models, w	ith Autom	atic Light	Control and	HC Objectives			]
odel	Head	Eyepiece(s)	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
151ALC	Monocular, 30° inclined	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x	Fixed, 130x120 mm, with sample clips	Coaxial coarse and fine, limit stop	N.A. 0.65 fixed, with diaphragm	1 W X-LED <sup>1</sup> , manual and automatic brightness control
-152ALC -153ALC*	Monocular, 30° inclined	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 60x*	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual and automatic brightness control
-155ALC	Monocular, 30° inclined	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual and automatic brightness control
-156ALC -157ALC*	Binocular, 30° inclined	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 60x*	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual and automatic brightness control
-159ALC	Binocular, 30° inclined	WF 10x/18	Quadruple	HC (high contrast) 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual and automatic brightness control
-150 - Core	dless Mode	ls, with N	-PLAN Obj	ectives and L	i-Ion Rechargeab	le Batteries		
odel	Head	Eyepiece(s)	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
·151R-PL	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x	Fixed, 130x120 mm, with sample clips	Coaxial coarse and fine, limit stop	N.A. 0.65 fixed, with diaphragm	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-152R-PL	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-153R-PL	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x, 60x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-155R-PL	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-156R-PL -157R-PL*	Binocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x, 60x*	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-159R-PL	Binocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x, 100x	Double layer, 125x116 mm, moving range 70x30 mm	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-150 - Pola	rized Liaht	Cordless N	Models. wit	th N-PLAN O	ojectives and Li-lo	n Rechargea	ble Batteries	]
odel	Head	Eyepiece(s)		Objectives	Stage	Focusing	Condenser	Illumination
150P-MRPL	Monocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x	Round, 360° rotating, 120 mm diameter, with sample clips	Coaxial coarse and fine, limit stop	N.A. 1.25, iris diaphragm, fixed	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte
-150P-BRPL	Binocular, 30° inclined, 360° rotating	WF 10x/18	Quadruple	N-PLAN 4x, 10x, 40x	Round, 360° rotating, 120 mm diameter, with sample clips	Coaxial coarse and fine, limit stop	N.A. 1.25, iris diaphragm, fixed	1 W X-LED <sup>1</sup> , manual brightness contro Li-lon rechargeable batte

**Focusing** 

Coaxial coarse

and fine, limit

Coaxial coarse

and fine, limit

stop

Condenser

Abbe N.A. 1.25, iris diaphragm,

Abbe N.A. 1.25,

iris diaphragm, focusable

Illumination

1 W X-LED<sup>1</sup>,

manual brightness control,

Li-lon rechargeable battery

manual brightness control, Li-lon rechargeable battery

Model

B-150D-MRPL

B-150D-BRPL

Head

Monocular,

30° inclined,

360° rotating

Binocular, 30° inclined, 360° rotating

Nosepiece

Quadruple

Quadruple

Eyepiece(s)

WF 10x/18

WF 10x/18

**Objectives** 

4x, 10x, 40x

4x, 10x, 40x, 100x

N-PLAN

Stage

Double layer, 125x116

Double layer, 125x116

mm, moving range 70x30 mm

mm, moving range

# **B-150** Series - Accessories

#### **Eyecups & Eyepieces** M-001 Huygens 5x eyepiece

WF10x/18 eyepiece, high eyepoint M-002.1

WF10x/18 micrometric eyepiece, high eyepoint M-004

M-008 WF10x/18 eyepiece, high eyepoint, with pointer

M-003 WF16x/12 eyepiece M-162 WF20x/10 eyepiece

#### **Objectives**

#### HC

M-137	HC (high contrast) objective 4x/0.10
M-138	HC (high contrast) objective 10x/0.25
M-139	HC (high contrast) objective 20x/0.40
M-141	HC (high contrast) objective 40x/0.65
M-142	HC (high contrast) objective 60x/0.85
M-143	HC (high contrast) objective 100x/1.25 (oil)

#### N-PLAN

M-164	N-PLAN objective 4x/0.10
M-165	N-PLAN objective 10x/0.25
M-166	N-PLAN objective 20x/0.40
NA 167	N DI ANI I' 1' 40 10 CE

M-167 N-PLAN objective 40x/0.65 M-168 N-PLAN objective 60x/0.85

M-169 N-PLAN objective 100x/1.25 (oil)

#### Stages

M-040 Attachable mechanical stage (only for B-151, B-151ALC and B-151R-PL)

#### **Condensers & Filters**

M-974	Blue filter, 32mm diameter
M-976	Green filter, 32mm diameter
M-978	Yellow filter, 32mm diameter
M-988	Frosted glass filter, 32mm diameter
M-155	Polarising set (filters only)
_	J , , , , , , , , , , , , , , , , , , ,

#### **Camera Adapters**

M-115	0.35x C-Mount projection lens
M-114	0.5x C-Mount projection lens
M-118	0.75x C-Mount projection lens

#### Miscellaneous

IVIISCEIIAIIC	eous e
<u>15104</u>	Cleaning kit
15008	Immersion oil, 10ml
15009	Immersion oil, 100ml
DC-002	Plastic dust cover, medium, 490(l)x490(h) mm
M-005	Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)
M-069	Solar charger
M-972	Plane-concave mirror, with base
AB-010	Antibacterial surface treatment, only for newly purchased microscope

#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. Output voltage: 5 Vdc. -

Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>).

Charging models: with solar panel (12h), with external USB power supply (2.5h)



#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.





#### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com




# B-190 Series




**Advanced Biological Microscopes For Students** 

# Fulfill the Next Generation Learning Challenges

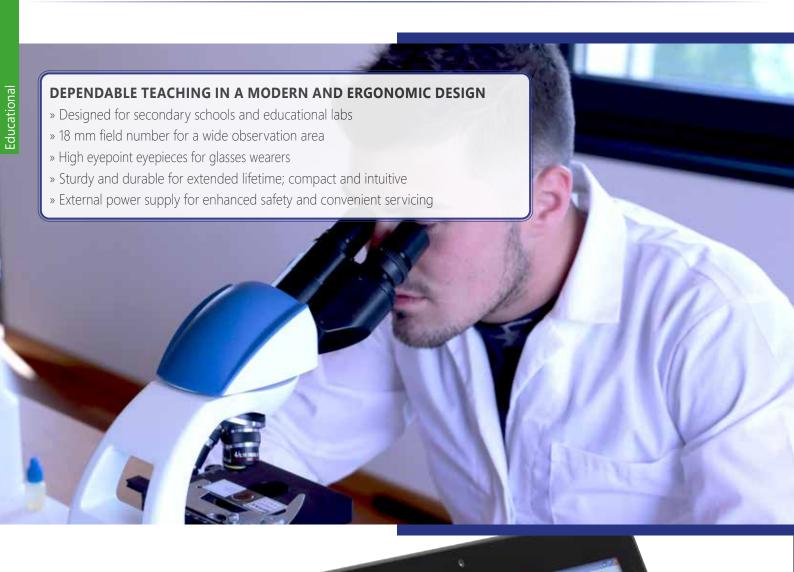


#### **DIGITAL IMAGING UNIT - ONLY AVAILABLE AT OPTIKA**

- » A completely new, revolutionary experience for unparalleled comfort
- » Responsive and smooth control for accurate results in few clicks
- » Large touch screen with 360° rotating and tilting holding solution
- » Simultaneous camera and power connection for long-term operation
- » Easily detachable to be used as a laptop (keyboard sold separately)

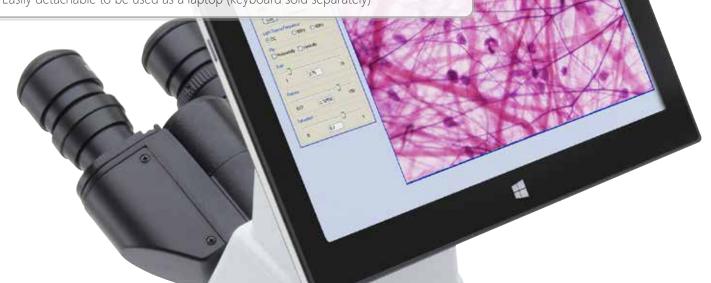

# Inventors Of A New Way To Teach Microscopy

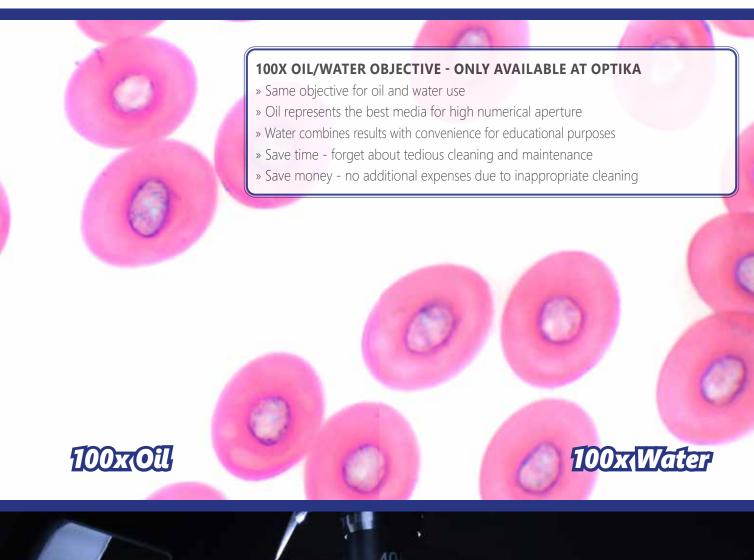



# **B-190** Series

Valuable configurations of modern microscopes ideal for teachers and secondary schools, with four achromatic lenses, FN 18 high eyepoint eyepieces, finite optical system, coaxial focusing, mechanical stage and powerful, uniform, white color temperature 3 W X-LED<sup>2</sup> illumination. Slim and easy to carry, all the models are equipped with long lasting LED illumination to provide over 20 years of use

#### High eyepoint eyepieces for glasses wearers

These eyepieces are designed in such a way that the exit pupil is further away from the eye lens than standard eyepieces, being are well suited for eyeglasses wearers

#### X-LED<sup>2</sup> State-of-the-art illumination system for incomparable light intensity

Provided with an exclusive lens & collector design, OPTIKA X-LED technology ensures unmatched uniformity & brightness (more than a 30 W halogen lamp) for excellent color fidelity with constant pure-white color temperature





#### 100x oil/water objective: same objective for dual use

This new, revolutionary objective is something you've never seen before! Oil ensures the best performance achievable; water represents the most convenient solution as eliminates tedious cleaning

#### A completely new, revolutionary experience for unparalleled comfort

B-190TBPL includes built-in camera and Digital Imaging Unit with large touch screen for a smooth and responsive control, with dependable results in few clicks, and providing an extremely comfortable solution for open discussions

#### External power supply for enhanced safety and convenient servicing


OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit

Advanced Biological Microscopes For Students And Teachers

#### **Optimum And Unparalleled Comfort In Use**

The B-190TBPL offers you a unique, incomparable solution. It includes a built-in camera and Digital Imaging Unit with large touch screen, for a responsive and smooth control. Simultaneous camera and power connection ensure long-term operation, with dependable results in one click. It provides a reliable and comfortable solution for open discussion: 360° rotating and tilting screen, easily detachable, that can be used as a laptop.



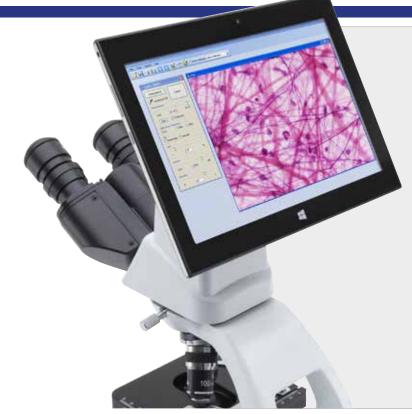

## Get the most out of our accessories



# M-174 - Polarizing set Set for simple polarization analysis. Upgrade your B-190 to a polarizing microscope and look at birefringent samples.

**ANALYZER** 



## **B-191sPL**















Monocular microscope ideal for teachers and secondary schools, with four N-PLAN objectives (600x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W **X-LED**<sup>2</sup> illumination. Slim and easy to carry, yet sturdy and resistant, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

Head: Monocular, 30° inclined; 360° rotating.

Eyepiece: WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 60x/0.85
- All with anti-fungus treatment.


**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-191PL**

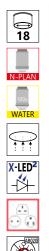



Monocular microscope ideal for teachers and secondary schools, with four N-PLAN objectives (1000x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W **X-LED**<sup>2</sup> illumination. Slim and easy to carry, yet sturdy and resistant, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

**Head:** Monocular, 30° inclined; 360° rotating.

Eyepiece: WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

Condenser: Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-192sPL**



Binocular microscope ideal for teachers and secondary schools, with four N-PLAN objectives (600x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W *X-LED*<sup>2</sup> illumination. Slim and easy to carry, yet sturdy and resistant, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

Head: Binocular, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 60x/0.85

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-192PL**



Binocular microscope ideal for teachers and secondary schools, with four N-PLAN objectives (1000x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W **X-LED**<sup>2</sup> illumination. Slim and easy to carry, yet sturdy and resistant, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

**Head:** Binocular, 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-193PL**

















**B-190TBPL** 



Trinocular microscope for camera connection ideal for teachers and secondary schools, with four N-PLAN objectives (1000x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W **X-LED**<sup>2</sup> illumination. Slim and easy to carry, yet sturdy and resistant, it is equipped with all the main controls to start learning how to use an advanced microscope and with long lasting LED illumination to provide over 20 years of use

Head: Trinocular (split ratio: 50/50), 30° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/18 mm, high eyepoint, secured by screw.

Nosepiece: Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

Digital binocular microscope ideal for teachers and secondary schools, with four N-PLAN objectives (1000x), FN 18 high eyepoint, finite optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 3 W **X-LED**<sup>2</sup> illumination. The 3.1 MP CMOS camera ensures excellent colour reproduction and is connected to Digital Imaging Unit with vivid color graphic display for unparalleled comfort and performance. The Digital Imaging Unit represents a 2-in-1 solution as it can be disconnected and used as a real PC, being Windows-based, with powerful Intel processor and large touch screen of 10.1" for fast, responsive and smooth control

**Head:** Binocular with integrated 3.1 MP camera, 30° inclined; 360° rotating. Digital Imaging Unit detachable included, rotating and tilting.

Dioptric adjustment: Left eyepiece.

Eyepieces: WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 100x/1.25 (Oil/Water). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

**Illumination:** X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply

10.1″

MP

## **B-190TBPL+**



#### **DIGITAL FOCUSING CONTROL**

B-190TBPL+ features a revolutionary focusing system based on unique digital control unit. There is a single adjustment knob located on the left side of the instrument.

The knob controls the focus, which has three different resolutions:

- Macro, 250μm per step (7500μm each knob revolution)
- Semi-Micro, 8μm per step (240μm each knob revolution)
- Micro, 1μm per step (30μm each knob revolution)

The selection of the three resolutions is done simply by pressing the knob itself. Repeated operations cyclically switch the type of focus in following order: Macro → Semi-Micro → Micro, then back to Macro and so on.

On the basis of chosen resolution, three specific LEDs located on the instrument's front side illumininate:





m - Micro





sm - Semi Micro

M - Macro



B-190TBPL+ offers the same features than B-190TBPL but with an upgraded 5Mp camera and a fully electronically controlled triple motorized focus adjustment: Macro, Semi-Micro and Micro. The new **X-OEF** system (Optika EasyFocus system) replaces the classic mechanical focusing system with an innovative solution that allows more functions.

EasyFocus allows the selection of three different resolutions (shift for each step): Macro (250 μm); Semi-Micro (8 μm); Micro (1 μm).

The new X-OEF system includes a stage retraction function that facilitates and speeds up immersion oil laying and slide changing operations. Each focus resolution (Macro, Semi-Micro and Micro) is clearly displayed to the user by three dedicated blue LEDs.

**Head:** Binocular with integrated 5 MP camera, 30° inclined; 360° rotating. Digital Imaging Unit detachable included, rotating and tilting.

Dioptric adjustment: Left eyepiece.

Eyepieces: WF10x/18 mm, high eyepoint, secured by screw.

**Nosepiece:** Quadruple ball bearings revolving nosepiece, reversed.

#### **Objectives:**

- N-PLAN objective 4x/0.10
- N-PLAN objective 10x/0.25
- N-PLAN objective 40x/0.65
- N-PLAN objective 100x/1.25 (Oil/Water). All with anti-fungus treatment.

**Specimen stage:** Mechanical stage, 125x115 mm, 70x30 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

Focusing: X-OEF system with three different resolutions (Macro; Semi-Micro; Micro) and stage retraction function.

**Condenser:** Abbe N.A. 1.25, pre-centered, focusable, with iris diaphragm.

Illumination: X-LED<sup>2</sup> with white 3 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

#### STAGE RETRACTION SYSTEM - ONE TOUCH MEMORY FUNCTIONS

This revolutionary system stores the last focus point to facilitate and speed up routine operations

normally cost the operator time and effort:

#### 1. Easy immersion oil laying operation

The last focusing point (height of stage) is automatically stored by the digital processor unit. This feature is really useful in case of oil immersion objectives use


The procedure is very simple and it permits to switch from 40x to 100 objective with an ease never achieved before:

- focus the 40x objective. Its focusing point is automatically stored by the system;
- push the focusing knob for 2 seconds; the stage reaches its lower position, allowing you to easily put a drop of oil on the sample,
- put a drop of oil on the sample and switch from 40x to 100x objective;
- push once again the dedicated button; the stage rises reaching the stored fucus position;
- adjust the micrometric focus.

#### 2. Ergonomic slide changing operation

By following above mentioned procedure, it is possible to significantly save time when you change

- focus the objective. Its focusing point is automatically stored by the system;
- push the dedicated botton located on the microscope body; the stage reaches its lower position, allowing you to easily change the sample;
- change slide;
- push once again the dedicated button; the stage rises reaching the stored fucus position;
- adjust the micrometric focus

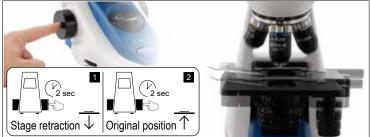



### 1

# **B-190TBPL & B-190TBPL+**

#### **Unique Features**

- > Simultaneous camera & power connection
- > Equipped with the latest **Windows OS** & **Intel** processor
- > Easily detachable, can be used as a laptop (keyboard sold separately only for B-190TBPL)









Simple and user-friendly, ideal for students and experienced users.

TB-KBD2 - Accessory keyboard


# Digital microscope with camera & Digital Imaging Unit

	B-190TBPL	B-190TBPL+
Camera resolution (n° of pixels: W x H)	2048x1536	2592x1944
Color / Monochrome	Color	Color
Sensor Size	1/2.5"	1/2.5"
Sensor Technology	CMOS	CMOS
Sensor Type	APTINA CMOS	APTINA CMOS
Image Format	4/3	4/3
Pixel Size (mm)	2.8x2.8	2.2x2.2
Acquisition technology	Rolling shutter	Rolling shutter
Sensitivity	0.53 V/lux-second	18.8ke/lus
Signal/Noise Ratio (DB)	40.5	40
Dynamic range (DB)	66.5	73
ADC conversion	12 Bit	12 Bit
Color Depth	8 Bit	8 Bit
Exposure time	0.1 - 1.5 sec	0.1 msec - 1 sec
Binning	1x1	1x1
IR filter	Yes	Yes
IR filter range (nm)	380-650 (IR CUT)	380-650 (IR CUT)
Camera resolution (MP)	3.1	5.1
USB standard	USB2.0	USB2.0
Frame rate full resolution (fps)	4@2048x1536	30@2592x1944
Frame rate other resolution (fps)	8@1280x1024; 30@640x480	30@1280x960; 30@640x480




**OPTIKA N-PLAN** objectives ensure bright, clear images with excellent flatness and compensation for chromatic aberration.


# **B-190** Series - **Zoom comparison**

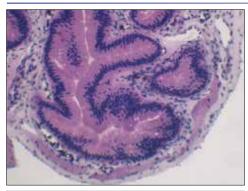

Frog small intestine - B-193PL - 10x objective

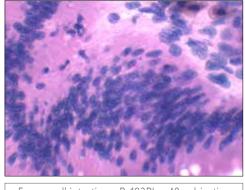

Frog small intestine - B-193PL - 40x objective

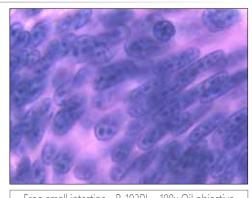

Frog small intestine - B-193PL - 100x Oil objective

## **B-190** Series - Accessories

#### **Eyecups & Eyepieces**

M-001 Huygens 5x eyepiece

M-002.1 WF10x/18 eyepiece, high eyepoint

M-004 WF10x/18 micrometric eyepiece, high eyepoint

M-008 WF10x/18 eyepiece, high eyepoint, with pointer

M-003 WF16x/12 eyepiece

M-162 WF20x/10 eyepiece

#### **Objectives**

M-164 N-PLAN objective 4x/0.10

M-165 N-PLAN objective 10x/0.25 M-166 N-PLAN objective 20x/0.40

M-167 N-PLAN objective 40x/0.65

M-168 N-PLAN objective 60x/0.85 M-169 N-PLAN objective 100x/1.25 (oil)

#### **Condensers & Filters**

M-174 Polarising set (filters only)

M-974 Blue filter, 32 mm diameter

M-976 Green filter, 32 mm diameter

M-978 Yellow filter, 32 mm diameter

M-988 Frosted glass filter, 32 mm diameter

#### **Camera Adapters**

M-115 0.35x C-Mount projection lens

M-114 0.5x C-Mount projection lens

M-118 0.75x C-Mount projection lens

M-173 C-Mount projection lens for APS-C/full frame reflex cameras (trino)

#### Miscellaneous

15104 Cleaning kit

15008 Immersion oil, 10ml

15009 Immersion oil, 100ml

DC-002 Plastic dust cover, medium, 490(l)x490(h) mm (except for B-190TB)

DC-003 TNT dust cover, medium, 600(l)x550(h) mm (only for B-190TB)

TB-KBD2 Keyboard for Digital Imaging Unit (only for B-190TB)

M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)

M-069 Solar charger

M-971 Plane-concave mirror, with base

VP-190 IQ/OQ/PQ manual for B-190 series

VP-TB IQ/OQ/PQ manual for TB series

AB-020 Antibacterial surface treatment, only for newly purchased microscope

#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. Output voltage: 5 Vdc. -

Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>).

Charging models: with solar panel (12h), with external USB power supply (2.5h)



#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.





#### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA Spain
OPTIKA China
OPTIKA India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com OPTIKA<sup>®</sup> North America OPTIKA<sup>®</sup> Central America OPTIKA<sup>®</sup> Africa namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# MS/SFX Series



**Entry-Level Monoscopes & Stereomicroscopes For Students** 

# Addressed For Simplicity & The Youngest Users



## TOUCH CONTROL, A NEW FRONTIER IN ILLUMINATION ADJUSTMENT » Light intensity settable via a simple click

- » Cordless use, totally independent from mains/batteries connection
- » Battery status indicator informs when the battery is ready to be used
- » Freely settable illumination for incident and transmitted light
- » External power supply for enhanced safety and convenient servicing



(1)



# **MS/SFX** Series

A great variety of mainly cordless binocular stereomicroscopes with turnable objectives, FN 20 eyepieces, 1 W LED transmitted/incident illumination and different stands to start exploring and discovering sciences and materials, including biology, entomology, rocks, plants, and many more specimens. Most of the models are equipped with premium features, such as the exclusive, comfortable touch control and the rechargeable batteries. Slim and easy to carry, all the models are equipped with long lasting LED illumination to provide over 20 years of use



#### Turnable objectives, up to 3 magnifications on 20 mm field number

SFX series is designed for basic operations, being provided with widefield eyepieces and objective turret magnification changer (up to 3 magnifications) for step magnification

#### Touch control - light intensity settable via a simple click

SFX Series represents a new frontier in illumination adjustment with an incredibly comfortable solution: the touch control, consisting in 10 pre-set intensity levels (5 for transmitted, 5 for incident)

#### Cordless use, totally independent from the mains/batteries connection

SFX Series works with or without the batteries in place and are provided with three NiMH rechargeable batteries for the longest autonomy in outdoor use (4-hour autonomy, at medium intensity)

1

# **MS/SFX** Series



#### Battery status indicator informs when batteries are ready to be used

SFX Series has a smart charging indicator which indicates current charging status at all times - even when not in charge or during storage: if it is on, it means it is immediately ready to work

#### Longlife LED illumination (providing over 20 years of use)

Money & energy saving thanks to LED long lifetime (65.000 hours, 22 years in case of 8 hours/day) which is more than 20 times compared to a standard halogen bulb

#### External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit







#### ST-50Led - When Long Working Distance Is Required

ST-50Led has a special objective with long working distance that allows you to inspect bulky samples, thanks also to its overhanging arm stand and LED flexible incident light

#### **MS-01 - Multifunctional Testing Equipment**

MS-01 is a portable monoscope ideal as versatile testing equipment especially for surface analysis and measurements to be used directly on the specimen with 10x objective and penlight for incident illumination

#### **GREENOUGH**









Portable monoscope ideal as multifunctional testing equipment especially for surface analysis and measurements to be used directly on the specimen. Equipped with fixed objective (10x), FN 15.5 eyepiece with crosshair and 0.3 W LED penlight for incident illumination powered by rechargeable batteries (not provided)

Eyepiece: WF10x/15.5 mm, micrometric with crosshair.

**Objective:** Achromatic 10x with anti-fungus treatment.

Working distance: 6 mm.

**Focusing:** Rack and pinion focusing mechanism.

Illumination: 0.3 W LED incident penlight with batteries (not provided).

## **SFX-31**



Cordless binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (2x-4x), FN 20 eyepieces, pillar stand and 1 W LED transmitted / incident illumination with rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

**Head:** Binocular, 45° inclined

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/20 mm, secured by screw and with rubber cups.

**Objective:** Achromatic 2x-4x with anti-fungus treatment.

Working distance: 57 mm.

Stand: Pillar with focus.

**Focusing:** Rack and pinion focusing mechanism.

**Illumination:** 1 W LED incident and transmitted, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K.

Multi-plug 100-240Vac/5Vdc external power supply.

## **SFX-33**



Cordless binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (2x-4x), FN 20 eyepieces, fixed arm with handle and 1 W LED transmitted / incident illumination with comfortable touch control and rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use

**Head:** Binocular, 45° inclined.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/20 mm, secured by screw and with rubber cups.

**Objective:** Achromatic 2x-4x with anti-fungus treatment.

Working distance: 57 mm

Stand: Fixed with handle and focus.

**Focusing:** Rack and pinion focusing mechanism.

**Illumination:** 1 W LED incident and transmitted, with touch brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

## **SFX-51**

**20** 

1W LED

1W LED

2x-4x

7



Cordless binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (2x-4x), FN 20 eyepieces, 360° rotating head, fixed arm with handle and 1 W LED transmitted / incident illumination with comfortable touch control and rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use

Head: Binocular, 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

Eyepieces: WF10x/20 mm, secured by screw and with rubber cups.

**Objective:** Achromatic 2x-4x with anti-fungus treatment.

Working distance: 76 mm

Stand: Fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** 1 W LED incident and transmitted, with touch brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply

## **SFX-91**



Cordless binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (1x-2x-4x), FN 20 eyepieces, precision fixed arm with handle and 1 W LED transmitted / incident illumination with comfortable touch control and rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Binocular, 45° inclined.

**Dioptric adjustment:** Both eyepieces.

Eyepieces: WF10x/20 mm, secured by screw and with rubber cups.

**Objective:** Achromatic 1x-2x-4x with anti-fungus treatment.

Working distance: 60 mm.

**Stand:** High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** 1 W LED swiveling incident and transmitted, with touch brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

# SFX-91D/SFX-91DW



Model SFX-91DW Connectivity Option: Wi-Fi 5 (AP Mode): Direct connection Wi-Fi 5 (STA Mode): LAN Wi-fi LAN (Ethernet)

Cordless digital binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (1x-2x-4x), FN 20 eyepieces, precision fixed arm with handle and 1 W LED transmitted / incident illumination with comfortable touch control and rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

**Head:** Binocular with integrated 5 MP camera, 45° inclined. **Model SFX-91DW:** Binocular head, 45° inclined, with built-in 4Mp USB/LAN/WiFi digital camera.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/20 mm, secured by screw and with rubber cups.

**Objective:** Achromatic 1x-2x-4x with anti-fungus treatment.

Working distance: 60 mm.

**Stand:** High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** 1 W LED swiveling incident and transmitted, with touch brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.



## ST-50Led













Binocular stereomicroscope ideal for large specimens, with long working distance fixed objective (2x), FN 20 eyepieces, overhanging stand and 1 W LED swiveling incident illumination. Additional objectives available for different magnifications

**Head:** Binocular, 45° inclined.

Dioptric adjustment: Left eyepiece.

Eyepieces: WF10x/20 mm, secured by screw.

**Objective:** Achromatic 2x with anti-fungus treatment.

Working distance: 110 mm.

Stand: Overhanging with focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** 1 W LED swiveling incident on flexible arm.

Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

## SFX-85V



Digital Zoom microscope with tilting/folding 9" monitor with camera and software integrated, LED Incident & transmitted light, precision fixed arm with handle and incredibly bright and uniform led illumination (transmitted and incident), settable through the exclusive touch control. Powered by external power supply.

Head: Monitor 9" screen tilting and folding.

**Zoom Optics:** 0.7X to 4.5X

**Total Magnification on Monitor:** 8x to 54x

Camera: 5MP, CMOS.

Working distance: 100 mm.

**Stand:** High-grade fixed stand with focus and handle.

Focusing: Rack and pinion focusing mechanism

**Illumination:** Incident: 21 LEDs, 1W in total. Transmitted: 1W LED. Touch brightness control. Multi-plug 100-240Vac, 50/60 Hz external

power supply.





## **MS/SFX** Series - Accessories

#### **ACCESSORIES FOR MS-1**

	•			
١n	isce	Пои	200	

DC-001 Plastic dust cover, small, 340(l)x400(h) mm

M-899 Pen illuminator <u>15104</u> Cleaning kit

#### **ACCESSORIES FOR SFX SERIES & ST-50Led**

**Eyecups & Eyepieces** 

<u>ST-001</u> <u>WF5x/22 eyepieces (pair), 30.5mm diameter (except for ST-50Led)</u>

ST-002 ST-003 WF10x/20 eyepieces (pair) WF15x/15 eyepieces (pair) WF20x/13 eyepieces (pair)

ST-004 ST-005 WF10x/20 micrometric eyepiece

ST-001.1 WF5x/22 eyepieces (pair), 30mm diameter (only for ST-50Led)

**Objectives** 

ST-025 1x objective (only for ST-50Led) ST-026 3.5x objective (only for ST-50Led)

Stage

Glass object-plate, 95mm diameter (only for ST-30FX)

ST-014 ST-015 Glass object-plate, 60mm diameter (except for ST-30FX & ST-50Led)

ST-011 White/black object-plate, 60mm diameter (except for ST-30FX & ST-50Led)

ST-012 White/black object-plate, 95mm diameter (only for ST-30FX)

Camera Adapters

0.5x C-Mount projection lens 0.35x C-Mount projection lens M-114 M-115 M-118 0.75x C-Mount projection lens

Miscellaneous

M-113.1 Ring adapter, 30mm (for monocular and binocular microscopes) (except ST-50Led)

M-113.2 Ring adapter, 30.5mm (for monocular and binocular microscopes) (only for ST-50Led)

DC-001 Plastic dust cover, small, 340(l)x400(h) mm (except for ST-50Led) DC-002 Plastic dust cover, medium, 490(l)x490(h) mm (only for ST-50Led) Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100) M-005

Cleaning kit <u>15104</u>

ST-041 Sample clip (only for ST-30FX)

AB-010 Antibacterial surface treatment, only for newly purchased microscope

#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.





How to connect the cameras to our microscopes. Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### Headquarters and Manufacturing Facilities

OPTIKA S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### Optika Sales branches

**OPTIKA**° Spain **OPTIKA**° China **OPTIKA**° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

**OPTIKA**<sup>®</sup> North America **OPTIKA**<sup>®</sup> Central America **OPTIKA**® Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# SLX Series



**Stereomicroscopes For Higher Education & Laboratory** 

# Extremely Versatile Cordless Stereo & Stereozoom Microscopes




#### THE LONGEST AUTONOMY ON THE MARKET

- » Longlife LED illumination (providing over 20 years of use)
- » Ultra-flat base with Ø 100 mm disc for diffused transmitted light
- » Cordless use, totally independent from mains/batteries connection
- » Freely settable illumination incident, oblique and transmitted light



1

# **SLX** Series

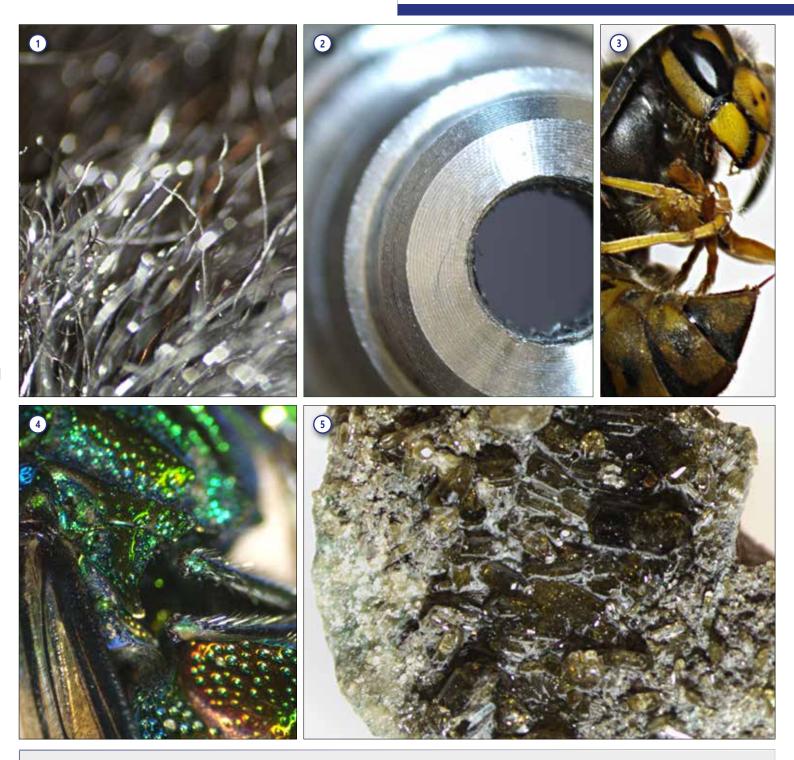

#### Legend

- Aluminum SLX-1 and 4x objective.
   Component worked on lathe SLX-2 and 3x zoom.
- 3. Wasp SLX-3 and 4x zoom.

- 4. Fly, detail SLX-2 and 4.5x zoom.
- 5. Rock SLX-2 1.5x zoom.

1

### **SLX** Series

Valuable configurations of cordless and modern stereo & stereozoom microscopes ideal for a variety of applications, including industrial purposes as well as dissection, biology, entomology, anatomy, chemistry and material science among the others.

Provided with dual magnification or 6.43:1 zoom ratio, FN 21 high eyepoint eyepieces, highgrade precise fixed arm with focus and handle with the latest technology of **EcoLED™** illumination plus rechargeable batteries. Slim and easy to carry, all the models are equipped with long lasting LED illumination to provide over 20 years of use

#### High eyepoint eyepieces for glasses wearers

These eyepieces are designed in such a way that the exit pupil is further away from the eye lens than standard eyepieces, being are well suited for eyeglasses wearers

#### The longest autonomy on the market ensured by EcoLED™

OPTIKA has re-designed illumination in microscopy, once again: a special coating process on optics combined with a new, higher ratio between low consumptions and ultra-efficiency has addressed us to top brightness levels

#### 6.43:1 zoom ratio - zoom magnification from 7x to 45x

Purposely designed for professional routine inspections, the total magnification can be even extended to 180x with 20x eyepieces and 2x additional lens, obtaining an excellent results in this class



#### Ultra-flat base with Ø 100 mm disc for diffused transmitted light

A new level of ergonomy and comfort is achieved during operations, with the ultra-flat base of only 3 cm height to ensure smooth specimen movement and the Ø 100 mm for top class diffusion of the transmitted ligh

### Stereomicroscopes For Higher Education & Laboratory





#### Longlife LED illumination (providing over 20 years of use)

Money & energy saving thanks to LED long lifetime (65.000 hours, 22 years in case of 8 hours/day) which is more than 20 times compared to a standard halogen bulb

#### Cordless use, totally independent from mains/batteries connection

All models work with or without the batteries in place and are provided with three NiMH rechargeable batteries for the longest autonomy in outdoor use (12-hour autonomy, at medium intensity)





#### External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit

### **SLX** Series - Get the most out of our accessories

#### **Additional Lenses**

Simply to be screwed into the threads below the objectives of SLX-2 and SLX-3 to either increase or decrease total magnification, or to increase the working distance when users need to work with hands under the microscope



#### ST-040.1 - Darkfield condenser

This is a darkfield condenser for stereo microscopes with bottom light and 100 mm round working plate to provide darkfield microscopy features, fitting all OPTIKA stereomicroscopes with 100 mm mounting size and transmitted light



### **SLX** Series - Range

### SLX-1

















Cordless binocular stereomicroscope ideal for industrial purposes and students/amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (2x-4x), FN 21 high eyepoint eyepieces, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use

Head: Binocular, 45° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups

**Objective:** Achromatic 2x-4x with anti-fungus treatment.

Working distance: 100 mm.

Stand: High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** EcoLED™ swiveling incident and transmitted, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

### SLX-2

















Cordless binocular stereozoom microscope ideal for industrial purposes and students/amateurs to dissect and discover mainly biology and materials science in 3D, with 0.7x...4.5x zoom, FN 21 high eyepoint eyepieces, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use

**Head:** Binocular, 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups

**Objective:** Parfocal achromatic zoom 0.7x...4.5x (6.43:1 ratio) with anti-fungus treatment.

Working distance: 100 mm.

**Stand:** High-grade, precision fixed with handle and focus.

**Focusing:** Rack and pinion focusing mechanism.

**Illumination:** EcoLED™ swiveling incident and transmitted, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

### **SLX** Series - Range

### SLX-3



Cordless trinocular stereozoom microscope ideal for industrial purposes and students/amateurs to dissect and discover mainly biology and materials science in 3D, with 0.7x...4.5x zoom, FN 21 high eyepoint eyepieces, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use

Head: Trinocular (split ratio: 50/50), 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups

**Objective:** Parfocal achromatic zoom 0.7x...4.5x (6.43:1 ratio) with anti-fungus treatment.

Working distance: 100 mm.

**Stand:** High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination:** EcoLED™ swiveling incident and transmitted, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.

### **SLX** Series - Comparison Chart

Model	Head	Eyepieces	Objective	Working Distance	Stand	Illumination
SLX-1	Binocular 45° inclined 360° rotating	WF 10x/21	2x – 4x selectable	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)
SLX-2	Binocular 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)
SLX-3	Trinocular (50/50) 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)

#### **Optical performance SLX-1**

Eyepiece	10x (ST-081)		15x (ST-082)		20x (ST-083)		10x (ST-084)	
Field number (mm)	21		15		10		21	
Additional lens	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
1x	20x - 40x	10.50 - 5.25	30x - 60x	7.50 - 3.75	40x - 80x	5.00 - 2.50	20x - 40x	10.50 - 5.25

#### **Optical performance SLX-2 - SLX-3**

Eyepiece 10x (ST-081)		15x (ST-082)		20x (S	T-083)	10x (ST-084) 21		
Field number (mm)	21		15		10			
Additional lens	Total magnification	Field of View (mm)						
0.5x	3.5x - 22.5x	60.00 - 9.33	5.25x - 33.75x	42.86 - 6.67	7x - 45x	28.57 - 4.44	3.5x - 22.5x	60.00 - 9.33
0.75x	5,25x - 33.75x	40.00 - 6.22	7.875x - 50.625x	28.57 - 4.44	10.5x - 67.5x	19.05 - 2.96	5.25x - 33.75x	40.00 - 6.22
1x	7x - 45x	30.00 - 4.67	10.5x - 67.5x	21.43 - 3.33	14x - 90x	14.29 - 2.22	7x - 45x	30.00 - 4.67
1.5x	10.5x - 67.5x	20.00 - 3.11	15.75x - 101.25x	14.29 - 2.22	21x - 135x	9.52 - 1.48	10.5x - 67.5x	20.00 - 3.11

### $(\mathbf{T})$

### **SLX** Series - Accessories

### Eyecups & Eyepieces

<u>S1-036</u>	<u>Eyecups (pair), flat</u>
ST-081	EW10x/21 eyepieces (pair), high eyepoint, with rubber cup
CRU-TS	WE15y/15 avaniages (nair) high avanaint

ST-082 WF15x/15 eyepieces (pair), high eyepoint ST-083 WF20x/10 eyepieces (pair), high eyepoint

ST-084 WF10x/21 micrometric eyepiece, high eyepoint, with rubber cup

#### **Additional Lenses**

ST-085.1	Additional	lens 0.5x	(w.d.	165mm	) with	SZ-EXT	(onl	v for SL	X-2 8	ι SLX-3	()
			1				1	/			+

ST-091 Additional lens 0.75x (w.d. 105mm) (only for SLX-2 & SLX-3)

ST-086.1 Additional lens 1.5x (w.d. 45mm) with compensating disc (only for SLX-2 & SLX-3)

#### **Stages**

ST-100.1 Hand moving stage, 100mm diameter

ST-110.1 Moving stage, coaxial knobs, 100mm diameter

ST-111.1 Moving stage, micrometric screws, 100mm diameter

#### **Condensers & Filters**

ST-040.1 Darkfield condenser, 100mm diameter

ST-088.1 Polarising set (filters and rotating stage), 100mm diameter

#### **Camera Adapters**

M-113.1 Ring adapter, 30mm (for monocular and binocular microscopes)

M-115 0.35x C-Mount projection lens

M-114 0.5x C-Mount projection lens M-118 0.75x C-Mount projection lens

M-173 C-Mount projection lens for APS-C/full frame reflex cameras (trino)

M-699 Universal adapter for C-Mount projection lens (trino)

M-620 0.35x focusable C-Mount adapter

M-620.1 0.5x focusable C-Mount adapter
M-620.2 0.65x focusable C-Mount adapter

M-620.3 1x focusable C-Mount adapter

#### Miscellaneous

15104	Cleaning	l/it
13107	Cicariiiiq	NIL

DC-002 Plastic dust cover, medium, 490(l)x490(h) mm

M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)

ST-041 Sample clip

ST-042 White/black object-plate, 100mm diameter

ST-043 Glass object-plate, 100mm diameter

ST-092 Protective glass for stereohead

VP-SLX IQ/OQ/PQ manual for SLX series

AB-020 Antibacterial surface treatment, only for newly purchased microscope



#### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA' S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA° Spain
OPTIKA° China
OPTIKA° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor.

Ideal for precision lens or prism cleaning.

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com






LABORATORY Microscopes

# **LABORATORY** Microscopes

B-290 SERIES - Entry-Level Lab Upright Microscopes	page 81
B-380 SERIES - Middle-Level Routine Lab Upright Microscopes	page 97
B-510 SERIES - Advanced Routine Lab Upright Microscopes	page 121
B-810/B-1000 SERIES - Research Lab Upright Microscopes	page 147
IM-300 SERIES - Routine Lab Inverted Microscopes	page 193
IM-5 SERIES - Routine & Research Lab Inverted Microscopes	page 211
IM-7 MODEL - Inverted Research Microscope	page 231
POL SERIES - Routine & Research Lab Polarizing Microscopes	page 243
FLUO SERIES - Routine & Research Lab Fluorescence Microscopes	page 253
SZR-180 MODEL - Research Stereomicroscope	page 275

### Icons

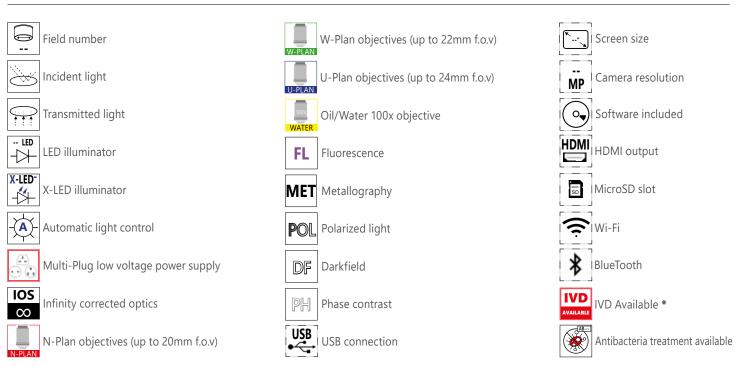

<sup>\*</sup> The IVD code must be requested at order



# B-290 Series



**Entry-Level Lab Upright Microscopes** 

# Best Value-for-Money Solutions & Versatile Use

#### **SUITABLE FOR UNIVERSITIES, EXPERTS & ROUTINE LABS**

- » The ideal choice for common lab requirements
- » Simply engineered for life-science
- » Ready for phase contrast and darkfield

#### **EXCELLENT PRICE/PERFORMANCE RATIO**

- » N-PLAN objectives (160 mm or IOS) for flat images on 20 mm FN
- » Fixed Koehler illumination for crisp and contrasted images
- » Rounded edge, rackless stage to prevent scratches




#### **B-290TB - BREAK NEW GROUND WITH DIGITAL IMAGING UNIT**

- » Large touch-screen of with fast, responsive and smooth control
- » 360° rotatable, tiltable and easily detachable
- » Simultaneous camera & power connection for long-term operation

#### **B-290LD SERIES - ON-FIELD TBC & MALARIA DIAGNOSIS**

- » Ultra-convenient LED fluorescence, blue filter
- » No waiting time, immediate operation
- » Cost-effective, money-saving technology

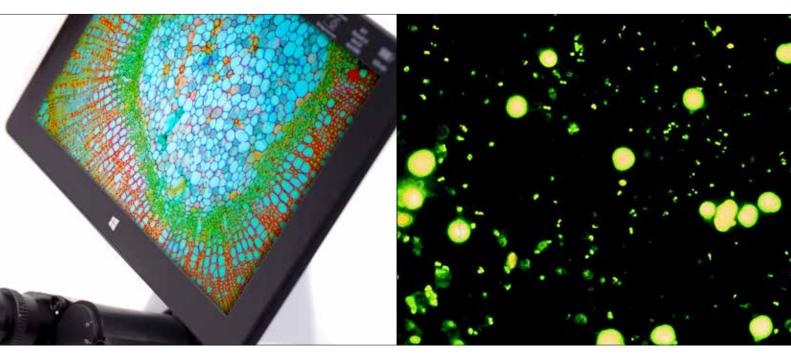

# 100x Oil/Water Objective — Only at OPTIKA

#### SAME OBJECTIVE FOR OIL AND WATER USE

- » Oil represents the best media for high numerical aperture
- » Water combines results with convenience
- » Water is recommended especially for educational purposes

#### **UNPARALLELED TIME & MONEY SAVING**

- » Save time by forgetting about tedious cleaning
- » No time-wasting procedures for maintainance
- » No additional expenses due to inappropriate cleaning



# **X-LED**<sup>3</sup> – Only Available at OPTIKA

#### STATE-OF-THE-ART ILLUMINATION SYSTEM

- » Uncomparable light intensity, exclusive lens & collector design
- » Constant pure-white colour temperature at all intensity levels
- » Unmatched color fidelity, uniformity and brightness

#### **CUT ELECTRICITY BILLS BY 90%**

- » Money & energy saving, only 3.6 W
- » More efficient brightness than a 50 W halogen lamp
- » LED long lifetime (65,000 hours = 22 years at 8 hours/day usage)



# **IOS** & Professional Features

#### LABORATORY GRADE OPTICS FOR OUTSTANDING IMAGES

- » Planachromatic optics
- » Designed to ensure field flatness on 20 mm (N-PLAN)
- » IOS Infinity corrected optical system

#### **FULL CONTROL OF YOUR IMAGES**

- » Fully settable condenser for perfect imaging
- » Easy to set objective-coded iris diaphragm, focusable & centrable
- » Phase contrast and darkfield slider available





# **B-290** Series

This series incorporates all the experience gathered by OPTIKA Microscopes in the field of light microscopy, adapted specifically for common laboratory applications. Suitable for routine microscopy with brightfield, darkfield, phase contrast and LED fluorescence, designed to last.

#### X-LED<sup>3</sup> Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature (6,300K colour temperature).

Relevant money and energy saving thanks to the incredibly low energy consumptions which allows you to cut the electricity bills by 90%!

The electric consumption (3.6 W only) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.

#### **Safe And Convenient Operations**


Rounded edge rackless stage has been designed with a belt-driven mechanism that allows a smooth movement without any protruding part.

This design gives you a more compact solution and lowers any risk of injury after accidentally hitting the rack with your hands.



#### **Clear Observation With 100x Objective**


Students and basic users will enjoy B-290 Series for the clear and sharp images they can get using 100x objective with water, thanks to the extremely bright X-LED<sup>3</sup> light source and the fully centerable Abbe condenser. Forget about the tedious lens cleaning you are used to when using 100x objective: dirt and dust will not affect your objective.

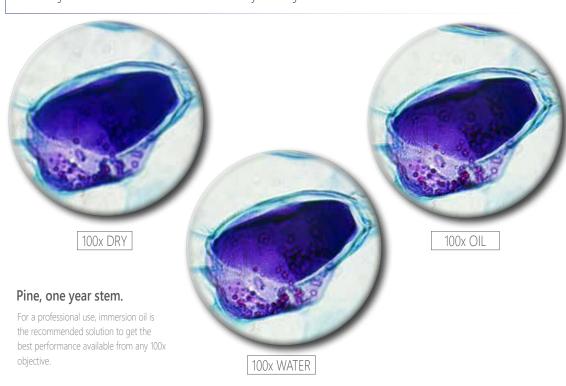



#### Laboratory Grade Optics, N-PLAN & IOS N-PLAN System

OPTIKA N-PLAN objectives ensure bright, clear images with excellent flatness and compensation for chromatic aberration. IOS Infinity-corrected optical system prevents image deterioration even if other optical components are added, such as polarizers, beamsplitters and so on.

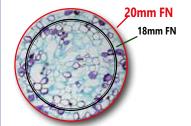

#### **Large Specimen View (20 mm Field Number)**

The **F.O.V.** (field of view) is based on a comfortable diameter of 20 mm. This means that a wide area of the sample can be inspected and allows a natural and easy view, particularly needed in a laboratory environment.

#### **Easy Transportation**

B-290 has a carefully shaped design that gives stability and ease in the transportation. Thanks to the handle on the back, it can be safely moved around the classroom or the laboratory.



### Get the most out of our accessories



#### M-666.290 Heating stage

Use the temperature controlled heating stage for specific applications where accurate temperature selection up to 50° C and programmable settings are needed. Requested in research arenas such as microbiology, biochemistry, material science, pharmacology, crystallographic characterization, surface quality control, and in all the temperature related inspections.

#### M-184 Darkfield condenser

With optional diaphragm M-184 you can easily obtain a darkfield view for dry objectives.
M-184 can be placed directly under the condenser without the need of additional

operations.





(/0.25


### Legend

- 1. B-290TB with 360° rotating Digital Imaging Unit for discussion.
- 2. Pine, one year stem, with B-292 and 10x objective.
- 3. B-290 Series exclusive X-LED<sup>3</sup> illumination system.
- 4. B-293 used by a teacher.
- 5. Lilly pollen, with B-292PLi and 20x objective.
- 6. Ascaris female, with B-293 and 40x objective.
- 7. Student setting the right focusing on a B-290TB.

# **B-290** Series

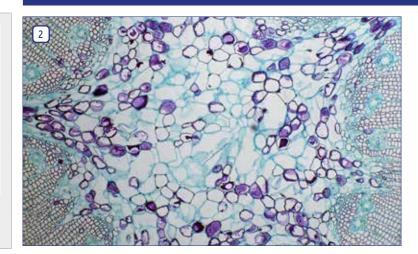

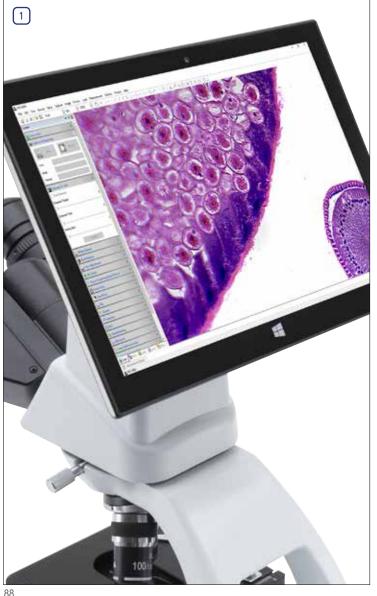


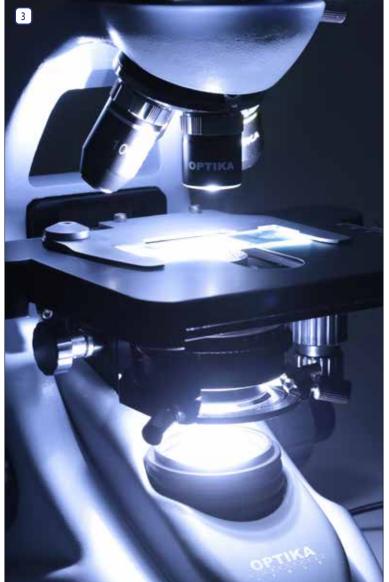

# Entry-Level Lab Upright Microscopes





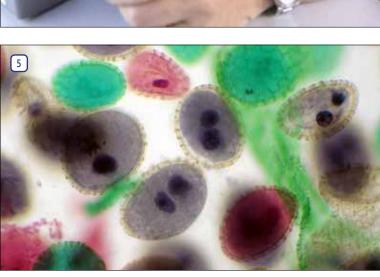



### **B-290** Series - Brightfield Models

### **B-292**



Binocular head with N-PLAN objectives, rackless stage and exclusive **X-LED**<sup>3</sup> for unmatchable performance, powerful and uniform illumination.

**Observation mode:** Brightfield.

**Head:** Binocular, 30° inclined and 360° rotating. Interpupillary distance 48-75 mm

**Dioptric adjustement:** On both eyepieces.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

Nosepiece: Quadruple revolving nosepiece, rotation on ball bearings.

#### **Objectives:**

N-PLAN 4x/0.10 N-PLAN 10x/0.25 N-PLAN 40x/0.65 N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300 K) and light intensity control.

Multi-plug 100-240Vac/6Vdc external power supply.

**B-293** 



Trinocular head with N-PLAN objectives, rackless stage and exclusive **X-LED**<sup>3</sup> for unmatchable performance, powerful and uniform illumination.

Observation mode: Brightfield.

**Head:** Trinocular (fixed, 50/50), 360° rotating and 30° inclined. Interpupillary distance 55-75 mm.

**Dioptric adjustement:** On both eyepieces.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quadruple revolving nosepiece, rotation on ball bearings.

#### **Objectives:**

N-PLAN 4x/0.10 N-PLAN 10x/0.25 N-PLAN 40x/0.65 N-PLAN 100x/1.25 (Oil/Water) All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.


**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> type with white 3.6 W LED (6,300 K) and light intensity control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-290** Series - Brightfield Models

### **B-292PLi**

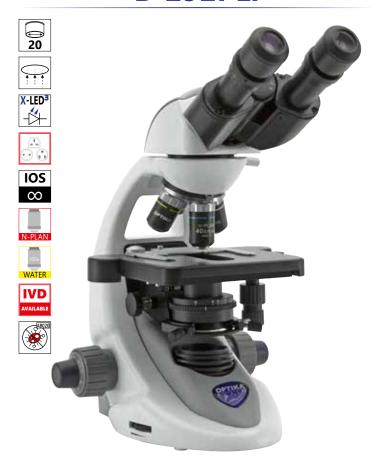

Binocular head with IOS N-PLAN (infinity corrected) objectives, rackless stage and exclusive **X-LED**<sup>3</sup> for incredibly bright illumination.

Observation mode: Brightfield.

**Head:** Binocular, 360° rotating and 30° inclined. Interpupillary distance 48-75 mm

**Dioptric adjustement:** On both eyepieces.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quadruple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25

IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> type with white 3.6 W LED (6,300 K) and light intensity control.

Multi-plug 100-240Vac/6Vdc external power supply.

### **B-293PLi**



Trinocular head with IOS N-PLAN (infinity corrected) objectives, rackless stage and exclusive **X-LED**<sup>3</sup> for incredibly bright illumination.

**Observation mode:** Brightfield.

**Head:** Trinocular (fixed, 50/50), 360° rotating and 30° inclined. Interpupillary distance 55-75 mm.

Dioptric adjustement: On both eyepieces.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quadruple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25

IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> type with white 3.6 W LED (6,300 K) and light intensity control.

Multi-plug 100-240Vac/6Vdc external power supply.

## **B-290TB** - Microscope with Camera & Digital Imaging Unit

### **B-290TB**

























Built-in camera and Digital Imaging Unit with N-PLAN objectives, rackless stage and exclusive **X-LED**<sup>3</sup> for unmatchable performance in illumination. Ideal for discussion group with 360° rotating Digital Imaging Unit.

Observation mode: Brightfield.

**Head:** Binocular, 360° rotating and 30° inclined. Interpupillary distance from 48 to 75 mm; dioptric adjuctment on left eyepiece. Built-in 3.1 MP camera.

**Dioptric adjustement:** On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

Nosepiece: Quadruple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

N-PLAN 4x/0.10 N-PLAN 10x/0.25

N-PLAN 40x/0.65 N-PLAN 100x/1.25 (Oil/Water)


All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> type with white 3.6 W LED (6,300 K) and light intensity control. Multi-plug 100-240Vac/6Vdc external power supply.

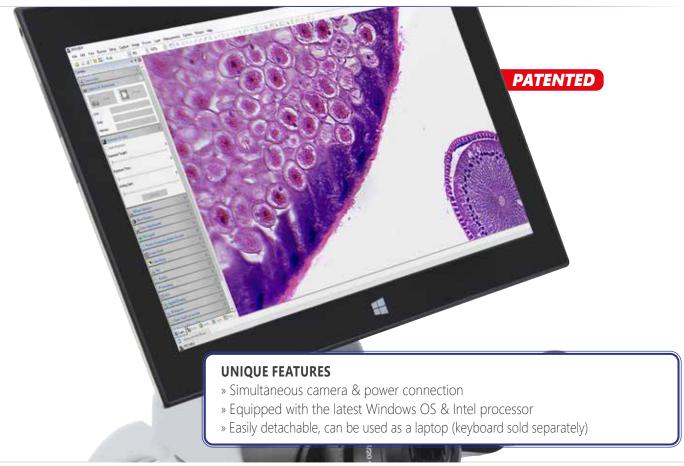

# **B-290TB** - Microscope with Camera & Digital Imaging Unit

The latest OPTIKA digital microscopes with Digital Imaging Unit, open new microscopy horizons, combining high-end optics with innovative digital technology for microscopic imaging. B-290TB includes a 3MP camera with a Digital Imaging Unit. View, capture, analyze and share your images with simplicity and reliability.

DIGITAL IMAGING UNIT TEC	CHNICAL SPECIFICATIONS
Operating system	Windows 11
CPU	Intel® Celeron
CPU speed	1.10 GHz
Graphic card	Intel® HD Graphics 500
RAM	Ram 6 GB DDR3
Display size	LED 10.1" IPS Multi Touch Screen
Display resolution	1920x1200
Storage	Hdd 128 GB (EMMC)
Network	802.11b/g/n/ac,2.4GHz
Input ports	USB 3.0, Type C, Micro SD card reader, Microphone, DC in
Output ports	Headphone - Micro HDMI
Battery Type	Lithium-ion
Battery capacity	3500 mAh
Power consumption	24.05W
Power supply	12V 2A EU / US
Dimensions (mm)	247 x 165 X 9
Weight (Kg)	0.57
Language	Multilanguage

Simple and user-friendly, ideal for students and experienced users.




### 60° Rotaling & Tilling



MP B 2.0 2.5" MOS
2.5"
MOS
4\3
x 1536
2 micron
nes\sec
- 30 fps (640x480)
- Man
- Man
- Man
- Man

### **Detachable**

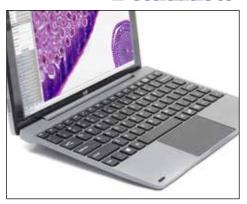

TB-KBD2 Accessory keyboard

### **B-290LD** - LED Fluorescence Microscopes



Fluorescence binocular and trinocular microscopes especially designed for tubercolosis and malaria analysis.

#### **Observation mode:** Brightfield.

**Head:** Binocular or trinocular, 360° rotating and 30° inclined. Interpupillary distance 48-75mm (bino version) 55-75mm (trino version).

**Dioptric adjustement:** On both eyepieces.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by a screw.

Nosepiece: Quadruple revolving nosepiece, rotation on ball bearings.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Brightfield Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300 K) and light intensity control.

Multi-plug 100-240Vac/6Vdc external power supply.

**Fluorescence Illumination:** Extra efficiency LED, with light intensity control. Peak wavelength: 465 nm, Power: 3.6W.

**Epi Fluorescence Attachment:** Slider with 3 positions (2 fluorescence, 1 brightfield), with 1 included filterset: Fluorescence B: EX 460-490, DM 505, EM 515LP: Acridine Yellow, Acridine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP, etc.

#### Part number: B-292LD1.50

Equipped with binocular head and following objectives:

IOS N-PLAN 10x/0.25, with anti-fungus treatment

IOS N-PLAN 20x/0.40, with anti-fungus treatment

IOS N-PLAN 40x/0.65, with anti-fungus treatment

IOS W-PLAN MET 50x/0.75 (No Cover), with anti-fungus treatment.

#### Part number: B-293LD1.50

Trinocular version of B-292LD1.50.

#### Part number: B-292LD1

Equipped with binocular head and following objectives:

IOS N-PLAN 10x/0.25, with anti-fungus treatment

IOS N-PLAN 20x/0.40, with anti-fungus treatment

IOS N-PLAN 40x/0.65, with anti-fungus treatment

IOS W-PLAN 100x/0.80 (No Cover, Dry), with anti-fungus treatment.

#### Part number: B-293LD1

Trinocular version of B-292LD1.

#### **Standard filterset**

Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)	
B (Blue)	460 - 490	505	515LP	

# **B-290** Series - Comparison Chart

Model	Head	Eyepieces	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
B-292	Binocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	N-PLAN 4x, 10x, 40x, 100x (oil/water)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-292PLi	Binocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 4x, 10x, 40x, 100x (oil/water)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-293	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	N-PLAN 4x, 10x, 40x, 100x (oil/water)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-293PLi	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 4x, 10x, 40x, 100x (oil/water)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-290TB	Binocular, 30° inclined, 360° rotating, Digital Imaging Unit	WF 10x/20	Quadruple, reversed	N-PLAN 4x, 10x, 40x, 100x (oil/water)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Coaxial coarse and fine, limit stop	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler

Model	Head	Eyepieces	Nosepiece	Objectives	Stage	Condenser	Fluorescence	Illumination
B-292LD1.50	Binocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 10x, 20x, 40x. W-PLAN MET 50x	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W Fluo LED; Blue Filterset	3.6 W X-LED <sup>3</sup> , brightness control
B-292LD1	Binocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 10x, 20x, 40x. W-PLAN 100x (dry)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W Fluo LED; Blue Filterset	3.6 W X-LED <sup>3</sup> , brightness control
B-293LD1.50	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 10x, 20x, 40x. W-PLAN MET 50x	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W Fluo LED; Blue Filterset	3.6 W X-LED <sup>3</sup> , brightness control
B-293LD1	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quadruple, reversed	IOS N-PLAN 10x, 20x, 40x. W-PLAN 100x (dry)	Double layer, 150x139 mm, moving range 75x33 mm, X-axis rackless	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W Fluo LED; Blue Filterset	3.6 W X-LED <sup>3</sup> , brightness control

#### Optical performance B-292, B-293 and B-290TB

Eyepiece			10x (M-160)		15x (M-161)		20x (M-162)	
Field number (mm)			20		16		10	
Objective	N.A.	W.D (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
4x (N-PLAN 160 mm)	0.10	15.20	40x	5	60x	4	80x	2.5
10x (N-PLAN 160 mm)	0.25	5.50	100x	2	150x	1.6	200x	1
20x (N-PLAN 160 mm)	0.40	5.00	200x	1	300x	0.8	400x	0.5
40x (N-PLAN 160 mm)	0.65	0.45	400x	0.5	300x	0.4	800x	0.25
60x (N-PLAN 160 mm)	0.80	0.13	600x	0.33	900x	0.26	1200x	0.16
100x (N-PLAN 160 mm)	1.25 (oil/water)	0.13	1000x	0.2	1500x	0.16	2000x	0.1

#### Optical performance B-292PLi, B-293PLi and LD models

Eyepiece			10x (M	10x (M-160)		15x (M-161)		20x (M-162)	
Field number (mm)			20	20		16		10	
Objective	N.A.	W.D (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	
4x (IOS N-PLAN)	0.10	16.80	40x	5	60x	4	80x	2.5	
10x (IOS N-PLAN)	0.25	5.80	100x	2	150x	1.6	200x	1	
20x (IOS N-PLAN)	0.40	5.10	200x	1	300x	0.8	400x	0.5	
40x (IOS N-PLAN)	0.65	0.43	400x	0.5	300x	0.4	800x	0.25	
50x (IOS W-PLAN MET, no cover)	0.75	0.32	500x	0.4	750x	0.32	1000x	0.2	
60x (IOS N-PLAN)	0.80	0.14	600x	0.33	900x	0.26	1200x	0.16	
100x (IOS N-PLAN)	1.25 (oil/water)	0.13	1000x	0.2	1500x	0.16	2000x	0.1	
100x (IOS W-PLAN, no cover)	0.80 (dry)	3.20	1000x	0.2	1500x	0.16	2000x	0.1	

### **B-290** Series - Accessories

Eyecups	& Ey	epieces
---------	------	---------

M-001	Huygens 5x eyepiece
M-008.1	WF10x/20 eyepiece, high eyepoint, with pointer, rubber cup
M-160	EW10x/20 eyepiece, high eyepoint, with rubber cup
M-161	EW15x/16 eyepiece, with rubber cup
M-162	WF20x/10 eyepiece
M-163	EW10x/20 micrometric eyepiece, high eyepoint, with rubber cup
Ohiectives	

### N-PLAN

M-164	N-PLAN objective 4x/0.10
M-165	N-PLAN objective 10x/0.25
M-166	N-PLAN objective 20x/0.40
M-167	N-PLAN objective 40x/0.65
M-168	N-PLAN objective 60x/0.85
M-169	N-PLAN objective 100x/1.25 (oil)

W-PLAN	
M-170	W-PLAN PH objective 10x/0.25
M-171	W-PLAN PH objective 20x/0.40
M-172	W-PLAN PH objective 40x/0.65
M-182	W-PLAN PH objective 100x/1.25 (oil)

#### **IOS N-PLAN**

M-144	IOS N-PLAN objective 4x/0.10
M-145	IOS N-PLAN objective 10x/0.25
M-146	IOS N-PLAN objective 20x/0.40
M-147	IOS N-PLAN objective 40x/0.65
M-149	IOS N-PLAN objective 60x/0.80
M-148	IOS N-PLAN objective 100x/1.25 (oil)

IOS W-PLA	AN
M-634.1	IOS W-PLAN objective 50x/0.95 (oil)
M-1120.N	IOS W-PLAN PH objective 10x/0.25
	IOS W-PLAN PH objective 20x/0.40
	IOS W-PLAN PH objective 40x/0.65
	IOS W-PLAN PH objective 100x/1.25 (oil)
	IOS W-PLAN MET objective 50x/0.75
M-698.2	IOS W-PLAN MET objective 100x/0.80 (dry)
CI	

#### Stages

M-175	Rotating stage for polarising set (for 150x139mm rackless stage)
M-666.290	Heating stage (on newly purchased microscopes, for 150x139mm stage), multiplug

#### **Condensers & Filters**

M-174	Polarising set (filters only)
M-184	Darkfield stop for condenser
M-971	Plane-concave mirror, with base
M-975	Blue filter, 45mm diameter
M-977	Green filter, 45mm diameter
M-979	Yellow filter, 45mm diameter
M-989	Frosted glass filter, 45mm diameter
M-1124.1	Brightfield condenser (with phase slider slot)
	(except for B-292, B-293 and B-290TB)

#### **Camera Adapters**

M-114	0.5x C-Mount projection lens
M-115	0.35x C-Mount projection lens
M-118	0.75x C-Mount projection lens
M-173	Photo adapter for APS-C and full frame reflex cameras (tr

### Please refer to the Adapter reference list on Digital section.

How to connect the cameras to our microscopes.

Miscellaneous 15104 Cleaning kit

· · · · · · · · · · · · · · · · · · ·
Immersion oil, 10ml
Immersion oil, 100ml
Plastic dust cover, medium, 490(l)x490(h) mm (except for B-290TB)
TNT dust cover, medium, 600(l)x550(h) mm (only for B-290TB)
Keyboard for Digital Imaging Unit (only for B-290TB)
Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)
Slide PH 10/20 160mm
<u>Slide PH 40 160mm</u>
<u>Slide PH 100 160mm</u>
Slide PH 10/20 IOS
Slide PH 40 IOS
Slide PH 100 IOS
Slide DK DRY
Solar charger
Centering telescope, 23mm diameter (except for B-292, B-293 and B-290TB)
IQ/OQ/PQ manual for B-290 series
IQ/OQ/PQ manual for TB series
Antibacterial surface treatment, only for newly purchased microscope

#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. - Output voltage: 5 Vdc. Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>). Charging models: with solar panel (12h), with external USB power supply (2.5h)





It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.



v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

OPTIKA S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### Optika Sales branches

<b>OPTIKA</b> ° Spain	spain@optikamicroscopes.com
<b>OPTIKA</b> ° China	china@optikamicroscopes.com
<b>OPTIKA</b> ° India	india@optikamicroscopes.com

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# B-380 Series



**Middle-Level Routine Lab Upright Microscopes** 

# Just What You Need. Right When Is Needed

#### **ROUTINE IN UNIVERSITIES, LABS & INDUSTRIES**


- » Wide range to fullfil specific lab requirements
- » Valuable solutions for life and material sciences
- » Compliant with several observation methods

#### THE PREFERRED PARTNER FOR ROUTINE TASKS

- » Full planarity optics on 20 mm (N-PLAN) according to ISO 19012-1
- » Fixed Koehler illumination for crisp and contrasted images
- » Rounded edge, rackless stage to prevent scratches



# Multiple Observation Methods

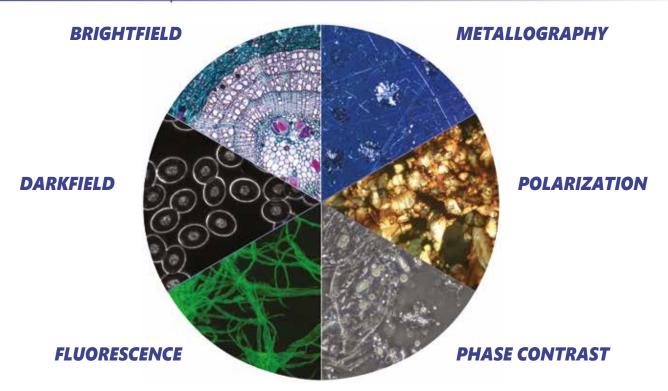

# 100x Oil/Water Objective – Only at OPTIKA

#### SAME OBJECTIVE FOR OIL AND WATER USE

- » Oil represents the best media for high numerical aperture
- » Water combines results with convenience
- » Water is recommended especially for educational purposes

#### **UNPARALLELED TIME & MONEY SAVING**

- » Save time by forgetting about tedious cleaning
- » No time-wasting procedures for maintainance
- » No additional expenses due to inappropriate cleaning



# **X-LED**<sup>3</sup> – Only Available at OPTIKA

#### STATE-OF-THE-ART ILLUMINATION SYSTEM

- » Uncomparable light intensity, exclusive lens & collector design
- » Constant pure-white colour temperature at all intensity levels
- » Unmatched color fidelity, uniformity and brightness

#### **CUT ELECTRICITY BILLS BY 90%**

- » Money & energy saving, only 3.6 W
- » More efficient brightness than a 50 W halogen lamp
- » LED long lifetime (65,000 hours = 22 years at 8 hours/day usage)



# **ALC** - Only Available At OPTIKA


#### **AUTOMATIC LIGHT CONTROL IN 3 STEPS**

- » When another objective is used
- » When the diphragm aperture changes
- » When processing samples with different opacity

#### FORGET ABOUT MANUAL LIGHT ADJUSTEMENT

- » Choose the light intensity according to your preference
- » Press the ALC button and the brightness is saved
- » The microscope will automatically regulate the light



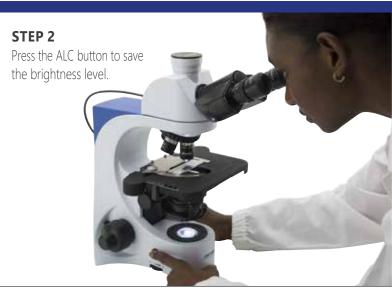

STEP 3

#### Forget about the illumination!

The microscope will automatically adjust the brightness for you, in case of:


- Another objective is used
- The diaphragm aperture is changed
- Another specimen with different opacity is processed







## **B-380** Series

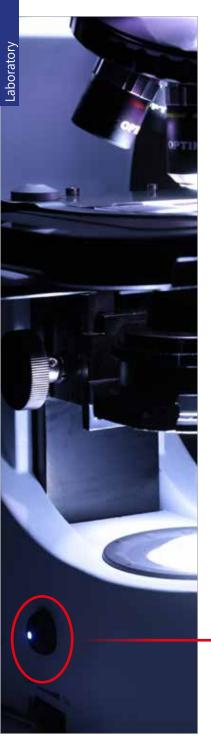

This series incorporates all the experience gathered by OPTIKA Microscopes in the field of light microscopy, adapted specifically for common laboratory applications.

Suitable for routine microscopy with brightfield, darkfield (oil and dry), phase contrast, fluorescence and polarized light, designed to be extremely stable on the bench and last long.

#### **Purposely Designed For Intense Use, Effortless**


Full of features that help being more comfortable especially in case of long-term use. All the main controls are located close to each other to enable minimal movements and reinforce the advantages that the ergonomy brings to this series.

#### X-LED<sup>3</sup> Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature (6,300K colour temperature).

Relevant money and energy saving thanks to the incredibly low energy consumptions which allows you to cut the electricity bills by 90%!

The electric consumption (3.6 W only) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.

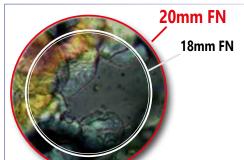

#### **Large Specimen View (20 mm Field Number)**

The **F.O.V.** (field of view) is based on a comfortable diameter of 20 mm.

This means that a wide area of the sample can be inspected and allows a natural and easy view, particularly needed in a laboratory environment.

#### ALC - Automatic Light Control, Only Available At OPTIKA

#### **Incomparable Comfort With The Exclusive Automatic Light Control (ALC)**

Light intensity is automatically adjusted by the microscope itself in order to maintain the same level as the one the user has previously chosen.

No matter if the aperture of the diaphragm changes, if another objective is used, and if the opacity of the sample is different...the microscope will set the light for you according to your preferences.

On ALC Models.

#### **Safe And Convenient Operations**

Rounded edge rackless stage has been designed with a belt-driven mechanism that allows a smooth movement without any protruding part. This design gives you a more compact solution and lowers any risk of injury after accidentally hitting the rack with your hands.



### Middle-Level Routine Lab Upright Microscopes

#### Universal Condenser For Brightfield, Darkfield & Phase

OPTIKA B-380 phase contrast microscopes are equipped with a 5-position dedicated rotating condenser for brightfield (standard use), phase contrast (10x/20x, 40x and 100x phase diaphragms), and a darkfield position for dry objectives.



#### **Exclusive X-LED<sup>3</sup> Darkfield Condenser**

The special condenser with integrated, exclusive X-LED<sup>3</sup> illuminator replaces any other external and expensive lighting source required for these applications and is ideal for great-looking, rich and high-quality specimen view.





#### In fluorescence we can offer several options.

According to your application and to the fluorochromes you are using, we can help you to identify the best light source.

#### **Traditional, HBO Fluorescence**

- » The most used and diffused method, worldwide
- » Wide spectrum range for future upgrades



#### **Innovative, LED Fluorescence**

- » Recommended for routine applications
- » Cost-effective, money saving technology
  - » Ready for immediate operation
  - » Eliminate warm-up/cool-down times
  - » Forget lamp replacement & centering



### Middle-Level Routine Lab Upright Microscopes

### **Get the most of our accessories**

M-181

Complete Phase Contrast Set with IOS W-PLAN PH obj. 10x, 20x, 40x, 100x, with Darkfield position

The B-380 series can be upgraded at any time with phase contrast kits (M-179 with W-PLAN PH objectives and M-181 with IOS W-PLAN PH objectives) including all the components you need to inspect transparent specimens such as microorganisms, thin tissue slices, lithographic patterns, fibers, glass, etc.

#### M-975.1 Ring with blue filter;

Increase the colour temperature of light (toward the blue).

#### M-977.1 Ring with green filter;

Optimize the resolution of phase contrast.

#### M-979.1 Ring with yellow filter;

Decrease the colour temperature of light (toward the red).

#### M-989.1 Ring with frosted glass filter;

Increase the uniformity of illumination, even further.







#### M-156 Koehler field diaphragm

Additional field diaphragm for upgrading the Fixed Koehler illumination system to a Full Koehler type. To be ordered on newly purchased B-380 microscope.

#### M-975 Blue filter;

Increase the colour temperature of light (toward the blue).

#### M-977 Green filter;


Optimize the resolution of phase contrast.

#### M-979 Yellow filter;

Decrease the colour temperature of light (toward the red).

#### M-989 Frosted glass filter;

Increase the uniformity of illumination, even further.

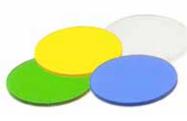


## **B-380** Series


#### Legend


- 1. Planachromatic Phase Contrast objectives.
- 2. Coded iris diaphragm for each objective.
- 3. B-383POL, tuff observed under polarized light.
- 4. Tilia three year stem at 4x magnification, B-383PL.
- 5. B-380 head with built-in Automatic Light Control system.
- 6. Head with Siedentopf adjustment system.
- 7. B-383POL attachment with Bertrand lens.
- 8. Handle for easy and comfortable transportation.
- 9. Coin at 4x magnification, B-383MET.
- 10. Coin at 50x magnification, B-383MET.
- 11. Innovative design of B-380 series.





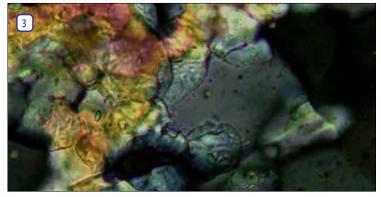

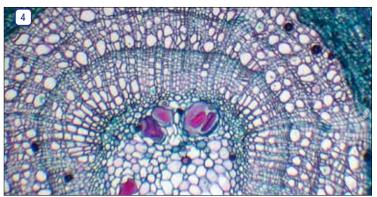


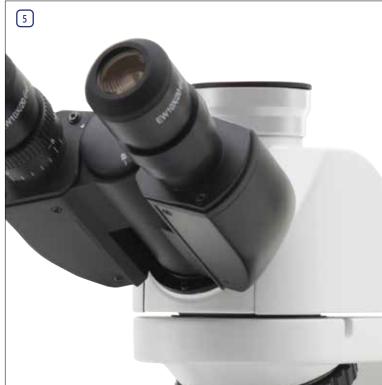







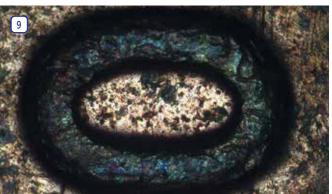





### **B-380** Series - Brightfield Models

### B-382PL-ALC



Brightfield binocular microscope with N-PLAN objectives, rackless stage and combining the exclusive **X-LED**<sup>3</sup> with **ALC** (Automatic Light Control) for great-looking, rich and high-quality view.

**Observation mode:** Brightfield.

**Head:** Binocular, 30° inclined, 360° rotating (when ALC cable is unplugged).

Interpupillary distance: Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

Nosepiece: Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

N-PLAN 4x/0.10 N-PLAN 10x/0.25

N-PLAN 40x/0.65 N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 150x139 mm, 75x33 mm X-Y range.


**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) and brightness control. **ALC** system.

Multi-plug 100-240Vac/6Vdc external power supply.

### **B-383PL**

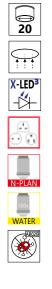



Brightfield trinocular microscope with N-PLAN objectives, rackless stage and the exclusive **X-LED**<sup>3</sup> for great-looking, rich and high-quality view.

Observation mode: Brightfield.

**Head:** Trinocular (fixed 50/50), 30° inclined, 360° rotating.

Interpupillary distance: Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

N-PLAN 4x/0.10 N-PLAN 10x/0.25

N-PLAN 40x/0.65 N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 150x139 mm, 75x33 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.

Multi-plug 100-240Vac/6Vdc external power supply.

### **B-380** Series - Brightfield Models

### B-382PLi-ALC





















Brightfield binocular microscope with IOS N-PLAN (Infinity Corrected) objectives, rackless stage and combining the exclusive **X-LED**<sup>3</sup> with **ALC** (Automatic Light Control) for great-looking, rich and high-quality view.

Observation mode: Brightfield.

**Head:** Binocular, 30° inclined, 360° rotating (when ALC cable is unplugged).

**Interpupillary distance:** Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

Nosepiece: Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25

IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) and brightness control. **ALC** system.

Multi-plug 100-240Vac/6Vdc external power supply.

**B-383PLi** 



Brightfield trinocular microscope with IOS N-PLAN (Infinity Corrected) objectives, rackless stage and the exclusive **X-LED**<sup>3</sup> for great-looking, rich and high-quality view.

Observation mode: Brightfield.

**Head:** Trinocular (fixed 50/50), 30° inclined, 360° rotating.

Interpupillary distance: Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25

IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

**⊕** 20

X-LED<sup>3</sup>

W-PLAN

DF

### **B-380** Series - Phase Contrast Models

### B-382PH-ALC



Phase contrast, darkfield and brightfield binocular microscope with W-PLAN objectives, rackless stage and combining the exclusive **X-LED³** with **ALC** (Automatic Light Control) for great-looking, rich and high-quality view.

**Observation mode:** Brightfield, phase contrast and darkfield (dry).

**Head:** Binocular, 30° inclined, 360° rotating (when ALC cable is unplugged).

**Interpupillary distance:** Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

N-PLAN 4x/0.10 W-PLAN PH 10x/0.25

W-PLAN PH 40x/0.65 W-PLAN PH 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 150x139 mm, 75x33 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) and brightness control. **ALC** system.

Multi-plug 100-240Vac/6Vdc external power supply.

### **B-383PH**



Phase contrast, darkfield and brightfield trinocular microscope with W-PLAN objectives, rackless stage and the exclusive **X-LED**<sup>3</sup> for great-looking, rich and high-quality view.

**Observation mode:** Brightfield, phase contrast and darkfield (dry).

**Head:** Trinocular (fixed 50/50), 30° inclined, 360° rotating.

Interpupillary distance: Adjustable between 48 and 75 mm.

Dioptric adjustment: On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

**Objectives:** 

N-PLAN 4x/0.10 W-PLAN PH 10x/0.25

W-PLAN PH 40x/0.65 W-PLAN PH 100x/1.25 (Oil/Water)

All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 150x139 mm, 75x33 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

### **B-380** Series - Phase Contrast Models

### B-382PHi-ALC























Phase contrast, darkfield and brightfield binocular microscope with IOS W-PLAN (Infinity Corrected) objectives, rackless stage and combining the exclusive **X-LED**<sup>3</sup> with **ALC** (Automatic Light Control) for great-looking, rich and high-quality view.

**Observation mode:** Brightfield, phase contrast and darkfield (dry).

Head: Binocular, 30° inclined, 360° rotating (when ALC cable is unplugged).

Interpupillary distance: Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

#### **Objectives:**

IOŚ W-PLAN PH 10x/0.25 IOS W-PLAN PH 20x/0.40 IOS W-PLAN PH 40x/0.65 IOS W-PLAN PH 100x/1.25 (Oil) All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) and brightness control. **ALC** system. Multi-plug 100-240Vac/6Vdc external power supply.

### **B-383PHi**























Phase contrast, darkfield and brightfield microscope with IOS W-PLAN (Infinity Corrected) objectives, rackless stage and the exclusive **X-LED**<sup>3</sup> for great-looking, rich and high-quality view.

**Observation mode:** Brightfield, phase contrast and darkfield (dry).

**Head:** (50/50), 30° inclined, 360° rotating.

Interpupillary distance: Adjustable between 48 and 75 mm.

**Dioptric adjustment:** On the left eyepiece tube.

Eyepieces: WF10x/20 mm, high eye-point and secured by screw.

**Nosepiece:** Quintuple revolving nosepiece, rotation on ball bearings.

#### **Objectives:**

IOS W-PLAN PH 10x/0.25 IOS W-PLAN PH 20x/0.40 IOS W-PLAN PH 40x/0.65 IOS W-PLAN PH 100x/1.25 (Oil) All with anti-fungus treatment.

**Specimen stage:** Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.

**Focusing:** Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

**Condenser:** Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield.

**Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

### **B-383DK** - Darkfield Microscope

Laboratory upright microscope for brightfield and darkfield observations with N-PLAN objectives (and W-PLAN 100x with iris diaphragm) for biology and especially darkfield fresh blood analysis and the exclusive **X-LED³** illumination system. The special condenser with integrated, exclusive X-LED3 illuminator replaces any other external and expensive lighting source required for these applications and is ideal for great-looking, rich and high-quality specimen view. Our immersion darkfield system provides the same result achieved by 150W external illuminators in combination with traditional cardioid darkfield condenser.


















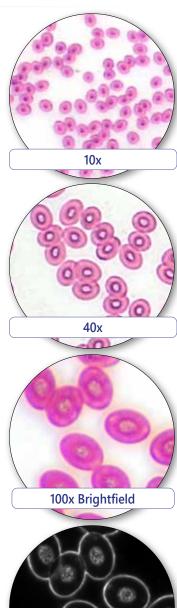

100x Darkfield

Part	Description					
<b>Observation mode:</b>	Brightfield, oil immersion darkfield.					
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.					
Interpupillary distance:	Adjustable between 48 and 75 mm.					
Dioptric adjustment:	On the left eyepiece tube.					
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.					
Nosepiece	Quintuple revolving nosepiece, rotation on ball bearings.					
Objectives:	N-PLAN 4x/0.10 N-PLAN 10x/0.25 N-PLAN 40x/0.65 W-PLAN 100x/1.25 (oil) with iris All with anti-fungus treatment.					

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 150x139 mm, 75x33 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Brightfield condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.
Darkfieldfield condenser:	Darkfield N.A. 1.36 (oil immersion) with built-in X-LED <sup>3</sup> .
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-383FL** - HBO Fluorescence Microscope

Laboratory upright microscope for brightfield and fluorescence observations with IOS N-PLAN objectives.

The HBO fluorescence illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Auramine O, Acridine Yellow, FITC and GFP (with blue filter set), Rhodamine, Nile Red, Cy5.5 and Propidium lodide (with the green one), yet giving the possibility to combine any other specific filter sets for future upgrade. Transmitted light through the exclusive **X-LED<sup>3</sup>** to ensure great-looking, rich and high-quality specimen view.




Part	Description				
<b>Observation mode:</b>	Brightfield, HBO fluorescence.				
<b>Epi-illumination and filters:</b>	HBO 100 W high pressure mercury lamp. 3-position filter holder; blue and green included.				
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.				
Interpupillary distance:	Adjustable between 48 and 75 mm.				
Dioptric adjustment:	On the left eyepiece tube.				
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25 IOS N-PLAN 20x/0.40 IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water) All with anti-fungus treatment.				

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-383LD** - LED Fluorescence Microscope

Entry-level laboratory upright microscope for brightfield and fluorescence observations with IOS N-PLAN objectives. The extremely powerful LED fluorescence illuminator is combined with blue excitation filter set for the visualization of the following fluorochromes: Acridine Yellow, Acridine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP, etc. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>3</sup> to ensure great-looking, rich and high-quality specimen view.

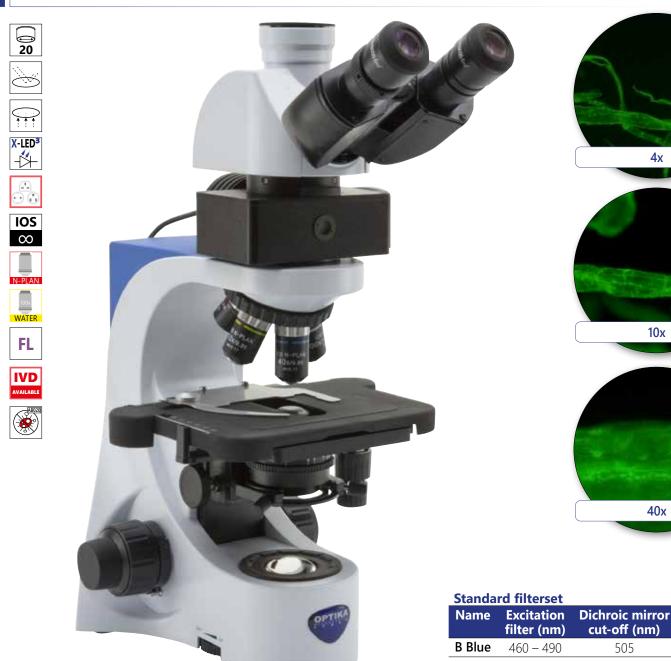

Part	Description					
<b>Observation mode:</b>	Brightfield, LED fluorescence.					
<b>Epi-illumination and filter:</b>	High-power blue LED with brightness control. 3-position filter holder; blue included.					
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.					
Interpupillary distance:	Adjustable between 48 and 75 mm.					
Dioptric adjustment:	On the left eyepiece tube.					
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.					
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.					
Objectives:	4x/0.10, W.D. 16.8 mm 20x/0.40, W.D. 5.1 mm 100x/1.25 (Oil/Water), W.D. 0.13 mm All with anti-fungus treatment.					

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

**4**x

40x

**Emission** 

filter (nm)

515LP

# **B-383POL** - Polarizing Microscope

Upright microscope for brightfield and polarizing light observations with strain-free IOS N-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED<sup>3</sup>** illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.























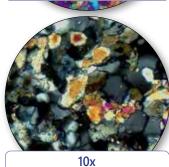



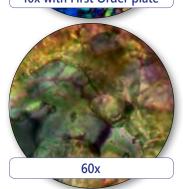

Part	Description					
Observation mode:	Brightfield, transmitted polarized light and conoscopy.					
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale.  Tint plates: 1° order red (λ); λ/4; Quartz wedge.					
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.					
Interpupillary distance:	Adjustable between 48 and 75 mm.					
Dioptric adjustment:	On the left eyepiece tube.					
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw. One with crosshair.					
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings Centering system for each objective.					

Part	Description			
Objectives (strain-free):	IOS N-PLAN POL 4x/0.10 IOS N-PLAN POL 10x/0.25 IOS N-PLAN POL 40x/0.65 IOS N-PLAN POL 60x/0.80 All with anti-fungus treatment.			
Specimen stage:	Rotatable stage with locking mechanism. Vernier scale with accuracy 0.1 mm. Diameter 160 mm. Specimen slide clamps.			
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.			
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable. With rotating polarizing filter.			
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.			

## **B-383MET** - Metallurgical Microscope

Brightfield upright microscope with IOS W-PLAN MET objectives and metallurgical attachment combining the exclusive **X-LED³** lighting source both for incident and transmitted illumination. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.
























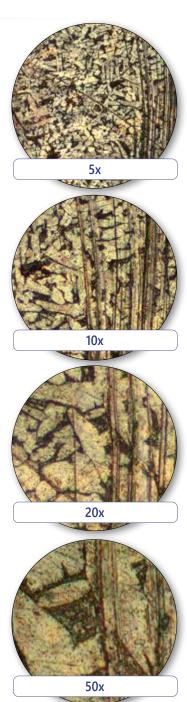

Part	Description					
Observation mode:	Brightfield, incident polarized light.					
Epi-illumination and polarizing filters:	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Field and aperture diaphragms, polarizer & analyzer filters.					
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.					
Interpupillary distance:	Adjustable between 48 and 75 mm.					
Dioptric adjustment:	On the left eyepiece tube.					
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.					
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.					

Part	Description			
Objectives (strain-free):	$\begin{array}{llllllllllllllllllllllllllllllllllll$			
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range. With tempered glass plate.			
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.			
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.			
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.			

# **B-380** Series - Comparison chart

Model	Head	Eyepieces	Nosepiece	Objectives	Stage	Focusing	Condenser	Illumination
B-382PL-ALC	Binocular, 30° inclined	WF 10x/20	Quintuple, reversed	N-PLAN 4x, 10x, 40x, 100x (oil/ water)	Rackless double layer, 150x139mm, moving range 75x33mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control, ALC control. Fixed Koehler
B-383PL	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	N-PLAN 4x, 10x, 40x, 100x (oil/ water)	Rackless double layer, 150x139mm, moving range 75x33mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-382PLi-ALC	Binocular, 30° inclined	WF 10x/20	Quintuple, reversed	IOS N-PLAN 4x, 10x, 40x, 100x (oil/ water)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control, ALC control. Fixed Koehler
B-383PLi	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	IOS N-PLAN 4x, 10x, 40x, 100x (oil/ water)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-382PH-ALC	Binocular, 30° inclined	WF 10x/20	Quintuple, reversed	W-PLAN 4x, 10xPH, 40xPH, 100xPH (oil)	Rackless double layer, 150x139mm, moving range 75x33mm	Coaxial coarse and fine, limit stop, adjustable tension	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield	3.6 W X-LED <sup>3</sup> , brightness control, ALC control. Fixed Koehler
B-383PH	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	W-PLAN 4x, 10xPH, 40xPH, 100xPH (oil)	Rackless double layer, 150x139mm, moving range 75x33mm	Coaxial coarse and fine, limit stop, adjustable tension	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-382PHi-ALC	Binocular, 30° inclined	WF 10x/20	Quintuple, reversed	IOS W-PLAN 10xPH, 20xPH, 40xPH, 100xPH (oil)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield	3.6 W X-LED <sup>3</sup> , brightness control, ALC control. Fixed Koehler
В-383РНі	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	IOS W-PLAN 10xPH, 20xPH, 40xPH, 100xPH (oil)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-383DK	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	N-PLAN 4x, 10x, 40x, W-PLAN 100x (oil, with iris diaphragm)	Rackless double layer, 150x139mm, moving range 75x33mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable. Additional darkfield condenser, N.A. 1.36, built-in X-LED <sup>3</sup>	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-383FL	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	IOS N-PLAN 4x, 10x, 20x, 40x, 100x (oil/water)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	Transmitted: 3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler Incident: HBO 100 W high- pressure mercury bulb
B-383LD	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	IOS N-PLAN 4x, 10x, 20x, 40x, 100x (oil/water)	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	Transmitted: 3.6 W X-LED³, brightness control. Fixed Koehler Incident: High-power blue LED
B-383POL	Trinocular, 30° inclined, 360° rotating	WF 10x/20 (one with crosshair reticle)	Quadruple, reversed	Strain-free IOS N-PLAN POL 4x, 10x, 40x, 60x	Round, 360° rotating, 160mm diameter, with sample clips and stop knob	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable. With rotating polarizer	3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler
B-383MET	Trinocular, 30° inclined, 360° rotating	WF 10x/20	Quintuple, reversed	IOS W-PLAN MET 5x, 10x, 20x, 50x	Rackless double layer, 233x147mm, moving range 78x54 mm	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, iris diaphragm, focusable and centrable	<b>Transmitted and incident:</b> 3.6 W X-LED <sup>3</sup> , brightness control. Fixed Koehler

### **B-380** Series - Accessories

Eyecups &	Eyepieces
M-001	<u>Huygens 5x eyepiece</u>
M-008.1	WF10x/20 eyepiece, high eyepoint, with pointer, rubber cup
M-160	EW10x/20 eyepiece, high eyepoint, with rubber cup
M-161	EW15x/16 eyepiece, with rubber cup
M-162	WF20x/10 eyepiece
M-163	EW10x/20 micrometric eyepiece, high eyepoint, with rubber cu
<b>Objectives</b>	
N-PLAN	
M-164	N-PLAN objective 4x/0.10
M-165	N-PLAN objective 10x/0.25
M-166	N-PLAN objective 20x/0.40
M-167	N-PLAN objective 40x/0.65
M-168	N-PLAN objective 60x/0.85
M-169	N-PLAN objective 100x/1.25 (oil)
IOS N-PLA	
M-144	IOS N-PLAN objective 4x/0.10
M-145	IOS N-PLAN objective 10x/0.25
M-146	IOS N-PLAN objective 20x/0.40
M-147	IOS N-PLAN objective 40x/0.65
M-149	IOS N-PLAN objective 60x/0.80
M-148	IOS N-PLAN objective 100x/1.25 (oil)
M-144P	IOS N-PLAN POL objective 4x/0.10
M-145P	IOS N-PLAN POL objective 10x/0.25
M-146P	IOS N-PLAN POL objective 20x/0.40
M-147P	IOS N-PLAN POL objective 40x/0.65
M-149P	IOS N-PLAN POL objective 60x/0.80
M-148P	IOS N-PLAN POL objective 100x/1.25 (oil)
W-PLAN	<b>,</b> , , , , , , , , , , , , , , , , , ,
M-059	W-PLAN objective 100x/1.25OI - (oil) objective with iris for DF
M-170	W-PLAN PH objective 10x/0.25
M-171	W-PLAN PH objective 20x/0.40
M-172	W-PLAN PH objective 40x/0.65
M-182	W-PLAN PH objective 100x/1.25 (oil)
IOS W-PL/	
M-634.1	IOS W-PLAN objective 50x/0.95 (oil)
M-336	IOS W-PLAN MET objective 5x/0.12
M-338	IOS W-PLAN MET objective 10x/0.25
M-339	IOS W-PLAN MET objective 20x/0.40
M-335	IOS W-PLAN MET objective 50x/0.75
M-698.2	IOS W-PLAN MET objective 100x/0.80 (dry)
M-1120.N	IOS W-PLAN PH objective 10x/0.25
	IOS W-PLAN PH objective 20x/0.40
M-1122.N	IOS W-PLAN PH objective 40x/0.65
M 1122 N	100 W DI AN DI -1:

Stages
M-175

M-1123.N IOS W-PLAN PH objective 100x/1.25 (oil)

Rotating stage for polarising set (for 150x139mm rackless stage) (except for B-382PH-ALC, B-383PH and B-383DK) Rotating stage for polarising set (for 233x147mm rackless stage) (except for B-382PHi-ALC and B-383PHi) M-175.1 M-635 Heating stage (on newly purchased microscopes, for 233x147mm stage), multiplug M-666.290 Heating stage (on newly purchased microscopes, for 150x139mm stage), multiplug





### **B-380** Series - Accessories

Condensers	& Filters	
M-174.1	Polarising set (filters only) (except for B-383POL)	
M-179	PH set - 10x, 40x, 100x W-PLAN PH obj. & BF/DF/PH condenser	
M-181	PH set - 10x, 20x, 40x, 100x IOS W-PLAN PH obj. & BF/DF/PH condenser	A second
M-185	Darkfield condenser for dry objectives	
M-975.1	Ring with blue filter, 45mm diameter	
M-977.1	Ring with green filter, 45mm diameter	
M-979.1	Ring with yellow filter, 45mm diameter	
M-989.1	Ring with frosted glass filter, 45mm diameter	
M-1124.1	Brightfield condenser (with phase slider slot)	
Camera Ada	pters	9
M-115	0.35x C-Mount projection lens	
M-114	0.5x C-Mount projection lens	
M-118	0.75x C-Mount projection lens	
M-173	C-Mount projection lens for APS-C/full frame reflex cameras (trino)	
M-620	0.35x focusable C-Mount adapter (biological microscopes)	
M-620.1	0.5x focusable C-Mount adapter (biological microscopes)	
M-620.2	0.65x focusable C-Mount adapter (biological microscopes)	M-181 - PH set - 10x, 20x, 40x, 100x IOS W-PLAN PH o
M-620.3	1x focusable C-Mount adapter (biological & stereomicroscopes)	
M-699	Universal adapter for C-Mount projection lens (trino)	
Miscellaneo		
15008	Immersion oil, 10ml	
15009	Immersion oil, 100ml	
15104	Cleaning kit	
DC-002	Plastic dust cover, medium, 490(l)x490(h) mm (except for B-383POL, B-38	
DC-003	TNT dust cover, medium, 600(l)x550(h) mm (only for B-383POL, B-383ME	T, B-383LD1, B-383LD2 and B-383FL)
M-005	Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)	
M-102	Slide PH 10/20 160mm	
M-103	<u>Slide PH 40 160mm</u>	
M-104	Slide PH 100 160mm	
M-105	Slide PH 10/20 IOS	
M-106	Slide PH 40 IOS	
M-107	Slide PH 100 IOS	
M-108	Slide DK DRY	
14000		

HBO 100W high-pressure mercury bulb for fluorescence (only for B-383FL)

Antibacterial surface treatment, only for newly purchased microscope

HBO 100W high-pressure mercury bulb for fluorescence (OSRAM) (only for B-383FL) Koehler field diaphragm (on newly purchased microscopes) (except for B-383POL)





M-069

M-151

M-151.1

M-156 M-1380

VP-380

AB-020

#### How to connect the cameras to our microscopes.

Solar charger

Please refer to the Adapter reference list on Digital section.

Centering telescope, 23mm diameter

IQ/OQ/PQ manual for B-380 series



#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. - Output voltage: 5 Vdc. Autonomy: over 6 hours at medium intensity (X-LED<sup>3</sup>). Charging models: with solar panel (12h), with external USB power supply (2.5h)



 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# B-510 Series



**Advanced Routine Lab Upright Microscopes** 

# Born To Be Professional

#### HIGH-GRADE CONFIGURATIONS FOR PROFESSIONALS

- » Wide range to fullfil specific lab requirements
- » Valuable solutions for life and material sciences
- » Compliant with several observation methods

#### AN AFFORDABLE PARTNER WITH HIGH-END FEATURES

- » IOS W-PLAN objectives for flat images on 22 mm FN
- » Full Koehler illumination for enhanced images
- » Rounded edge, rackless stage to prevent scratches



# Optically Impressive

#### **MAINTAINING GOOD EYESIGHT**

- » 10x/22 eyepieces for large specimen view
- » Comfortable rubber cup to get rid of annoying external light
- » High eye-point for glasses wearers, dioptric adjustment (left eyepiece)

#### **B-510 & IOS W-PLAN: THE PERFECT COMBINATION**

- » IOS Infinity corrected optical system
- » Full planarity optics on 22 mm (W-PLAN) according to ISO 19012-1
- » High-grade Semi-Apo lens available ideal for fluorescence






# An Extensive Range of Different Configurations


### THE MODELS FOR ASBESTOS TESTING IN ENVIRONMENTAL MONITORING

- » B-510ASB measures fiber concentrations in air and includes 40xPH lens and 12.5x Walton & Beckett eyepieces
- » B-510POL & B-510POL-I for bulk/fiber class identification


### B-510DK - THE DEDICATED MODEL FOR FRESH BLOOD ANALYSIS

- » Exclusive X-LED³ darkfield cardioid condenser with high N.A. 1.36 and the new IOS W-PLAN 100x oil iris objective
- » Brightfield condenser also supplied

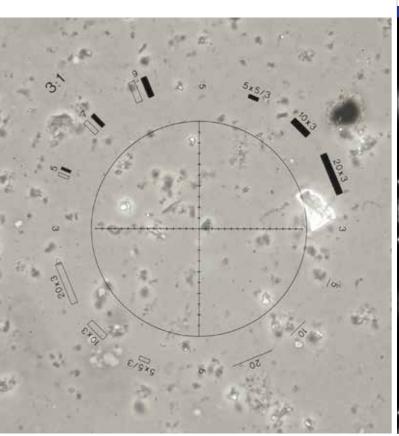

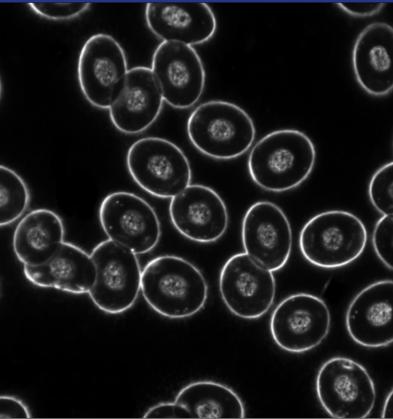

# Multiple Observation Methods

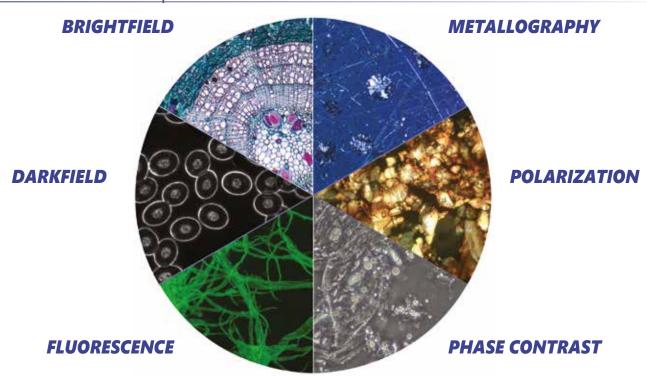

# Many Specimens, Many Observers - Intense Productivity

#### **INCREASE YOUR SAMPLE THROUGHPUT**

- » Large, resistant stage to easily and quickly process 2 samples
- » Ergonomic design and controls for extended operation
- » Convenient handle for easy transportation

### DISCUSSION BRIDGES FOR SIMULTANEOUS OBSERVATIONS AND TEACHING

- » RGB pointer with brightness adjustment for the main observer
- » Face-to-face attachment with 1 extra viewing head, 20mm FN
- » Side attachment with 1, 2 & 4 extra viewing heads, 20mm FN




# **X-LED<sup>3</sup>** - Only Available at OPTIKA

#### STATE-OF-THE-ART ILLUMINATION SYSTEM

- » Uncomparable light intensity, exclusive lens & collector design
- » Constant pure-white color temperature at all intensity levels
- » Unmatched color fidelity, uniformity and brightness

#### **CUT ELECTRICITY BILLS BY 90%**

- » Money & energy saving, only 3.6 W
- » More efficient brightness than a 50 W halogen lamp
- » LED long lifetime (65,000 hours = 22 years at 8 hours/day usage)

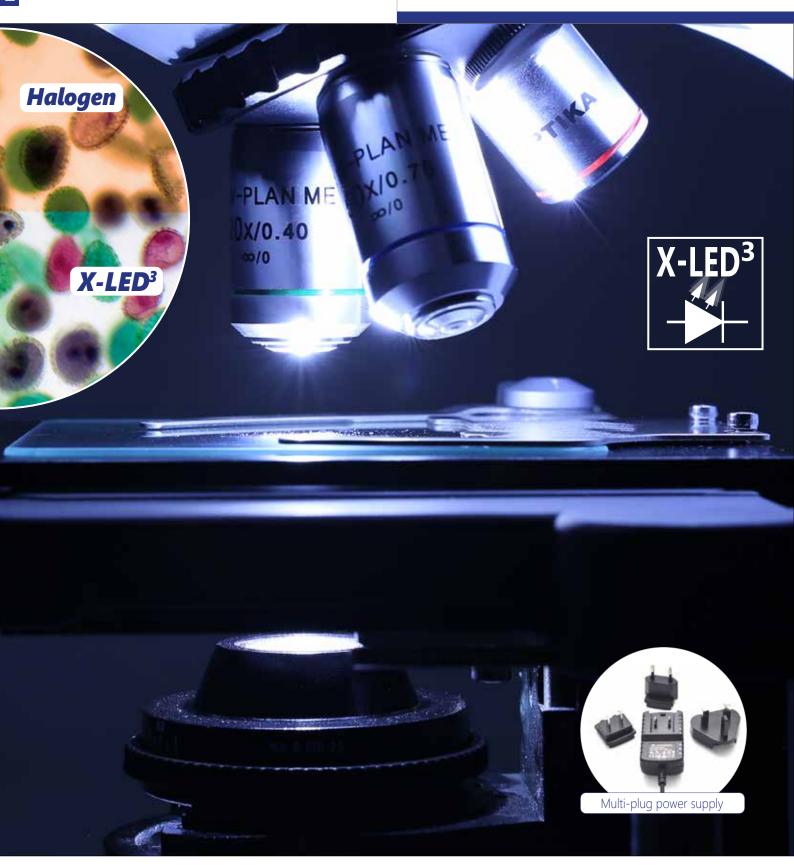

# Go Digital - Vivid Colors & Contrast For Stunning Images

#### STAY CONNECTED WITH YOUR SPECIMEN, EASILY

- » Trincular port to be always updated with the latest technology cameras, even in the future
- » Wide range of cameras matching all the needs, including the more specific ones
- » Modern C-mount focusable professional adapters for all kinds of cameras

#### **PROFESSIONAL IMAGE ANALYSIS**

- » Multi-language software for live-view, picture and video in different file formats
- » Advanced functions for pictures processing (EDF, stitching, multi-fluorescence combine)
- » Powerful tools to perform measurements and generate custom reports



### **B-510** Series

OPTIKA B-510 Series meets a wide variety of analysis applications, thanks to the comprehensive range of microscope models equipped with enhanced and impressive optics, a wide field of view of 22 mm, the state-of-the-art, exclusive X-LED lighting source and Koehler illuminator to produce high sample contrast and homogeneous bright light.

#### **A Perfect Downgrading of Top-Level Series**


Many components of B-510 come from the B-810/1000 Series, the top-level in OPTIKA range, to ensure the state-of-the-art performance and at the same time an incredible level of reliability and durability. Its excellent quality/price ratio is achieved through an intelligent rationalization of production costs and choice of materials.

#### X-LED<sup>3</sup> Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature (6,300K colour temperature).

Relevant money and energy saving thanks to the incredibly low energy consumptions which allows you to cut the electricity bills by 90%!

The electric consumption (3.6 W only) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.

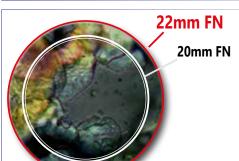

#### **Large Specimen View (22 mm Field Number)**

The **F.O.V.** (field of view) is based on a comfortable diameter of 22 mm.

This means that an extra wide area of the sample can be inspected and allows a natural and easy view, particularly needed in a laboratory environment.

#### **Safe And Convenient Operations**

Rounded edge rackless stage has been designed with a belt-driven mechanism that allows a smooth movement without any protruding part.

This design gives you a more compact solution and lowers any risk of injury after accidentally hitting the rack with your hands.

The electric consumption (3.6 W only) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.



### Advanced Routine Lab Upright Microscopes

### Universal Condenser For Brightfield, Darkfield & Phase

OPTIKA B-510 phase contrast microscopes are equipped with a 5-position dedicated rotating condenser for brightfield (standard use), phase contrast (10x/20x, 40x and 100x phase diaphragms), and a darkfield position for dry objectives.



#### **Exclusive X-LED<sup>3</sup> Darkfield Condenser**

The special condenser with integrated, exclusive X-LED<sup>3</sup> illuminator replaces any other external and expensive lighting source required for these applications and is ideal for great-looking, rich and high-quality specimen view.



#### In fluorescence we can offer several options.

According to your application and to the fluorochromes you are using, we can help you to identify the best light source.

#### **Traditional, HBO Fluorescence**

- » The most used and diffused method, worldwide
- » Wide spectrum range for future upgrades

#### **Innovative, LED Fluorescence**

- » Recommended for routine applications
- » Cost-effective, money saving technology
  - » Ready for immediate operation
  - » Eliminate warm-up/cool-down times
  - » Forget lamp replacement & centering



# **B-510BF / B-510ERGO** - Brightfield Microscope

Advanced routine laboratory microscope for brightfield observations with IOS W-PLAN objectives and rackless stage.

The high-efficiency **X-LED³** makes it reliable for all transmitted light observations for great-looking, rich and high-quality view.










Part	Description
<b>Observation mode:</b>	Brightfield.
	Trinocular (fixed 50/50), 30° inclined, 360° rotating. Binocular ergonomical head, 30°-60° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Dioptric adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS W-PLAN 4x/0.10 IOS W-PLAN 10x/0.25 IOS W-PLAN 40x/0.65 IOS W-PLAN 100x/1.25 (Oil) All with anti-fungus treatment.

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

# **B-510PH** - Phase Contrast Microscope

Advanced routine laboratory microscope for brightfield, darkfield and phase contrast observations with IOS W-PLAN PH objectives and rackless stage. Especially dedicated to phase contrast observation, the microscope ensures a high image sharpness even with complex specimens. The high-efficiency **X-LED<sup>3</sup>** makes it reliable for all transmitted light observations.



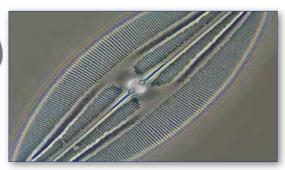








Part	Description
Observation mode:	Brightfield, phase contrast and darkfield (dry).
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS W-PLAN PH 10x/0.25 IOS W-PLAN PH 20x/0.40 IOS W-PLAN PH 40x/0.65 IOS W-PLAN PH 100x/1.25 (Oil) All with anti-fungus treatment.

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510ASB** - Asbestos Analysis Microscope

Advanced routine laboratory microscope for brightfield and phase contrast observations with IOS W-PLAN objectives and rackless stage. Ideal for Asbestos analysis in accordance to international rules with 12.5x eyepieces and Walton & Becket graticule to perform perfect asbestos fibers analysis at a glance. The high-efficiency **X-LED³** makes it reliable for all transmitted light observations.



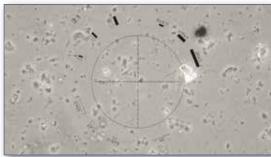






Part	Description	
<b>Observation mode:</b>	Brightfield, phase contrast.	
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Dioptric adjustment:	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and WF12.5x/15 mm with dioptric adjustment, one with Walton & Beckett graticule.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS W-PLAN 4x/0.10 IOS W-PLAN 10x/0.25 IOS W-PLAN PH 40x/0.65 IOS W-PLAN 100x/1.25 (Oil) All with anti-fungus treatment.	

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable. With 40x phase contrast slider.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510DK** - Immersion Darkfield Microscope

Advanced routine laboratory microscope for brightfield and darkfield observations with IOS W-PLAN objectives (including 100x with iris diaphragm) and rackless stage for biology and especially darkfield fresh blood analysis and the exclusive **X-LED³** illumination system. The special condenser with integrated, exclusive X-LED³ illuminator replaces any other external and expensive lighting source required for these applications and is ideal for great-looking, rich and high-quality specimen view. Our immersion darkfield system provides the same result achieved by 150W external illuminators in combination with traditional cardioid darkfield condenser.



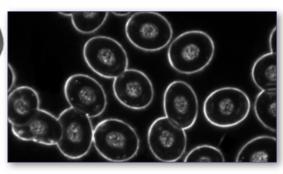








Part	Description
Observation mode:	Brightfield, oil immersion darkfield.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS W-PLAN 4x/0.10 IOS W-PLAN 10x/0.25 IOS W-PLAN 40x/0.65 IOS W-PLAN 100x/1.25 (oil) with iris All with anti-fungus treatment.

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Brightfield condenser:	Darkfield N.A. 1.36 (oil immersion) with built-in X-LED <sup>3</sup> .
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510FL** - HBO Fluorescence Microscope

Advanced routine laboratory microscope for brightfield and fluorescence observations with Semi-Apo IOS W-PLAN F objectives to enhance the visibility of the sample and increase the overall contrast. The **HBO fluorescence** illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Auramine, FITC, GFP and YFP (with blue filter set) plus Rhodamine, Texas Red and TRITC (with the green one), yet giving the possibility to combine any other specific filter sets for future upgrade. Transmitted light through the exclusive **X-LED<sup>3</sup>** to ensure great-looking, rich and high-quality specimen view.

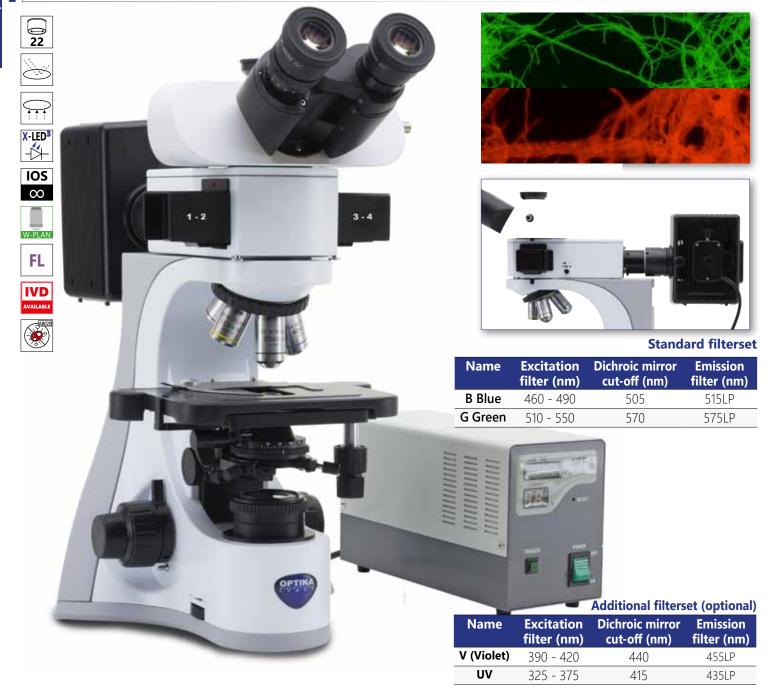

Part	Description	
<b>Observation mode:</b>	Brightfield, HBO fluorescence.	
<b>Epi-illumination and filter:</b>	HBO 100 W high pressure mercury lamp. 4-position filter holder; blue & green included.	
Head:	Trinocular (3-position 100/0, 50/50, 0/100), 30° inclined, 360° rotating.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Dioptric adjustment:	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS W-PLAN F 4x/0.13 IOS W-PLAN F 10x/0.30 IOS W-PLAN F 20x/0.50 IOS W-PLAN F 40x/0.75 All with anti-fungus treatment.	

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510LD4** - LED Fluorescence Microscope

Advanced routine fluorescence microscope for transmitted brightfield and fluorescence observations with IOS W-PLAN objectives. The extremely powerful LED Fluorescence Illuminators are combined with corresponding excitation filter sets for the visualization of most fluorochromes. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED<sup>3</sup>** to ensure great-looking, rich and high-quality specimen view.












### A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

Part	Description	
<b>Observation mode:</b>	Brightfield, LED fluorescence.	
<b>Epi-illumination and</b>	High-power LED with brightness control.	
filter:	4-position filter holder; none included.	
Head:	Trinocular (3-position 100/0, 50/50, 0/100), 30° inclined,	
	360° rotating.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Dioptric adjustment:	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Object.: B-510LD4:	IOS W-PLAN 4x/0.10 IOS W-PLAN 10x/0.25	
B-510LD4-SA:	IOS W-PLAN 40x/0.65 IOS W-PLAN 100x/1.25 (Oil) IOS W-PLAN F 4x/0.13 IOS W-PLAN F 10x/0.30 IOS W-PLAN F 40x/0.75 IOS W-PLAN F 100x/1.3 (Oil) All with anti-fungus treatment.	

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510LD4** - LED Fluorescence Microscope

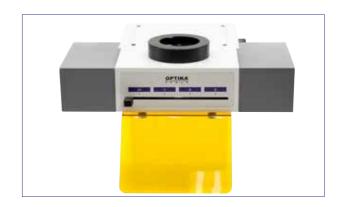

OPTIKA LED Fluorescence attachment is a revolutionary solution

It consists of a 4-position selector for the use of 4 fluorescent illuminators, called LED Fluorescence Cubes.

Each Cube is composed of a filterset mounted on a filterblock and a high power LED with emission corresponding to the filters installed. In this way the selection of each filter controls the lighting up of the corresponding LED.

The microscope is supplied without any LED Fluorescence Cube. A selection of 9 types is available, as shown in the table below.

**B-510LD4**Trinocular brightfield microscope, 1000x, IOS W-PLAN objectives



B-510LD4-SA
Trinocular brightfield microscope, 1000x,
IOS W-PLAN F (Semi-Apo) objectives. For UV Fluorescence



M-1220 ... M-1228 LED Fluorescence Cubes (LED + Filterset)


























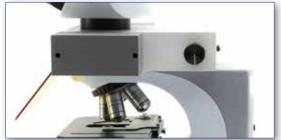



#### **LED Fluorescence Cubes available (LED + Filterset)**


Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1220 - Blue	460	455 - 495	500	510LP
M-1220.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1221 - Green	523	510 - 550	570	575LP
M-1221.1 - Green (pass band)	523	510 - 550	570	585-625
M-1222 - Violet	405	390 - 420	440	450LP
M-1223 - UV	365	325 - 375	415	435LP
M-1223.1 - UV (pass band)	365	340 - 390	405	420-470
M-1224 - Red 1 *	623	590 - 650	660	665LP
M-1225 - Red 2 *	623	595 - 645	655	665 - 715
M-1226 - Deep Red *	660	623 - 678	685	690 - 750
M-1227 - Far Red *	740	720 - 760	770	780LP
M-1228 - Amber	590	582 - 603	610	615 - 645


<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.


# **B-510POL** - Polarizing Microscope

Advanced routine laboratory microscope for transmitted light in brightfield and polarized light observations with strain-free IOS W-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED<sup>3</sup>** illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.



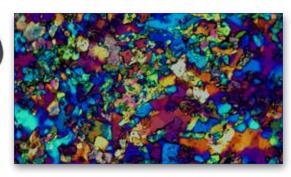

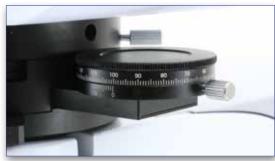

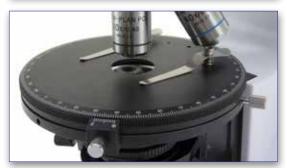

Part	Description
<b>Observation mode:</b>	Brightfield, transmitted polarized light and conoscopy.
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale. Tint plates: 1° order red (\(\lambda\); \(\lambda/4\); Quartz wedge.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups. One with crosshair.
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings. Centering system for each objective.

Part	Description	
Objectives (strain-free):	IOS W-PLAN POL 4x/0.10 IOS W-PLAN POL 10x/0.25 IOS W-PLAN POL 20x/0.45 IOS W-PLAN POL 40x/0.65 All with anti-fungus treatment.	
Specimen stage:	Rotatable stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm. Specimen slide clamps.	
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.	
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable. With rotating polarizing filter.	
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.	

### **B-510POL-I** - Polarizing Microscope

Advanced routine laboratory microscope for brightfield and polarized light observations in transmitted and incident light with strain-free IOS LWD W-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED**<sup>3</sup> illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.






















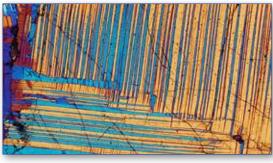





#### **Incident/transmitted light Objectives included**

### **Description**

IOS LWD W-PLAN POL 5x/0.12, W.D. 15.5 mm
IOS LWD W-PLAN POL 10x/0.25, W.D. 10.0 mm
IOS LWD W-PLAN POL 20x/0.40, W.D. 5.8 mm

IOS LWD W-PLAN POL 50x/0.75, W.D. 0.32 mm


Part	Description
<b>Observation mode:</b>	Brightfield, transmitted/incident polarized light and conoscopy.
Epi-illumination and filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With polarizer and rotating analyzer for incident illumination, aperture and field diaphragm. With additional filter holder.
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale. Tint plates: 1° order red (λ); λ/4; Quartz wedge.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Dioptric adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups. One with crosshair.
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings. Centering system for each objective.

Part	Description
Objectives (strain-free):	IOS LWD W-PLAN POL 5x/0.12 IOS LWD W-PLAN POL 10x/0.25 IOS LWD W-PLAN POL 20x/0.40 IOS LWD W-PLAN POL 50x/0.75 All with anti-fungus treatment.
Specimen stage:	Rotatable stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm. Specimen slide clamps.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable. With rotating polarizing filter.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.


## **B-510MET** - Metallurgical Microscope

Advanced routine laboratory microscope with IOS W-PLAN MET objectives and metallurgical attachment with the exclusive **X-LED³** lighting source for incident illumination only. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.



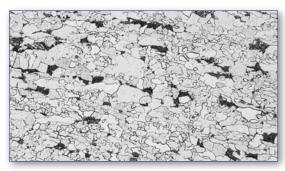



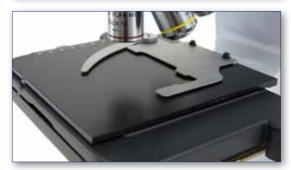




Part	Description
Observation mode:	Brightfield, simple polarized light, oblique illumination on incident light.
Epi-illumination and polarizing filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With aperture and field diaphragms, and oblique illumination system. With polarizer and analyzer. Multi-plug 100-240Vac/6Vdc external power supply.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.


Part	Description
Objectives (strain-free):	IOS W-PLAN MET 5x/0.12 IOS W-PLAN MET 10x/0.25 IOS W-PLAN MET 20x/0.40 IOS W-PLAN MET 50x/0.75 All with anti-fungus treatment.
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.


## **B-510METR** - Metallurgical Microscope

Advanced routine laboratory microscope with IOS W-PLAN MET objectives and metallurgical attachment with the exclusive **X-LED³** lighting source for both transmitted and incident illumination. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.



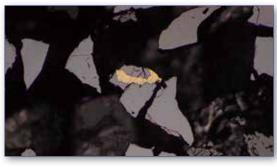

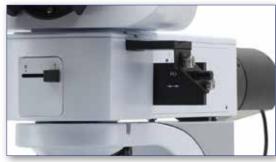

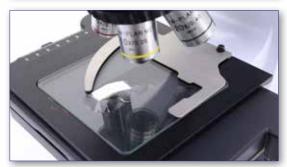



Part	Description
Observation mode:	Brightfield on transmitted light. Brightfield, simple polarized light, oblique illumination on incident light.
Epi-illumination and polarizing filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With aperture and field diaphragms, and oblique illumination system. With polarizer and analyzer.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.


Part	Description
Objectives (strain-free):	IOS W-PLAN MET 5x/0.12 IOS W-PLAN MET 10x/0.25 IOS W-PLAN MET 20x/0.40 IOS W-PLAN MET 50x/0.75 All with anti-fungus treatment.
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range. With tempered glass plate.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

### **B-510** - Discussion Microscopes

Advanced routine laboratory microscope for brightfield observations with IOS W-PLAN objectives and rackless stage. Ideal for discussion groups and teaching purpose for multiple observers, up to five users simultaneously. A three-color LED pointer facilitates the indication and identification of the object observed.

The high-efficiency **X-LED<sup>3</sup>** makes it reliable for all transmitted light observations for great-looking, rich and high-quality view.

B-510-2F

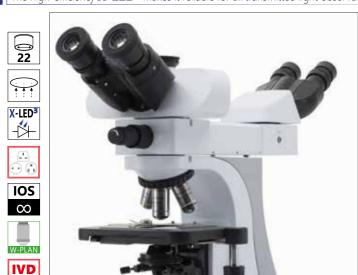







Part	Description	
<b>Observation mode:</b>	Brightfield	
Head:	Trinocular (fixed photo port 50/50), 30° inclined, 360° rotating.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Dioptric adjustment:	On the left eyepiece tube.	
Eyepieces:	Main head: WF10x/22 mm, high eye-point and with rubber cups. Additional head(s): WF10x/20 mm, high eye-point.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS W-PLAN 4x/0.10 IOS W-PLAN 10x/0.25 IOS W-PLAN 40x/0.65 IOS W-PLAN 100x/1.25 (Oil) All with anti-fungus treatment.	

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

# B-510 Series - Comparison chart

Model	Head	Eyepieces	Nosepiece	Attachment	Objectives
B-510BF	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	-	IOS W-PLAN 4x, 10x, 40x, 100x (oil)
B-510ERGO	Binocular ergonomical, 30°-60° inclined 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	-	IOS W-PLAN 4x, 10x, 40x, 100x (oil)
B-510PH	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	-	IOS W-PLAN PH 10x, 20x, 40x, 100x (oil)
B-510ASB	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point WF12.5x/15mm w/W&B reticle	Quintuple, reversed	-	IOS W-PLAN 4x, 10x, 40xPH, 100x (oil)
B-510DK	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	-	IOS W-PLAN 4x, 10x, 40x, 100x with iris diaphragm (oil)
B-510FL	Trinocular (100/0, 50/50, 0/100), 30° in- clined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Epi-Florescence attachment, with 4-position filterset slider. Equipped with Blue (FITC) and Green (TRITC) filtersets	IOS W-PLAN F 4x, 10x, 20x, 40x (oil)
B-510LD4	Trinocular (3-position 100/0, 50/50, 0/100), 30° inclined, 360° rotating.	WF10x/22mm, high eye-point	Quintuple, reversed	Epi-Florescence attachment, with 4-position filterset slider.	<b>B-510LD4:</b> IOS W-PLAN 4x, 10x, 40x, 100x (oil) <b>B-510LD4-SA:</b> IOS W-PLAN F 4x, 10x, 40x, 100x (oil)
B-510POL	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quadruple, reversed. Objective positions centrable.	Swing-out Bertrand lens with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale. Tint plates: 1° order red ( $\lambda$ ); $\lambda/4$ ; Quartz wedge	IOS W-PLAN POL 4x, 10x, 20x, 40x
B-510POL-I	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quadruple, reversed. Objective positions centrable.	Incident light attachment with Polarizer for incident illumination, with Aperture & Field diaphragms and additional filter holder. Swing-out Bertrand lens with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale. Tint plates: 1° order red (λ); λ/4; Quartz wedge	IOS LWD W-PLAN POL 5x, 10x, 20x, 50x
B-510MET	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Incident light attachment with Polarizer and rotating Analyzer for incident illumination, with Aperture & Field diaphragms and 2 additional filter holders. Oblique illumination	IOS W-PLAN MET 5x, 10x, 20x, 50x
B-510METR	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Incident light attachment with Polarizer and rotating Analyzer for incident illumination, with Aperture & Field diaphragms and 2 additional filter holders; Epi/Transmitted light selector. Oblique illumination	IOS W-PLAN MET 5x, 10x, 20x, 50x
B-510-2F	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Discussion bridge for 2 observers , Face-to-Face type. 2nd binocular head with WF10x/20mm eyepieces. Equipped with RGB discussion pointer	IOS W-PLAN 4x, 10x, 40x, 100x (oil)
B-510-2	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Discussion bridge for 2 observers , Side-by-Side type. 2nd binocular head with WF10x/20mm eyepieces. Equipped with RGB discussion pointer	IOS W-PLAN 4x, 10x, 40x, 100x (oil)
B-510-3	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Discussion bridge for 3 observers. Additional binocular heads with WF10x/20mm eyepieces. Equipped with RGB discussion pointer	IOS W-PLAN 4x, 10x, 40x, 100x (oil)
B-510-5	Trinocular (50/50), 30° inclined, 360° rotating	WF10x/22mm, high eye-point	Quintuple, reversed	Discussion bridge for 5 observers. Additional binocular heads with WF10x/20mm eyepieces. Equipped with RGB discussion pointer	IOS W-PLAN 4x, 10x, 40x, 100x (oil)

# **B-510** Series - Comparison chart

Stage	Focusing	Condenser	Incident Illumination	Transmitted Illumination
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Phase condenser (10x/20x, 40x, 100x) with darkfield (dry) and brightfield	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Abbe N.A. 1.25, with objective- coded iris diaphragm, focusable and centerable. With 40x phase contrast slider.	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system. Additional darkfield condenser, N.A. 1.36, built-in X-LED <sup>3</sup>	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	High Pressure HBO 100 W mercury bulb	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	High-power blue LED with brightness control	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Rotating stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm; Specimen slide clamps	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system. With rotating polarizing filter	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Rotating stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm; Specimen slide clamps	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system. With rotating polarizing filter	X-LED <sup>8</sup> with white 8 W LED (6,300K), brightness control	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-a- xis rackless; with metal plate for metallurgical samples	Coaxial coarse and fine, limit stop, adjustable tension	-	X-LED <sup>8</sup> with white 8 W LED (6,300K), brightness control	-
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless; with glass plate for metallurgical samples	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	X-LED <sup>8</sup> with white 8 W LED (6,300K), brightness control	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.
Double layer, 233x147 mm, moving range 78x54 mm, X-axis rackless	Coaxial coarse and fine, limit stop, adjustable tension	Swing-out type, sliding-in, N.A. 0.2/0.9, with centering system	-	X-LED <sup>3</sup> with white 3.6 W LED (6,300K), brightness control. Full Koehler type.

### **B-510** Series - Accessories

Eyecups &	Eyepieces
M-601	WF15x/16 eyepiece, high eyepoint
M-690	Eyecups (pair)
M-780	PL10x/22 eyepiece, high eyepoint, with rubber cup
M-781	PL10x/22 micrometric eyepiece, high eyepoint, with rubber cup
AA-01	12.5x/18 eyepieces (pair), high eyepoint, focusable, W&B reticle
Objectives	
IOS W-PLA	
M-1049	IOS W-PLAN objective 2x/0.05
M-1125	IOS W-PLAN objective 4x/0.10
M-1126	IOS W-PLAN objective 10x/0.25
M-1127	IOS W-PLAN objective 20x/0.40
<u>M-1128</u>	IOS W-PLAN objective 40x/0.65
M-634.1	IOS W-PLAN objective 50x/0.95 (oil)
M-1129	IOS W-PLAN objective 60x/0.80
M-1130	IOS W-PLAN objective 100x/1.25 (oil)
M-1130.1	IOS W-PLAN objective 100x/1.25OI - (oil) with iris for DF
IOS W-PLA	
M-1060	IOS W-PLAN F objective 4x/0.13
<u>M-1061</u>	IOS W-PLAN F objective 10x/0.30
<u>M-1062</u>	IOS W-PLAN F objective 20x/0.50
M-1063	IOS W-PLAN F objective 40x/0.75
<u>M-1064</u>	IOS W-PLAN F objective 100x/1.30 (oil)
IOS W-PLA	
M-336	IOS W-PLAN MET objective 5x/0.12
M-338	IOS W-PLAN MET objective 10x/0.25
M-339	IOS W-PLAN MET objective 20x/0.40
M-335	IOS W-PLAN MET objective 50x/0.75
	IOS W-PLAN MET objective 100x/0.80 (dry)
IOS W-PLA	
	IOS W-PLAN PH objective 10x/0.25
	IOS W-PLAN PH objective 20x/0.40
	IOS W-PLAN PH objective 40x/0.65
	IOS W-PLAN PH objective 100x/1.25 (oil)
IOS W-PLA M-1131	
	IOS W-PLAN POL objective 4x/0.10
M-1132 M-1133	IOS W-PLAN POL objective 10x/0.25 IOS W-PLAN POL objective 20x/0.45
M-1134	
M-1134 M-1135	IOS W-PLAN POL objective 40x/0.65 IOS W-PLAN POL objective 60x/0.80
M-1136	IOS LWD W-PLAN POL objective 60x/0.80
M-1137	IOS LWD W-PLAN POL objective 3x/0.12
M-1138	IOS LWD W-PLAN POL objective 10x/0.25
M-1139	IOS LWD W-PLAN POL objective 20x/0.40
141-1133	103 EVVD VV-1 ENIV I OF ODJECTIVE 30X/0.13





### **B-510** Series - Accessories

**Stages** M-175.1 Rotating stage for polarising set (for 233x147mm rackless stage) (except for B-510PH, B-510ASB and B-510DK) M-635 Heating stage (on newly purchased microscopes, for 233x147mm), multiplug **Condensers & Filters** M-181 PH set - 10x, 20x, 40x, 100x IOS W-PLAN PH obj. & BF/DF/PH condenser Darkfield condenser for dry objectives (except for B-510DK, B-510MET and B-510METR) M-185 M-636 Polarising set (filters only) (except for B-510POL & B-510POL-I) Fluorescence filter set V (filterblock not needed) (only for B-510FL) M-637 M-638 Fluorescence filter set UV-DAPI (filterblock not needed) (only for B-510FL) M-975 Blue filter, 45mm diameter M-977 Green filter, 45mm diameter M-979 Yellow filter, 45mm diameter M-989 Frosted glass filter, 45mm diameter M-1220 Blue LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) Blue (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) M-1220.1 M-1221 Green LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) M-1221.1 Green (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) M-1222 Violet LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) UV LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) M-1223 M-1223.1 UV (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) M-1224 Red 1 LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) 7 Red 2 LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) \* M-1225 Deep Red LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) \* M-1226 M-1227 Far Red LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) \* M-1228 Amber LED Fluorescence Cube (LED+Filterset), (only for B-510LD4) **Camera Adapters** M-113.1 Ring adapter, 30mm (for monocular and binocular microscopes) M-115 0.35x C-Mount projection lens 0.5x C-Mount projection lens M-114 M-118 0.75x C-Mount projection lens C-Mount projection lens for APS-C/full frame reflex cameras (trino) M-173 M-620 0.35x focusable C-Mount adapter (biological microscopes) M-620.1 0.5x focusable C-Mount adapter (biological microscopes) 0.65x focusable C-Mount adapter (biological microscopes) M-620.2 1x focusable C-Mount adapter (biological & stereomicroscopes) M-620.3 M-699 Universal adapter for C-Mount projection lens (trino) Miscellaneous 15008 Immersion oil, 10ml 15009 Immersion oil, 100ml 15104 Cleaning kit DC-003 TNT dust cover, medium, 600(l)x550(h) mm M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100) M-069 Solar charger (except for B-510MET, B-510METR and B-510POL-I) M-151 HBO 100W high-pressure mercury bulb for fluorescence (only for B-510FL) HBO 100W high-pressure mercury bulb for fluorescence (OSRAM) (only for B-510FL) M-151.1

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order it by specifying with "AR GLASS"



M-1004.N

M-1037 VP-510

AB-020

#### How to connect the cameras to our microscopes.

IQ/OQ/PQ manual for B-510 series

Please refer to the Adapter reference list on Digital section.

Gout analysis kit (only for B-510BF, B-510ERGO and B-510PH)

Antibacterial surface treatment, only for newly purchased microscope

#### 15104 - Cleaning kit

It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.



Centering telescope, 30mm diameter (except for B-510MET, B-510METR, B-510POL and B-510POL-I)

#### M-069 - Solar charger

Included battery: rechargeable – Lithium-Poly. Capacity: 2500 mAh. - Output voltage: 5 Vdc. Autonomy: over 6 hours at medium intensity (X-LED³). Charging models: with solar panel (12h), with external USB power supply (2.5h)



v7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# B-810 & B-1000



**Research Lab Upright Microscopes** 



OPTIKA Microscopes, thanks to the long experience achieved in microscopy development, has conceived the new B-1000: a major leap in our technological offer. As a flagship instrument, B-1000 originates from customer most demanding feedbacks and needs. Its modularity and versatility will allow to find the perfect place in any clinical or basic reasearch laboratory. All controls are easily accessible and comfortable also for extended periods of observation.

B-1000 is built on IOS Infinity Corrected optical system, which gives both top-notch optical performances, and the possibility to extend your instrument with the broad range of accessories and modules. X-LED illumination is the best solution to have pure white light, very intense even at higher magnification, and optimum power efficiency given by solid state source.

If you are a looking for our best solution to your present and future professional demands, B-1000 is the answer.





Highest category of optical equipment among our product range guarantees a sharp and clear view in any situation, while top level mechanical design offers sturdiness and long lifetime.





#### **Solid Stand – Extra Stability**

A completely new design and a die-cast aluminum stand offer solidity and durability, even for the most demanding laboratory use.

This new microscope can seamlessly be upgraded with many attachments that extend its field of use.

#### **X-LED White Illumination**

X-LED illumination system is based on a pure white high-efficiency LED and a special optics. It guarantees constant color temperature, no heat, and an extreme electrical consumption efficiency. The whole system is pre-aligned and boasts a lifetime of 50,000 hours.

#### **X-LED** benefits

Powerful pure white LED illumination, ideal for brightfield, darkfield and phase contrast applications. Constant color temperature through all the intensity levels. No heat generation, avoiding damage of the specimen.

Factory pre-centering assures uniform illumination over the field of view, yet providing perfect Koehler alignment. Very long lifetime and high power efficiency.

#### **Light under control**

Intelligent control of the microscope illumination: the "AUTO-OFF" function automatically switches the light off after a user-selectable time period. "BOOST" gives an extra high level of illumination for light-demanding applications. "AUTO" allows to store an illumination level, and to maintain it throughout the inspection.





Low position focus and stage controls allow a fast and comfortable operation. Frequently used controls as light intensity adjustment and diaphragm are also placed in the lower part of the stand and enable operation without having to take the eyes off the specimen. All optical heads are equipped with high-point eyepieces and dioptric adjustment, for the best viewing experience.



#### **Modularity – Build your own solution**

Many worlds in one instrument. Modularity allows to build the desired solution (brightfield, darkfield, phase contrast, material science, fluorescence, motorized automation and so on). B-1000 has the flexibility to help your work the best way.

#### **Comfortable Stage**

Refined belt-driven stage, with a wide working surface and a highly precise XY movement.

### **High Quality IOS Optical System**

Infinity corrected optical system, based on planachromatic, fluorite and semi-apochromatic objectives, designed to give sharp and clear images, both for the user and the digital camera. Quintuple and sextuple nosepieces give the flexibility to build the objectives that best suits your needs. The system is completed with wide field, high-point eyepieces, with a field number of 24mm.

#### **Ready for Digital Imaging**

Range of adapters can accommodate for C-mount digital cameras, as well as reflex cameras. Focus adjustment gives perfectly clear digital images. Our cameras include specific software for capturing, measuring, marking and storing your images. Pro View software allows to perform image acquisition, post-processing, measurements and storage of your images. User can save a preset for later work, or even create a multi-focus composition.





### Remote-controlled microscope


The stage can be remote-controlled through a dedicated software: X, Y, Z axes, as well as nosepiece, can be moved with a single click. Communication protocol is available for interfacing with custom software, such as automated analysis or autofocus.

### X-Y-Z motorized stage Motorized nosepiece






## B-810/B-1000 - Observation Methods

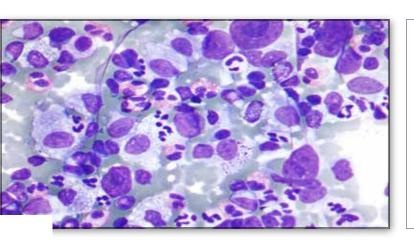

#### Pathology / Cytology


Since B-810 / B-1000 use white LED illumination, they can maintain the same color temperature even if the brightness is changed. "AUTO" function automatically adjusts the light intensity when the objective is changed or the aperture diaphragm is set to a different value.

These feautures, along with motorized stage and ergonomic controls, make your workflow easier.

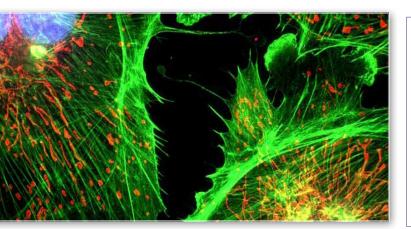

#### **Fluorescence Microscopy**

A new attachment for epi-fluorescence provides the ultimate solution in the field of fluorescence diagnostic. Vibrationfree six positions filter wheel with shutter, field and aperture diaphragms offer all you need for a complete analysis. Custom filtersets are available and mounted on request. For application where efficiency, rapidity and ease of use are crucial, this model offers also a LED epi-fluorescence attachment, with very high power standard illuminators.

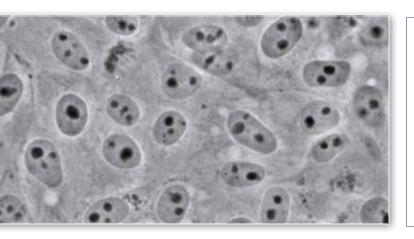

#### **Phase Contrast Microscopy**


The bright LED illuminator brings a comfortable view in phase contrast with all magnifications. Universal wheel condenser allows to quickly switch between brightfield, darkfield and phase contrast.

Ideal for clinical laboratories or fibers (e.g. asbestos) analysis.

## **B-810/B-1000** - Observation Methods


#### **Darkfield Microscopy**

Ideal for observing blood cells, diatoms, small insects, bone, fibers, unstained bacteria, yeast, protozoa, mineral and chemical crystals, colloidal particles, dust-count specimens, and thin sections of polymers and ceramics.

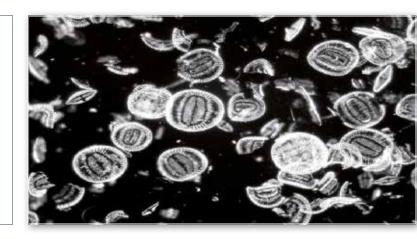

#### **Material Science**


A new attachment designed specifically for metallographic inspection, with dedicated objectives set, for the most complete epi-illumination analysis: brightfield, darkfield and polarizing view observations.

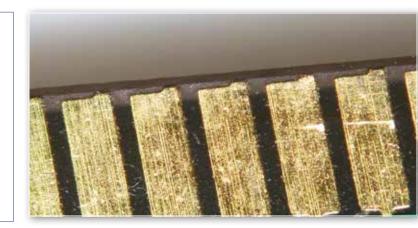


#### **Polarizing Microscopy**

Polarized light microscopy is used in geological applications or also for both natural and industrial minerals, composites such as concretes, ceramics, mineral fibers and polymers, and crystalline or biological molecules such as DNA, starch, wood and urea.

Attachments for a full polarization analysis are available (both for transmitted and incident light), so it's possible to look at color fringes right away.

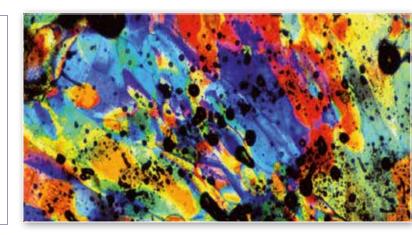

# B-810/B-1000 - Design and Production

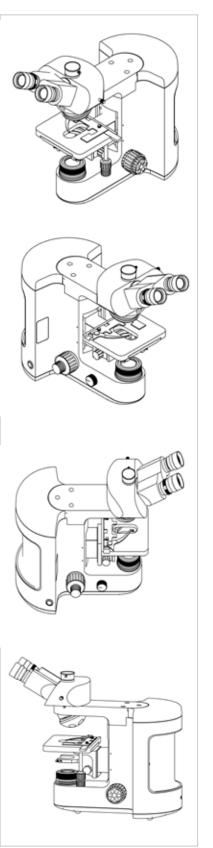



Laboratory

## B-810/B-1000 - Design and Production

**OPTIKA** workshop provides the facilities for precise and reliable optomechanical manufacturing, essential for this kind of instruments.

CNC machining department, equipped with 5-axis milling machine and lathe.



Work in progress inside the milling machine.



Die-cast stands ready to be processed.



Microscope stands exiting from the internal varnishing facility.





All processes are carefully monitored through the application of ISO 9001 Quality System standards.

### **Eyepieces**



#### M-1001

PL10x/22 eyepieces (pair), high eyepoint, with rubber cup (retractable)



#### M-1002

PL10x/24 eyepieces (pair), high eyepoint, with rubber cup (retractable); for B-1000 only



#### M-1003

PL15x/16 eyepieces (pair), high eyepoint



#### M-781

PL10x/22 eyepiece, high eyepoint, with micrometric scale (10mm/100um) & rubber cup (retractable)



#### **AA-01**

12.5x/15 eyepieces (pair), high eyepoint, with dioptric adjustment, one with W&B graticule for Asbestos



#### M-1004.N

Centering telescope



### Heads













### **Bodies**



#### M-1187

Main body with focus system and X-LED<sup>3</sup> illumination for B-810.



#### M-1021B

Main body with focus system and X-LED<sup>8</sup> illumination, for general purposes for B-1000.



#### M-1021M

Main body with focus system and X-LED8 illumination, for brightfield metallurgical B-1000.



#### M-1022M

Main body with focus system (incident light only) for brightfield metallurgical B-1000.



#### M-1021MD

Main body with focus system and X-LED8 illumination, for brightfield/darkfield metallurgical B-1000B-1000.



#### M-1022MD

Main body with focus system (incident light only), for brightfield/darkfield metallurgical B-1000.



#### M-1156

Upgrade (controller) for any kind of motorization (stage, Z-axis, nosepiece, or all of them together); for B-1000

#### M-1149

Motorization of Z-axis for B-1000

### **Nosepieces**



#### M-1040

Quintuple reversed nosepiece, for RMS objectives; for B-810



#### M-1041

Sextuple reversed nosepiece, for RMS objectives; for B-810



#### M-1042

Sextuple reversed nosepiece, for RMS objectives with DIC slot; for B-1000



#### M-1043

Sextuple motorized reversed nosepiece, for RMS objectives with DIC slot; for B-1000



#### M-1044

Quintuple reversed nosepiece, with centrable positions for POL objectives; for B-1000



### **W-PLAN Objectives**



#### Plan Objectives - IOS W-PLAN Series



M-1049 IOS W-PLAN 2x/0.05

**M-1125** IOS W-PLAN 4x/0.10

M-1126 IOS W-PLAN 10x/0.25

M-1127 IOS W-PLAN 20x/0.40

M-1128 IOS W-PLAN 40x/0.65

M-1129 IOS W-PLAN 60x/0.80

M-1130 IOS W-PLAN objective 100x/1.25 (oil)

**M-1130.1** IOS W-PLAN objective 100x/1.25 (Iris)

#### Plan Objectives - IOS W-PLAN F PH Series



M-1325 IOS W-PLAN F PH 10x/0.25

**M-1326** IOS W-PLAN F PH 20x/0.40

M-1327 IOS W-PLAN F PH 40x/0.65

M-1328 IOS W-PLAN F PH 100x/1.25 (oil)

#### Plan Semi-APO Objectives - IOS W-PLAN F Series



M-1060 IOS W-PLAN F 4x/0.13

M-1061 IOS W-PLAN F 10x/0.30

M-1062 IOS W-PLAN F 20x/0.50

M-1063 IOS W-PLAN F 40x/0.75

M-1064 IOS W-PLAN F 100x/1.30 (oil)

#### Plan Objectives - IOS U-PLAN POL Series



**M-1080** IOS U-PLAN POL 4x/0.10 **M-1081** IOS U-PLAN POL 10x/0.25

M-1081.5 IOS U-PLAN POL 20x/0.45

**M-1082** IOS U-PLAN POL 40x/0.65 **M-1083** IOS U-PLAN POL 60x/0.85

#### Plan Objectives - IOS LWD W-PLAN MET Series



M-1099 IOS LWD W-PLAN MET 2.5x/0.08 (with depolarizer)

### **U-PLAN Objectives**



#### Plan Semi-APO Objectives - IOS U-PLAN F Series



M-1075 IOS U-PLAN F 4x/0.13

**M-1076** IOS U-PLAN F 10x/0.30

**M-1077** IOS U-PLAN F 20x/0.50

M-1078 IOS U-PLAN F 40x/0.75

M-1079 IOS U-PLAN F 100x/1.30 (oil)

#### Plan APO Objectives - IOS U-PLAN APO Series



M-1301 IOS U-PLAN APO 2x/0.08

**M-1302** IOS U-PLAN APO 4x/0.13

**M-1303** IOS U-PLAN APO 10x/0.40

**M-1304** IOS U-PLAN APO 20x/0.75

**M-1305** IOS U-PLAN APO 40x/0.95

**M-1306** IOS U-PLAN APO 60x/0.90 **M-1307** IOS U-PLAN APO 100x/1.35

#### Plan Semi-APO Objectives - IOS U-PLAN F PH Series



M-1310 IOS U-PLAN F PH 4x/0.13

M-1311 IOS U-PLAN F PH 10x/0.40

M-1312 IOS U-PLAN F PH 20x/0.75

M-1313 IOS U-PLAN F PH 40x/0.95

**M-1314** IOS U-PLAN F PH 60x/0.90

M-1315 IOS U-PLAN F PH 100x/1.35

M-1157 universal condenser required

#### Plan Objectives - IOS LWD U-PLAN POL Series



M-1090 IOS LWD U-PLAN POL 5x/0.15

M-1091 IOS LWD U-PLAN POL 10x/0.30

M-1092 IOS LWD U-PLAN POL 20x/0.45

M-1093 IOS LWD U-PLAN POL 50x/0.55

#### Plan Objectives - IOS LWD U-PLAN MET Series



M-1100 IOS LWD U-PLAN MET 5x/0.15

M-1101 IOS LWD U-PLAN MET 10x/0.30

M-1102 IOS LWD U-PLAN MET 20x/0.45

**M-1103** IOS LWD U-PLAN MET 50x/0.55

M-1104 IOS LWD U-PLAN MET 100x/0.80 (dry)

#### Plan Semi-APO Objectives - IOS LWD U-PLAN F MET Series



**M-1171** IOS LWD U-PLAN F MET 5x/0.15

**M-1172** IOS LWD U-PLAN F MET 10x/0.30

M-1173 IOS LWD U-PLAN F MET 20x/0.50

**M-1174** IOS LWD U-PLAN F MET 50x/0.80

M-1175 IOS LWD U-PLAN F MET 100x/0.90 (dry)

#### Plan Objectives - IOS LWD U-PLAN MET BD (Darkfield) Series



M-1094 IOS LWD U-PLAN MET BD 5x/0.15

M-1095 IOS LWD U-PLAN MET BD 10x/0.30

M-1096 IOS LWD U-PLAN MET BD 20x/0.45

M-1097 IOS LWD U-PLAN MET BD 50x/0.55

**M-1098** IOS LWD U-PLAN MET BD 100x/0.80 (dry)

#### Plan Semi-APO Objectives - IOS LWD U-PLAN F MET BD Series



**M-1180** IOS LWD U-PLAN F MET BD 5x/0.15

M-1181 IOS LWD U-PLAN F MET BD 10x/0.30

M-1182 IOS LWD U-PLAN F MET BD 20x/0.50

**M-1183** IOS LWD U-PLAN F MET BD 50x/0.80

M-1184 IOS LWD U-PLAN F MET BD 100x/0.90 (dry)

### **Stages**



#### M-1140

Standard mechanical stage; 175x145mm for B-1000



#### M-1141

Rackless mechanical stage; movement knobs with friction adjustment control; 242x157mm; for B-1000.



#### M-1143

MPC (Mineral Solid Surface) rackless mechanical stage; movement knobs with friction adjustment control; 242x157mm; for B-1000.



#### M-1148

Metallurgical stage with glass, for metallurgical model; 175x145mm; for B-1000



#### M-1190

Rackless mechanical stage; 220x149mm; for B-810



## **Stages**



#### M-1145 + M-1146

Rotating Stage, centrable + attachable XY stage; dia. 172mm; for B-1000



#### M-1147

Motorized mechanical stage; 242x157mm; for B-1000



#### M-1144

Heating stage, with digital temperature controller; 175x145mm; for B-1000.



#### M-1190H

Heating stage, with digital temperature controller; 220x149mm; for B-810.



### **Condensers**



#### M-1189

0.90 N.A. swing-out condenser for B-810



#### M-1150

0.90 N.A. swing-out condenser for B-1000



1.2 N.A. swing-out condenser for B-810



#### M-1151

1.2 N.A. swing-out condenser for B-1000



#### M-1155

0.9/0.25 N.A. swing-out condenser (to be used with objective M-1049) for B-1000



#### M-1154

0.70 N.A. swing-out MET condenser for B-1000



### **Condensers**

#### M-1157 8-Position universal condenser for B-1000



#### Parts for the universal condenser M-1157

M-1205 Top lens 0.2 N.A.
M-1206 Top lens 0.9 N.A.
M-1207 Top lens 1.4 N.A.
M-1208 DIC 10x prism
M-1209 DIC 20x prism
M-1210 DIC 40x/60x prism
M-1212 10x/20x phase ring
M-1213 40x/60x phase ring
M-1214 100x phase ring
M-1215 4x phase ring
M-1216 Darkfield stop (dry)
M-1217 Darkfield stop (oil)

M-1153

0.90 N.A. swing-out POL condenser for B-1000



M-618

Darkfield condenser for dry objectives for B-1000



#### M-185

Darkfield condenser for dry objectives for B-810



#### M-1253.NO

Phase contrast condenser with insert slide 10x/20x-40x for B-810



#### M-1253

Phase contrast condenser with insert slide 10x/20x-40x for B-1000



#### M-1255.NO

Phase contrast condenser 10x, 20x, 40x, 100x, BF, DF for B-810




#### M-1255

Phase contrast condenser 10x, 20x, 40x, 100x, BF, DF for B-1000




### **Attachments**

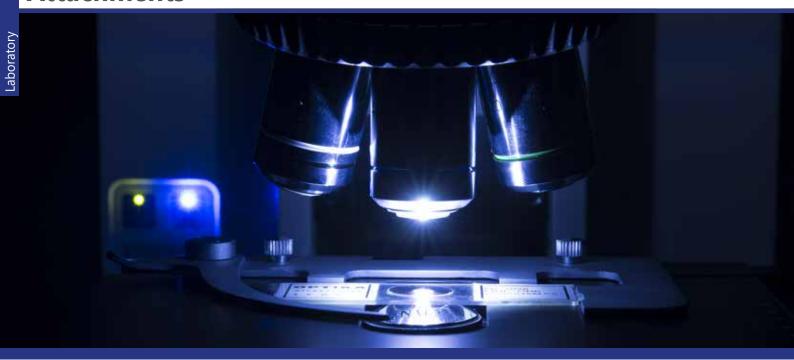

### **Fluorescence attachments**

#### M-1031M

4-position LED Fluorescence attachment. Without LED-Filtersets (see B-1000 4FLD page for choice); for B-1000

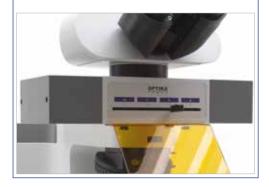

#### M-1032

6-position HBO Fluorescence attachment, with standard Blue and Green filtersets (FITC & TRITC). With Aperture & Field Diaphragms; for B-1000.



### **Polarizing attachments**

#### M-1033

Bertrand lens with analyzer and slot for sliders (with Lambda, 1/4 Lambda and Quartz Edge); for B-1000



#### M-1034


8W X-LED incident polarizing attachment. With Aperture & Field Diaphragms; for B-1000



### **ALC** attachments

#### M-1030

Automatic Light Control (ALC) system for B-1000

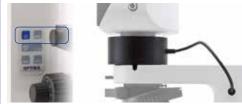

### **Attachments**



### **Metallurgical attachments**

Metallurgical attachment for incident brightfield observation, 18 W LED.

Equipped with Polarizer and rotating Analyzer.

With Aperture & Field Diaphragms.

For B-1000.



#### M-1039MD

Metallurgical attachment for incident brightfield and darkfield observation, 18 W LED. Equipped with Polarizer and rotating Analyzer; with Aperture & Field Diaphragms.

Built-in MET nosepiece with 6-positions for 26 mm-Darkfield Metallurgical Objectives.

For B-1000.



### **Discussion Bridges for B-1000**









M-1162

### M-1159

Discussion Bridge with 1 extra head. Face-to-Face type

#### M-1160

Discussion Bridge with 1 extra head. Side-by-Side type

#### M-1161

Discussion Bridge with 2 extra heads

#### M-1162

Discussion Bridge with 4 extra heads

Discussion Bridge with 9 extra heads

All Discussion Bridges are equipped with extra heads with WF10x/20mm eyepieces

### **Application Sets**

### **Koehler DIC transmitted**

M-550 - Interferential green filter IF550.

M-190P - Polarizer for transmitted light.

M-1201 - Analyzer for transmitted light.

M-1202 - DIC prism for transmitted light.



Usable objective series: IOS W-PLAN F, IOS U-PLAN F. (See page 164)

### Transmitted Koehler DIC combined with Fluo HBO

M-550 - Interferential green filter IF550.

M-190P - Polarizer for transmitted light.

M-1203 - Analyzer for reflected light.

M-1202 - DIC prism for transmitted light.



Usable objective series: IOS W-PLAN F, IOS U-PLAN F. (See page 164)

### Nomarski DIC transmitted

M-1157 - 8-Position universal condenser.

M-1206 - Top lens 0.9 N.A.

M-1208 - DIC 10x prism for universal condenser.

M-1209 - DIC 20x prism for universal condenser.

M-1210 - DIC 40x/60x prism for universal condenser.

M-1201 - Analyzer for transmitted light.

M-1202 - DIC prism for transmitted light.



Usable objective series: IOS U-PLAN F, IOS U-PLAN APO, IOS U-PLAN F PH. (See page 162)

### Transmitted Nomarski DIC combined with Fluo HBO

M-1157 - 8-Position universal condenser.

M-1206 - Top lens 0.9 N.A.

M-1208 - DIC 10x prism for universal condenser.

M-1209 - DIC 20x prism for universal condenser.

M-1210 - DIC 40x/60x prism for universal condenser.

M-1203 - Analyzer for reflected light.

M-1202 - DIC prism for transmitted light.



Usable objective series: IOS U-PLAN F, IOS U-PLAN APO, IOS U-PLAN F PH. (See page 162)

### Nomarski DIC reflected for metallurgical appl.

M-870 - DIC prism for metallurgical reflected light.



Usable objective series: IOS LWD U-PLAN F MET Series IOS LWD U-PLAN F MET BD Series.

(See page 162)

### **GOUT** analysis kit

M-1047 - GOUT analysis kit





IOS ∞

W-PLAN

DF

## B-810 - Brightfield & Phase Contrast Microscope

B-810 is the result of the long experience gathered by OPTIKA Microscopes in the field of light microscopy, offering an extremely valuable product for routine and research laboratory brightfield & phase contrast applications.

Ergonomic design for comfortable long-term use and smooth operation with minimal movements meets the concept of modularity, to offer a the possibility to create a tailored version and match all the personal requirements.

It is qualified as an particularly performing and robust solution, considering the field of view of 22 mm, the state-of-the-art, exclusive **X-LED³** lighting source (3.6 W) and the sturdy dye-cast frame for high stability combined with a wide variety of heads, objectives and condensers to get the most out of a microscope.



## **B-810** - Configuration Chart



## B-1000BF - Brightfield Microscope

The modular OPTIKA B-1000 helps you working in a comfortable way during extended periods of use and let you perform reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Semi-Apo objectives), the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W)

illumination system, designed by OPTIKA and the Koehler diaphragm. B-1000 gives multiple options as manual or motorized configuration, with the possibility of

**ALC** (Automatic Light Control) and a variety of objectives, stages and condensers.





## **B-1000BF** - Configuration Chart




## B-1000PH - Phase Contrast Microscope

The modular OPTIKA B-1000 is available in phase contrast and helps you working in a comfortable way during extended periods of use performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W)illumination system, designed by OPTIKA and the Koehler diaphragm. B-1000 gives multiple options as manual or motorized configuration, with the possibility of ALC (Automatic Light Control) and a variety of objectives, stages and condensers.



## **B-1000PH** - Configuration Chart

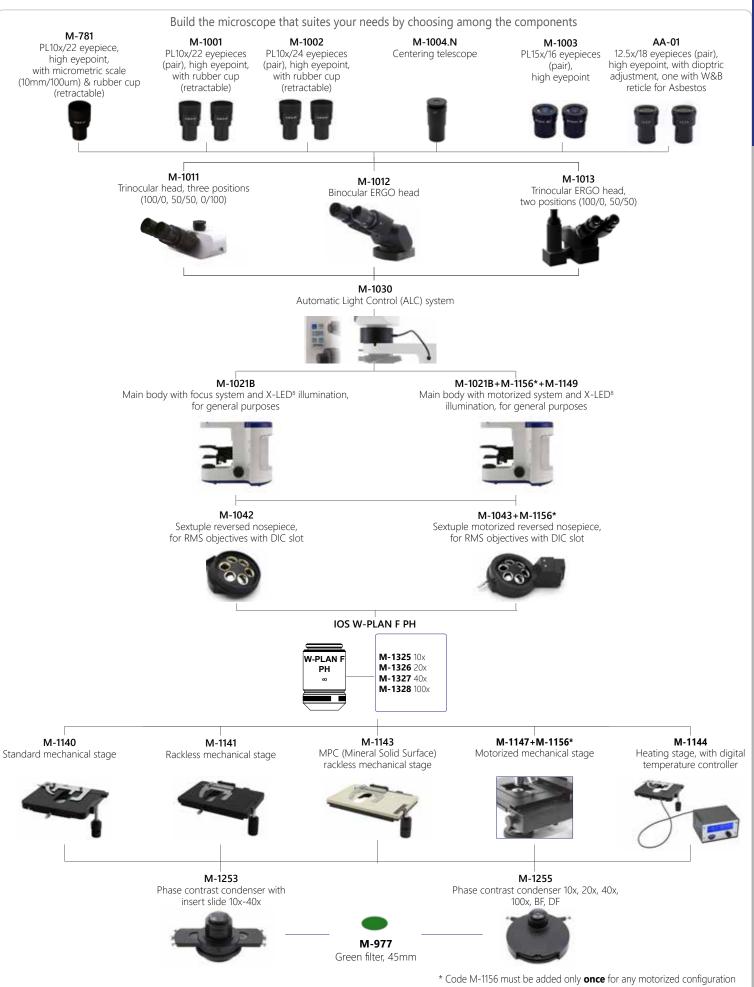

# B-1000FL-HBO - HBO Fluorescence Microscope

The modular OPTIKA B-1000 can stand a HBO fluorescence attachment, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Semi-Apo objectives), the state-of-the-art, exclusive X-LED<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration, with a variety of objectives, stages and condensers.

Standard	filterset		
Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
B Blue	460 - 490	500	520LP
G Green	510 - 550	570	590LP
	l filterset (opt		
Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
V (Violet)	400 - 410	455	455LP
22	330 - 385	400	420LP
X-LED <sup>®</sup> O O O			
W-PLAN U-PLAN FL IVD AVAILABLE			
OPTIKA			ОРТІКА

## **B-1000FL-HBO** - Configuration Chart



#### (2)

### **B-1000LD4** - LED Fluorescence Microscope

The modular OPTIKA B-1000 can stand a LED fluorescence attachment, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Apo and Semi-Apo objectives), the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration, with a variety of objectives, stages and condensers.



## A new milestone achieved in Fluorescence Microscopy


- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

#### **LED Fluorescence Cubes available (LED + Filterset)**

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1220 - Blue	460	455 - 495	500	510LP
M-1220.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1221 - Green	523	510 - 550	570	575LP
M-1221.1 - Green (pass band)	523	510 - 550	570	585-625
M-1222 - Violet	405	390 - 420	440	450LP
M-1223 - UV	365	325 - 375	415	435LP
M-1223.1 - UV (pass band)	365	340 - 390	405	420-470
M-1224 - Red 1 *	623	590 - 650	660	665LP
M-1225 - Red 2 *	623	595 - 645	655	665 - 715
M-1226 - Deep Red *	660	623 - 678	685	690 - 750
M-1227 - Far Red *	740	720 - 760	770	780LP
M-1228 - Amber	590	582 - 603	610	615 - 645

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.

### **B-1000LD4** - Configuration Chart

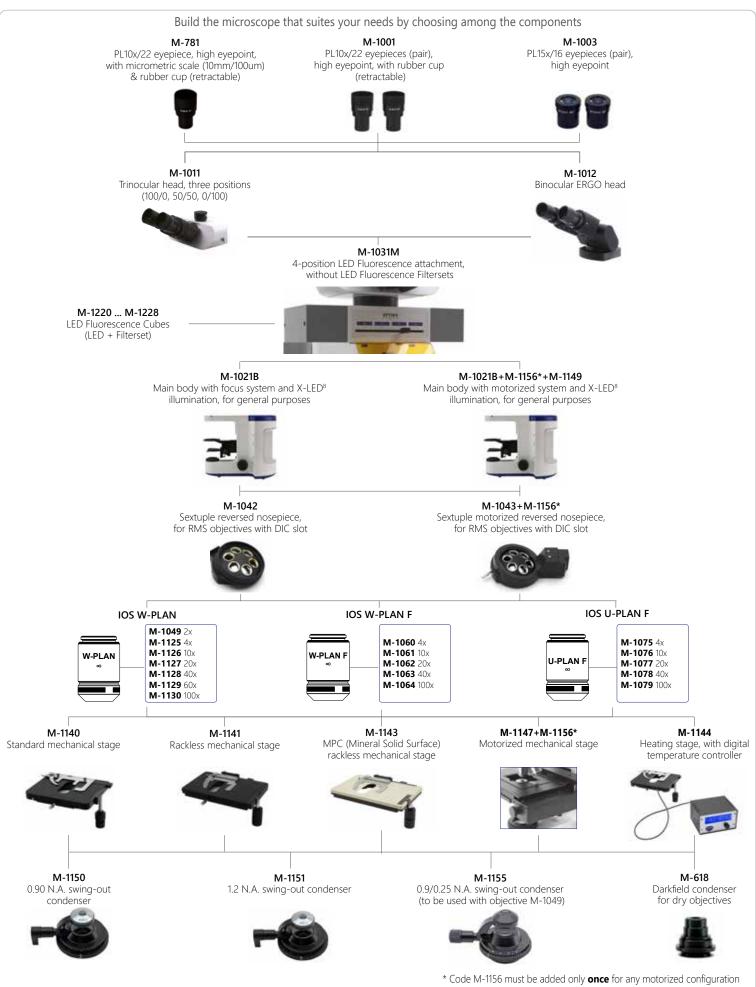


### **B-1000POL** - Polarizing Microscope

The modular OPTIKA B-1000 is available with transmitted polarized light, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration.




### **B-1000POL** - Configuration Chart

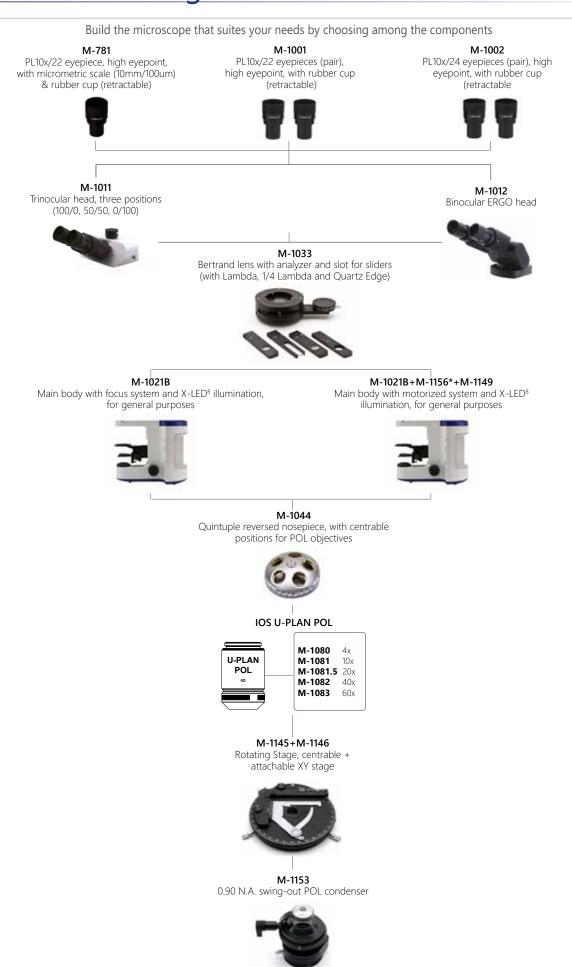

### B-1000POL-I - Polarizing Microscope

The modular OPTIKA B-1000 is available with transmitted and incident polarized light, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive X-LED8 (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

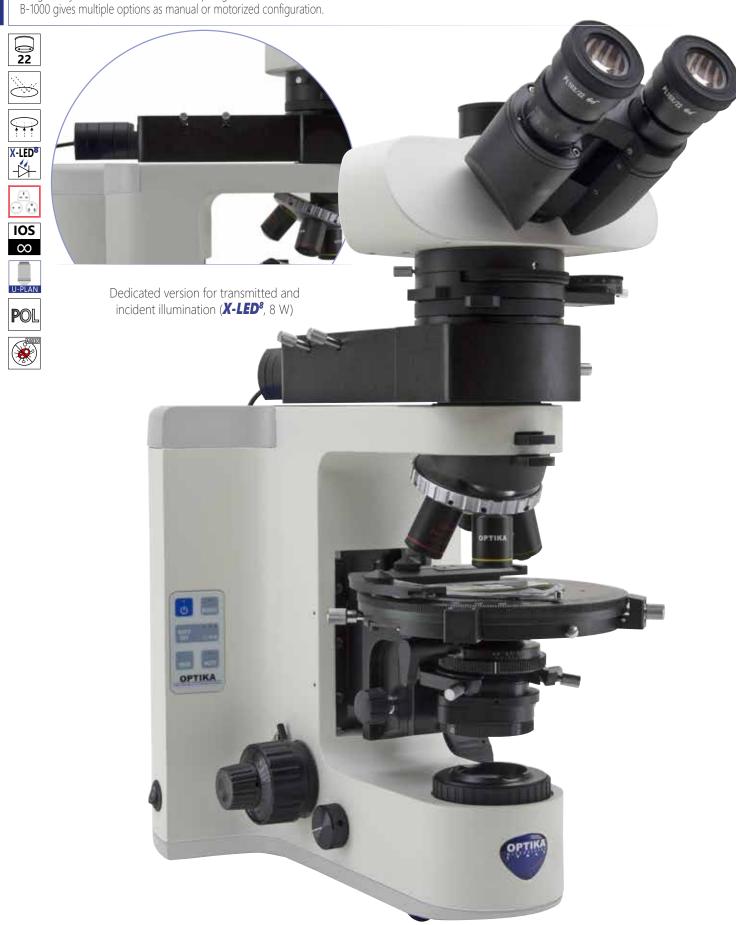

### **B-1000POL-I** - Configuration Chart



### **B-1000METBF** - Brightfield Metallurgical Microscope

The modular OPTIKA B-1000METBF offers superior quality **brightfield incident light**, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) transmitted illumination (Koehler system). The incident light relies on an incredibly bright 18 W LED illumination, designed by OPTIKA.

B-1000 gives multiple options as manual or motorized configuration.



### **B-1000METBF** - Configuration Chart





\* Code M-1156 must be added only **once** for any motorized configuration

### **B-1000METDK** - Darkfield Metallurgical Microscope

The modular OPTIKA B-1000METDK offers superior quality **brightfield and darkfield incident light**, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) transmitted illumination (Koehler system). The incident light relies on an incredibly bright 18 W LED illumination, designed by OPTIKA. B-1000 gives multiple options as manual or motorized configuration.



### **B-1000METDK** - Configuration Chart

Build the microscope that suites your needs by choosing among the components




### **B-1000 Multi-Head** - Discussion Microscopes


The modular OPTIKA B-1000 helps you working in a comfortable way during extended periods of use and let you perform reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Semi-Apo objectives), the state-of-the-art, exclusive **X-LED** illumination

system, designed by OPTIKA and the Koehler diaphragm. B-1000 gives multiple options as manual or motorized configuration, with a variety of objectives, stages and condensers. Ideal for discussion groups and teaching purpose for multiple observers, up to ten users simultaneously.



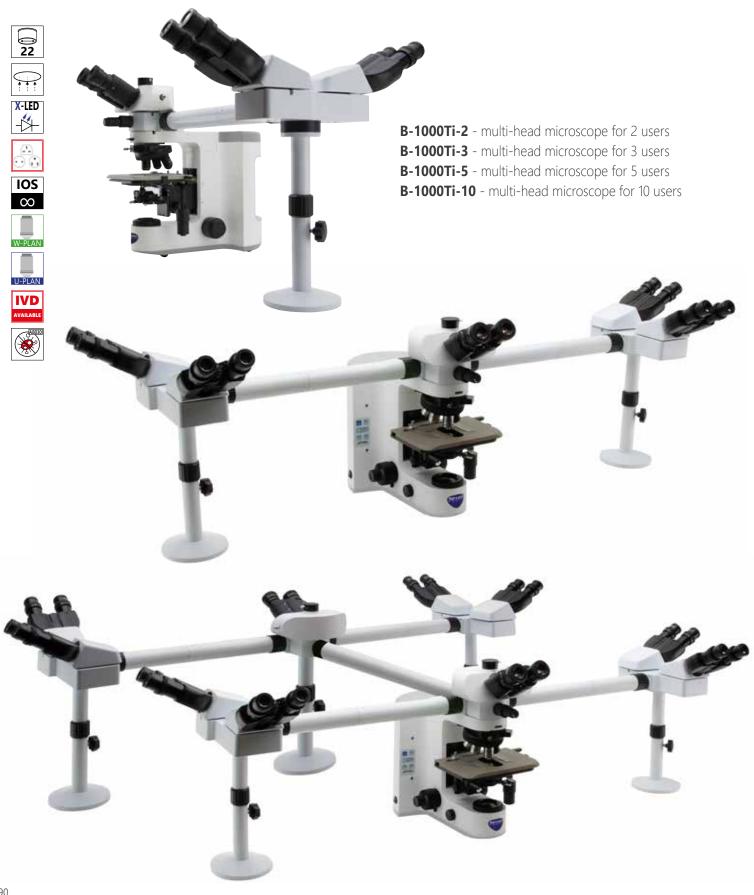

### **B-1000 Multi-Head** - Configuration Chart

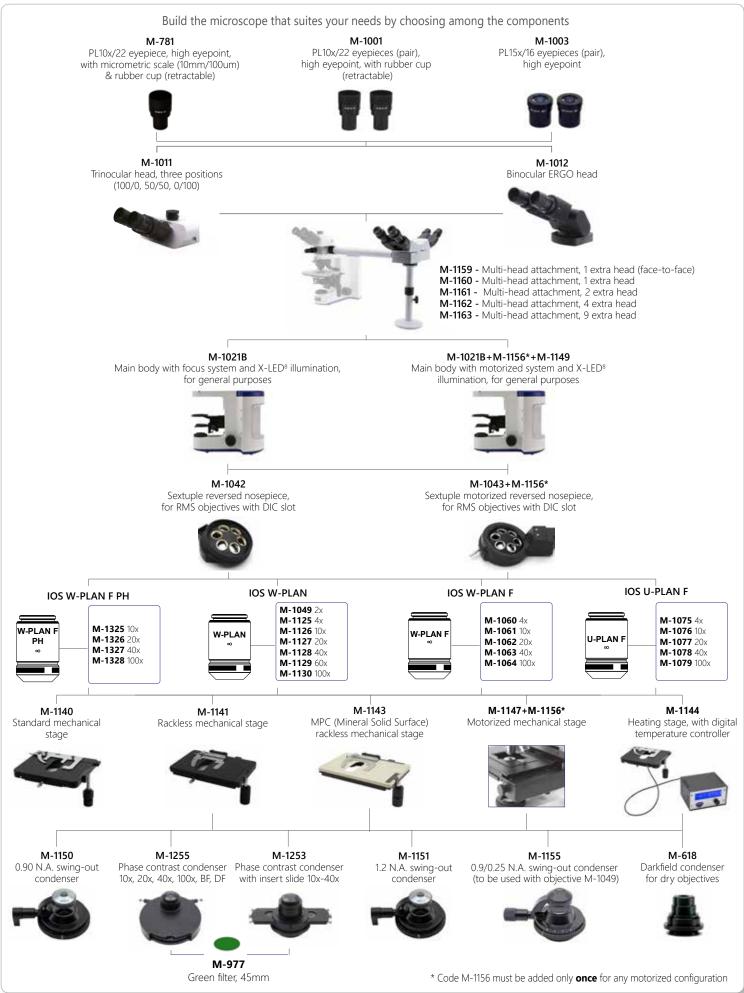

### **B-810/B-1000** - Accessories

#### **Eyecups & Eyepieces**

M-690 Eyecups (pair)

#### **Condensers & Filters**

- M-613 Polarizing set (filters only) (except for B-1000FL-LED, B-1000MET, B-1000POL, and B-1000POL-I)
- M-615 Lambda filter for polarizing set (except for B-1000FL-LED, B-1000MET, B-1000POL, and B-1000POL-I)
- M-975 Blue filter, 45mm diameter (except for B-1000MET)
- M-977 Green filter, 45mm diameter (except for B-1000MET)
- M-979 Yellow filter, 45mm diameter (except for B-1000MET)
- M-989 Frosted glass filter, 45mm diameter (except for B-1000MET)
- M-1047 Gout analisys kit
- M-1164 Empty fluorescence filterblock (only for B-1000FL-HBO)
- M-1165 Fluorescence filter set V (filterblock included) (only for B-1000FL-HBO)
- M-1166 Fluorescence filter set UV-DAPI (filterblock included) (only for B-1000FL-HBO)
- M-1220 Blue LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1220.1 Blue (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1221 Green LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1221.1 Green (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1222 Violet LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1223 UV LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1223.1 UV (pass band) LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-1224 Red 1 LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4) \*
- M-1225 Red 2 LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4) \*
- M-1226 Deep Red LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4) \*
- M-1227 Far Red LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4) \*
- M-1228 Amber LED Fluorescence Cube (LED+Filterset), (only for B-1000LD4)
- M-ND25 Neutral density filter, 25% transmission (only for B-1000FL-HBO)

#### **Camera Adapters**

- M-113.1 Ring adapter, 30mm (for monocular and binocular microscopes)
- M-115 0.35x C-Mount projection lens
- M-114 0.5x C-Mount projection lens
- M-118 0.75x C-Mount projection lens
- M-173 C-Mount projection lens for APS-C/full frame reflex cameras (trino)
- M-620 0.35x focusable C-Mount adapter (biological microscopes)
- M-620.1 0.5x focusable C-Mount adapter (biological microscopes)
- M-620.2 0.65x focusable C-Mount adapter (biological microscopes)
- M-620.3 1x focusable C-Mount adapter (biological & stereomicroscopes)
- M-699 Universal adapter for C-Mount projection lens (trino)

#### Miscellaneous

- 15008 Immersion oil, 10ml
- 15009 Immersion oil, 100ml
- 15104 Cleaning kit
- DC-005 TNT dust cover, extra large, 820(l)x550(h) mm
- M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)
- M-151 HBO 100W high-pressure mercury bulb for fluorescence
- M-151.1 HBO 100W high-pressure mercury bulb for fluorescence (OSRAM)
- VP-1000 IQ/OQ/PQ manual for B-1000 series (Brightfield)
- VP-1000MET IQ/OQ/PQ manual for B-1000 series (Metallographic)
- VP-1000POL IQ/OQ/PQ manual for B-1000 series (Polarizing)
- VP-810 IQ/OQ/PQ manual for B-810 series
- VP-1000PH IQ/OQ/PQ manual for B-1000 series (Phase Contrast)
- VP-1000FL IQ/OQ/PQ manual for B-1000 series (Fluorescence)
- AB-030 Antibacterial surface treatment, only for newly purchased microscope



It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.





How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

\* If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order it by specifying with "AR GLASS"

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA° Spain
OPTIKA° China
OPTIKA° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com OPTIKA<sup>®</sup> North America OPTIKA<sup>®</sup> Central America OPTIKA<sup>®</sup> Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# IM-300 Series



**Routine Lab Inverted Microscopes** 


### Your Preferred Inverted Microscope for Routine

#### **ROUTINE IN UNIVERSITIES, LABS & INDUSTRIES**

- » Wide range to fullfil specific lab requirements
- » Valuable solutions for life and material sciences
- » Compliant with several observation methods

#### AN AFFORDABLE PARTNER WITH HIGH-END FEATURES

- » IOS LWD W-PLAN objectives for flat images on 22 mm FN
- » Fast, efficient investigation with no particular sample prep
- » Trinocular port with beam splitter for most light-demanding needs

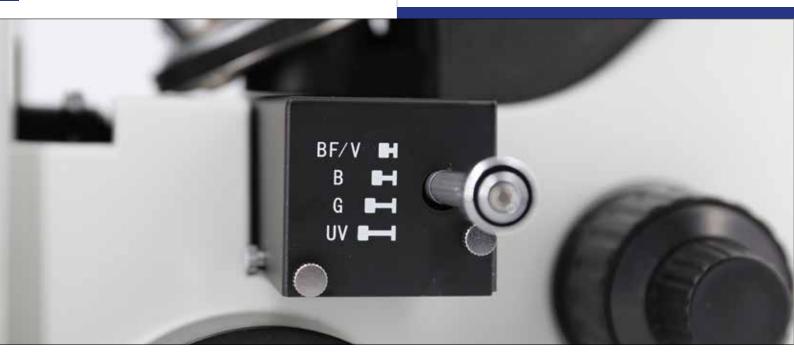

## Optically Impressive

#### **MAINTAINING GOOD EYESIGHT**

- » 10x/22 eyepieces for large specimen view
- » Comfortable rubber cup to get rid of annoying external light
- » High eye-point for glasses wearers, diopter adjustment (left eyepiece)

#### IM-300 & IOS W-PLAN: THE PERFECT COMBINATION

- » IOS Infinity corrected optical system
- » Full planarity optics on 22 mm (W-PLAN) according to ISO 19012-1
- » High-grade Semi-Apo lens available ideal for fluorescence




### An Extensive Range of Different Configurations

#### **OBSERVE EVEN THE MOST COMPLEX SAMPLES**

- » Phase contrast objective for transparent sample examination
- » LED and HBO fluorescence available for specific purposes
- » High quality no cover glass objectives for material science

#### CONCEPTUAL INNOVATION IN LED FLUORESCENCE

- » Choose the lowest operational cost, LED lifetime of 65,000 hours
- » Immediate operation, eliminating warm-up/cool-down times
- » Forget about lamp centering, adjustment and maintenance

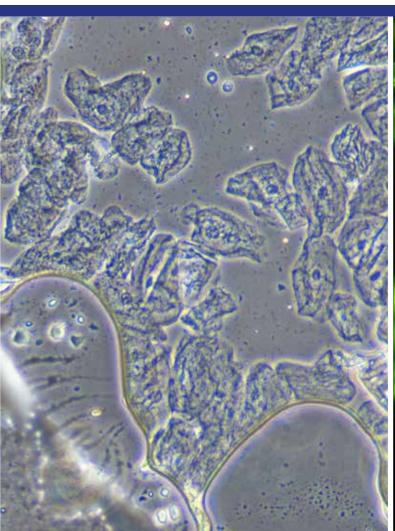




RPC - Relief Phase Contrast

**RPC** system is designed to increase visibility and contrast in unstained and living material by detecting optical gradients (or slopes) and converting them into variations of light intensity. Typical applications are transparent specimens, bacteria, tissue culture work, spermatozoa, cells in glass containers, protozoa, mites, fibers, etc.

When observed under modulation contrast optics (RPC), transparent objects that are essentially not visible in ordinary brightfield microscopy take on an apparent three-dimensional appearance dictated by phase gradients in the specimen. There are also no halos exhibited in the image, unlike the images produced with phase contrast optics.

**RPC** is recommended over **DIC** technique in case of specimens like crystals (with effects upon polarized light), or contained in specimen carriers such as plastic culture vessels, Petri dishes, etc.

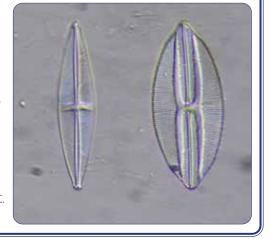

# **X-LED**<sup>8</sup> - Only Available at OPTIKA

#### STATE-OF-THE-ART ILLUMINATION SYSTEM

- » Uncomparable light intensity, exclusive lens & collector design
- » Constant pure-white color temperature at all intensity levels
- » Unmatched color fidelity, uniformity and brightness

#### **CUT ELECTRICITY BILLS BY 90%**

- » Money & energy saving, 8 W
- » More efficient brightness than a 100 W halogen lamp
- » LED long lifetime (65,000 hours = 22 years at 8 hours/day usage)



## Go Digital - Vivid Colors & Contrast For Stunning Images

#### STAY CONNECTED WITH YOUR SPECIMEN, EASILY


- » Trincular port to be always updated with the latest technology cameras, even in the future
- » Wide range of cameras matching all the needs, including the more specific ones
- » Modern C-mount focusable professional adapters for all kinds of cameras

#### **PROFESSIONAL IMAGE ANALYSIS**

- » Multi-language software for live-view, picture and video in different file formats
- » Advanced functions for pictures processing (EDF, stitching, multi-fluorescence combine)
- » Powerful tools to perform measurements and generate custom reports




# Multiple Observation Methods

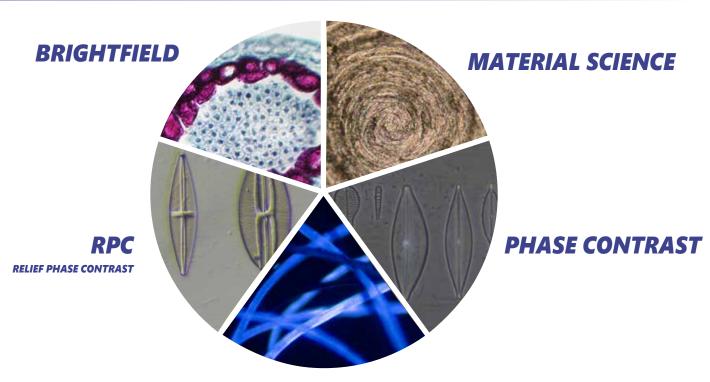

**FLUORESCENCE** 

### **IM-300** Series

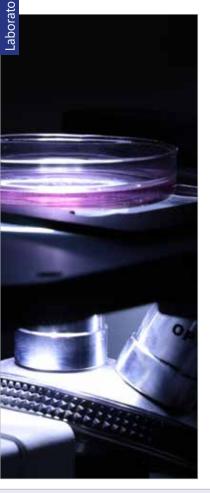


Inverted microscopes are useful for observing living cells or organisms at the bottom of a large container (e.g., a tissue culture flask) under more natural conditions than on a glass slide, as it occurs with a conventional microscope. IM-300 Series is engineered and designed to be your ideal solution for fast and reliable routine inspections, with the exclusive, state-of-the-art X-LED<sup>8</sup> illumination system. A particularly simple and ingenious optical design allows stable alignments and smooth and accurate movements. OPTIKA provides different configurations, including the innovative LED fluorescence technology for a new, enhanced experience.

#### X-LED<sup>8</sup> Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature (6,300K colour temperature).

Relevant money and energy saving thanks to the incredibly low energy consumptions which allows you to cut the electricity bills by 90%!

The electric consumption (8 W only) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.

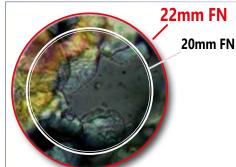

### Large Specimen View (22 mm Field Number)

The **F.O.V.** (field of view) is based on a comfortable diameter of 22 mm.

This means that an extra wide area of the sample can be inspected and allows a natural and easy view, particularly needed in a laboratory environment.



According to your application and to the fluorochromes you are using, we can help you to identify the best light source.



#### **Traditional, HBO Fluorescence**

- » The most used and diffused method, worldwide
- » Wide spectrum range for future upgrades
- » HBO Fluorescence power supply



### Innovative, LED Fluorescence

- » Recommended for routine applications
- » Cost-effective, money saving technology
  - » Ready for immediate operation
  - » Eliminate warm-up/cool-down times
  - » Forget lamp replacement & centering

### Routine Lab Inverted Microscopes

### Get the most out of our accessories



#### **DESIGNED TO FACILITATE YOUR DAILY ROUTINE**

- » Removable condenser to increase the working distance
- » Mechanical stage and side extensions for great comfort (as optional)
- » Different inserts available according to the container used (as optional)



M-793.1 Holder for Petri 38mm diameter (M-793.2 needed)



**M-793.2** Holder for Terasaki and Petri 65mm diameter



**M-793.3** Holder for slides and Petri 54mm diameter



M-793.4 Holder for 2+2 slides



**M-793.5**Holder for metallurgical samples



**M-793.6** Holder for Utermohl-Chamber (M-793.3 needed)



M-793.7 Load bearing side extension



M-792
Mechanical stage for IM-300
(except IM-300LD2 and IM-300LD4)

M-792.2 Mechanical stage for IM-300 (IM-300LD2 only)

### **IM-300** - Brightfield & Phase Contrast Microscope

IM-300 looks at the challenge of the future with confidence, offering first-class optical quality and mechanical versatility, to extend its use with several accessories. Ensuring top-level brightfield and phase contrast observation, as it comes with a set of 3 IOS LWD W-PLAN PH objectives (10x, 20x and 40x). The high-efficiency **X-LED<sup>8</sup>** makes it reliable for all transmitted light observations.

For a more complete solution, choose among the several accessories available (objectives, mechanical stage, side extensions, holders and stage inserts).





















Part	Description				
Observation mode:	Brightfield, phase contrast.				
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.				
Interpupillary distance:	Adjustable between 50 and 75 mm.				
Diopter adjustment:	On the left eyepiece tube.				
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40x/0.65 All with anti-fungus treatment.				

Part	Description			
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.			
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.			
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield. Supplied with green filter IF550.			
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/6Vdc external power supply.			

### IM-300D- Brightfield & Phase Contrast Microscope




















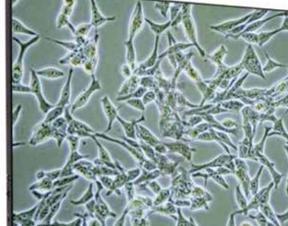

IM-300D looks at the challenge of the future with confidence, offering first-class optical quality and mechanical versatility, to extend its use with several accessories. Ensuring top-level brightfield and phase contrast observation, as it comes with a set of 3 IOS LWD W-PLAN PH objectives (10x, 20x and 40x).
The high-efficiency **X-LED**<sup>8</sup> makes it reliable for all transmitted

light observations.

For a more complete solution, choose among the several accessories available (objectives, mechanical stage, side extensions, holders and stage inserts). This model is equipped with an Intel micro PC, a 15.6" 4K touch screen, 6Mpx high-sensitivity color camera, Optika ProView image analysis software.







Part	Description				
<b>Observation mode:</b>	Brightfield, phase contrast.				
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.				
Interpupillary distance:	Adjustable between 50 and 75 mm.				
<b>Diopter adjustment:</b>	On the left eyepiece tube.				
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40x/0.65 All with anti-fungus treatment.				

Part	Description			
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.			
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.			
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield. Supplied with green filter IF550.			
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/6Vdc external power supply.			
Digital equipment:	Intel micro PC with Image analysis software for imaging. 15.6" 4K touch screen; 6MP high-sensitivity color camera.			

### IM-300F - HBO Fluorescence Microscope

Routine inverted fluorescence microscope for transmitted brightfield, phase contrast and fluorescence observations with IOS LWD W-PLAN objectives. The HBO fluorescence illuminator is combined with blue and green excitation filter set for the visualization of the following fluorochromes: Ethidium Bromide, Acridine Orange, Alexa Fluor 488, Fluo-4, FITC Plus, Spectrum Green, Ethidium Homodimer I, Propidium Iodide (PI) and Spectrum Gold. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



Part	Description				
<b>Observation mode:</b>	Brightfield, phase contrast, HBO fluorescence.				
<b>Epi-illumination and filter:</b>	HBO 100 W high pressure mercury lamp. 3-position filter holder; blue & green included.				
Head:	Trinocular (2-position 100/0, 0/100), 45° inclined.				
Interpupillary distance:	Adjustable between 50 and 75 mm.				
<b>Diopter adjustment:</b>	On the left eyepiece tube.				
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40X/0.65 All with anti-fungus treatment.				

Part	Description				
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.				
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.				
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield. Supplied with green filter IF550.				
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/6Vdc external power supply.				

### IM-300FL4 - HBO Fluorescence Microscope

Advanced inverted microscope for brightfield and fluorescence observations with Semi-Apo IOS LWD U-PLAN F objectives to enhance the visibility of the sample and increase the overall contrast. The HBO fluorescence illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Ethidium Bromide, Acridine Orange, Alexa Fluor 488, Fluo-4, FITC Plus, Spectrum Green, Ethidium Homodimer I, Propidium Iodide (PI) and Spectrum Gold, with additional filterset: Alexa Fluor 350, Aminocoumarin, Hoechst, DAPI, BFP (Blue Fluorescent Protein), Atto 390, Tetracycline, Pacific Orange and Spectrum Blue. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



Part	Description				
Observation mode:	Brightfield, HBO fluorescence.				
<b>Epi-illumination and filter:</b>	HBO 100 W high pressure mercury lamp. 4-position filter holder; blue & green included.				
Head:	Trinocular (2-position 100/0, 0/100), 45° inclined.				
Interpupillary distance:	Adjustable between 50 and 75 mm.				
<b>Diopter adjustment:</b>	On the left eyepiece tube.				
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.				


rait	Description			
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.			
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.			
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm.			
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/6Vdc external power supply.			

### **IM-300LD2** - LED Fluorescence Microscope

Routine inverted fluorescence microscope for transmitted brightfield, phase contrast and fluorescence observations with IOS LWD W-PLAN PH objectives.

The LED fluorescence illuminators are combined with blue and green excitation filter set for the visualization of the following fluorochromes: GFP, Alexa Fluor 488, Calcein, SYBR Green, FITC, Fluo-4, MitoTracker Green, Spectrum Gold, Propidium Iodide, Ethidium Homodimer I. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



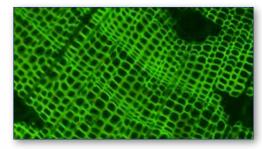



### A new milestone achieved in Fluorescence Microscopy

- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for routine applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

#### **LED Fluorescence Cubes (LED + Filterset) included**

Name	LED emission (nm)	Excita- tion filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
Blue	460	455 - 495	500	510LP
Green	523	510 - 550	570	575LP

Part	Description				
Observation mode:	Brightfield, phase contrast, LED fluorescence.				
<b>Epi-illumination and filter:</b>	High-power LED with brightness control. 3-position filter holder; blue and green filtesets included.				
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.				
Interpupillary distance:	Adjustable between 50 and 75 mm.				
<b>Diopter adjustment:</b>	On the left eyepiece tube.				
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.				
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.				
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40x/0.65 All with anti-fungus treatment.				

Part	Description				
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.				
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.				
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield Supplied with green filter IF550.				
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. 100-240Vac/12Vdc external power supply.				

### IM-300LD4 - LED Fluorescence Microscope

Advanced fluorescence inverted microscope for transmitted brightfield and fluorescence observations with IOS LWD U-PLAN F objectives.

The extremely powerful LED Fluorescence Illuminators are combined with corresponding excitation filter sets for the visualization of most fluorochromes. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



### A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

#### LED Fluorescence Cubes available (LED + Filterset)

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1230 - Blue	460	455 - 495	500	510LP
M-1230.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1231 - Green	523	510 - 550	570	575LP
M-1231.1 - Green (pass band)	523	510 - 550	570	585-625
M-1232 - Violet	405	390 - 420	440	450LP
M-1233 - UV	365	325 - 375	415	435LP
M-1233.1 - UV (pass band)	365	340 - 390	405	420-470
M-1234 - Red 1 *	623	590 - 650	660	665LP
M-1235 - Red 2 *	623	595 - 645	655	665 - 715
M-1236 - Deep Red *	660	623 - 678	685	690 - 750
M-1237 - Far Red *	740	720 - 760	770	780LP
M-1238 - Amber	590	582 - 603	610	615 - 645

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.

Description

Part	Description
Observation mode:	Brightfield, LED fluorescence.
<b>Epi-illumination and filter:</b>	High -Power LED with brightness control. 4-position filter holder; none included.
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Diopter adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.

rait	Description
Specimen stage:	Mechanical stage, 250x290 mm, with round glass and metal stage, inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/12Vdc external power supply.

### IM-300LD4D - LED Fluorescence Microscope























Advanced fluorescence inverted microscope for transmitted brightfield and fluorescence observations with IOS LWD U-PLAN F objectives. The extremely powerful LED Fluorescence Illuminators are combined with corresponding excitation filter sets for the visualization of most fluorochromes. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view. This model is equipped with an Intel micro PC, a 15.6" 4K touch screen, 6Mpx high-sensitivity color camera and Optika ProView image analysis software for fluorescence.

### A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1230 - Blue	460	455 - 495	500	510LP
M-1230.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1231 - Green	523	510 - 550	570	575LP
M-1231.1 - Green (pass band)	523	510 - 550	570	585-625
M-1232 - Violet	405	390 - 420	440	450LP
M-1233 - UV	365	325 - 375	415	435LP
M-1233.1 - UV (pass band)	365	340 - 390	405	420-470
M-1234 - Red 1	623	590 - 650	660	665LP
M-1235 - Red 2	623	595 - 645	655	665 - 715
M-1236 - Deep Red	660	623 - 678	685	690 - 750
M-1237 - Far Red	740	720 - 760	770	780LP
M-1238 - Amber	590	582 - 603	610	615 - 645

Part	Description
Observation mode:	Brightfield, LED fluorescence.
<b>Epi-illumination and filter:</b>	High -Power LED with brightness control. 4-position filter holder; none included.
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Diopter adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.

Part	Description
Specimen stage:	Mechanical stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/12Vdc external power supply.
Digital equipment:	Intel micro PC with Image analysis software for Fluorescence. 15.6" 4K touch screen; 6MP IR sensitive high-sensitivity color camera.

### IM-300METLD- Metallurgical Microscope

LED routine inverted microscope with IOS LWD U-PLAN MET objectives for material science and metallographic applications, combining a sturdy yet compact structure with dedicated components required in this field, like the NCG (no cover glass) objectives working without cover slide ideal for metallographic samples and other opaque specimens. A particularly simple and ingenious optical design allows stable alignments and smooth and accurate movements. This model is equipped with an 18W LED lighting system.



Part	Description
<b>Observation mode:</b>	Brightfield, simple polarized light.
Epi-illumination and polarizing filters:	LED 18 W with brightness control. With centrable aperture and field diaphragms. With polarizer and 360° analyzer. Supplied with blue (LBD) filter.
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Diopter adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.

Part	Description
Objectives:	IOS LWD U-PLAN MET 5x/0.15 IOS LWD U-PLAN MET 10x/0.30 IOS LWD U-PLAN MET 20x/0.45 IOS LWD U-PLAN MET 50x/0.55 All with anti-fungus treatment.
Specimen stage:	Fixed stage, 250x160 mm, with round metal stage insert.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.

## IM-300 Series - Comparison chart

Model	Туре	Objectives	Condenser	Incident illumination	Fluorescence slider	Transmitted illumination	Digital equipment
IM-300	BF, PH	IOS LWD W-PLAN PH 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	-	-	8 W X-LED <sup>8</sup> , brightness control	-
IM-300D	BF, PH	IOS LWD W-PLAN PH 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	-	-	8 W X-LED <sup>8</sup> , brightness control	Intel micro PC 15.6" 4K touch screen; 6MP highsensitivity color camera.
IM-300F	BF, FL, PH	IOS LWD W-PLAN PH 10x, 20, 40x	LWD, N.A. 0.30, iris diaphragm	FL HBO with blue and green filtersets	2-position +BF	8 W X-LED <sup>8</sup> , brightness control	-
IM-300FL4	BF, FL	IOS LWD U-PLAN F 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	FL HBO with blue and green filtersets	3-position +BF	8 W X-LED <sup>8</sup> , brightness control	-
IM-300LD2	BF, FL, PH	IOS LWD W-PLAN PH 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	FL LED with blue and green filtersets	2-position +BF	8 W X-LED <sup>8</sup> , brightness control	-
IM-300LD4	BF, FL	IOS LWD U-PLAN F 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	LED Fluorescence Cubes as optional	4-position	8 W X-LED <sup>8</sup> , brightness control	-
IM-300LD4D	BF, FL	IOS LWD U-PLAN F 10x, 20x, 40x	LWD, N.A. 0.30, iris diaphragm	LED Fluorescence Cubes as optional	2-position	8 W X-LED <sup>8</sup> , brightness control	Intel micro PC 15.6" 4K touch screen, 6MP IR sensitive high-sensitivity color camera.
IM-300METLD	BF, PO	IOS LWD U-PLAN MET 5x, 10x, 20x, 50x	-	LED 18 W, brightness control	-	-	-

### IM-300 Series - Optical performance

#### IM-300 / IM-300LD2 / IM-300F

Eyepiece			10x (M-780)		
Field number (mm)			2	22	
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
2x	0.08	19.40	20x	11	
4x	0.13	10.40	40x	5.50	
10x PH	0.25	7.30	100x	2.20	
20x PH	0.40	6.80	200x	1.10	
40x PH	0.60	3.00	400x	0.55	
40x	0.60	3.00	400x	0.55	
60x	0.70	1.70	600x	0.37	

#### IM-300FL4 / IM-300LD4 / IM-300LD4D

Eyepiece			10x (M-780)			
Field number (mm)			2	22		
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)		
4x	0.13	18.52	40x	5.50		
10x	0.30	7.11	100x	2.20		
20x	0.45	5.91	200x	1.10		
40x	0.65	1.61	400x	0.55		
60x	0.75	1.04	600x	0.37		

#### IM-300MFTI D

IIVI-300IVIETED						
Eyepiece			10x (l	M-780)	15x (I	VI-601)
Field number (mm)				22	1	16
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	Total magnification	Field of view (mm)
5x	0.15	10.80	50x	4.40	75x	3.20
10x	0.30	10	100x	2.20	150x	1.60
20x	0.45	4	200x	1.10	300x	0.80
50x	0.55	7.90	500x	0.44	750x	0.32
100x	0.80	2.10	1000x	0.22	1500x	0.16

### **IM-300** Series - Accessories

	£ Eyepieces	M-7
M-601		<u>M-7</u>
M-780		<u>M-7</u>
M-781		M-7
Objectives		<u>M-7</u>
IOS W-PL		<u>M-7</u>
	IOS W-PLAN objective 2x/0.08	Con
M-782	IOS LWD W-PLAN objective 4x/0.13	<u>M-1</u>
M-773		<u>M-7</u>
M-786	IOS LWD W-PLAN objective 60x/0.70	<u>M-8</u>
IOS W-PL		
	IOS LWD W-PLAN PH objective 4x/0.13	<u>M-8</u>
	IOS LWD W-PLAN PH objective 10x/0.25	
<u>M-784N</u>	,	Can
<u>M-785</u>	,	M-1
IOS U-PLA		M-1 M-1
M-800		M-1
M-801		M-1
	IOS LWD U-PLAN F objective 20x/0.45	141 1
M-803		M-6
M-804	IOS LWD U-PLAN F objective 60x/0.75	M-6
IOS U-PLA		M-6
	IOS LWD U-PLAN F PH objective 20x/0.45	M-6
	IOS LWD U-PLAN F PH objective 40x/0.65	<u>M-6</u>
IOS U-PLA		١.
	IOS LWD U-PLAN MET objective 5x/0.15	Acc
	IOS LWD U-PLAN MET objective 10x/0.30	<u>AB-(</u>
	IOS LWD U-PLAN MET objective 20x/0.45	<u>AB-</u>
	IOS LWD U-PLAN MET objective 50x/0.55	DC-
M-1104	IOS LWD U-PLAN MET objective 100x/0.80 (dry)	<u>M-0</u>
ı	J-PLAN RPC	M-1
M-861		M-7
M-862	IOS LWD U-PLAN RPC Objective 10x/0.25	<u>M-9</u>
M-863		VP-I
M-864	IOS LWD U-PLAN RPC Objective 40x/0.65	+ IC
FLuo Atta		* If t
M-676	Empty fluorescence filterblock	leng
M-677	Fluorescence filter set V (filterblock included)	
M-677.1	Fluorescence filter set V (filterblock NOT included)	
<u>M-677ND</u>	Neutral density ND filter, 25% transmission	
<u>M-678</u>	Fluorescence filter set UV-DAPI (filterblock included)	
<u>M-678.1</u>	Fluorescence filter set UV-DAPI (filterblock NOT included)	
<u>M-678ND</u>	Neutral density ND filter, 50% transmission	
M-797	HBO fluo attachment, 2-pos. (B&G filter set)	
M-798	HBO fluo attachment, 4-pos. (B&G filter set)	
M-1230	Blue LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1230.1	Blue (pass band)LED Fluorescence Cube (LED+Filterset), for IM-300	<u>LD4</u>
M-1231	Green LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1231.1	Green (pass band) LED Fluorescence Cube (LED+Filterset), for IM-3	<u>00LD4</u>
M-1232	Violet LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1233	UV LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1233.1	UV (pass band) LED Fluorescence Cube (LED+Filterset), for IM-3001	<u>_D4</u>
M-1234	Red 1 LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1235	Red 2 LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
<u>M-1236</u>	Deep Red LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1237	Far Red LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
M-1238	Amber LED Fluorescence Cube (LED+Filterset), for IM-300LD4	
Stages		
M-792	Mechanical stage for IM-300 (except IM-300LD2 and IM-300LD4)	

M-793.2	Holder for Terasaki and Petri 65mm diameter	
M-793.3	Holder for slides and Petri 54mm diameter	
M-793.4	Holder for 2+2 slides	
M-793.5	Holder for metallurgical samples	
M-793.6	Holder for Utermohl-Chamber (M-793.3 needed)	
M-793.7	Load bearing side extension	
Condensers & Filters		
M-1004.N	Centering telescope (30mm Ø)	
M-785.2N	Slider with phase rings (4x/10x, 20x/40x, BF)	
M-860	Slider with rotating 10x-20x-40x OPTIKA	
	Modulation Contrast (OMC) (IM-300/IM-5)	
M-860.1	Slider with rotating 4x	
	OPTIKA Modulation Contrast (OMC) (IM-300/IM-5)	
Camera Adapters		
M-113.1	Ring adapter, 30mm (for monocular and binocular microscopes)	
M-114	0.5x C-Mount projection lens	
M-115	0.35x C-Mount projection lens	
M-118	0.75x C-Mount projection lens	
M-173	T2-Mount projection lens for APS-C/full frame reflex cameras	
	(trino) (T2 ring not included)	
M-620	0.35x focusable C-Mount adapter	
M-620.1	0.5x focusable C-Mount adapter	
M-620.2	0.65x focusable C-Mount adapter	
M-620.3	1x focusable C-Mount adapter	
M-699	Universal adapter for projection lens (trino)	
Accessories		

710005501105	
AB-030	Antibacterial surface treatment, for B-800/1000, IM-300/5/7, IS
AB-040	Antibacterial surface treatment, for B-800/1000, IM-300/5/7, IS
DC-004	TNT dust cover, large, 700(l)x550(h) mm
M-005	Micrometric slide, 26x76mm, 2 scales (1mm/100 & 10mm/100)
M-151.1	HBO 100W high-pressure mercury bulb for fluorescence
M-787	Cut-off filter (infrared)
M-977	Green filter, 45mm diameter
VP-IM300	IQ/OQ/PQ manual for IM-300

the use of a camera is needed, when used for red emission fluorescence waveths above 650nm, please order it by specifying with "AR GLASS"





M-792.2

Mechanical stage for IM-300 (IM-300LD2 only) Holder for Petri 38mm diameter (M-793.2 needed)

 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# IM-5 Series



**Routine & Research Lab Inverted Microscopes** 

### The Best Option for Routine & Research

#### INTUITIVE YET SUPERIOR CONFIGURATIONS FOR PROFESSIONALS

- » Wide range to fullfil specific lab requirements
- » Valuable solutions for life and material sciences
- » Compliant with several observation methods

#### AN AFFORDABLE PARTNER WITH UNIQUE HIGH-END FEATURES

- » IOS LWD U-PLAN objectives for flat images on 24 mm FN
- » Fast, efficient investigation with no particular sample prep
- » Trinocular port with beam splitter for most light-demanding needs




## Optically Impressive

#### **MAINTAINING GOOD EYESIGHT**

- » 10x/24 eyepieces for the highest F.O.V. on an inverted microscope
- » Comfortable rubber cup to get rid of annoying external light
- » High eye-point for glasses wearers and dioptric adjustment

#### IM-5 & IOS U-PLAN: THE PERFECT COMBINATION

- » IOS Infinity corrected optical system
- » Full planarity optics on 24 mm (U-PLAN) according to ISO 19012-1
- » High-grade Semi-Apo lens available ideal for fluorescence

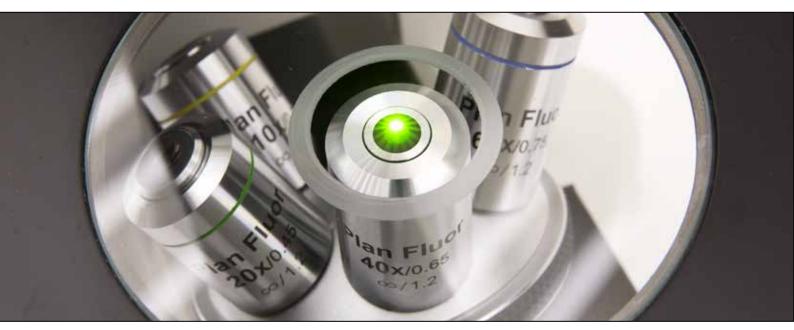


### An Extensive Range of Different Configurations

#### **OBSERVE EVEN THE MOST COMPLEX SAMPLES**

- » Phase contrast lens for transparent sample examination
- » Motorized LED fluorescence available for specific purposes
- » High quality no cover glass objectives for material science

### AUTOMATIC LED SELECTION & CONCEPTUAL INNOVATION IN LED FLUORESCENCE

- » Choose the fluorescence filter for motorized LED selection
- » Immediate operation, eliminating warm-up/cool-down times
- » Forget about lamp centering, adjustment and maintenance

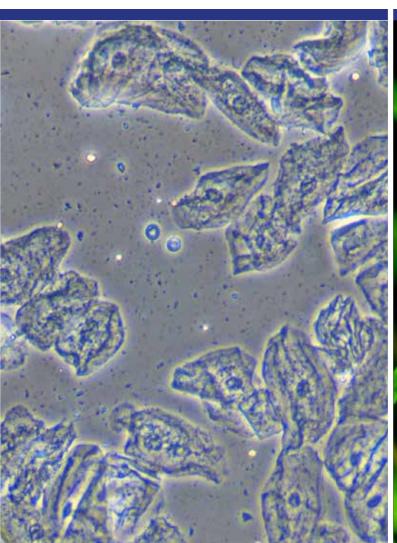




### RPC Modulation Contrast

**RPC** system is designed to increase visibility and contrast in unstained and living material by detecting optical gradients (or slopes) and converting them into variations of light intensity. Typical applications are transparent specimens, bacteria, tissue culture work, spermatozoa, cells in glass containers, protozoa, mites, fibers, etc.

When viewed under modulation contrast optics, transparent objects that are essentially invisible in ordinary brightfield microscopy take on an apparent three-dimensional appearance dictated by phase gradients in the specimen. There are also no halos exhibited in the image, unlike the images produced with phase contrast optics.

**RPC** is recommended over **DIC** technique in case of specimens like crystals (with effects upon polarized light), or contained in specimen carriers such as plastic culture vessels, Petri dishes, etc.

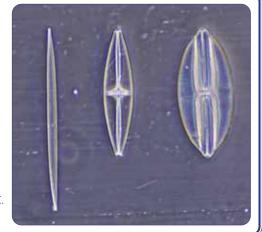

### Born To Be Professional

#### **DESIGNED TO FACILITATE YOUR DAILY ACTIVITIES**

- » Mechanical stage and side extensions for great comfort
- » Large, resistant stage to easily and quickly process samples
- » Different inserts available according to the container used

### CREATE YOUR COMPLETE, FLEXIBLE WORKING STATION

- » Integrable micromanipulation system available
- » Hoffman® modulation contrast available
- » Stage top incubation system available



# Go Digital - Vivid Colors & Contrast For Stunning Images

### STAY CONNECTED WITH YOUR SPECIMEN, EASILY

- » Trincular port to be always updated with the latest technology cameras, even in the future
- » Wide range of cameras matching all the needs, including the more specific ones
- » Modern C-mount focusable professional adapters for all kinds of cameras

#### **PROFESSIONAL IMAGE ANALYSIS**

- » Multi-language software for live-view, picture and video in different file formats
- » Advanced functions for pictures processing (EDF, stitching, multi-fluorescence combine)
- » Powerful tools to perform measurements and generate custom reports



# **X-LED**<sup>8</sup> - Only Available at OPTIKA

#### STATE-OF-THE-ART ILLUMINATION SYSTEM

- » Uncomparable light intensity, exclusive lens & collector design
- » Constant pure-white color temperature at all intensity levels
- » Unmatched color fidelity, uniformity and brightness

#### **CUT ELECTRICITY BILLS BY 90%**

- » Money & energy saving, 8 W (on X-LED8)
- » More efficient brightness than a 100 W (for X-LED8) halogen lamp
- » LED long lifetime (65,000 hours = 22 years at 8 hours/day usage)



# Adjust It To Your Individual Needs

### FULLY SETTABLE, ADJUSTABLE IN HEIGHT CONDENSER FOR PERFECT IMAGING


- » Full Koehler illumination for enhanced images
- » Field & aperture diaphragms, centrable; N.A. 0.50 condenser
- » Removable/rotatable condenser to increase the working distance

### IMPROVED OPTICAL PERFORMANCE, LONG WORKING DISTANCE

- » Superior image quality, crisp and bright details
- » Excellent contrast and resolution due to high numerical apertures
- » Comprehensive range of objectives for extended versatility



### **IM-5** Series

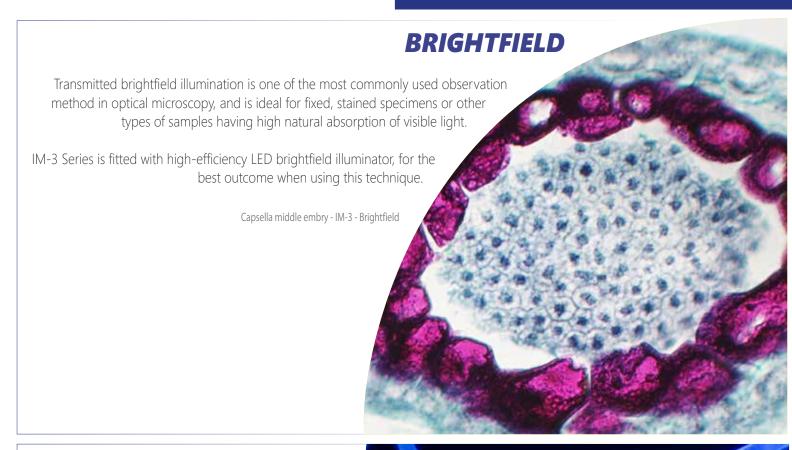

### **FLUORESCENCE**

The fluorescence microscopy is the most demanding technique in biology and biomedical sciences, as well as in materials science.

This method is capable to study organic and inorganic samples thanks to primary fluorescence (auto-fluorescence) or secondary (staining and labelling with fluorochromes)


IM-Series is tailored for applications in research, clinical and pharmaceutical diagnostic field.


Fluorescence illuminators available as mercury lamp (IM-3F & IM-3FL4) and also as LED (IM-3LD).

Cotton fibers - IM-3FL4 - UV Fluorescence



## Multiple Observation Methods

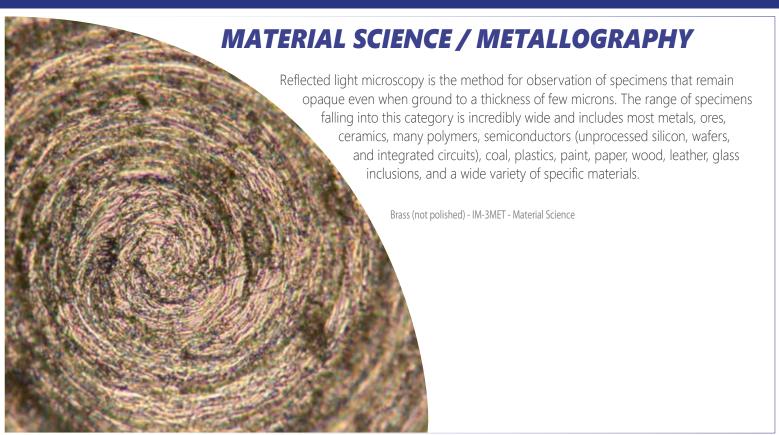

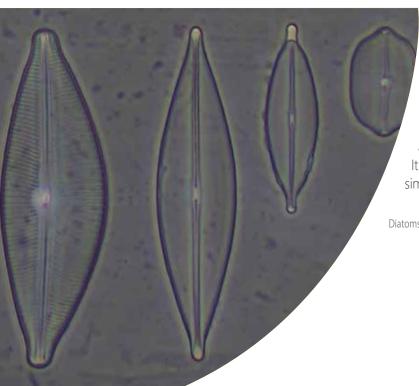

### **PHASE CONTRAST**

Phase-contrast microscopy is a particular technique applied in transparent, non-stainable, samples like culture of living cells, microorganisms, lithographic patterns, latex dispersions, fibers, asbestos and subcellular particles.

It reveals many cellular structures that are not visible with a simple brightfield microscope.

Diatoms - IM-3 - Phase contrast

# **IM-5** Series

### Significant Time And Money Saving

The Series has been designed to increase comfort and achieve significant benefits, especially in terms of time saving with quick and intuitive installation, pre-aligned phase contrast system and pre-aligned LED light source.


As time is money, these features bring to a drastic impact on cost reduction, even more evident thanks to the exclusive illumination system provided by OPTIKA.

### X-LED<sup>8</sup> Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature (6,300K colour temperature).

Relevant money and energy saving thanks to the incredibly low energy consumptions which allows you to cut the electricity bills by 90%!

The electric consumption (8 W) proves the high efficiency of this system: incredibly high light intensity combined with low consumption.

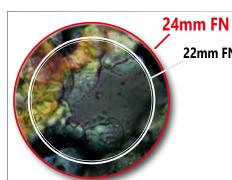

### **The Widest Specimen Area Available** (24mm Field Number) 22mm FN


The **F.O.V.** (field of view) is based on a very comfortable diameter of 24 mm.

This means that an extra wide area of the sample can be inspected and allows a natural and easy view, particularly needed in a laboratory environment.



IM-5 ensures significant repeatibility since the level of light intensity can be seen at any time from the frontal panel in order to reproduce the same conditions. "ECO" button makes the microscope more environmentally sensitive, with automatic switch-off after 20 minutes of inactivity.



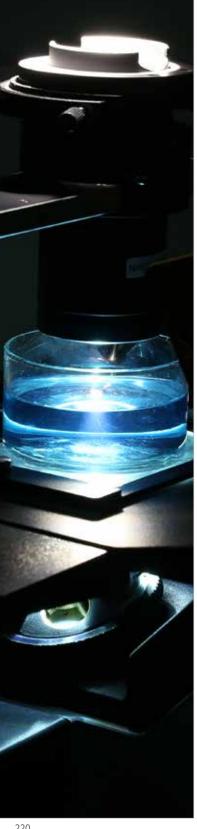

### Routine & Research Lab Inverted Microscopes

### In fluorescence we offer the latest technology.

IM-5FLD is a state-of-the art LED fluorescence microscope, equipped with motorized selection of the best LED according to the filter selected (blue, green, UV and an empty position fo optional filter) by using the filter holder slide.

### **Innovative, LED Fluorescence**

- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering



### Get the most out of our accessories



#### M-793.4

Holder for 2+2 slides. (Only for IM-5 and IM-5FLD)



#### M-793.5

Holder for small metallurgical samples. (Only for IM-5MET)



#### M-793.6

Holder for Utermöhl-Chamber. (Only for IM-5 and IM-5FLD)

### **Accessories included**



Holder for Petri diameter 38mm. (Included with IM-5 and IM-5FLD)



Holder for Terasaki and Petri diameter 65mm. (Included with IM-5 and IM-5FLD)



Holder for slide and Petri diameter 54mm. (Included with IM-5 and IM-5FLD)

# IM-5 - Brightfield & Phase Contrast Microscope

Phase contrast, brightfield and darkfield (dry) trinocular inverted microscope ideal for laboratory requirements (especially cell culture), with freely configurable lenses according to customer's preferences, FN 24 high eyepoint, infinity corrected optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 8 W X-LED\*. Sturdy and incredibly reliable, it is equipped with all the main controls in ergonomic position and with long lasting LED illumination to provide over 20 years of use.



## **IM-5** - Specifications

**Condenser:** 

(Full Koehler):

**Transmitted illumination** 



Abbe N.A. 0.50, removable, with iris diaphragm and slider for phase contrast.

X-LED8 with white 8 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power

### IM-5 is freely configurable in terms of objectives, by choosing among:

Included ■ Optional □

Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:		
M-782	IOS LWD W-PLAN objective 4x/0.13	
M-773	IOS LWD W-PLAN objective 40x/0.60	
M-786	IOS LWD W-PLAN objective 60x/0.70	

Positive Phase Contrast Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:		
M-782.1	IOS LWD W-PLAN PH objective 4x/0.13	
M-783N	IOS LWD W-PLAN PH objective 10x/0.25	
M-784N	IOS LWD W-PLAN PH objective 20x/0.40	
M-785	IOS LWD W-PLAN PH objective 40x/0.65	

Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-800	IOS LWD U-PLAN F objective 4x/0.13	
M-801	IOS LWD U-PLAN F objective 10x/0.30	
M-802	IOS LWD U-PLAN F objective 20x/0.45	
M-803	IOS LWD U-PLAN F objective 40x/0.65	
M-804	IOS LWD U-PLAN F objective 60x/0.75	

Positive Phase Contrast Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1177	IOS LWD U-PLAN F PH objective 20x/0.45	
M-1178	IOS LWD U-PLAN F PHobjective 40x/0.65	

### IM-5FLD - LED Fluorescence Microscope

Phase contrast, brightfield and darkfield (dry) LED fluorescence trinocular inverted microscope, with freely configurable lenses according to customer's preferences, FN 24 high eyepoint, infinity corrected optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 8 W *X-LED8*. The 4-position epi-fluorescence attachment is powered by extremely powerful 5 W LEDs fluorescence illuminator and combined with blue, green and UV excitation filters for the visualization of the following fluorochromes: Acridine Yellow, Acridine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP (blue filter) plus Rhodamine, Texas Red and TRITC (green filter) plus Alexa Fluor® 350, 7- Amino-4-methylcoumarin, 6-Aminoquinoline, Calcofluor® White, Dansyl cadaverine, DAPI, Dapoxyl, DIDS, Europium (III) Chloride, Fluoro-Gold™, Fura-2, Hoechst 33342 & 33258, 1,5 IAEDANS, Indo-1, Marina Blue®, 4-Methylumbelliferone, PBF1, Pyrene, SBFI, Y66F, Y66H (UV filter) among the others. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Sturdy and incredibly reliable, it is equipped with all the main controls in ergonomic position and with long lasting, efficient LED illumination to provide over 20 years of use.



# **IM-5FLD** - Specifications



Part	Description
Head:	Trinocular (split ratio: 100/0, 0/100), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	Both eyepieces.
Eyepieces:	WF10x/24 mm, high eyepoint and with retractable rubber cups.
Epi-fluorescence illumination & filters:	High-power 5 W LEDs with brightness control, motorized LED selection with centrable field diaphragm, 4-position filter holder; blue (EX 450-490, DM 495, EM 500-550), green (EX 540-580, DM 585, EM 608-682) and UV (EX 340-390, DM 400, EM 420LP) excitation filters included.
Nosepiece:	Quintuple ball bearings revolving nosepiece, reversed.
Objectives:	Selectable according to customer's preferences. All with anti-fungus treatment.
Specimen stage:	Fixed stage, 215x250 mm and attachable mechanical stage, 290x250 mm, 120x80 mm X-Y movement range.
Focusing:	Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.
Condenser:	Abbe N.A. 0.50, removable, with iris diaphragm and slider for phase contrast.
Transmitted illumination (Full Koehler):	X-LED <sup>8</sup> with white 8 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power supply.

#### **Fluorescence filtersets**

Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
B (Blue)	450 – 490	495	500 - 550
G (Green)	540 – 580	585	607 - 682
UV (Ultraviolet)	340 -390	400	420LP

### IM-5FLD is freely configurable in terms of objectives, by choosing among:

Included ■ Optional □

Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:		
M-782	IOS LWD W-PLAN objective 4x/0.13	
M-773	IOS LWD W-PLAN objective 40x/0.60	
M-786	IOS LWD W-PLAN objective 60x/0.70	

Positive Phase Contrast Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:			
M-782.1	IOS LWD W-PLAN PH objective 4x/0.13		
M-783N	IOS LWD W-PLAN PH objective 10x/0.25		
M-784N	IOS LWD W-PLAN PH objective 20x/0.40		
M-785	IOS LWD W-PLAN PH objective 40x/0.65		

	rected Semi-Apochromatic, Long Working Distance field flatness up to F.N. 25:	
M-800	IOS LWD U-PLAN F objective 4x/0.13	
M-801	IOS LWD U-PLAN F objective 10x/0.30	
M-802	IOS LWD U-PLAN F objective 20x/0.45	
M-803	IOS LWD U-PLAN F objective 40x/0.65	
M-804	IOS LWD U-PLAN F objective 60x/0.75	

Positive Phase Contrast Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1177	IOS LWD U-PLAN F PH objective 20x/0.45	
M-1178	IOS LWD U-PLAN F PHobjective 40x/0.65	

## **IM-5MET** - Metallurgical Microscope

Industrial and materials science inverted microscope especially designed for opaque specimens (including metals microstructure investigation and studies such as grain size, grain boundaries, phases, transformation, inclusions, and non-metals, as well as sample preparation and treatment) in metallography labs. Freely configurable lenses according to customer's preferences, FN 24 high eyepoint, infinity corrected optical system, coaxial focusing, mechanical stage, and epi-illumination attachment powered by halogen 12 V/100 W with brightness control. Sturdy and incredibly reliable, it is equipped with all the main controls in ergonomic position and with long lasting, efficient LED illumination to provide over 20 years of use.



# **IM-5MET** - Specifications



Part	Description
Head:	Trinocular (split ratio: 100/0, 50/50), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	Both eyepieces.
Eyepieces:	WF10x/24 mm, high eyepoint, secured by screw and with retractable rubber cups.
Epi-illumination & filters:	Halogen 12 V/100 W with brightness control. With field and aperture diaphragms, polarizer and analyzer filters.
Nosepiece:	Quintuple ball bearings revolving nosepiece, reversed.
Objectives:	Selectable according to customer's preferences. All with anti-fungus treatment.
Specimen stage:	Mechanical stage, 240x250 mm.
Focusing:	Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

### IM-5MET is freely configurable in terms of objectives, by choosing among:

Included ■ Optional □

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1100	IOS LWD U-PLAN MET objective 5x/0.15	
M-1101	IOS LWD U-PLAN MET objective 10x/0.30	
M-1102	IOS LWD U-PLAN MET objective 20x/0.45	
M-1103	IOS LWD U-PLAN MET objective 50x/0.55	
M-1104	IOS LWD U-PLAN MET objective 100x/0.80 (dry)	

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:				
M-1171	IOS LWD U-PLAN F MET objective 5x/0.15			
M-1172	IOS LWD U-PLAN F MET objective 10x/0.30			
M-1173	IOS LWD U-PLAN F MET objective 20x/0.50			
M-1174	IOS LWD U-PLAN F MET objective 50x/0.80			
M-1175	IOS LWD U-PLAN F MET objective 100x/0.90 (dry)			

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:							
M-1094	M-1094 IOS LWD U-PLAN MET BD objective 5x/0.15						
M-1095 IOS LWD U-PLAN MET BD objective 10x/0.30							
M-1096 IOS LWD U-PLAN MET BD objective 20x/0.45							
M-1097 IOS LWD U-PLAN MET BD objective 50x/0.55							
M-1098	IOS LWD U-PLAN MET BD objective 100x/0.80 (dry)						

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:						
M-1180	M-1180 IOS LWD U-PLAN F MET BD objective 5x/0.15					
M-1181	M-1181 IOS LWD U-PLAN F MET BD objective 10x/0.30					
M-1182 IOS LWD U-PLAN F MET BD objective 20x/0.50						
M-1183	M-1183 IOS LWD U-PLAN F MET BD objective 50x/0.80					
M-1184	IOS LWD U-PLAN F MET BD objective 100x/0.90 (dry)					

# IM-5 Series - Comparison Chart

#### **Common features:**

- Head: Trinocular (2-position 100/0, 0/100), 45° inclined.
  Eyepieces: PL10x/24 mm, with dioptric adjustment, high eye-point and rubber cups. Dioptric adjustment on both eyepieces.
  Focusing mechanism: Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.

Model	Туре	Nosepiece	Stage	Condenser	Incident illumination	Fluorescence slider	Transmitted illumination
IM-5	BF, PH	Quintuple revolving nosepiece, rotation on ball bearings	Fixed, 215x250 mm can be equipped with mechanical (included), 290x250 mm, 120x80 mm movement range	N.A. 0.50 Köhler, W.D. 28 mm, rotatable to extend the W.D.	-	-	8 W X-LED <sup>8</sup> , brightness control and ECO function
IM-5FLD	BF, FL, PH	Quintuple revolving nosepiece, rotation on ball bearings	Fixed, 215x250 mm can be equipped with mechanical (included), 290x250 mm, 120x80 mm movement range	N.A. 0.50 Köhler, W.D. 28 mm, rotatable to extend the W.D.	FL LED with Blue, Green and UV filtersets	4-position	8 W X-LED <sup>8</sup> , brightness control and ECO function
IM-5MET	BF MET, DF MET	Quintuple revolving nosepiece, rotation on ball bearings. With 26 mm thread holes, 5 adapter rings (for RMS objectives) and DIC slot	Rackless, mechanical, 240x250 mm, 50x50 mm movement range	-	Halogen bulb, 12 V/100 W, brightness control and ECO function	-	-



### IM-5 Series - Optical Performance

Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:

Eyepiece			10x (M-880)		
Field number			24 (mm)		
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
4x	0.13	10.40	40x	6.0	
40x	0.60	3.10	400x	0.60	
60x	0.70	1.70	600x	0.40	

Positive Phase Contrast Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:

Eyepiece			10x (M-880)		
Field number			24 (r	nm)	
Objective	N.A.	W.D. (mm)	Total magnifi- cation	Field of view (mm)	
4x	0.13	10.40	40x	6.0	
10x	0.25	7.30	100x	2.4	
20x	0.40	6.80	200x	1.2	
40x	0.60	3.00	400x	0.60	

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 25:

Eyepiece			10x (M-880)		
Field number			24 (r	nm)	
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
5x	0.15	10.80	50x	4.8	
10x	0.30	10.00	100x	2.40	
20x	0.45	4.00	200x	1.20	
50x	0.55	7.90	500x	0.48	
100x	0.80	2.10	1000x	0.24	

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:

Eyepiece			10x (M-880)		
Field number			24 (m	nm)	
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
5x	0.15	9.00	50x	4.8	
10x	0.30	9.00	100x	2.40	
20x	0.45	3.40	200x	1.20	
50x	0.55	7.50	500x	0.48	
100x	0.80	2.00	1000x	0.24	

Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:

Eyepiece			10x (M-880)	
Field number			24 (ı	mm)
Objective	N.A.	W.D. (mm)	Total magnifi- cation	Field of view (mm)
4x	0.13	18.52	40x	6.0
10x	0.30	7.11	100x	2.4
20x	0.45	5.91	200x	1.2
40x	0.65	1.61	400x	0.60
60x	0.75	1.04	600x	0.40

Positive Phase Contrast Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:

Eyepiece			10x (M-880)		
Field number			24 (mm)		
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
20x	0.45	5.91	20x	1.2	
40x	0.65	1.61	400x	0.60	

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:

Eyepiece			10x (M-880)	
Field number			24 (mm)	
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)
5x	0.15	19.50	50x	4.8
10x	0.30	10.9	100x	2.40
20x	0.50	3.20	200x	1.20
50x	0.80	1.2	500x	0.48
1000x	0.90	1.00	1000x	0.24

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:

Eyepiece			10x (M-880)		
Field number			24 (r	mm)	
Objective	N.A.	W.D. (mm)	Total magnification	Field of view (mm)	
5x	0.15	13.50	50x	4.8	
10x	0.30	9.00	100x	2.40	
20x	0.50	2.50	200x	1.20	
50x	0.80	1.00	500x	0.48	
100x	0.90	1.00	1000x	0.24	



### **IM-5** Series - Accessories

M-880	
M-881	
M-882	
Objective IOS W-PL	
	IOS LWD W-PLAN objective 4x/0.10
M-773	
M-786 IOS W-PL	
	IOS LWD W-PLAN PH objective 4x/0.13
	IOS LWD W-PLAN PH objective 4x/0.15
	IOS LWD W-PLAN PH objective 10x/0.25
M-785	
IOS U-PL	
M-800	IOS LWD U-PLAN F objective 4x/0.13
M-801	
M-802	IOS LWD U-PLAN F objective 20x/0.45
M-803	IOS LWD U-PLAN F objective 40x/0.65
M-804	•
IOS U-PL	
	IOS LWD U-PLAN F PH objective 20x/0.45
M-1178	
	AN MET (Brightfield)
M-1100	
M-1101	
M-1102	IOS LWD U-PLAN MET objective 10x/0.30
M-1103	
M-1104	
	AN MET (Brightfield & Darkfield)
M-1094	
M-1095	
M-1096	
M-1097	IOS LWD U-PLAN MET BD objective 50x/0.55
M-1098	
	AN F MET (Brightfield)
M-1171	
M-1172	IOS LWD U-PLAN F MET objective 10x/0.30
M-1173	
M-1174	
	IOS LWD U-PLAN F MET objective 100x/0.90 (dry)
	AN F MET (Brightfield & Darkfield)
M-1180	IOS LWD U-PLAN F MET BD objective 5x/0.15
M-1181	IOS LWD U-PLAN F MET BD objective 10x/0.30
M-1182	IOS LWD U-PLAN F MET BD objective 20x/0.50
M-1183	IOS LWD U-PLAN F MET BD objective 50x/0.80
M-1184	IOS LWD U-PLAN F MET BD objective 100x/0.90 (dry)
	U-PLAN RPC
M-861	IOS LWD U-PLAN RPC objective 4x/0.13
M-862	IOS LWD U-PLAN RPC objective 10x/0.25
M-863	IOS LWD U-PLAN RPC objective 20x/0.40
M-864	IOS LWD U-PLAN RPC objective 40x/0.65

Condensers & Filters					
M-550	Interferential green filter IF550 (except for IM-5MET)				
<u>M-677ND</u>	Neutral density filter, 25% transmission (only for IM-5MET)				
M-678ND Neutral density filter, 50% transmission (only for IM-5MET)					
	Camera Adapters				
M-113.1	Ring adapter, 30mm (for monocular and binocular microscopes)				
M-115	0.35x C-Mount projection lens				
<u>M-114</u>	0.5x C-Mount projection lens				
M-118	0.75x C-Mount projection lens				
M-173	C-Mount projection lens for APS-C/full frame reflex cameras (trino)				
M-620	0.35x focusable C-Mount adapter (biological microscopes)				
M-620.1	0.5x focusable C-Mount adapter (biological microscopes)				
M-620.2	0.65x focusable C-Mount adapter (biological microscopes)				
M-620.3	1x focusable C-Mount adapter (biological & stereomicroscopes)				
M-699	Universal adapter for C-Mount projection lens (trino)				
Miscellane					
15104	Cleaning kit				
CL-36	Halogen bulb 12V/100W (only for IM-5MET)				
DC-005	TNT dust cover, extra large, 820(l)x550(h) mm				
M-005					
	(1mm/100 & 10mm/100)				
M-641	Adapter for micromanipulator plate (only for IM-5)				
M-793.1	Holder for Petri 38mm diameter (M-793.2 needed)				
=	(except for IM-5MET)				
M-793.2	Holder for Terasaki and Petri 65mm diameter				
	(except for IM-5MET)				
M-793.3	Holder for slides and Petri 54mm diameter				
	(except for IM-5MET)				
M-793.4	Holder for 2+2 slides (except for IM-5MET)				
M-793.5	Holder for metallurgical samples (only for IM-5MET)				
M-793.6	Holder for Utermohl-Chamber (M-793.3 needed)				
14 700 7	(except for IM-5MET)				
M-793.7	Load bearing side extension (except for IM-5MET)				
M-860	Slider with rotating 10x-20x-40x OPTIKA Modulation Contrast slit				
M-860.1	Slider with rotating 4x OPTIKA Modulation Contrast slit				
M-870	DIC slider with Nomarski prism for reflected light				
\/D I\/E	(only for IM-5MET)				
VP-IM5	IQ/ÓQ/PQ manual for IM-5 series				
<u>AB-030</u>	Antibacterial surface treatment, only for newly purchased microscope				

**Condensers & Filters** 



**How to connect the cameras to our microscopes.**Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA® S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA\* Spain spain@optikamicroscopes.com
OPTIKA\* China china@optikamicroscopes.com
OPTIKA\* India india@optikamicroscopes.com

**OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



IM-7



**Inverted Research Microscope** 


IM-7 represents the best of what Optika dedicates to the world of research. This model was created to meet all the needs related to research in life science and designed to be complemented by a series of packages dedicated to more advanced individual applications. For all intents and purposes, IM-7 is to be considered as an inverted imaging platform, due to its high expandability and state-of-the-art quality.



Top-level of optical equipment among our product range provides a sharp and clear view in any situation, while top-level mechanical design offers sturdiness and long lifetime.



### **IM-7**

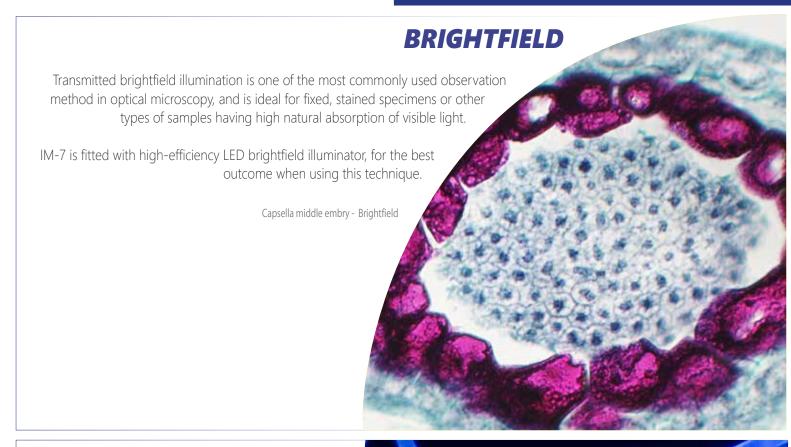

### **FLUORESCENCE**

The fluorescence microscopy is the most demanding technique in biology and biomedical sciences, as well as in materials science.

This method is capable to study organic and inorganic samples thanks to primary fluorescence (auto-fluorescence) or secondary (staining and labelling with fluorochromes)


IM-7 is tailored for applications in research, clinical and pharmaceutical diagnostic field.


Fluorescence illuminators available as mercury lamp.

Cotton fibers - UV Fluorescence



## Multiple Observation Methods

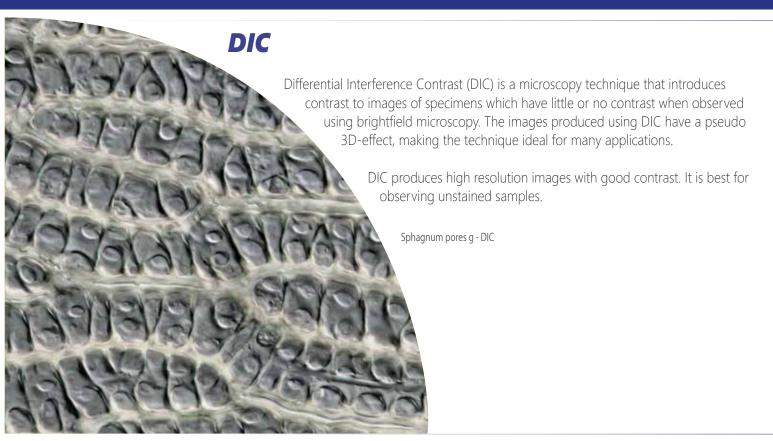

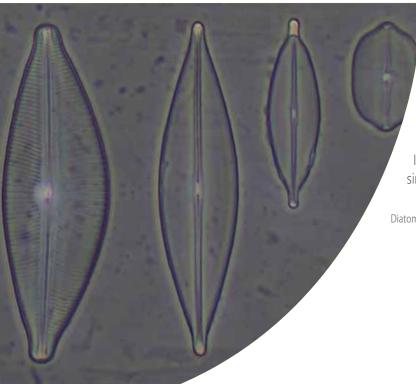

### **PHASE CONTRAST**

Phase-contrast microscopy is a particular technique applied in transparent, non-stainable, samples like culture of living cells, microorganisms, lithographic patterns, latex dispersions, fibers, asbestos and subcellular particles.

It reveals many cellular structures that are not visible with a simple brightfield microscope.

Diatoms - Phase contrast

### **IM-7** - Features



### **Revolving nosepiece for DIC and Fluorescence filter turret**


The six-position nosepiece has a slot (in each of the six positions) for inserting DIC prisms.

The filter turret can hold up to six fluorescence filterblocks. It is easily extractable, in order to facilitate the operation of inserting or replacing the filterblocks.

### Mechanical stage and universal condenser

The wide 3-layer mechanical stage comes with several interchangeable plates for the use of Petri dishes, flasks and slides. The movement of the stage is controlled by a long tilting handle equipped with a pair of knobs for X/Y axes.

The universal condenser is a 6-position type, designed for brightfield, phase contrast and DIC.

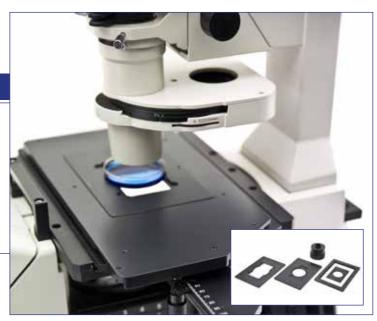



#### **Illuminator arm**

The arm of the transmitted illuminator is backward tilting up to 30 degrees and it allows to use flasks and big bottles.

### **IM-7** - Features



### Main photo tube and Bertrand lens

The main photo tube located on the binocular head is easily controllable by using its control knob. 3 positions selectable: 100/0, 50/50, 0/100.

For phase contrast centering operations a Bertrand lens is available and it can be easily inserted by means of a dedicated knob.

### Fluorescence equipment

A complete package of accessory dedicatd to Fluorescence technique is available as option.







#### Fluorescence filterblocks

In addition to the four standard fluorescence B-G-V-UV filterblocks, many others are available upon request to satisfy every kind of need.

# IM-7 - Features / Controls

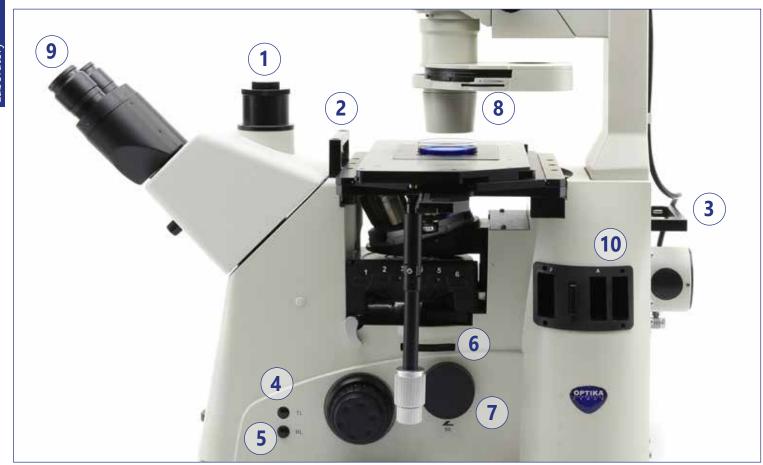

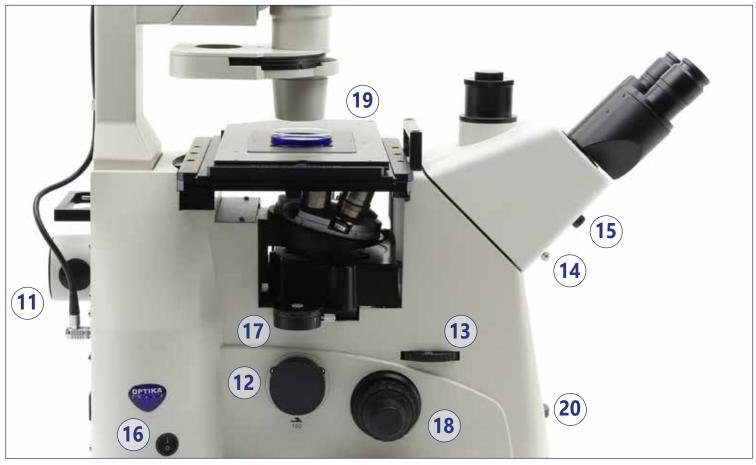

### IM-7 - Features / Controls

- Main photo tube, on binocular head
- **2 -** Handle, for transportation
- **3 -** Handle, for transportation
- **4 -** Transmitted illuminator switch
- **5 -** Fluorescence illuminator switch
- **6 -** Magnification changer, 1x-1.5x
- **7 -** Right side photo tube
- 8 Universal condenser, for brightfield, PH and DIC
- **9 -** Binocular head, with standard WF10x/22mm eyepieces or WF10x/25mm (optional)
- **10 -** Slots for field diaphragm, aperture diaphragm, ND filters

- 11 Input port for fluorescence illuminator
- **12 -** Left side photo tube
- **13 -** Side photo tubes control (100/0; 20/80; 0/100)
- **14 -** Main photo tube control (100/0; 50/50; 0/100)
- 15 Bertrand lens insertion control
- **16** Main switch
- 17 Slot for analyzer
- **18 -** Focusing knobs
- **19 -** 3-layer mechanical stage
- 20 Transmitted illuminator brightness control



# IM-7 - Standard Specifications

Part	Description
Optical system:	Infinity corrected.
Head:	Type: Trinocular (Siedentopf); Inclination 45°; Interpupillary distance 47-78 mm; Tube inner diameter 30 mm; Built-in Bertrand lens.
Photo tubes:	Trino Port/Binocular: 100/0 ; 50/50 ; 0/100; Left Side Port/Binocular: 100/0 ; 0/100.; Right Side Port/Binocular: 100/0; 20/80.
Eyepieces:	PLAN WF10x/25 mm, high eyepoint and with built-in rubber cups; focusable .
Nosepiece:	Sextuple ball bearings revolving nosepiece, reversed; DIC slots fo all positions.
Objectives:	Infinity corrected; 60 mm parfocal distance; All with anti-fungus treatment; Selectable according to customer's preferences (see objective table).
Magnification changer:	1x, 1.5x
Specimen stage:	340x230 mm, three-layers mechanical stage; 130x85 mm moving range; Moving mechanism: Rack and pinion (flexible knob); Holder for Petri dish, 160x110 mm; Holder for Terasaki plate (96 well); Holder for 1 slide. Anti-scratch painting.
Focusing:	Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen; Adjustable tension of coarse focusing knob; 10 mm coarse total travel; Coarse travel per single rotation: 2 mm; Fine travel per single rotation: 0.2 mm; Fine graduations: 100; Fine resolution: 2 µm.
Condenser:	Koehler type, 6-position, focusable, centrable; N.A. 0.55; Iris diaphragm; Phase contrast positions: 10x PH, 20x PH, 40x PH (rings included); DIC positions: DIC1 (10x), DIC2 (20x, 40x) (prisms not included); 26 mm working distance.
Transmitted illumination (Full Koehler):	X-LED <sup>10</sup> ; 10W LED, high efficiency, 6000K; Brightness control; 65.000-hour lifetime.
Incident illumination:	Upgradable to Incident Fluorescence illumination through the purchase of the components indicated in the fluoscence package table.
Differential Interference Contrast (DIC)	Upgradable to DIC through the purchase of the components indicated in the DIC package table.

# IM-7 - Objectives/Application Packages

IM-7 is freely configurable in terms of objectives, by choosing am	ong	g
--	-----	---

Included ■ Optional □

	Objectives 20x, 40x, and 60x feature a correction collar in order to compensate for various thicknesses of cover glasses or different containers.			
M-1320	IOS LWD U-PLAN F (SEMI-APO) PH 4x/0.13, W.D. 16.5 mm, Cover glass -			
M-1321	IOS LWD U-PLAN F (SEMI-APO) PH 10x/0.3, W.D. 7.4 mm, Cover glass 1.2 mm			
M-1322	IOS LWD U-PLAN F (SEMI-APO) PH 20x/0.45, W.D. 7.5 - 8.8 mm, Cover glass 0 - 2 mm			
M-1323	IOS LWD U-PLAN F (SEMI-APO) PH 40x/0.60, W.D. 3.0 – 4.4 mm, Cover glass 0 - 2 mm			
M-1324	IOS LWD U-PLAN F (SEMI-APO) PH 60x/0.70, W.D. 1.8 – 2.6 mm, Cover glass 0.1 - 1.3 mm			

#### **HBO Fluorescence package:**

rescence package.	
EPI Fluorescence internal attachment	
OSRAM 100W HBO high pressure mercury bulb	
HBO Lamp house	
Optika 100W HBO power supply	
6-position fluorescence filterbox turret	
UV protector orange shield	
B filterblock, filters included	
G filterblock, filters included	
V filterblock, filters included	
UV filterblock, filters included	
Aperture diaphragm slider	
Field diaphragm slider	
Slider with neutral density filter for HBO illumination	
Empty fluorescence filterblock	
	EPI Fluorescence internal attachment  OSRAM 100W HBO high pressure mercury bulb  HBO Lamp house  Optika 100W HBO power supply 6-position fluorescence filterbox turret  UV protector orange shield  B filterblock, filters included  G filterblock, filters included  V filterblock, filters included  UV filterblock, filters included  Aperture diaphragm slider  Field diaphragm slider  Slider with neutral density filter for HBO illumination

Name	Excitation filter (nm)	Dichroic cut-off mirror (nm)	Emission filter (nm)
B (Blue)	460 – 490	500	510LP
G (Green)	510 - 550	570	590LP
V (Violet)	400 – 410	455	455LP
UV (Ultraviolet)	330 – 385	400	420LP

#### **DIC - Differential Interference Contrast package:**

	<b>—</b> • · · ·	
Included	<ul><li>Optional</li></ul>	
Included		ш

M-1350	DIC Prism for 10x	
M-1351	DIC Prism for 20x	
M-1352	DIC Prism for 40x	
M-1354	DIC 1 prism 10x for condenser	
M-1355	DIC 2 prism 20x-40x for condenser	
M-1356	Slider with rotating analyzer	

### **PH - Phase Contrast package:**

IIICIUUEU	_	Optional	ш

M-1348	4X Phase ring	
M-1349	60X Phase ring	

### **IM-7** - Accessories

#### Eyecups & Eyepieces M-1360 PL10x/25, high eyepoint, focusable, rubber cup M-1361 PL10x/25, micrometric eyepiece, high eyepoint, focusable, rubber cup **Objectives** IOS U-PLAN F (SEMI-APO) PH 4x/0.13, W.D. 16.5 mm, Cover glass -M-1320 IOS U-PLAN F (SEMI-APO) PH 10x/0.3, W.D. 7.4 mm, Cover glass 1.2 mm IOS U-PLAN F (SEMI-APO) PH 20x/0.45, W.D. 7.5 - 8.8 mm, Cover glass 0 - 2 mm M-1321 M-1322 IOS U-PLAN F (SEMI-APO) PH 40x/0.60, W.D. 3.0 – 4.4 mm, Cover glass 0 - 2 mm M-1323 IOS U-PLAN F (SEMI-APO) PH 60x/0.70, W.D. 1.8 – 2.6 mm, Cover glass 0.1 - 1.3 mm M-1324 Fluorescence package M-1330 EPI Fluorescence internal attachment M-151.1 OSRAM 100W HBO high pressure mercury bulb M-1332 **HBO** lamp house PS-HBO 100W HBO power supply) M-1334 6-position fluorescence filterbox turret M-1335 UV protector orange shield M-1336 B filterblock, filters included M-1337 G filterblock, filters included M-1338 V filterblock, filters included M-1339 UV filterblock, filters included M-1340 Aperture diaphragm slider M-1341 Field diaphragm slider PH package M-1348 4X Phase ring M-1349 60X Phase ring DIC package M-1350 DIC Prism for 10x M-1351 DIC Prism for 20x M-1352 DIC Prism for 40x M-1354 DIC 1 prism 10x for condenser M-1355 DIC 2 prism 20x-40x for condenser M-1356 Slider with rotating analyzer **Camera Adapters** M-620 0.35x focusable C-Mount adapter (for main photo tube only) 0.5x focusable C-Mount adapter (for main photo tube only) M-620.1 M-620.2 0.65x focusable C-Mount adapter (for main photo tube only) M-620.3 1x focusable C-Mount adapter (for main photo tube only) M-1365 0.5x focusable C-Mount adapter (for left/right side photo tube only) 1x focusable C-Mount adapter (for left/right side photo tube only)

#### M-1366 Miscellaneous

Color temperature filter, 38mm M-1370

Green filter, 38mm <u>M-1371</u>

Yellow filter, 38mm M-1372

<u>M-1373</u> Frosted filter, 38mm

<u> 15104</u> <u>Cleaning kit</u>

DC-005 TNT dust cover, extra large, 820(l)x550(h) mm

Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100) M-005

VP-IM7 IQ/OQ/PQ manual for IM-7

<u>AB-040</u> Antibacterial surface treatment, only for newly purchased microscope





How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### Headquarters and Manufacturing Facilities

OPTIKA S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### Optika Sales branches

**OPTIKA**° Spain **OPTIKA**° China **OPTIKA**° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

**OPTIKA**<sup>®</sup> North America **OPTIKA**<sup>®</sup> Central America **OPTIKA**® Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# POL Series



**Routine & Research Lab Polarizing Microscopes** 

### **Polarized Light** Microscopy

Polarized light microscopy is an optical microscopy technique involving polarized light. Simple techniques include illumination of the sample with polarized light. Directly transmitted or incident light can, optionally, be blocked with a polariser orientated at 90 degrees to the illumination.

These illumination techniques are most commonly used on birefringent samples where the polarized light interacts strongly with the sample and so generating contrast with the background. Polarized light microscopy is used extensively in optical mineralogy.

As polarised light passes through a birefringent sample, the phase difference between the fast and slow directions varies with the thickness, and wavelength of light used. The optical path difference (o.p.d.) is defined as

$$o.p.d. = \Delta n x t$$

where t is the thickness of the sample.

This then leads to a phase difference between the light passing in the two vibration directions of

$$\delta = 2 \pi (\Delta n \times t / \lambda)$$

For example, if the optical path difference is  $\lambda$  / 2 , then the phase difference will be  $\pi$  , and so the polarisation will be perpendicular to the original, resulting in all of the light passing through the analyser for crossed polars. If the optical path difference is n x  $\lambda$ , then the phase difference will be 2 n x  $\pi$  , and so the polarisation will be parallel to the original. This means that no light will be able to pass through the analyser which it is now perpendicular to. The Michel-Levy Chart arises when polarised white light is passed through a birefringent sample. If the sample is of uniform thickness, then only one specific wavelength will meet the above condition described above, and be perpendicular to the direction of the analyser. This means that instead of polychromatic light being viewed at the analyser, one specific wavelength will have been removed. This information can be used in a number of ways:

- If the birefringence is known, then the thickness, t, of the sample can be determined
- If the thickness is known, then the birefringence of the sample can be determined


As the order of the optical path difference increases, then it is more likely that more wavelengths of light will be removed from the spectrum. This results in the appearance of the colour being "washed out", and it becomes more difficult to determine the properties of the sample. This, however, only occurs when the sample is relatively thick when compared to the wavelength of light.


# lcons

W-PLAN Objectives

(22mm field of view)
U-PLAN Objectives

(25mm field of view)

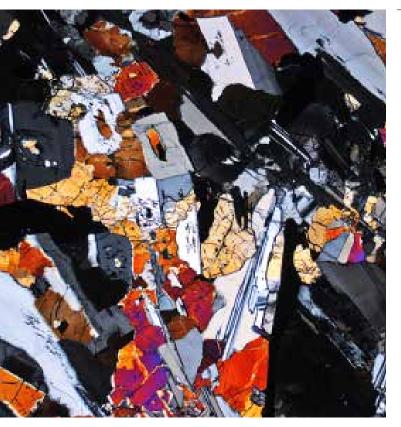


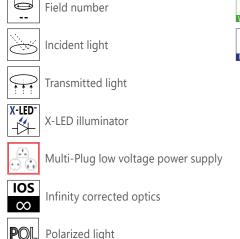

## **B-383POL** - Polarizing Microscope

Upright microscope for brightfield and polarizing light observations with strain-free IOS N-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED<sup>3</sup>** illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.



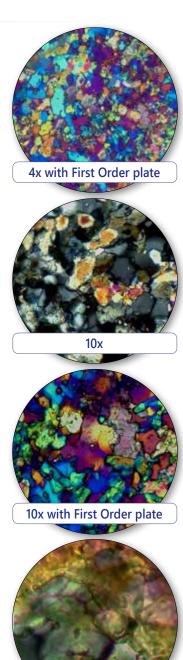

60x


Part	Description
Observation mode:	Brightfield, transmitted polarized light and conoscopy.
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale.  Tint plates: 1° order red (λ); λ/4; Quartz wedge.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 48 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw. One with crosshair.
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings. Centering system for each objective.


Part	Description
Objectives (strain-free):	IOS N-PLAN POL 4x/0.10 IOS N-PLAN POL 10x/0.25 IOS N-PLAN POL 40x/0.65 IOS N-PLAN POL 60x/0.80 All with anti-fungus treatment.
Specimen stage:	Rotatable stage with locking mechanism. Vernier scale with accuracy 0.1 mm. Diameter 160 mm. Specimen slide clamps.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable. With rotating polarizing filter.
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## B-510POL - Polarizing Microscope

Advanced routine laboratory microscope for transmitted light in brightfield and polarized light observations with strain-free IOS W-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED<sup>3</sup>** illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.



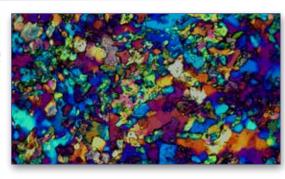

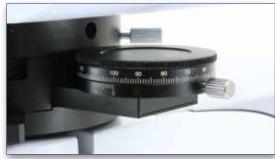



Part	Description
<b>Observation mode:</b>	Brightfield, transmitted polarized light and conoscopy.
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale.  Tint plates: 1° order red (λ); λ/4; Quartz wedge.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups. One with crosshair.
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings. Centering system for each objective.

Part	Description
Objectives (strain-free):	IOS W-PLAN POL 4x/0.10 IOS W-PLAN POL 10x/0.25 IOS W-PLAN POL 20x/0.45 IOS W-PLAN POL 40x/0.65 All with anti-fungus treatment.
Specimen stage:	Rotatable stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm. Specimen slide clamps.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable. With rotating polarizing filter.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-510POL-I** - Polarizing Microscope

Advanced routine laboratory microscope for brightfield and polarized light observations in transmitted and incident light with strain-free IOS LWD W-PLAN POL objectives. Complete of polarizer and analyzer filters, Bertrand lens for conoscopic observation, compensator plates and high-precision rotatable stages. It comes with the exclusive **X-LED**<sup>3</sup> illumination system to deliver bright and clear images, along with all the accessories to perform accurate polarization analysis in biology and materials science.






















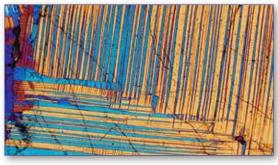





### **Incident/transmitted light Objectives included**

Description

IOS LWD W-PLAN POL 5x/0.12, W.D. 15.5 mm

IOS LWD W-PLAN POL 10x/0.25, W.D. 10.0 mm

IOS LWD W-PLAN POL 20x/0.40, W.D. 5.8 mm

IOS LWD W-PLAN POL 50x/0.75, W.D. 0.32 mm

Part	Description
<b>Observation mode:</b>	Brightfield, transmitted/incident polarized light and conoscopy.
Epi-illumination and filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With polarizer and rotating analyzer for incident illumination, aperture and field diaphragm. With additional filter holder.
Bertrand lens and polarizing attachment:	Swing-out type with centering mechanism for observation in conoscopy/orthoscopy. Rotatable analyzer from 0° to 90° with graduated scale. Tint plates: 1° order red (λ); λ/4; Quartz wedge.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Dioptric adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups. One with crosshair.
Nosepiece:	Quadruple revolving nosepiece, rotation on ball bearings. Centering system for each objective.

Part	Description
Objectives (strain-free):	IOS LWD W-PLAN POL 5x/0.12 IOS LWD W-PLAN POL 10x/0.25 IOS LWD W-PLAN POL 20x/0.40 IOS LWD W-PLAN POL 50x/0.75 All with anti-fungus treatment.
Specimen stage:	Rotatable stage with locking mechanism and centering knobs. Vernier scale with accuracy 0.1 mm. Diameter 145 mm. Specimen slide clamps.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable. With rotating polarizing filter.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-1000POL** - Polarizing Microscope


The modular OPTIKA B-1000 is available with transmitted polarized light, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration.



## **B-1000POL** - Configuration Chart


Build the microscope that suites your needs by choosing among the components

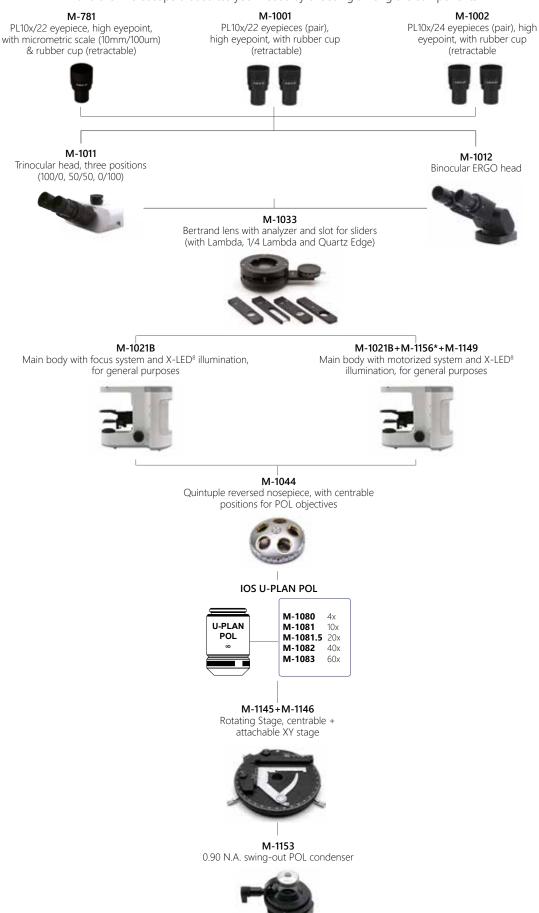

<sup>\*</sup> Code M-1156 must be added only **once** for any motorized configuration


## B-1000POL-I - Polarizing Microscope

The modular OPTIKA B-1000 is available with transmitted and incident polarized light, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive X-LED8 (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

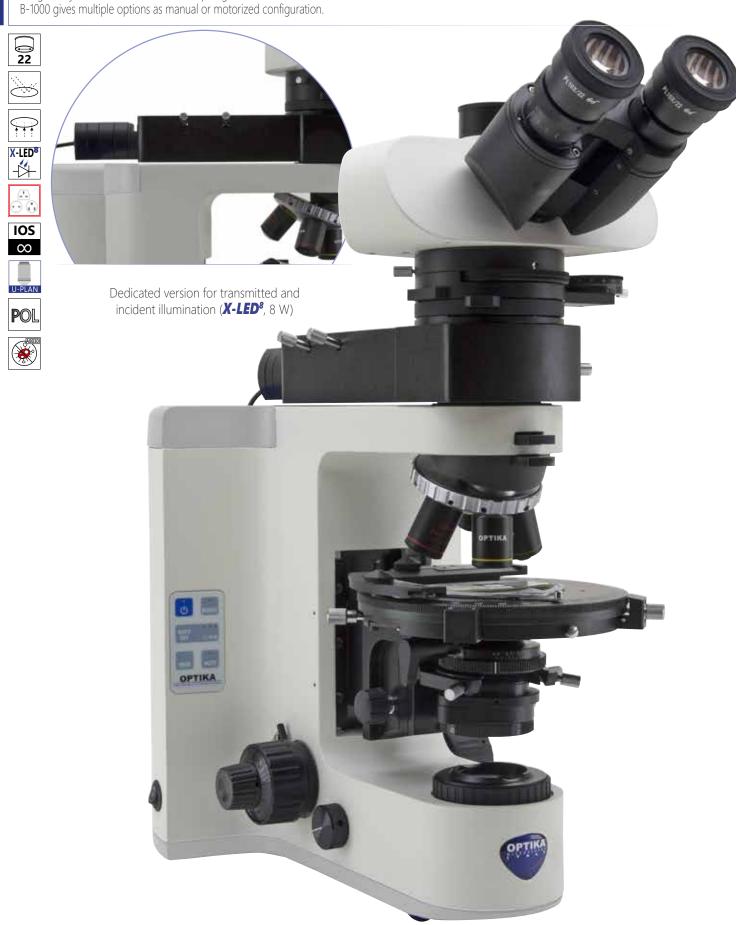

## **B-1000POL-I** - Configuration Chart

Build the microscope that suites your needs by choosing among the components

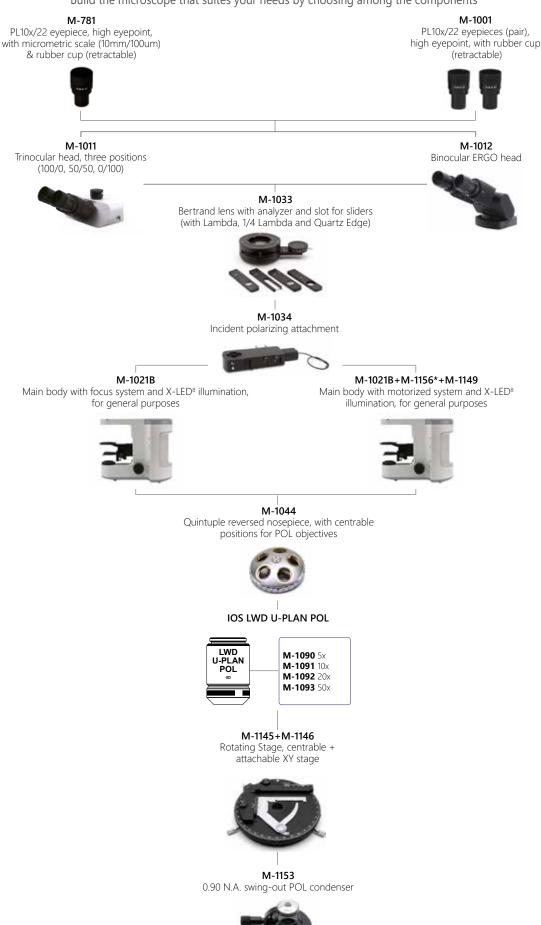

 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com


**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# FLUO Series

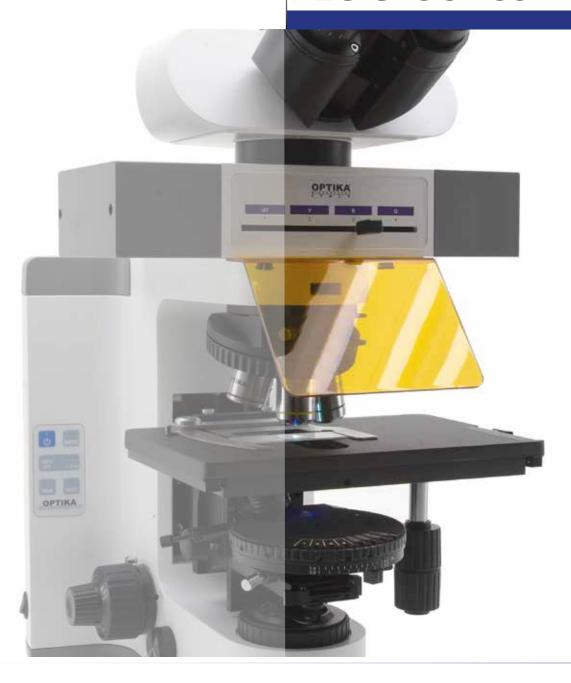

**Routine & Research Lab Fluorescence Microscopes** 

## Fluorescence Microscopy

### **Epi Fluorescence microscopes**


A fluorescence microscope is an optical microscope that uses fluorescence and phosphorescence instead of, or in addition to, reflection and absorption to study properties of organic or inorganic substances. The "fluorescence microscope" refers to any microscope that uses fluorescence to generate an image. The Epi Fluorescence microscope is equipped with a fluorescence illuminator wich generates incident fluorescence light.

### **Principle**

The specimen is illuminated with light of a specific wavelength (or wavelengths) which is absorbed by the fluorophores, causing them to emit light of longer wavelengths (i.e., of a different color than the absorbed light). The illumination light is separated from the much weaker emitted fluorescence through the use of a spectral emission filter. Typical components of a fluorescence microscope are a light source (HBO mercury-vapor lamps are common; more advanced forms are high-power LEDs), the excitation filter, the dichroic mirror, and the emission filter. The filters and the dichroic mirror are chosen to match the spectral excitation and emission characteristics of the fluorophore used to label the specimen. In this manner, the distribution of a single fluorophore (color) is imaged at a time. Multi-color images of several types of fluorophores must be composed by combining several single-color images.

Most fluorescence microscopes in use are epifluorescence microscopes, where excitation of the fluorophore and detection of the fluorescence are done through the same light path (through the objective). These microscopes are widely used in biology and are the basis for more advanced microscope designs.


### lcons

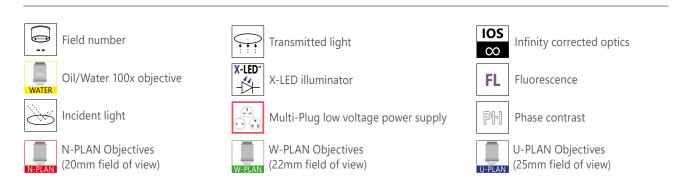

## Fluorescence Microscopy

### **Epifluorescence microscopy**

The majority of fluorescence microscopes, especially those used in the life sciences, are of the epifluorescence design. Light of the excitation wavelength illuminates the specimen through the objective lens. The fluorescence emitted by the specimen is focused to the detector by the same objective that is used for the excitation which for greater resolution will need objective lens with higher numerical aperture. Since most of the excitation light is transmitted through the specimen, only reflected excitatory light reaches the objective together with the emitted light and the epifluorescence method therefore gives a high signal-to-noise ratio. The dichroic beamsplitter acts as a wavelength specific filter, transmitting fluoresced light through to the eyepiece or detector, but reflecting any remaining excitation light back towards the source.

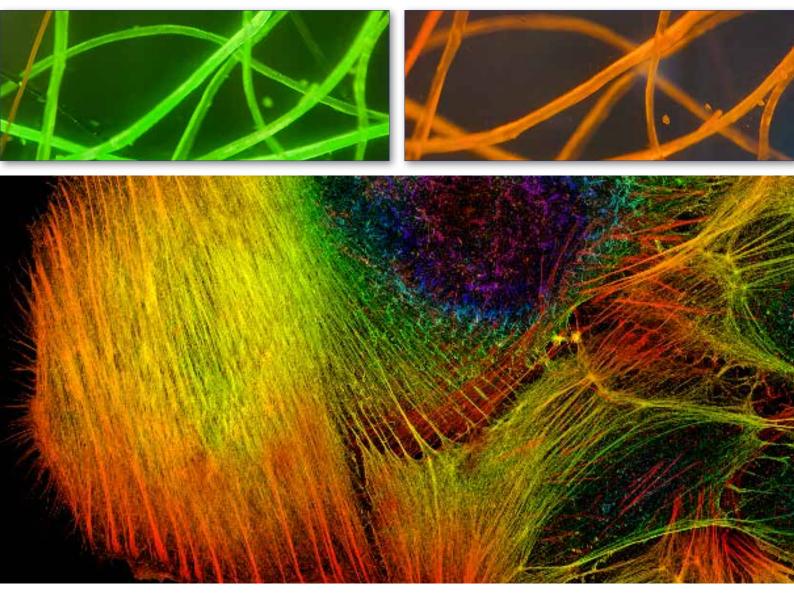

## **B-290LD** - LED Fluorescence Microscopes



Fluorescence binocular and trinocular microscopes especially designed for tubercolosis and malaria analysis.

### Observation mode: Brightfield.

**Head:** Binocular or trinocular, 360° rotating and 30° inclined. Interpupillary distance 48-75mm (bino version) 55-75mm (trino version).

**Dioptric adjustement:** On the left eyepiece tube.

**Eyepieces:** WF10x/20 mm, high eye-point and secured by a screw.

Nosepiece: Quadruple revolving nosepiece, rotation on ball bearings.

**Specimen stage:** Double layer rackless mechanical sliding stage, 150x139 mm, 75x33 mm X-Y movement range. Vernier scale on the two axes, accuracy: 0.1 mm.

**Focusing:** Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.

**Condenser:** Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.

**Brightfield Illumination (Fixed Koehler type):** X-LED<sup>3</sup> with white 3.6 W LED (6,300 K) and light intensity control.

Multi-plug 100-240Vac/6Vdc external power supply.

**Fluorescence Illumination:** Extra efficiency LED, with light intensity control. Peak wavelength: 465 nm, Power: 3.6W.

**Epi Fluorescence Attachment:** Slider with 3 positions (2 fluorescence, 1 brightfield), with 1 included filterset: Fluorescence B: EX 460-490, DM 505, EM 515LP: Acridine Yellow, Acridine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP, etc.

### Part number: B-292LD1.50

Equipped with binocular head and following objectives:

IOS N-PLAN 10x/0.25, with anti-fungus treatment

IOS N-PLAN 20x/0.40, with anti-fungus treatment

IOS N-PLAN 40x/0.65, with anti-fungus treatment

IOS W-PLAN MET 50x/0.75 (No Cover), with anti-fungus treatment.

### Part number: B-293LD1.50

Trinocular version of B-292LD1.50.

### Part number: B-292LD1

Equipped with binocular head and following objectives:

IOS N-PLAN 10x/0.25, with anti-fungus treatment

IOS N-PLAN 20x/0.40, with anti-fungus treatment

IOS N-PLAN 40x/0.65, with anti-fungus treatment

IOS W-PLAN 100x/0.80 (No Cover, Dry), with anti-fungus treatment.

### Part number: B-293LD1

Trinocular version of B-292LD1.

### **Standard filterset**

Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
B (Blue)	460 - 490	505	515LP

## **B-383FL** - HBO Fluorescence Microscope

Laboratory upright microscope for brightfield and fluorescence observations with IOS N-PLAN objectives.

The HBO fluorescence illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Auramine, FITC, GFP and YFP (with blue filter set) plus Rhodamine, Texas Red and TRITC (with the green one), yet giving the possibility to combine any other specific filter sets for future upgrade. Transmitted light through the exclusive **X-LED**<sup>3</sup> to ensure great-looking, rich and high-quality specimen view.




Part	Description		
<b>Observation mode:</b>	Brightfield, HBO fluorescence.		
<b>Epi-illumination and filters:</b>	HBO 100 W high pressure mercury lamp. 3-position filter holder; blue and green included.		
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.		
Interpupillary distance:	Adjustable between 48 and 75 mm.		
<b>Dioptric adjustment:</b>	On the left eyepiece tube.		
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.		
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.		
Objectives:	IOS N-PLAN 4x/0.10 IOS N-PLAN 10x/0.25 IOS N-PLAN 20x/0.40 IOS N-PLAN 40x/0.65 IOS N-PLAN 100x/1.25 (Oil/Water) All with anti-fungus treatment.		

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-383LD** - LED Fluorescence Microscope

Entry-level laboratory upright microscope for brightfield and fluorescence observations with IOS N-PLAN objectives. The extremely powerful LED fluorescence illuminator is combined with blue excitation filter set for the visualization of the following fluorochromes: Accidine Yellow, Accidine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP, etc. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED<sup>3</sup>** to ensure great-looking, rich and high-quality specimen view.

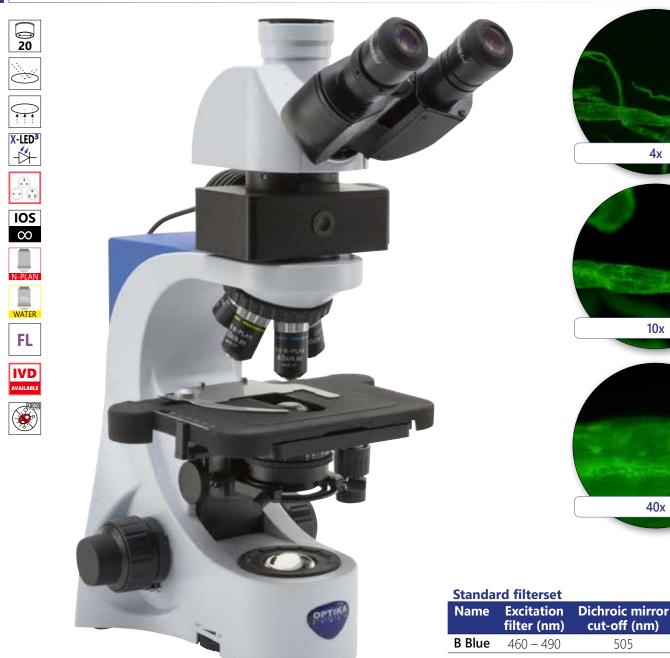

4x
10x
40x
TUA

**Emission** 

filter (nm)


515LP

Part	Description		
<b>Observation mode:</b>	Brightfield, LED fluorescence.		
<b>Epi-illumination and filter:</b>	High-power blue LED with brightness control. 3-position filter holder; blue included.		
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.		
Interpupillary distance:	Adjustable between 48 and 75 mm.		
Dioptric adjustment:	On the left eyepiece tube.		
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.		
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.		
Objectives:	4x/0.10, W.D. 16.8 mm 10x/0.25, W.D. 5.8 mm 20x/0.40, W.D. 5.1 mm 40x/0.65, W.D. 0.43 mm 100x/1.25 (Oil/Water), W.D. 0.13 mm All with anti-fungus treatment.		

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.

## **B-510FL** - HBO Fluorescence Microscope

Advanced routine laboratory microscope for brightfield and fluorescence observations with Semi-Apo IOS W-PLAN F objectives to enhance the visibility of the sample and increase the overall contrast. The **HBO fluorescence** illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Auramine, FITC, GFP and YFP (with blue filter set) plus Rhodamine, Texas Red and TRITC (with the green one), yet giving the possibility to combine any other specific filter sets for future upgrade. Transmitted light through the exclusive **X-LED**<sup>3</sup> to ensure great-looking, rich and high-quality specimen view.

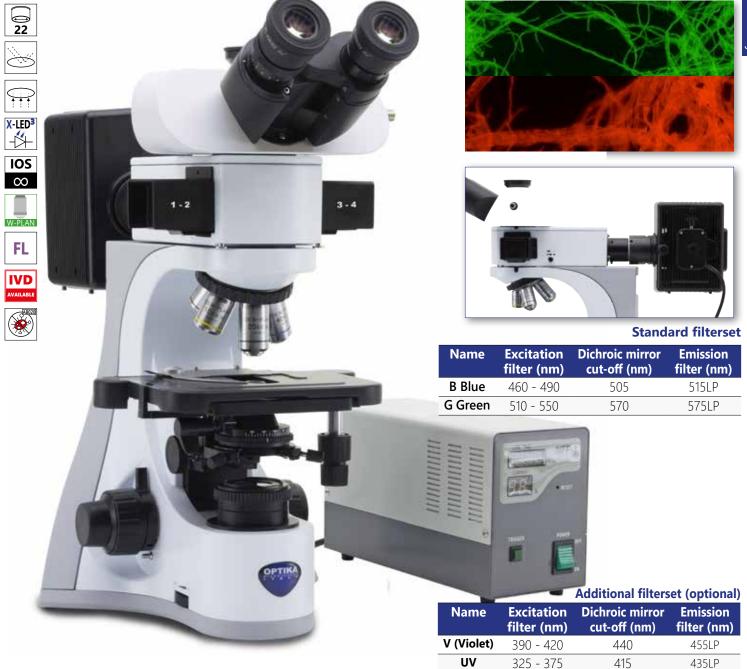

Part	Description			
<b>Observation mode:</b>	Brightfield, HBO fluorescence.			
<b>Epi-illumination and filter:</b>	HBO 100 W high pressure mercury lamp. 4-position filter holder; blue & green included.			
Head:	Trinocular (3-position 100/0, 50/50, 0/100), 30° inclined, 360° rotating.			
Interpupillary distance:	Adjustable between 50 and 75 mm.			
Dioptric adjustment:	On the left eyepiece tube.			
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.			
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.			
Objectives:	IOS W-PLAN F 4x/0.13 IOS W-PLAN F 20x/0.50 All with anti-fungus treatment.			

Part	Description
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.

## **B-510LD4** - LED Fluorescence Microscope

Advanced routine fluorescence microscope for transmitted brightfield and fluorescence observations with IOS W-PLAN objectives. The extremely powerful LED Fluorescence Illuminators are combined with corresponding excitation filter sets for the visualization of most fluorochromes. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED<sup>3</sup>** to ensure great-looking, rich and high-quality specimen view.
































## A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

Part	Description		
<b>Observation mode:</b>	Brightfield, LED fluorescence.		
<b>Epi-illumination and</b>	High-power LED with brightness control.		
filter:	4-position filter holder; none included.		
Head:	Trinocular (3-position 100/0, 50/50, 0/100), 30° inclined,		
	360° rotating.		
Interpupillary distance:	Adjustable between 50 and 75 mm.		
Dioptric adjustment:	On the left eyepiece tube.		
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.		
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.		
Object.: B-510LD4:	IOS W-PLAN 4x/0.10		
B-510LD4-SA:	IOS W-PLAN 40x/0.65 IOS W-PLAN 100x/1.25 (Oil) IOS W-PLAN F 4x/0.13 IOS W-PLAN F 10x/0.30 IOS W-PLAN F 40x/0.75 IOS W-PLAN F 100x/1.3 (Oil) All with anti-fungus treatment.		

Part	Description		
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.		
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.		
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.		
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6,300K) with brightness control.  Multi-plug 100-240Vac/6Vdc external power supply.		

## **B-510LD4** - LED Fluorescence Microscope

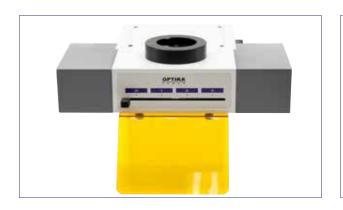

OPTIKA LED Fluorescence attachment is a revolutionary solution

It consists of a 4-position selector for the use of 4 fluorescent illuminators, called LED Fluorescence Cubes.

Each Cube is composed of a filterset mounted on a filterblock and a high power LED with emission corresponding to the filters installed. In this way the selection of each filter controls the lighting up of the corresponding LED.

The microscope is supplied without any LED Fluorescence Cube. A selection of 9 types is available, as shown in the table below.

**B-510LD4**Trinocular brightfield microscope, 1000x, IOS W-PLAN objectives



B-510LD4-SA
Trinocular brightfield microscope, 1000x,
IOS W-PLAN F (Semi-Apo) objectives. For UV Fluorescence



M-1220 ... M-1228 LED Fluorescence Qubes (LED + Filterset)



























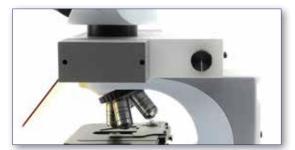

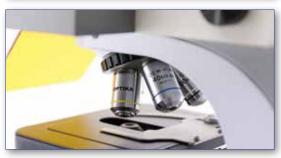

### LED Fluorescence Cubes available (LED + Filterset)

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1220 - Blue	460	455 - 495	500	510LP
M-1220.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1221 - Green	523	510 - 550	570	575LP
M-1221.1 - Green (pass band)	523	510 - 550	570	585-625
M-1222 - Violet	405	390 - 420	440	450LP
M-1223 - UV	365	325 - 375	415	435LP
M-1223.1 - UV (pass band)	365	340 - 390	405	420-470
M-1224 - Red 1 *	623	590 - 650	660	665LP
M-1225 - Red 2 *	623	595 - 645	655	665 - 715
M-1226 - Deep Red *	660	623 - 678	685	690 - 750
M-1227 - Far Red *	740	720 - 760	770	780LP
M-1228 - Amber	590	582 - 603	610	615 - 645

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.

## **B-1000FL-HBO** - HBO Fluorescence Microscope

The modular OPTIKA B-1000 can stand a HBO fluorescence attachment, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Semi-Apo objectives), the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration, with a variety of objectives, stages and condensers.

## **B-1000FL-HBO** - Configuration Chart

Build the microscope that suites your needs by choosing among the components



## **B-1000LD4** - LED Fluorescence Microscope

The modular OPTIKA B-1000 can stand a LED fluorescence attachment, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics (including Apo and Semi-Apo objectives), the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) illumination system, designed by OPTIKA and the Koehler diaphragm.

B-1000 gives multiple options as manual or motorized configuration, with a variety of objectives, stages and condensers.



## A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

### LED Fluorescence Cubes available (LED + Filterset)

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
M-1220 - Blue	460	455 - 495	500	510LP
M-1220.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1221 - Green	523	510 - 550	570	575LP
M-1221.1 - Green (pass band)	523	510 - 550	570	585-625
M-1222 - Violet	405	390 - 420	440	450LP
M-1223 - UV	365	325 - 375	415	435LP
M-1223.1 - UV (pass band)	365	340 - 390	405	420-470
M-1224 - Red 1 *	623	590 - 650	660	665LP
M-1225 - Red 2 *	623	595 - 645	655	665 - 715
M-1226 - Deep Red *	660	623 - 678	685	690 - 750
M-1227 - Far Red *	740	720 - 760	770	780LP
M-1228 - Amber	590	582 - 603	610	615 - 645

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.

## **B-1000LD4** - Configuration Chart

M-1140

M-1150

0.90 N.A. swing-out

condenser

Build the microscope that suites your needs by choosing among the components M-1001 PL10x/22 eyepiece, high eyepoint, PL10x/22 eyepieces (pair), PL15x/16 eyepieces (pair), with micrometric scale (10mm/100um) high eyepoint, with rubber cup high eyepoint & rubber cup (retractable) (retractable) M-1011 M-1012 Trinocular head, three positions Binocular ERGO head (100/0, 50/50, 0/100) M-1031M 4-position LED Fluorescence attachment, without LED Fluorescence Filtersets M-1220 ... M-1228 LED Fluorescence Cubes (LED + Filterset) M-1021B M-1021B+M-1156\*+M-1149 Main body with focus system and X-LED<sup>8</sup> Main body with motorized system and X-LED<sup>8</sup> illumination, for general purposes illumination, for general purposes M-1042 M-1043+M-1156\* Sextuple reversed nosepiece, Sextuple motorized reversed nosepiece, for RMS objectives with DIC slot for RMS objectives with DIC slot IOS U-PLAN F IOS W-PLAN IOS W-PLAN F M-1049 2x M-1125 4x M-1060 4x M-1075 4x M-1126 10x M-1061 10x M-1076 10x W-PLAN W-PLAN F U-PLAN F M-1127 20x M-1062 20x M-1077 20x M-1128 40x M-1063 40x M-1078 40x M-1129 60x M-1064 100x M-1079 100x M-1130 100x M-1141 M-1143 M-1147+M-1156\* M-1144 Motorized mechanical stage MPC (Mineral Solid Surface) Heating stage, with digital Standard mechanical stage Rackless mechanical stage rackless mechanical stage temperature controller M-1151 M-1155

0.9/0.25 N.A. swing-out condenser

(to be used with objective M-1049)

1.2 N.A. swing-out condenser


Darkfield condenser

for dry objectives

<sup>\*</sup> Code M-1156 must be added only **once** for any motorized configuration

## Fluorescence Filtersets

Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
B (Blue)	460 – 490	505	515LP
G (Green)	510 – 550	570	575LP
V (Violet) optional	385 – 425	440	455LP
UV (Ultraviolet) <b>optional</b>	325 – 375	415	435LP

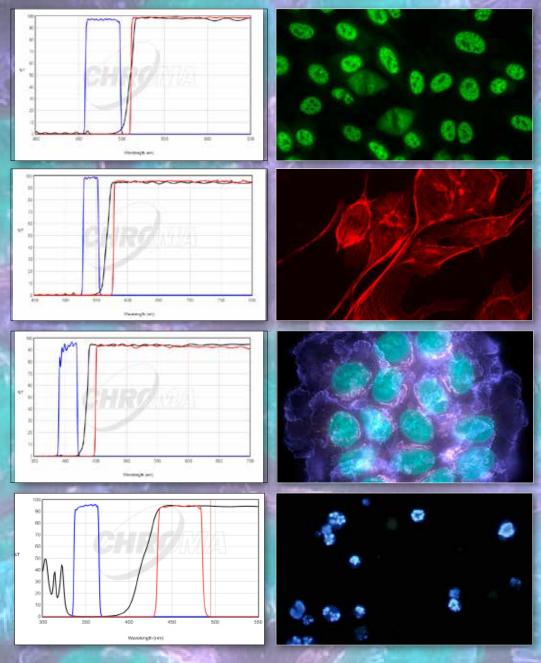

MANY MORE FILTERSETS AVAILABLE ON REQUEST

CHROMA TECHNOLOGY CORP

Specs are of B-510FL Filtersets

## **IM-300F** - HBO Fluorescence Microscope

Routine inverted fluorescence microscope for transmitted brightfield, phase contrast and fluorescence observations with IOS LWD W-PLAN objectives. The HBO fluorescence illuminator is combined with blue and green excitation filter set for the visualization of the following fluorochromes: Ethidium Bromide, Acridine Orange, Alexa Fluor 488, Fluo-4, FITC Plus, Spectrum Green, Ethidium Homodimer I, Propidium Iodide (PI) and Spectrum Gold. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



Part	Description
<b>Observation mode:</b>	Brightfield, phase contrast, HBO fluorescence.
<b>Epi-illumination and filter:</b>	HBO 100 W high pressure mercury lamp. 3-position filter holder; blue & green included.
Head:	Trinocular (2-position 100/0, 0/100), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Diopter adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40X/0.65 All with anti-fungus treatment.

Part	Description
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield. Supplied with blue (LBD) and green filter.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/6Vdc external power supply.

## IM-300FL4 - HBO Fluorescence Microscope

Advanced inverted microscope for brightfield and fluorescence observations with Semi-Apo IOS LWD U-PLAN F objectives to enhance the visibility of the sample and increase the overall contrast. The HBO fluorescence illuminator provides an outstanding flexibility of use, standing the blue and green filter sets (supplied as standard) for Ethidium Bromide, Acridine Orange, Alexa Fluor 488, Fluo-4, FITC Plus, Spectrum Green, Ethidium Homodimer I, Propidium lodide (PI) and Spectrum Gold, with additional filterset: Alexa Fluor 350, Aminocoumarin, Hoechst, DAPI, BFP (Blue Fluorescent Protein), Atto 390, Tetracycline, Pacific Orange and Spectrum Blue. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



Part	Description	
Observation mode:	Brightfield, HBO fluorescence.	
Epi-illumination and filter:	HBO 100 W high pressure mercury lamp. 4-position filter holder; blue & green included.	
Head:	Trinocular (2-position 100/0, 0/100), 45° inclined.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
<b>Diopter adjustment:</b>	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.	


Part	Description
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. Supplied with blue (LBD) filter.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm.  Multi-plug 100-240Vac/6Vdc external power supply.

## IM-300LD2 - LED Fluorescence Microscope

Routine inverted fluorescence microscope for transmitted brightfield, phase contrast and fluorescence observations with IOS LWD W-PLAN PH objectives.

The LED fluorescence illuminators are combined with blue and green excitation filter set for the visualization of the following fluorochromes: GFP, Alexa Fluor 488, Calcein, SYBR Green, FITC, Fluo-4, MitoTracker Green, Spectrum Gold, Propidium Iodide, Ethidium Homodimer I. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



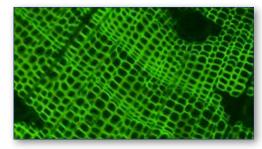



## A new milestone achieved in Fluorescence Microscopy

- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for routine applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

### **LED Fluorescence Cubes (LED + Filterset) included**

Name	LED emission (nm)	Excita- tion filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
Blue	460	455 - 495	500	510LP
Green	523	510 - 550	570	575LP

Part	Description	
<b>Observation mode:</b>	Brightfield, phase contrast, LED fluorescence.	
<b>Epi-illumination and filter:</b>	High-power LED with brightness control. 3-position filter holder; blue and green filtesets included.	
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
<b>Diopter adjustment:</b>	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS LWD W-PLAN PH 10x/0.25 IOS LWD W-PLAN PH 20x/0.40 IOS LWD W-PLAN PH 40x/0.65 All with anti-fungus treatment.	

Part	Description
Specimen stage:	Fixed stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. With 4x/10x, 20x/40x phase contrast slider and brightfield.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. 100-240Vac/12Vdc external power supply.

## IM-300LD4 - LED Fluorescence Microscope

Advanced fluorescence inverted microscope for transmitted brightfield and fluorescence observations with IOS LWD U-PLAN F objectives.

The extremely powerful LED Fluorescence Illuminators are combined with corresponding excitation filter sets for the visualization of most fluorochromes. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Transmitted light through the exclusive **X-LED**<sup>8</sup> to ensure great-looking, rich and high-quality specimen view.



## A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

LED Fluorescence Cubes available (LED + Filterset)

Name	LED emission (nm)	Excitation filter (nm)	Dichroic mirror cut-off (nm)	filter (nm)
M-1230 - Blue	460	455 - 495	500	510LP
M-1230.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1231 - Green	523	510 - 550	570	575LP
M-1231.1 - Green (pass band)	523	510 - 550	570	585-625
M-1232 - Violet	405	390 - 420	440	450LP
M-1233 - UV	365	325 - 375	415	435LP
M-1233.1 - UV (pass band)	365	340 - 390	405	420-470
M-1234 - Red 1 *	623	590 - 650	660	665LP
M-1235 - Red 2 *	623	595 - 645	655	665 - 715
M-1236 - Deep Red *	660	623 - 678	685	690 - 750
M-1237 - Far Red *	740	720 - 760	770	780LP
M-1238 - Amber	590	582 - 603	610	615 - 645

<sup>\*</sup> If the use of a camera is needed, when used for red emission fluorescence wavelengths above 650nm, please order OPTIKA camera code C-P6AR.

Part	Description	
Observation mode:	Brightfield, LED fluorescence.	
<b>Epi-illumination and filter:</b>	High -Power LED with brightness control. 4-position filter holder; none included.	
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
<b>Diopter adjustment:</b>	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.	

Part	Description
Specimen stage:	Mechanical stage, 250x290 mm, with round glass and metal stage, inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. Supplied with blue (LBD) filter.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/12Vdc external power supply.

## IM-300LD4D - LED Fluorescence Microscope


























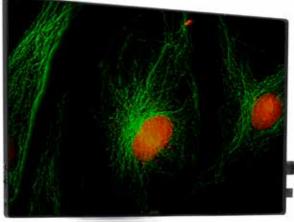

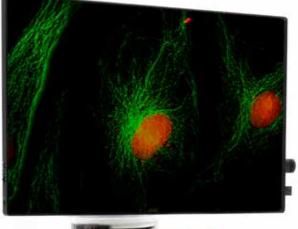

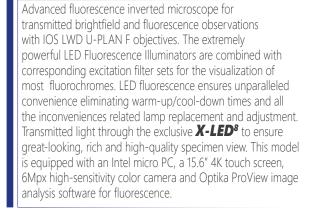

### A new milestone achieved in Fluorescence Microscopy

- » Full-modular Fluorescence System
- » Interchangeable LED-Filtersets
- » 4 LED-Filtersets slots
- » 10% higher light intensity than HBO
- » 35% higher light intensity than Metal-Halide
- » Adjustable light intensity
- » The selection of filtersets automatically involves the switching on of the corresponding LEDs
- » Recommended for research applications
- » Cost-effective, money saving technology
- » Ready for immediate operation
- » Eliminate warm-up/cool-down times
- » Forget lamp replacement & centering

Name	emission (nm)	filter (nm)	cut-off (nm)	filter (nm)
M-1230 - Blue	460	455 - 495	500	510LP
M-1230.1 - Blue (pass band)	460	455 - 495	500	518-542
M-1231 - Green	523	510 - 550	570	575LP
M-1231.1 - Green (pass band)	523	510 - 550	570	585-625
M-1232 - Violet	405	390 - 420	440	450LP
M-1233 - UV	365	325 - 375	415	435LP
M-1233.1 - UV (pass band)	365	340 - 390	405	420-470
M-1234 - Red 1	623	590 - 650	660	665LP
M-1235 - Red 2	623	595 - 645	655	665 - 715
M-1236 - Deep Red	660	623 - 678	685	690 - 750
M-1237 - Far Red	740	720 - 760	770	780LP
M-1238 - Amber	590	582 - 603	610	615 - 645

Part	Description
Observation mode:	Brightfield, LED fluorescence.
Epi-illumination and filter:	High -Power LED with brightness control. 4-position filter holder; none included.
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
<b>Diopter adjustment:</b>	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.
Objectives:	IOS LWD U-PLAN F 10x/0.30 IOS LWD U-PLAN F 20x/0.45 IOS LWD U-PLAN F 40x/0.65 All with anti-fungus treatment.

Part	Description
Specimen stage:	Mechanical stage, 250x290 mm, with round glass and metal stage inserts for slides and 54mm dia. Petri dishes.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.
Condenser:	LWD pre-centered condenser, N.A. 0.30, W.D. 72 mm. Supplied with blue (LBD) filter.
Transmitted illumination:	X-LED <sup>8</sup> with white 8 W LED (6,300K) with brightness control. With aperture diaphragm. Multi-plug 100-240Vac/12Vdc external power supply.
Digital equipment:	Intel micro PC with Image analysis software for Fluorescence. 15.6" 4K touch screen; 6MP IR sensitive high-sensitivity color camera.

## IM-5FLD - LED Fluorescence Microscope

Phase contrast, brightfield and darkfield (dry) LED fluorescence trinocular inverted microscope, with freely configurable lenses according to customer's preferences, FN 24 high eyepoint, infinity corrected optical system, coaxial focusing, mechanical stage, Abbe condenser and powerful, uniform, white color temperature 8 W *X-LED8*. The 4-position epi-fluorescence attachment is powered by extremely powerful 5 W LEDs fluorescence illuminator and combined with blue, green and UV excitation filters for the visualization of the following fluorochromes: Acridine Yellow, Acridine Orange, Auramine, DiO, DTAF, FITC, GFP, YFP (blue filter) plus Rhodamine, Texas Red and TRITC (green filter) plus Alexa Fluor® 350, 7- Amino-4-methylcoumarin, 6-Aminoquinoline, Calcofluor® White, Dansyl cadaverine, DAPI, Dapoxyl, DIDS, Europium (III) Chloride, Fluoro-Gold™, Fura-2, Hoechst 33342 & 33258, 1,5 IAEDANS, Indo-1, Marina Blue®, 4-Methylumbelliferone, PBF1, Pyrene, SBFI, Y66F, Y66H (UV filter) among the others. LED fluorescence ensures unparalleled convenience eliminating warm-up/cool-down times and all the inconveniences related lamp replacement and adjustment. Sturdy and incredibly reliable, it is equipped with all the main controls in ergonomic position and with long lasting, efficient LED illumination to provide over 20 years of use.



# **IM-5FLD** - Specifications



Part	Description
Head:	Trinocular (split ratio: 100/0, 0/100), 45° inclined.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	Both eyepieces.
Eyepieces:	WF10x/24 mm, high eyepoint and with retractable rubber cups.
Epi-fluorescence illumination & filters:	High-power 5 W LEDs with brightness control, motorized LED selection with centrable field diaphragm, 4-position filter holder; blue (EX 450-490, DM 495, EM 500-550), green (EX 540-580, DM 585, EM 608-682) and UV (EX 340-390, DM 400, EM 420LP) excitation filters included.
Nosepiece:	Quintuple ball bearings revolving nosepiece, reversed.
Objectives:	Selectable according to customer's preferences. All with anti-fungus treatment.
Specimen stage:	Fixed stage, 215x250 mm and attachable mechanical stage, 290x250 mm, 120x80 mm X-Y movement range.
Focusing:	Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.
Condenser:	Abbe N.A. 0.50, removable, with iris diaphragm and slider for phase contrast.
Transmitted illumination (Full Koehler):	X-LED <sup>8</sup> with white 8 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power supply.

### **Fluorescence filtersets**

Name	Excitation filter (nm)	Dichroic mirror cut-off (nm)	Emission filter (nm)
B (Blue)	450 – 490	495	500 - 550
G (Green)	540 – 580	585	607 - 682
UV (Ultraviolet)	340 -390	400	420LP

### IM-5FLD is freely configurable in terms of objectives, by choosing among:

Included ■ Optional □

Infinity-corrected Plan-Achromatic, Long Working Distance objectives field flatness up to F.N. 22:		ectives,
M-782	IOS LWD W-PLAN objective 4x/0.13	
M-773	IOS LWD W-PLAN objective 40x/0.60	
M-786	IOS LWD W-PLAN objective 60x/0.70	

Positive Phase Contrast Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 22:		
M-782.1	IOS LWD W-PLAN PH objective 4x/0.13	
M-783N	IOS LWD W-PLAN PH objective 10x/0.25	
M-784N	IOS LWD W-PLAN PH objective 20x/0.40	
M-785	IOS LWD W-PLAN PH objective 40x/0.65	

	rected Semi-Apochromatic, Long Working Distance field flatness up to F.N. 25:	
M-800	IOS LWD U-PLAN F objective 4x/0.13	
M-801	IOS LWD U-PLAN F objective 10x/0.30	
M-802	IOS LWD U-PLAN F objective 20x/0.45	
M-803	IOS LWD U-PLAN F objective 40x/0.65	
M-804	IOS LWD U-PLAN F objective 60x/0.75	

Positive Phase Contrast Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1177	IOS LWD U-PLAN F PH objective 20x/0.45	
M-1178	IOS LWD U-PLAN F PHobjective 40x/0.65	

 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# SZR-180



**Research Stereomicroscope** 

### SZR-180- Features



### 18:1 Zoom Ratio

The SZR-180 with its 18:1 zoom ratio puts it at the top of its class. The wide range of magnifications available (7.5x-135x), coupled with a precise click-stop mechanism for working with reproducible zoom settings, makes it ideal for both low-magnification observation and high-magnification screening and observation of small cell structures.

### **PLAN APO 1x objective**

The Plan APO 1x apochromatic main objective provides high-resolution images at both low and high magnifications, making it ideal for observing and capturing images of the smallest details free of aberrations and halos.

Combined with the zoom and the generous 23mm field number eyepieces, this provides a large FOV even at high magnifications (1.7mm).



### **Ergonomic stand with diascopic LED lighting**

Ergonomically designed, the SZR-180's wide base offers a huge space for sample holders, intuitive controls and an extremely low profile for easy access and reduced fatigue-free body movement. The LED diascopic illumination is designed to achieve high overall contrast and includes the innovative OIC (Oblique Illumination Contrast) system, ideal for transparent structures, and dark-field illumination.


## **SZR-180**- Specifications


The SZR-180 is the top-of-the-range stereo zoom microscope dedicated to the world of research.

With a wide magnification range and high optical resolution combined with effective ergonomics, the SZR-180 is the perfect companion in advanced research. The 18:1 zoom allows you to observe a wide range of samples, from single cells to entire microorganisms.





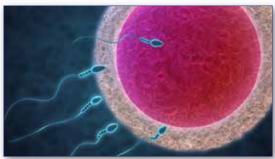

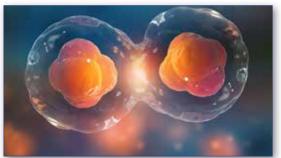

Part	Description
Head:	Trinocular 100/0 - 0/100 (on right tube), 20° inclined.
Interpupillary distance:	52-112 mm.
Dioptric adjustment:	Both eyepieces.
Eyepieces:	WF10x/23 mm, high eyepoint with rubber cups.
Objective:	PLAN Apochromatic 1x, N.A. 0.15, 7.5x-135x (zoom ratio 18:1), zoom click stops at 0.75x, 1x, 2x, 3x, 6x, 10x, 13.5x with aperture diaphragm.
Working distance:	60 mm.
Stand:	High-grade, precision fixed with handle and focus.
Focusing:	Coaxial coarse & fine, coarse total travel 125 mm, fine single rotation 2 mm. Rack and pinion focusing mechanism.
Stage:	Glass transparent plate for transmitted illumination diameter 180 mm.
Illumination:	LED transmitted, with OIC (Oblique Illumination Contrast) and darkfield illumination. Multi-plug 100-240Vac/12Vdc external power supply.

## **SZR-180** - Optical Performance

Eyepiece	10x	
Field number (mm)	23	
Objective	Total magnification	Field of View (mm)
1x (W.D: 60 mm)	7.5x-135x	30.60 - 1.70







### **SZR-180** - Accessories

### Miscellaneous

15104 Cleaning kit

DC-004 TNT dust cover, extra large, 700(l)x550(h) mm

M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)

AB-030 Antibacterial surface treatment, only for newly purchased microscope





### **How to connect the cameras to our microscopes.**Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com


### **Optika Sales branches**

OPTIKA Spain
OPTIKA China
OPTIKA India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



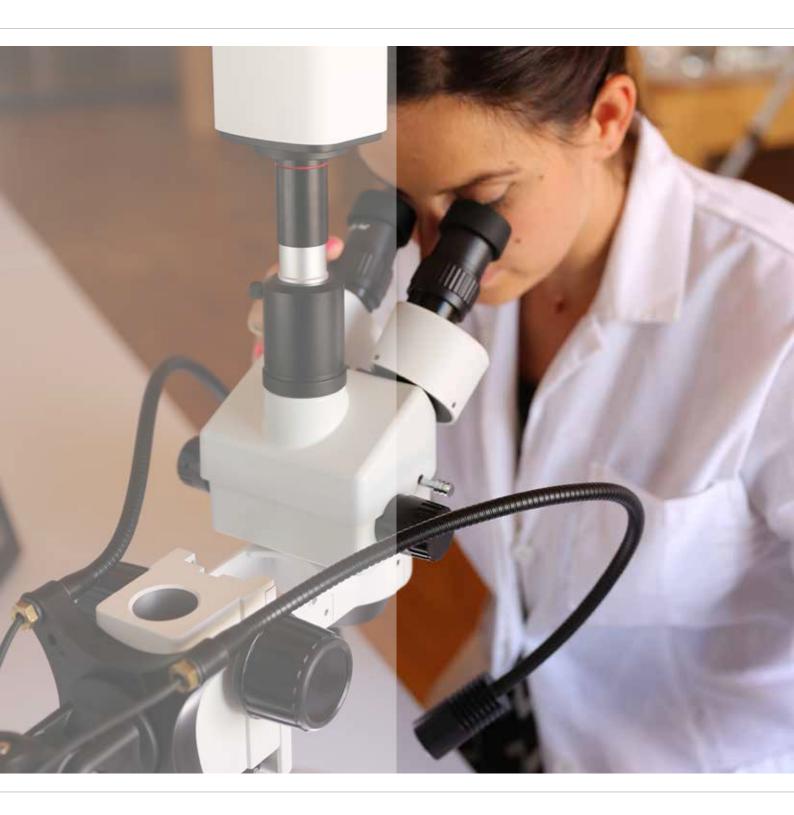


INSPECTION & INDUSTRIAL Microscopes

## **INSPECTION & INDUSTRIAL** Microscopes

Professional Stereo Microscopes	
<b>SLX SERIES</b> - Stereomicroscopes For Higher Education & Laboratory	page 281
SZ SERIES - Stereozoom Microscopes For Laboratory & Industry	page 291
SZP SERIES - Advanced CMO Stereozoom Microscopes	page 307
CL SERIES - Cold Light Illuminators	page 315
GEM SERIES - Stereozoom Microscopes For Gemology	page 319
Industrial Microscopes	
MET SERIES - Metallurgical Microscopes	page 329
IS SERIES - Inspection Video Microscopes	page 343

## Icons

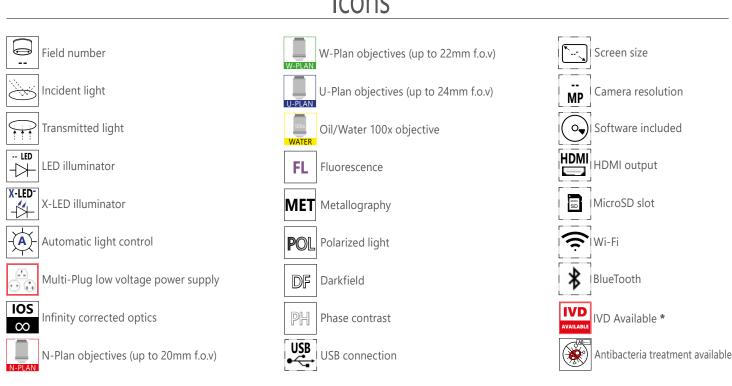



# SLX Series



**Stereomicroscopes For Higher Education & Laboratory** 

## Extremely Versatile Cordless Stereo & Stereozoom Microscopes

### PROFESSIONAL FEATURES FOR... WELL, EVERYONE


- » Level up skills and become a professional user
- » 3D Greenough view for high resoluted images & large field depth
- » 6.43:1 ratio 7x ... 45x or turnable objective 2x, 4x on 21 mm
- » Compact, practical and intuitive to use
- » Sturdy and durable for extended lifetime

### THE LONGEST AUTONOMY ON THE MARKET

- » Longlife LED illumination (providing over 20 years of use)
- » Ultra-flat base with Ø 100 mm disc for diffused transmitted light
- » Cordless use, totally independent from mains/batteries connection
- » Freely settable illumination incident, oblique and transmitted light
- » External power supply for enhanced safety and convenient servicing



## **SLX** Series

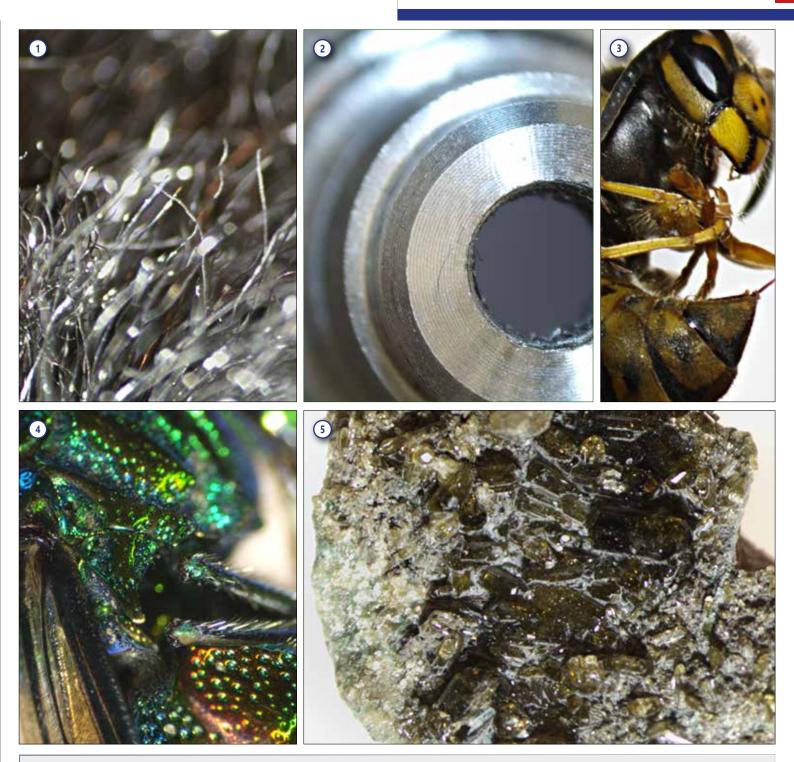

### Legend

- Aluminum SLX-1 and 4x objective.
   Component worked on lathe SLX-2 and 3x zoom.
- 3. Wasp SLX-3 and 4x zoom.

- 4. Fly, detail SLX-2 and 4.5x zoom.
- 5. Rock SLX-2 1.5x zoom.

3

## **SLX** Series

Valuable configurations of cordless and modern stereo & stereozoom microscopes ideal for a variety of applications, including dissection, biology, entomology, anatomy, chemistry, material science among the others and even industrial purposes.

Provided with dual magnification or 6.43:1 zoom ratio, FN 21 high eyepoint eyepieces, highgrade precise fixed arm with focus and handle with the latest technology of **EcoLED™** illumination plus rechargeable batteries. Slim and easy to carry, all the models with high-grade precise fixed arm are equipped with long lasting **LED** illumination to provide over 20 years of use.

### High eyepoint eyepieces for glasses wearers

These eyepieces are designed in such a way that the exit pupil is further away from the eye lens than standard eyepieces, being are well suited for eyeglasses wearers

### The longest autonomy on the market ensured by EcoLED™

OPTIKA has re-designed illumination in microscopy, once again: a special coating process on optics combined with a new, higher ratio between low consumptions and ultra-efficiency has addressed us to top brightness levels

### 6.43:1 zoom ratio - zoom magnification from 7x to 45x

Purposely designed for professional routine inspections, the total magnification can be even extended to 135x with 20x eyepieces and 1.5x additional lens, obtaining an excellent results in this class



### Ultra-flat base with Ø 100 mm disc for diffused transmitted light

A new level of ergonomy and comfort is achieved during operations, with the ultra-flat base of only 3 cm height to ensure smooth specimen movement and the Ø 100 mm for top class diffusion of the transmitted light

## Stereomicroscopes For Higher Education & Laboratory





### Longlife LED illumination (providing over 20 years of use)

Money & energy saving thanks to LED long lifetime (65.000 hours, 22 years in case of 8 hours/day) which is more than 20 times compared to a standard halogen bulb.

### Cordless use, totally independent from mains/batteries connection

All models work with or without the batteries in place with three NiMH rechargeable batteries (not included) for the longest autonomy in outdoor use (12-hour autonomy, at medium intensity).





### External power supply for enhanced safety and convenient servicing

OPTIKA's safety first approach drives to the use of a low voltage, multi-plug, external power supply in order to prevent any risk of electric shock and heatflow inside the unit.

### **SLX** Series - Get the most out of our accessories

### **Additional Lenses**

Simply to be screwed into the threads below the objectives of SLX-2 and SLX-3 to either increase or decrease total magnification, or to increase the working distance when users need to work with hands under the microscope



### ST-040.1 - Darkfield condenser

This is a darkfield condenser for stereo microscopes with bottom light and 100 mm round working plate to provide darkfield microscopy features, fitting all OPTIKA stereomicroscopes with 100 mm mounting size and transmitted light



## **SLX** Series - Range

### SLX-1

















Cordless binocular stereomicroscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with turnable objective (2x-4x), FN 21 high eyepoint, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Binocular, 45° inclined; 360° rotating.

Dioptric adjustment: Left eyepiece.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups

**Objective:** Achromatic 2x-4x with anti-fungus treatment.

Working distance: 100 mm.

Stand: High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination: EcoLED™** swiveling incident and transmitted, with

brightness control, rechargeable batteries (not provided).

Color temperature: 6,300 K.

Multi-plug 100-240Vac/5Vdc external power supply.

### SLX-2

















Cordless binocular stereozoom microscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with 0.7x...4.5x zoom, FN 21 high eyepoint, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Binocular, 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups

**Objective:** Parfocal achromatic zoom 0.7x...4.5x (6.43:1 ratio) with anti-fungus treatment.

Working distance: 100 mm.

**Stand:** High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.

**Illumination: EcoLED™** swiveling incident and transmitted, with brightness control, rechargeable batteries (not provided).

Color temperature: 6,300 K.

Multi-plug 100-240Vac/5Vdc external power supply.

## **SLX** Series - Range

### SLX-3



Cordless trinocular stereozoom microscope ideal for students, schools and amateurs to dissect and discover mainly biology and materials science in 3D, with 0.7x...4.5x zoom, FN 21 high eyepoint, precision fixed arm with handle and the latest technology of EcoLED™ illumination plus rechargeable batteries (not provided). Slim and easy to carry, it is equipped with long lasting LED illumination to provide over 20 years of use.

Head: Trinocular (split ratio: 50/50), 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups.

**Objective:** Parfocal achromatic zoom 0.7x...4.5x (6.43:1 ratio) with anti-fungus treatment.

Working distance: 100 mm

Stand: High-grade, precision fixed with handle and focus.

Focusing: Rack and pinion focusing mechanism.


**Illumination: EcoLED™** swiveling incident and transmitted, with brightness control, rechargeable batteries (not provided). Color temperature: 6,300 K. Multi-plug 100-240Vac/5Vdc external power supply.





## **SLX** Series - Range

## **SLX-4 / SLX-5**

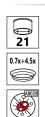





SLX-5

Binocular (SLX-4) or trinocular (SLX-5) stereomicroscope with Extremely stable, long overhanging stand for observation of particularly large specimens to be matched with a focusing system and stereozoom head. Smooth horizontal and vertical movements are ensured and the head can be easily swivelled for inspection at oblique angles. In case illumination is needed, choose from the wide choice of external illuminators available.

### Head:

SLX-4: Binocular, 45° inclined; 360° rotating.

SLX-5: Trinocular (split ratio: 50/50), 45° inclined; 360° rotating.

**Dioptric adjustment:** Both eyepieces.

**Eyepieces:** WF10x/21 mm, high eyepoint, secured by screw and with rubber cups.

**Objective:** Parfocal achromatic zoom 0.7x...4.5x (6.43:1 ratio) with anti-fungus treatment.

Working distance: 100 mm.

**Focusing:** Rack and pinion mechanism controlled by a pair of knobs placed on both sides of the stand.

**Stand:** Overhanging stand 420 mm high, 430 mm horizontal arm length, base 230x230 mm with following adjustments: heigh, longitudinal extension and head rotation angle (left-right).

# **SLX** Series - Comparison Chart

Model	Head	Eyepieces	Objective	Working Distance	Stand	Illumination
SLX-1	Binocular 45° inclined 360° rotating	WF 10x/21	2x – 4x selectable	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)
SLX-2	Binocular 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)
SLX-3	Trinocular (50/50) 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	High-grade, precision fixed with handle and focus	<b>EcoLED™</b> swiveling incident and transmitted with brightness control, rechargeable batteries (not provided)
SLX-4	Binocular 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	Overhanging stand	Without illumination; External light source needed
SLX-5	Trinocular (50/50) 45° inclined 360° rotating	WF 10x/21	0.7x4.5x zoom	100 mm	Overhanging stand	Without illumination; External light source needed

### **Optical performance SLX-1**

Eyepiece	10x (S	10x (ST-081)		15x (ST-082)		20x (ST-083)		10x (ST-084)	
Field number (mm)	2	21		15		10		21	
Additional lens	Total magnification	Field of View (mm)							
1x	20x - 40x	10.50 - 5.25	30x - 60x	7.50 - 3.75	40x - 80x	5.00 - 2.50	20x - 40x	10.50 - 5.25	

### Optical performance SLX-2 - SLX-3 - SLX-4 - SLX-5

Eyepiece	10x (ST-081)		15x (ST-082)		20x (ST-083)		10x (ST-084)	
Field number (mm)	21		15		10		21	
Additional lens	Total magnification	Field of View (mm)						
0.5x	3.5x - 22.5x	60.00 - 9.33	5.25x - 33.75x	42.86 - 6.67	7x - 45x	28.57 - 4.44	3.5x - 22.5x	60.00 - 9.33
0.75x	5,25x - 33.75x	40.00 - 6.22	7.875x - 50.625x	28.57 - 4.44	10.5x - 67.5x	19.05 - 2.96	5.25x - 33.75x	40.00 - 6.22
1x	7x - 45x	30.00 - 4.67	10.5x - 67.5x	21.43 - 3.33	14x - 90x	14.29 - 2.22	7x - 45x	30.00 - 4.67
1.5x	10.5x - 67.5x	20.00 - 3.11	15.75x - 101.25x	14.29 - 2.22	21x - 135x	9.52 - 1.48	10.5x - 67.5x	20.00 - 3.11



# **SLX** Series - Accessories

Eyecups &	Eyepieces
ST-036	Eyecups (pair), flat
ST-081	EW10x/21 eyepieces (pair), high eyepoint, with rubber cup
ST-082	WF15x/15 eyepieces (pair), high eyepoint
ST-083	WF20x/10 eyepieces (pair), high eyepoint
ST-084	WF10x/21 micrometric eyepiece, high eyepoint, with rubber cup
Additional	
ST-085.1	Additional lens 0.5x (w.d. 165mm) with SZ-EXT (except for SLX-1)
ST-091	Additional lens 0.75x (w.d. 105mm) (except for SLX-1)
ST-086.1	Additional lens 1.5x (w.d. 45mm) with compensating disc (except for SLX-1)
ST-087	Additional lens 2x (w.d. 33mm) (only for SLX-4 & SLX-5)
Stages	
ST-100.1	Hand moving stage, 100mm diameter
ST-110.1	Moving stage, coaxial knobs, 100mm diameter
<u>ST-111.1</u>	Moving stage, micrometric screws, 100mm diameter
Condensers	
ST-040.1	Darkfield condenser, 100mm diameter
ST-088.1	Polarising set (filters and rotating stage), 100mm diameter
Camera Ada	
<u>M-113.1</u>	Ring adapter, 30mm (for monocular and binocular microscopes)
<u>M-115</u>	0.35x C-Mount projection lens
M-114	0.5x C-Mount projection lens
M-118	0.75x C-Mount projection lens
M-173	C-Mount projection lens for APS-C/full frame reflex cameras (trino)
M-699	Universal adapter for C-Mount projection lens (trino)
M-620	0.35x focusable C-Mount adapter
M-620.1	0.5x focusable C-Mount adapter
M-620.2	0.65x focusable C-Mount adapter
<u>M-620.3</u>	1x focusable C-Mount adapter
Miscellaneo	
15104 DC 002	Cleaning kit
DC-002	Plastic dust cover, medium, 490(l)x490(h) mm (only for SLX-1, SLX-2 and SLX-3)
DC-004	TNT dust cover, large, 700(l)x550(h) mm (only for SLX-4 & SLX-5) Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)
M-005	
<u>ST-041</u> ST-042	Sample clip White / black phiest plate 100mm diameter
	White/black object-plate, 100mm diameter
ST-043 ST-092	Glass object-plate, 100mm diameter  Protective glass for stereohead
VP-SLX	IQ/OQ/PQ manual for SLX series
AB-020	
AD-UZU	Antibacterial surface treatment, only for newly purchased microscope



### How to connect the cameras to our microscopes. Please refer to the Adapter reference list on Digital section.

15104 - Cleaning kit It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.



v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

### **Optika Sales branches**

**OPTIKA**° Spain **OPTIKA**° China **OPTIKA**° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

**OPTIKA**<sup>®</sup> North America **OPTIKA**° Central America **OPTIKA**<sup>®</sup> Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# SZ Series



**Stereozoom Microscopes For Laboratory & Industry** 

# FIELD OF VIEW - 23 mm Full Plan Field of View

### 10x/22mm or 10x/23mm PLAN EYEPIECES

Plan eyepieces with 22 or 23 mm field of view; high eye-point type, also suitable for the use of eyeglasses.

### **ZOOM OBJECTIVE**

The High-Grade Zoom Objective ensures a sharp and clear vision. With a zoom ratio of 6.72:1 or 8.46:1 it makes this series a perfect instrument for any application.









# **X-LED** - Two Times Brighter Than Any Other


### **X-LED**: A SEA OF LIGHT

Where present, stands' illuminators are equipped with X-LED systems:

- the incident illuminator with **X-LED³** (single LED, 3.6W)
- the transmitted illuminator with **X-LED**<sup>T1-T3</sup> (12 or 60 LEDs)

### **X-LED**: RESPECTING COLORS

With 6,300K color temperature the specimens are illuminated with the most natural light. It allows to respect their colors, without altering the nuances.

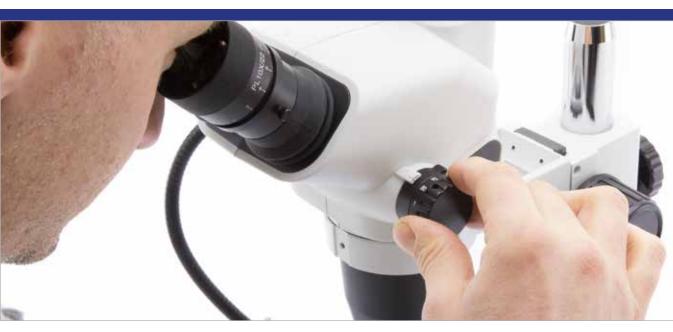



# **SZ** Series




### **Greenough Optical System**

The V-shape optical path of Greenough allows us to design a very compact and a slim unit, highly versatile and appreciated for the 3D viewing. Samples with significant depth can be quickly inspected. Binocular and trinocular heads are 45° inclined to ensure comfortable posture to the user even after several hours of operation.

### **Zoom Ratio**

SZ Series has 0.67x-4.5x or 0.65x-5.5x zoom range (6.72:1 or 8.46:1 zoom ratio, depending on the head), being purposely designed for routine and advanced inspections. These zoom ratios enable most samples to be observed at the appropriate magnifications. When combined with proper accessories (2x additional lens and 25x eyepieces), SZX-A delivers excellent image quality up to 275x.

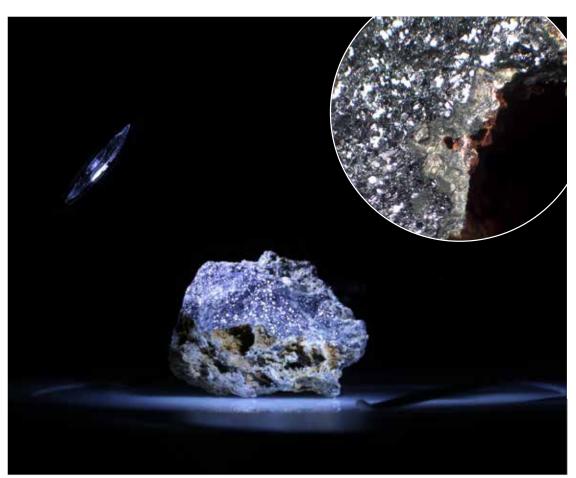

### **X-LED** Exclusive Lighting Source

Special technology able to double the light intensity for incomparable performance, ensuring constant pure-white colour temperature.


Relevant money & energy saving thanks to the incredibly low energy consumptions allow you to cut the electricity bills by 90%!

of great help.

### Get the most out of our accessories



# SZ Series - Modular Chart

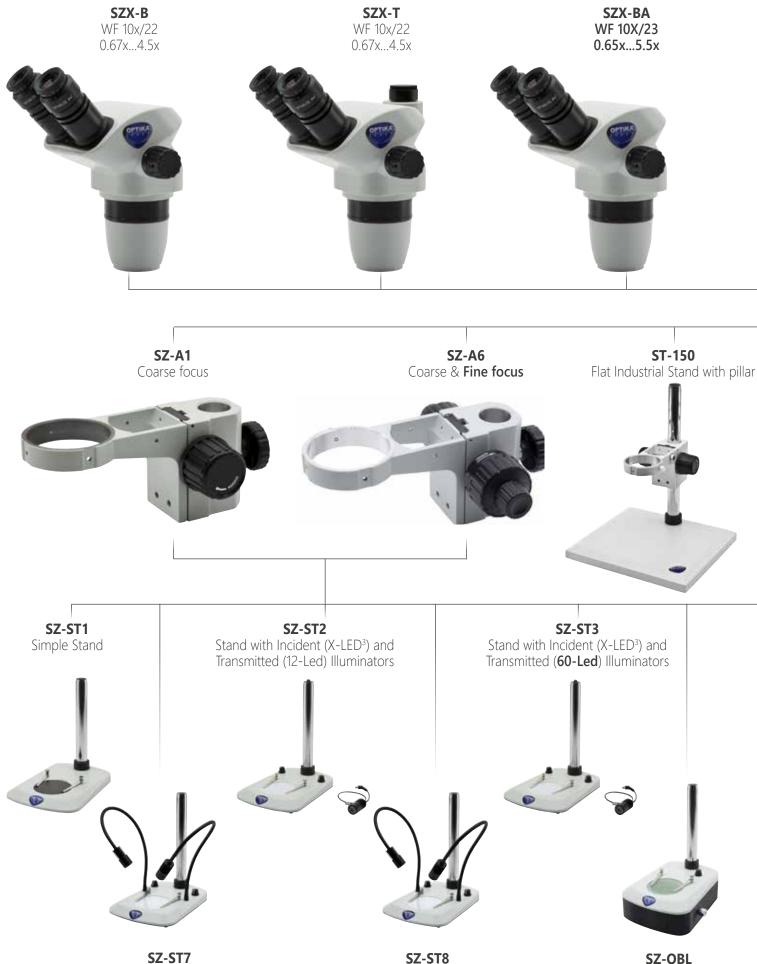

Stand with **Double Incident** (X-LED<sup>3</sup>) and Transmitted (12-Led) Illuminators

Stand with **Double Incident** (X-LED<sup>3</sup>) and Transmitted (**60-Led**) Illuminators

**SZ-OBL**Base with **oblique illumination**(X-LED<sup>3</sup>)

# SZ Series - Modular Chart

**SZX-TA** WF 10X23 0.65x - 5.5x



**SZO-B WF 10X/23**0.67x...4.5x **Click-stops** 



**SZO-T WF 10X/23**0.67x...4.5x **Click-stops** 



**ST-152**Flat Industrial Stand with fixed arm



SZ-STL1

Overhanging Stand

**SZ-STL3X**Advanced Precision Overhanging Stand with Coarse Focus Mechanism



**SZ-STL2**Advanced Overhanging Stand

**SZ-STL5**Flexible arm stand, 360° rotating



Stand Flexible arm stand, 360° rotating with 2x incident (X-LED³) illuminators







# **SZ** Series - Heads

# SZX-B / SZX-T



### SZX Heads - Excellent Price/Performance Ratio For Any Laboratory

- » Ideal for universities, experts & common routine lab requirements
- » 3D Greenough view for high resoluted images & large field depth
- » 22 mm field number and large working distance (up 110 mm)
- » 6.72:1 zoom ratio zoom magnification from 6.7x to 45x
- » Simultaneous eyepiece & camera observation (on SZX-T)
- » Cost-effective solution for diversified applications

### Head:

SZX-B: Binocular, 360° rotating on all stands and 45° inclined.

**SZX-T:** Trinocular (split ratio: 50/50), 360° rotating on all stands and 45° inclined

Interpupillary distance: Adjustable between 51 and 75 mm.

**Dioptric adjustment:** On both eyepieces.

**Eyepieces:** WF10x/22 mm, high eyepoint, secured by screw and with integrated rubber cups.

**Objective:** Parfocal achromatic zoom 0.67x...4.5x (zoom factor 6.72:1)

Working distance: 100 mm.

# SZX-BA / SZX-TA



### SZX-A Heads - Advanced Analysis With High Power Magnification

- » Purposely designed for particularly performing zoom conditions
- » 3D Greenough view for high resoluted images & large field depth
- » 23 mm field number and large working distance (up 110 mm)
- » 8.46:1 zoom ratio zoom magnification from 6.5x to 55x
- » Simultaneous eyepiece & camera observation (on SZX-TA)
- » High magnification change from overview to tiny details

### Head:

SZX-BA: Binocular, 360° rotating on all stands and 45° inclined.

**SZX-TA:** Trinocular (split ratio: 50/50), 360° rotating on all stands and 45° inclined.

**Interpupillary distance:** Adjustable between 51 and 75 mm.

**Dioptric adjustment:** On both eyepieces.

**Eyepieces:** WF10x/23 mm, high eyepoint, secured by screw and with integrated rubber cups.

**Objective:** Parfocal achromatic zoom 0.65x...5.5x (zoom factor 8.46:1).

Working distance: 102 mm.

# **SZ** Series - Heads

# SZO-B / SZO-T



### SZO Heads - Addressed For Extreme Reliability & Repeatability

- » Ensuring the sharpest vision, high productivity, repetitive analysis
- » 3D Greenough view for high resoluted images & large field depth
- » 23 mm field number and large working distance (up 110 mm)
- » 6.72:1 zoom ratio zoom magnification from 6.7x to 45x
- » Simultaneous eyepiece & camera observation (on SZO-T)
- » Multi-position click-stop no need to move eyes from eyepieces

### Head:

**SZO-B:** Binocular, 360° rotating on all stands and 45° inclined.

**SZO-T:** Trinocular (split ratio: 70/30), 360° rotating on all stands and 45° inclined

**Interpupillary distance:** Adjustable between 51 and 75 mm.

**Dioptric adjustment:** On both eyepieces.

**Eyepieces:** WF10x/23 mm, high eyepoint, secured by screw and with retractable rubber cups.

**Objective:** Parfocal achromatic zoom 0.67x...4.5x (zoom factor 6.72:1) with click-stop.

Working distance: 110 mm.



### **CLICK-STOP DETENTS**

On SZO Heads the zoom objective is equipped with a precise click-stop mechanism which makes easy to quickly find the default zoom positions: 0.67x, 0.80x, 1x, 1.2x, 1.5x, 2x, 2.5x, 3x, 3.5x, 4x, and 4.5x.



# **SZ** Series - Zoom Heads Comparison Chart

Model	Head	Eyepieces	Zoom Objective	Zoom Ratio	Click-Stops Detents	Working Distance	Photo Tube
SZX-B	Binocular 45° inclined	WF 10x/22	0.67x4.5x	6.72:1	-	100 mm	-
SZX-T	Trinocular 45° inclined	WF 10x/22	0.67x4.5x	6.72:1	-	100 mm	Fixed 50%-50%
SZX-BA	Binocular 45° inclined	WF 10x/23	0.65x5.5x	8.46:1	-	102 mm	-
SZX-TA	Trinocular 45° inclined	WF 10x/23	0.65x5.5x	8.46:1	-	102 mm	Fixed 50%-50%
SZO-B	Binocular 45° inclined	WF 10x/23	0.67x4.5x	6.72:1	0.67x, 0.80x, 1x, 1.2x, 1.5x, 2x, 2.5x, 3x, 3.5x, 4x, and 4.5x	110 mm	-
SZO-T	Trinocular 45° inclined	WF 10x/23	0.67x4.5x	6.72:1	0.67x, 0.80x, 1x, 1.2x, 1.5x, 2x, 2.5x, 3x, 3.5x, 4x, and 4.5x	110 mm	Fixed 70%-30%

# **SZ** Series - Pillar Stands



To be used in combination with a focusing mechanism (SZ-A1/A6)

Extremely stable, ultra-slim and modern stand to be matched with a focusing system and stereozoom head. The 25 mm height makes this stand very effective in preventing fatigue during operation, increasing the ergonomy and the performance as a result.

In case illumination is needed, choose from the wide choice of external illuminators available.

Width: 210 mm

Length: 270 mm

Thickness: 25 mm

Pillar diameter: 32 mm

Pillar height: 275 mm

Stage: 95 mm dia., black and white disc with sample clips

**Incident illumination:** None

**Transmitted illumination:** None



To be used in combination with a focusing mechanism (SZ-A1/A6)

# Top-class X-LED<sup>T</sup> transmitted light, with geometrically arranged LEDs and X-LED<sup>3</sup> incident light

The impressive incident illumination generated by the exclusive X-LED<sup>3</sup> swivelling illuminator (3.6 W) is combined with an extralarge, settable 12-LED X-LED<sup>T1</sup> disc (2 W, on SZ-ST2) or 60-LED X-LED<sup>T3</sup> disc (4 W, on SZ-ST3) for transmitted illumination.

**Incident illumination:** X-LED<sup>3</sup> with white 3.6 W LED and brightness control. Color temperature: 6,300 K.

### **Transmitted illumination:**


**SZ-ST2:** 12-LED X-LED<sup>T1</sup> disc with white 2 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply. **SZ-ST3:** 60-LED X-LED<sup>T3</sup> disc with white 4 W LED and brightness control. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power supply.

Base dimensions: 270×210 mm.

Pillar dimensions: Ø32×280 mm.

**Weight:** 1.90 kg.

# SZ-ST7 / SZ-ST8

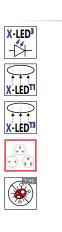



Top-class X-LED<sup>™</sup> transmitted light, with geometrically arranged LEDs and freely orientable, flexible double gooseneck X-LED<sup>®</sup> incident light

The impressive incident illumination generated by the exclusive X-LED<sup>3</sup> lighting system located in two flexible gooseneck arms (3.6 W each, for a total of 7.2 W) is combined with an extra-large, settable 12-LED X-LED<sup>T1</sup> disc (2 W, on SZ-ST7) or 60-LED X-LED<sup>T3</sup> disc (4 W, on SZ-ST8) for transmitted illumination.

**Incident illumination:** Two flexible X-LED<sup>3</sup> gooseneck arms with white 3.6 W LED/each and brightness control. Illuminance: 170,000 lux (at 10 cm distance). Color temperature: 6,300 K.

### **Transmitted illumination:**

**SZ-ST7:** 12-LED X-LED<sup>T1</sup> disc with white 2 W LED and brightness control. Illuminance: 4,000 lux. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power supply.

**SZ-ST8:** 60-LED X-LED<sup>T3</sup> disc with white 4 W LED and brightness control. Illuminance: 8,000 lux. Color temperature: 6,300 K. Multi-plug 100-240Vac/12Vdc external power supply.

Base dimensions: 270×210 mm.

Pillar dimensions: Ø32×280 mm.

**Stage:** 95 mm dia., translucent and black/white discs with sample clips

Weight: 2.40 kg.

# **SZ** Series - Pillar Stands

# SZ-OBL



To be used in combination with a focusing mechanism (SZ-A1/A6)

Extremely stable and modern stand, equipped with rotatable mirror and X-LED<sup>3</sup> transmitted light. Ideal for observing low-contrast samples that cannot be stained. The X-LED<sup>3</sup> light source with its high efficiency and no indirect heating is ideal for observing "live" specimens. Tiltable mirror with 2 faces, one smooth and one knurled, to get maximum contrast from your samples.

### **Transmitted illumination:**

X-LED<sup>3</sup> disc frosted white and brightness control. Illuminance: 8,000 lux. Color temperature: 6,300 K. Multi-plug 100-240Vac/6Vdc external power supply.

**Base dimensions:** 270×210 mm. **Pillar dimensions:** Ø32×280 mm. **Weight:** 2.4 kg.

# **SZ** Series - X-LED Comparison



SZ-ST2: 12-Led *X-LED<sup>TT</sup>*Transmitted Illuminator



SZ-ST3: 60-Led X-LED<sup>T3</sup>
Transmitted Illuminator



**SZ-ST7:** 12-Led **X-LED**<sup>TI</sup>
Transmitted Illuminator




**SZ-ST8:** 60-Led **X-LED**<sup>73</sup>
Transmitted Illuminator

# **SZ** Series - Pillar Stands Comparison Chart

Model	Base	Pillar	Stage	Incident Illumination	Transmitted Illumination
SZ-ST1	270x210x25h mm	32dia.x275h mm	95 mm dia., black and white disc with sample clips	None	None
SZ-ST2	270x210x30h mm	32dia.x280h mm	95 mm dia., white disc with sample clips	<b>X-LED³</b> (single Led, 3.6 W); brightness control	<b>X-LED</b> <sup>77</sup> (12-Led); brightness control
SZ-ST3	270x210x30h mm	32dia.x280h mm	95 mm dia., white disc with sample clips	<b>X-LED³</b> (single Led, 3.6 W); brightness control	<b>X-LED</b> <sup>T3</sup> (60-Led); brightness control
SZ-ST7	270x210x30h mm	32dia.x280h mm	95 mm dia., white disc with sample clips	Double <b>X-LED³</b> (total 7.2 W) on flexible goosencks; brightness control	<b>X-LED</b> <sup>TT</sup> (12-Led); brightness control
SZ-ST8	270x210x30h mm	32dia.x280h mm	95 mm dia., white disc with sample clips	Double <b>X-LED³</b> (total 7.2 W) on flexible goosencks; brightness control	<b>X-LED</b> <sup>73</sup> (60-Led); brightness control
SZ-OBL	270x210x85h mm	32dia.x280h mm	95 mm dia., white disc with sample clips	None	<b>X-LED</b> <sup>73</sup> brightness control

# **SZ** Series - Boom Stands

# SZ-STL3X

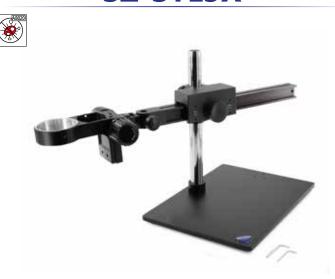

Extremely stable, long ball bearing stand complete of head holder and coarse focusing system for observation of particularly large specimens to be paired with a stereozoom head. Smooth horizontal and vertical movements are ensured and the head can be easily swivelled and tilted for inspection at oblique angles. In case illumination is needed, choose from the wide choice of external illuminators available.

**Base dimensions:** : 380x280 mm **Pillar dimensions:** Ø32×360 mm **Horizontal arm:** 460 mm

Maximum sample height: 270 mm

Swivelling movement: 360°
Tilting movement: 180°

Weight: 7 kg

# SZ-STL1



# SZ-STL2



Extremely stable, long overhanging stand for observation of particularly large specimens to be matched with a focusing system and stereozoom head. Smooth horizontal and vertical movements are ensured and the head can be easily swivelled for inspection at oblique angles. In case illumination is needed, choose from the wide choice of external illuminators available.

**Base dimensions:** 230×230 mm **Pillar dimensions:** Ø32×435 mm.

Horizontal arm: 415 mm

Maximum sample height: 400 mm

**Swivelling movement:** 360°

**Weight:** 14.00 kg

To be used in combination with a focusing mechanism (SZ-A1 or SZ-A6)

Extremely stable, hinged and long overhanging stand for observation of particularly large specimens to be matched with a focusing system and stereozoom head. Smooth horizontal and vertical movements are ensured and the head can be easily swivelled and tilted for inspection at oblique angles.

In case illumination is needed, choose from the wide choice of external illuminators available.

**Base dimensions:** 260×210 mm **Pillar dimensions:** Ø32×425 mm

Horizontal arm: 515 mm

Maximum sample height: 440 mm

**Swivelling movement:** 360°

Tilting movement: 180°

**Weight:** 18.8 kg

To be used in combination with a focusing mechanism (SZ-A1 or SZ-A6)

# **SZ** Series - Boom Stands / Flat Industrial Stands

# SZ-STL5



Highly versatile flexible arm stand, 360° rotating, ideal for industrial applications. It comes complete of head holder with focusing system and all the supports for table clamp. Small footprint is ensured when not in use, saving valuable space on the bench.

Table clamp included.

Total height: 330 mm

Total depth: 880 mm

Horizontal arm: 820 mm

Weight: 4.8 kg



Highly versatile flexible arm stand, 360° rotating, ideal for industrial applications. It comes complete of the settable and exclusive X-LED3 incident illumination system, head holder with focusing system and all the supports for table clamp.

Small footprint is ensured when not in use, saving valuable space on the bench. Multi-plug 100-240Vac/6Vdc external power supply. Head shown in picture not included

Total height: 330 mm
Total depth: 880 mm
Horizontal arm: 820 mm

**Illuminators:** 2x X-LED3 on gooseneck, with brightness control unit **Power supply:** Multi-plug 100/240Vac, 50-60Hz; output 12Vdc, 1.5A

Weight: 5.7 kg

ST-150 / ST-152



Large, simple plain stand with head holder and focusing mechanism.

ST-150: with pillar; ST-152: with fixed arm.

Base dimensions: 320x290 mm

**Pillar dimensions (ST-150):** Ø32×410 mm

Fixed arm dimensions (ST-152): 300 mm

303

### <u>3</u>

# **SZ** Series - Focusing Mechanisms

# SZ-A1



Entry-level coarse focusing to raise or lower the head to focus Coarse focusing system (76 mm head holder, 32 mm pillar) with adjustable tension for standard stereomicroscope requirements.

Type: Coarse

Coarse total travel: 50 mm Adjustable tension: Yes

**Head holder internal diameter:** 76 mm

Pillar diameter: 32 mm

# SZ-A6



Coaxial coarse and fine focusing results in easy, precise conditions

Coaxial coarse and fine focusing system (76 mm head holder, 32 mm pillar) with 2  $\mu$ m fine resolution and adjustable tension ideal for extremely precise fine adjustments and smooth movements.

**Type:** Coaxial coarse and fine

Coarse total travel: 50 mm

Fine total travel (per single rotation): 0.2 mm.

Fine graduations: 100
Fine resolution: 2 µm
Adjustable tension: Yes

Head holder internal diameter: 76 mm

Pillar diameter: 32 mm

# **SZ** Series - Best Selling Configurations



SZO-B + SZ-A1 + SZ-ST1



SZX-B + SZ-A1 + SZ-ST3



SZX-BA + SZ-A1 + SZ-ST7

# **SZ** Series - Optical Performance

# **SZX-B and SZX-T** Heads

Eyepiece	Eyepiece 10x (ST-141)		15x (ST-302)		20x (ST-303)		25x (ST-144)	
Field number (mm)	2	2	1	6	1	2	9	
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
0.3x (W.D: 241 mm)	2.01x-13.5x	109.45-16.30	3.02x-20.25x	79.60-11.85	4.02x-27x	59.70-8.89	5.02x-33.75x	44.78-6.67
0.5x (W.D: 160 mm)	3.35x-22.5x	65.67-9.78	5.02x-33.75x	47.76-7.11	6.7x-45x	35.82-5.33	8.37x-56.25x	26.87-4.0
0.75x (W.D: 116 mm)	5.02x-33.75x	43.78-6.52	7.54x-50.62x	31.84-4.74	10.05x-67.5x	23.88-3.56	12.56x-84.37x	17.91-2.67
1x (W.D: 100 mm)	6.7x-45x	32.84-4.89	10.05x-67.5x	23.88-3.56	13.4x-90x	17.91-2.67	16.75x-112.5x	13.43-2.0
1.5x (W.D: 44 mm)	10.05x-67.5x	21.89-3.26	15.07x-101.25	15.92-2.37	20.1x-135x	11.94-1.78	25.12x-168.75x	8.96-1.33
2x (W.D: 26 mm)	13.4x-90.0x	16.42-2.44	20.1x-135x	11.94-1.78	26.8x-180x	8.96-1.33	33.5x-225x	6.72-1.0

# **SZX-BA and SZX-TA** Heads

Eyepiece 10x (ST-306)		T-306)	15x (S	15x (ST-302)		20x (ST-303)		25x (ST-144)	
Field number (mm)	23		16		12		9		
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	
0.3x (W.D: 228 mm)	1.95x-15.5x	117.95-13.94	2.93x-23.25x	82.05-9.70	3.9x-31.0x	61.54-7.27	4.88x-38.75x	46.15-5.45	
0.5x (W.D: 164 mm)	3.25x-27.5x	70.77-8.36	4.88x-41.25x	49.23-5.82	6.5x-55x	36.92-4.36	8.13x-68.75x	27.69-3.27	
0.75x (W.D: 109 mm)	4.87x-41.25x	47.18-5.58	7.31x-61.88x	32.82-3.88	9.74x-82.5x	24.62-2.91	12.18x-103.13x	18.46-2.18	
1x (W.D: 102 mm)	6.5x-55x	35.38-4.18	9.75x-82.5x	24.62-2.91	13.0x-110x	18.46-2.18	16.25x-137.5x	13.85-1.64	
1.5x (W.D: 46 mm)	9.75x-82.5x	23.59-2.79	14.63x-123.75	16.41-1.94	19.5x-165x	12.31-1.45	24.38x-206.25x	9.23-1.09	
2x (W.D: 26 mm)	13.0x-110.0x	17.69-2.09	19.5x-165x	12.31-1.45	26.0x-210x	9.23-1.09	32.5x-262.5x	6.92-0.82	

# **SZO-B and SZO-T** Heads

Eyepiece	Eyepiece 10x (ST-301)		15x (S	15x (ST-302)		20x (ST-303)		25x (ST-144)	
Field number (mm)	Field number (mm) 23		16		12		9		
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	
0.3x (W.D: 287 mm)	2.01x-13.5x	114.43-17.04	3.02x-20.25x	79.60-11.85	4.02x-27x	59.70-8.89	5.02x-33.75x	44.78-6.67	
0.5x (W.D: 177 mm)	3.35x-22.5x	68.66-10.22	5.02x-33.75x	47.76-7.11	6.7x-45x	35.82-5.33	8.37x-56.25x	26.87-4.0	
0.75x (W.D: 120 mm)	5.02x-33.75x	45.77-6.81	7.54x-50.62x	31.84-4.74	10.05x-67.5x	23.88-3.56	12.56x-84.37x	17.91-2.67	
1x (W.D: 110 mm)	6.7x-45x	34.33-5.11	10.05x-67.5x	23.88-3.56	13.4x-90x	17.91-2.67	16.75x-112.5x	13.43-2.0	
1.5x (W.D: 47 mm)	10.05x-67.5x	22.89-3.41	15.07x-101.25	15.92-2.37	20.1x-135x	11.94-1.78	25.12x-168.75x	8.96-1.33	
2x (W.D: 26 mm)	13.4x-90.0x	17.16-2.56	20.1x-135x	11.94-1.78	26.8x-180x	8.96-1.33	33.5x-225x	6.72-1.0	

# **SZ** Series - Best Selling Configurations



SZO-T + SZ-A1 + SZ-STL1



SZO-T + SZ-A1 + SZ-STL2

# **SZ** Series - Accessories

### Evecups & Evepieces

ST-141	WF10x/22 eyepieces (pair), high eyepoint, focusable, rubber cup
31 171	
	(only for SZX-B & SZX-T)

WF10x/22 micrometric eyepiece, high eyepoint, focusable, rubber cup ST-145

(only for SZX-B & SZX-T)

ST-144 WF25x/9 eyepieces (pair), high eyepoint, focusable, with rubber cup ST-301 WF10x/23 eyepieces (pair), high eyepoint, focusable, with rubber cup

(only for SZO-B & SZO-T)

ST-302 WF15x/16 eyepieces (pair), high eyepoint, focusable, with rubber cup

ST-303 WF20x/12 eyepieces (pair), high eyepoint, focusable, with rubber cup ST-305 WF10x/23 micrometric eyepiece, high eyepoint, focusable, rubber cup (only for SZO-B & SZO-T)

ST-306 WF10x/23 eyepieces (pair), high eyepoint, focusable, rubber cup (only for SZX-BA & SZX-TA)

ST-310 WF10x/23 micrometric eyepiece, high eyepoint, focusable, rubber cup (only for SZX-BA & SZX-TA)

### **Additional Lenses**

ST-102 Additional lens 0.3x ST-103 Additional lens 0.5x ST-104 Additional lens 0.75x

ST-105 Additional lens 1.5x

ST-106 Additional lens 2x

### **Stages**

ST-100 Hand moving stage, for insert diameter of 95 mm

ST-110 Moving stage, coaxial knobs, for insert diameter of 95 mm

ST-111 Moving stage, micrometric screws, fort insert diameter of 95 mm

Applicable heating stage, for instert diameter of 95 mm, multiplug ST-666

### Condenser & Filters

ST-040 Darkfield condenser, 95mm diameter

ST-088 Polarising set (filters and rotating stage), 95mm diameter

### Camera Adapters

Ring adapter, 30mm (for monocular and binocular microscopes) <u>M-113.1</u>

M-115 0.35x C-Mount projection lens

0.5x C-Mount projection lens M-114

M-118 0.75x C-Mount projection lens

C-Mount projection lens for APS-C/full frame reflex cameras (trino) M-173

M-699 Universal adapter for C-Mount projection lens (trino)

0.35x focusable C-Mount adapter M-620

M-620.1 0.5x focusable C-Mount adapter

M-620.2 0.65x focusable C-Mount adapter

M-620.3 1x focusable C-Mount adapter

### Miscellaneous

15104 Cleaning kit

DC-002 Plastic dust cover, medium, 490(l)x490(h) mm

DC-003 TNT dust cover, medium, 600(l)x550(h) mm

DC-004 TNT dust cover, large, 700(l)x550(h) mm

M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)

M-069 Solar charger

ST-014 Glass object-plate, 95mm diameter

ST-041 Sample clip

ST-092 Protective glass for stereohead

VP-SZ IQ/OQ/PQ manual for SZ series

AB-010 Antibacterial surface treatment, only for newly purchased microscope (Only for SZO-B & SZO-T)

AB-020 Antibacterial surface treatment, only for newly purchased microscope



### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.









v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

OPTIKA S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

### **Optika Sales branches**

**OPTIKA**° Spain **OPTIKA**° China **OPTIKA**° India

spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

**OPTIKA**<sup>®</sup> North America **OPTIKA**<sup>®</sup> Central America **OPTIKA**<sup>®</sup> Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# SZP Series



**Advanced CMO Stereozoom Microscopes** 

OPTIKA

# **SZP** Series

### **Greenough Vs. CMO Stereo Microscopes**

The Greenough design, introduced in the early twentieth century, consists of two identical, symmetrical optical systems each containing a separate eyepiece and objective, arranged in careful alignment with a single housing. The two optical paths are arranged at a precise angle, due to which they converge to the point that determines the working distance and thus the focal length. This design allows high numerical apertures because the objectives are very similar to those used in compound microscopes. The lower section of the microscope contains the objectives, while the upper end of the body tubes projects an image to the eyes.

CMO (Common Main Objective) stereo microscopes have



3

# CMO Stereozoom Microscopes


the ability to collect more light than the Greenough design and feature greater optical aberration correction.

CMO stereo microscopes, unlike Greenough stereo microscopes, are by definition instruments whose optics are infinity corrected. This allows you to easily introduce accessories such as beamsplitters or aperture diaphragms thanks to which you can get a greater depth **ST-172** of field.

### Get the most out of our accessories

### ST-172 - Iris diaphragm module


With this accessory, the depth of field of the final image can be adjusted, for applications where different planes all in focus are needed.

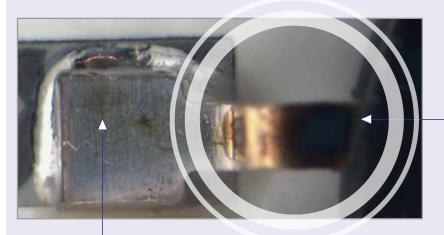

### **OUT-OF-FOCUS AREA**

With iris diaphragm open

### **FOCUS PLANE**

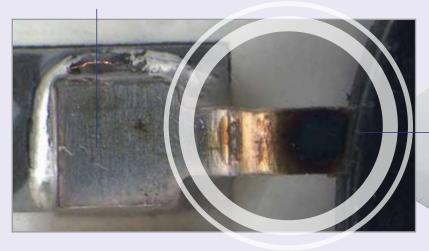

### **IN-FOCUS AREA**

With iris diaphragm close

OPTIKA

Soldered Led Pin - SZP-8 with iris diaphragm module and ST-156 stand.

# **SZP** Series - Heads

# **SZP-6, SZP-8, SZP-10**












Part	Description
Head:	Binocular, 360° rotating on all stands and 30° inclined.
Interpup. distance:	Adjustable between 52 and 75 mm.
Dioptric adjustment:	On both eyepieces.
Eyepieces:	WF10x/24 mm, high eye-point.
SZP-6 zoom body:	Parfocal achromatic zoom 0.8x5.0x (zoom factor 6.25:1).
SZP-8 zoom body:	Parfocal achromatic zoom 0.8x6.4x (zoom factor 8:1).
SZP-10 zoom body:	Parfocal achromatic zoom 0.8x8.0x (zoom factor 10:1).
Zoom positioning:	By click-stop mechanism
Objective lens:	Plan Achromatic 1x.
Working distance:	80 mm.
Optical system:	Galilean (Parallel, infinity corrected).

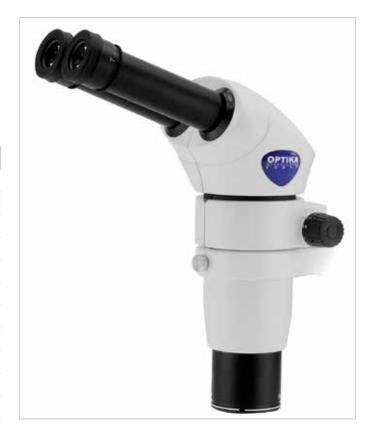

# SZP-6e, SZP-8e, SZP-10e









Part	Description
Head:	Ergonomical binocular, 360° rotating on all stands and freely inclinable from 0° to 35°.
Interpup. distance:	Adjustable between 55 and 80 mm.
Dioptric adjustment:	On both eyepieces.
Eyepieces:	WF10x/24 mm, high eye-point.
SZP-6e zoom body:	Parfocal achromatic zoom 0.8x5.0x (zoom factor 6.25:1).
SZP-8e zoom body:	Parfocal achromatic zoom 0.8x6.4x (zoom factor 8:1).
SZP-10e zoom body:	Parfocal achromatic zoom 0.8x8.0x (zoom factor 10:1).
Zoom positioning:	By click-stop mechanism
Objective lens:	Plan Achromatic 1x.
Working distance:	80 mm.
Optical system:	Galilean (Parallel, infinity corrected).



# **SZP-FL** - Epi-Fluorescence Attachment

### Attachment for fluorescence applications for SZP stereomicroscopes only.

Used in many applications like biology, botany, electronics, materials, forensics. Equipped with HBO 100W mercury lamp illuminator. *To be combined with any SZP model.* 





SZP-FL	HBO fluorescence attachment for SZP heads
Description:	SZP fluorescence attachment for biology, industrial inspection, forensics, etc. Essential tool for security printing and mineral research.
Illumination:	100W HBO high-pressure mercury bulb. Average lamp lifetime: 400 hours. Input voltage: 110/240Vac, 50/60Hz, 1A; Fuse: F8AL 250V. Maximum input power: 125W. Current and time counter LED displays.
Photo Attachment:	Trinocular output Photo/Video port.

# SZP Series - Stands

# **ST-155**











Modern, large plain stand equipped with **LED transmitted and incident** illumination, both with intensity control. It comes complete of head holder and focusing mechanism.

Base size: 330x290 mm. Height: 40 mm. Fixed arm height: 340 mm.

Head not included.



# **ST-156**











Modern, large plain stand equipped with LED transmitted and incident illumination, both with intensity control. It comes complete of head holder and coaxial coarse and fine focusing system. Base size: 330x290 mm. Height: 40 mm. Fixed arm height: 340 mm. Head not included.



# **SZP** Series - Optical Performance

### **SZP** Heads

### **Optical performance SZP-6 / SZP-6e**

Eyepiece	10x (ST-160)		15x (ST-161)		20x (ST-162)	
Field number (mm)	24		15		10	
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
0.3x (W.D: 280 mm)	2.4x-15x	100.00-16.00	3.6x-22.5x	62.50-10.00	4.8x-30x	41.67-6.67
0.5x (W.D: 165 mm)	4x-25x	60.00-9.60	6x-37.5x	37.50-6.00	8x-50x	25.00-4.00
1x (W.D: 80 mm)	8x-50x	30-4.80	12x-75x	18.75-3.00	16x-100x	12.50-2.00
2x (W.D: 32.5 mm)	16x-100x	15-2.40	24x-150x	9.38-1.50	32x-200x	6.25-1.00

### **Optical performance SZP-8 / SZP-8e**

Eyepiece	10x (ST-160)		15x (ST-161)		20x (ST-162)	
Field number (mm)	24		15		10	
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
0.3x (W.D: 280 mm)	2.4x-19.2x	100.00-12.50	3.6x-28.8x	62.50-7.81	4.8x-38.4x	41.67-5.21
0.5x (W.D: 165 mm)	4x-32x	60.00-7.50	6x-48x	37.50-4.69	8x-64x	25.00-3.13
1x (W.D: 80 mm)	8x-64x	30.00-3.75	12x-96x	18.75-2.34	16x-128x	12.50-1.56
2x (W.D: 32.5 mm)	16x-128x	15.00-1.88	24x-192x	9.38-1.17	32x-256x	6.25-0.78

### **Optical performance SZP-10 / SZP-10e**

	-					
Eyepiece	10x (ST-160)		15x (ST-161)		20x (ST-162)	
Field number (mm)	24		15		10	
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)
0.3x (W.D: 280 mm)	2.4x-24x	100.00-10.00	3.6x-36x	62.50-6.25	4.8x-48x	41.67-4.17
0.5x (W.D: 165 mm)	4x-40x	60.00-6.00	6x-60x	37.50-3.75	8x-80x	25.00-2.50
1x (W.D: 80 mm)	8x-80x	30.00-3.00	12x-120x	18.75-1.88	16x-160x	12.50-1.25
2x (W.D: 32.5 mm)	16x-160x	15.00-1.50	24x-240x	9.38-0.94	32x-320x	6.25-0.63

# **SZP** Series - Accessories

### **Eyecups & Eyepieces**

CT 4 CO	1 11/210 /04			( 11	*41 1.1
VI - I 6 1	1 \//\=       \/ / / / /	avaniaca hian	AVADAIDT	tocucanio	With runnar cun
ST-160	)	EVENIECE, IIIUII	EVEDUIII	, iucusanie,	with rubber cup

- ST-161 WF15x/16 eyepieces (pair), focusable, with rubber cup
- ST-162 WF20x/12 eyepieces (pair), focusable, with rubber cup
- ST-163 WF10x/24 micrometric eyepiece, high eyepoint, focusable, rubber cup

### **Additional Lenses**

ST-165	0.3x obje	ective (w.d	<u>d. 280mm)</u>

- ST-166 0.5x objective (w.d. 118mm)
- ST-167 2x objective (w.d. 32.5mm)

### **Camera Adapters**

- M-113.1 Ring adapter, 30mm (for monocular and binocular microscopes)
- M-115 0.35x C-Mount projection lens
- M-114 0.5x C-Mount projection lens
- M-118 0.75x C-Mount projection lens
- M-173 C-Mount projection lens for APS-C/full frame reflex cameras (trino)
- M-699 Universal adapter for C-Mount projection lens (trino)
- M-620 0.35x focusable C-Mount adapter
- M-620.1 0.5x focusable C-Mount adapter
- M-620.2 0.65x focusable C-Mount adapter
- M-620.3 1x focusable C-Mount adapter

### Miscellaneous

- 15104 Cleaning kit
- DC-002 Plastic dust cover, medium, 490(l)x490(h) mm
- DC-003 TNT dust cover, medium, 600(l)x550(h) mm
- DC-004 TNT dust cover, large, 700(l)x550(h) mm
- M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)
- M-151 HBO 100W high-pressure mercury bulb for fluorescence
- M-151.1 HBO 100W high-pressure mercury bulb for fluorescence (OSRAM)
- ST-098 Polarising set (filters and rotating stage)
- ST-170 Photo/Video beam splitter 1 port
- ST-171 Photo/Video beam splitter 2 port
- ST-172 Iris diaphgram module
- ST-176 Protective glass for stereohead (only works with standard 1x lens)
- SZP-FL HBO fluo attachment, 3-pos. (B & G filter set), multi-plug
- VP-SZP IQ/OQ/PQ manual for SZP series
- AB-020 Antibacterial surface treatment, only for newly purchased microscope

# 15104 - Cleaning kit It cleans glass quickly and effectively, without leaving residue or odor. Ideal for precision lens or prism cleaning.







### How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

### **Optika Sales branches**

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com

**OPTIKA**<sup>®</sup> North America **OPTIKA**<sup>®</sup> Central America **OPTIKA**<sup>®</sup> Africa namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# CL Series



**Cold Light Illuminators** 

# **CL** Series - Illuminators

## **CLD-01**



**CL-17.1:** Additional polarizing filters for CL-11.1

Polarizing filters (pair) for CL-11.1. CL-17.1: Ring optical fiber guide (length 700 mm). Guide diameter 16 mm. CL-12:

Diameter of the fixing ring: 55mm.

The circular end is suitable for all stereo series.

# **CL-11.1**





### Double-arm guide for CLD-01, with focusing lenses.

Arm length 500 mm. Each arm fitted with focusable lens adapter. Illuminance (at 10 cm distance): 100.000 lux.

Optional accessory:

CL-17.1: Polarizing filters (pair) for CL-11.1.





### Ring optical fiber guide for CLD-01.

Lenght 700 mm, diameter 16 mm.

The circular end is suitable for all series, by using the three locking screws. Diameter of the fixing ring: 55 mm.

Illuminance (at 10 cm distance): 40,000 lux.

# **CL** Series - Illuminators



Arm illuminator with magnetic fixing plate, USB powered.

**Light sources:** LED; With 2W high efficiency LED. **Color temperature:** pure white 6,300 K;

Luminous flux: 200 lm;

Illuminance: 20,000 (on 5 cm diameter).

Pair of arm illuminators with magnetic fixing plates, metal base, USB. Light sources: LED; With 2W high efficiency LEDs. Color temperature: pure white 6,300 K; Luminous flux: 200 lm (x2);

Illuminance: 20,000 (on 5 cm diameter) (x2).

# CL-32/CL-33 quick connection





# **CL-41**



Double-arm X-LED<sup>3</sup> illuminator, with brightness control.

**Light sources:** With 3,6 W high efficiency X-LED<sup>3</sup>.

Color temperature: pure white 6,300 K; Luminous flux: 400 lm each arm;

**Illuminance:** 170,000 lux (at 10 cm distance). Multi-plug 100-240Vac/12Vdc external power supply.

### (3)

# **CL** Series - Illuminators

## **CL-14**

# **CL-16.1**



**56-LED ring light illuminator, with brightness control.** Compact, with 360° rotating ring connector.

Illuminance (at 10 cm distance): >8,000 lux.



**Professional lighting system** including 144 LEDs (2W total power) for enhanced light uniformity and brightness. The ring light illumination is divided into 4 different zones, individually adjustable for **selectable light zones**. The **separated, external control panel** prevents interferences during use, whilst the sturdy metal structure makes it **more durable and resistant**.

**Diameter of the fixing ring:** 60mm. **Color temperature:** pure white 6,300 K. **Illuminance:** 6,000 lux (at 10 cm distance).

### **CL-18**





**Professional lighting system** including 96 LEDs for enhanced light uniformity and brightness, and **built-in rotating polarizing filter** and analyzer filters for complete polarization, ideal to reduce glares and reflections when viewing metallic parts.

Provided with on-board intensity level adjustement.

**Diameter of the fixing ring:** 60mm. **Color temperature:** pure white 6,300 K. **Illuminance:** 9,000 lux (at 10 cm distance).

v7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

### **Optika Sales branches**

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com OPTIKA® North America
OPTIKA® Central America
OPTIKA® Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# GEM Series



**Stereozoom Microscopes For Gemology** 

# **GEM** Series

As a leading company in the supply of gemological microscopes, OPTIKA offers a series of microscopes purposely designed for this sector by using both brightfield and darkfield methods OPTIGEM-10 and OPTIGEM-20. This series has been designed and manufactured in order to satisfy the requests of a very demanding industry; brightfield/darkfield, immersion analysis, light color temperature: no detail has been left to chance.

### **Specifically Designed for Specialists**

Gemological stereomicroscopes are meant to help with stone inspection. Jewels and gems have a variety of grades (or quality levels), which ultimately influence their value and cost on the market; therefore it is important to have solutions that are purposely designed for gemology. These stereomicroscopes are equipped with iris and darkfield condenser on the bottom light source, and with a set of on stage tweezers to hold the stone in place.

### **Much More Than Gemological Stereomicroscopes**

**OPTIGEM-10 & OPTIGEM-20** are two-in-one gemology instruments that can be used both in vertical and horizontal position in a very easy way, just by turning one knob (no disassembling and re-assembling operations are required). The horizontal position extends the use of a gemological microscope, giving the possibility to perform immersion analysis by submerging a sample in liquid. If the stone's refractive index is close to the liquid's one, immersion makes the interior more visible by reducing the effects of refraction and surface reflection. This enables you to see a gem's inclusions or color distribution more easily.

Immersion is also necessary to see crystal growth structures, which might help you separate natural from synthetic corundum. Features like curved growth striae in flame-fusion synthetics, or separation planes in assembled stones, are often far easier to see when the stone is immersed.



Vertical position for standard gem analysis with darkfield illumination and polarizing tecnique



Horizontal position for immersion gem analysis

### **Incredibly Versatile Operations**

**OPTIGEM Series** offers multiple options for illumination and contrast techniques, such as incident, transmitted and oblique brightfield darkfield, polarization and immersion analysis only on Optigem 10 & 20. They come with a special optical condenser configuration to ensure real, perfect darkfield application (see the dedicated chapter for further information).

# Stereozoom Microscopes For Gemology

### **Ultrabright LED Condenser for Optimized Illumination**

An ultrabright LED-based electronic condenser with intensity control allows to switch from brightfield to darkfield; it also produces perfectly the colour of daylight.

The condenser uses a new optical configuration especially created in order to obtain a perfect darkfield application.

With darkfield observation, the unscattered beams from the image are excluded: as a result, the field around the specimen is generally dark.

An additional flexible arm and velvet-field slider produce extra contrast for crisp and vibrant images. The illumination of OPTIGEM microscopes is greatly performing and this brings this series to be ideal for precious stones and jewels evaluation.



### Get the most out of our accessories



ST-202 - Polarizing analysis kit

Polarization technique allows to quickly determine if the stone at hand is isotropic or anisotropic or, at best, to determine the optic character of gemstones (twin planes, strain, pleochroism, etc.). It is also the preferred tool for separating synthetic Quartz from its natural counterparts. In addition, the polarizing microscope may be very useful for distinguishing solid inclusions from negative inclusions as well as for spotting polysynthetic twinning.

### **Applications**


Some application examples demonstrating the performance of OPTIGEM Series, especially designed to observe samples of precious stones and jewels and provided with specific features for gemological needs.


### Legend


- 1. Inspection of stones with OPTIGEM-10.
- 2. Inspection of stones with pure white darkfield illumination.
- 3. Immersion cell (ST-203) on a translating support (ST-204).
- 4. Sample of Citrine.
- 5. ST-201 accessory creates a soft darkfield illumination ideal for diamond analysis.
- 6. Optigem can be easily rotated to a horizontal working position.
- 7. ST-201 accessory for analysis under polarized light.
- 8. Working with Optigem and its accessories (they can be stacked for increased functionality).

# **GEM** Series



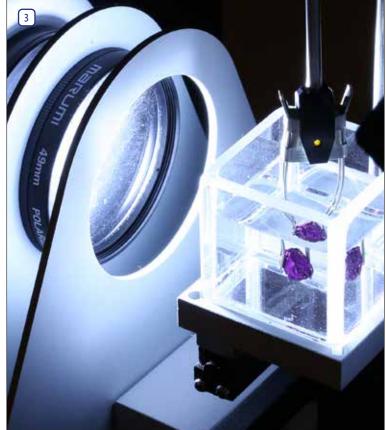

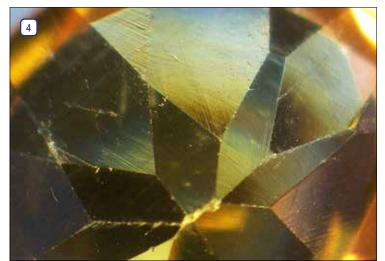

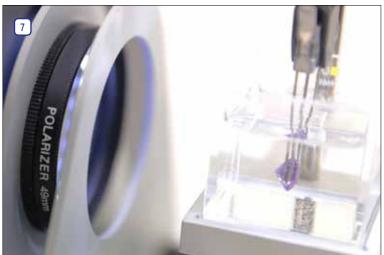







# **GEM** Series - **OPTIGEM10** model

Binocular gemological stereomicroscopes for brightfield and darkfield applications with special side-emitting **LED** illumination ring for true darkfield illumination. Equipped also with **4 incident LED** flexible arm and a diffusive **LED** disc for transmitted illumination. The instrument can be easily tilted horizontally for immersion analysis.



















Part	Description
<b>Observation mode:</b>	Brightfield, darkfield.
Heads	Binocular, 45° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 51 and 75 mm
Dioptric adjustment:	On both eyepiece tubes.
Eyepieces:	WF10x/21 mm, high eye-point.
Objective:	Parfocal achromatic zoom 0.7x4.5x (zoom factor 6.43:1).
Working distance:	100 mm
Specimen stage:	Gemological stage, with clamp for holding gems.
Focusing:	Rack and pinion mechanism controlled by a pair of knobs.
Stand:	Fixed arm stand with tilting system (with position lock control). It can be totally overturned and rotated allowing in this way to obtain a 2-in-one instrument: an instrument for standard observations and one for observation of water-immersed gems.

Part	Description
Darkfield illumination:	Equipped with a state-of-the-art illuminator for darkfield observation. It consist of an innovative side-emmiting LEDs ring with an emission angle of 38°. With brightness control.
Transmitted light llumination:	Equipped with a LED illuminator, located under the stage. With brightness control.
Incident illumination:	Equipped with a flexible gooseneck-arm 4-LED illuminator. With brightness control.
<b>Color temperature:</b>	Pure white 6,300 K

## **GEM** Series - **OPTIGEM20** model

Trinocular gemological stereomicroscopes for brightfield and darkfield applications with special side-emitting **LED** illumination ring for true darkfield illumination. Equipped also with **4 incident LED** flexible arm and a diffusive **LED** disc for transmitted illumination. The instrument can be easily tilted horizontally for immersion analysis.



















Part	Description
Observation mode:	Brightfield, darkfield.
Heads	Trnocular, 45° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 51 and 75 mm
Dioptric adjustment:	On both eyepiece tubes.
Eyepieces:	WF10x/21 mm, high eye-point.
Objective:	Parfocal achromatic zoom 0.7x4.5x (zoom factor 6.43:1).
Working distance:	100 mm
Specimen stage:	Gemological stage, with clamp for holding gems.
Focusing:	Rack and pinion mechanism controlled by a pair of knobs.
Stand:	Fixed arm stand with tilting system (with position lock control). It can be totally overturned and rotated allowing in this way to obtain a 2-in-one instrument: an instrument for standard observations and one for observation of water-immersed gems.

Part	Description
Darkfield illumination:	Equipped with a state-of-the-art illuminator for darkfield observation. It consist of an innovative side-emmiting LEDs ring with an emission angle of 38°. With brightness control.
Transmitted light llumination:	Equipped with a LED illuminator, located under the stage. With brightness control.
Incident illumination:	Equipped with a flexible gooseneck-arm 4-LED illuminator. With brightness control.
Color temperature:	Pure white 6,300 K

## **GEM** Series - Comparison chart


Model	Head	Eyepieces	Objective	Stand	Illumination
OPTIGEM-10	Binocular, 360° rotating, 45° inclined	Wide Field 10x/21mm	0.7 4.5x Zoom	Gemological stand	Incident illumination: 4 LED flexible arm with brightness adjustment.  Transmitted illumination: Diffusive LED disc for observation in brightfield and side-emitting LED ring for observation in darkfield.
OPTIGEM-20	Trinocular, 360° rotating, 45° inclined	Wide Field 10x/21mm	0.7 4.5x Zoom	Gemological stand	Incident illumination: 4 LED flexible arm with brightness adjustment.  Transmitted illumination: Diffusive LED disc for observation in brightfield and side-emitting LED ring for observation in darkfield.

## **GEM** Series - Optical Performance


### **OPTIGEM-10 / OPTIGEM-20 - Optical performance**

Eyepiece	10x (ST	-081)	15x (5	ST-082)	20x (ST-083)			
Field number (mm)	21		1	15	10			
Objective	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)	Total magnification	Field of View (mm)		
0.5x (W.D. 165 mm)	3.5x-22.5x	60.00-9.23	5.25x-33.75x	42.86-6.67	7x-45x	28.57-4.44		
0.75x (W.D. 117 mm)	5.25x-33.75x	40.00-6.22	7.875x-50.625x	28.57-4.44	10.50x-67.5x	19.05-2.96		
1x (W.D. 100 mm)	7x-45x	30.00-4.66	10.5x-67.5x	21.43-3.33	14x-90x	14.29-2.22		
1.5x (W.D. 47 mm)	10.50x-67.5x	20.00-3.11	15.75x-101.25x	14.29-2.22	21x-135x	9.52-1.48		
2x (W.D. 33 mm)	14x-90x	15.00-2.33	21x-135x	10.71-1.67	28x-180x	7.14-1.11		

## **GEM** Series - Contrast method comparison

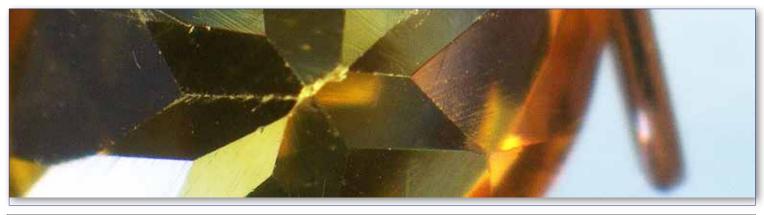

Citrine - OPTIGEM-20 - 0.7x zoom - Brightfield

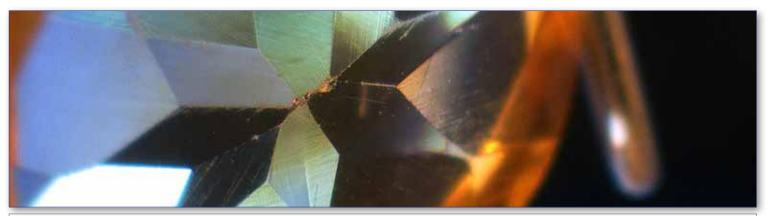

Citrine - OPTIGEM-20 - 0.7x zoom - Darkfield

## **GEM** Series - Accessories

### **Eyecups & Eyepieces**

ST-081	AVADIACAC I	nair) hi	th avanaint	with rubber cup

ST-082 WF15x/15 eyepieces (pair), high eyepoint

ST-083 WF20x/10 eyepieces (pair), high eyepoint

ST-084 WF10x/21 micrometric eyepiece, high eyepoint, with rubber cup

### **Additional Lenses**

ST-086	Additional lens 1.5x (w.d. 45mm)
ST-087	Additional lens 2x (w.d. 33mm)

### **Condenser & Filters**

ST-202 Polarizing analisys kit

### **Camera Adapters**

scor	micr	ar	TH.	()(	ım	l D	M	ar	ar	Ш	$\alpha$	วท	m	()r	(T)	ım	Jn	51	ter	ınar	1(1	KII	. I	1.5	-1	IVI	ш

M-115	0.35x	C-Mount	pro	jection	lens

0.5x C-Mount projection lens M-114

#### M-118 0.75x C-Mount projection lens

C-Mount projection lens for APS-C/full frame reflex cameras (trino) Universal adapter for C-Mount projection lens (trino) M-173

M-699

0.35x focusable C-Mount adapter M-620

0.5x focusable C-Mount adapter M-620.1

M-620.2 0.65x focusable C-Mount adapter

M-620.3 1x focusable C-Mount adapter

### Miscellaneous

4 = 4 0 4			1.4.
15104	(   🗅 :	anino	ı bıt

DC-002 Plastic dust cover, medium, 490(l)x490(h) mm

M-005 Micrometric slide, 26x76mm, with 2 scales (1mm/100 & 10mm/100)

ST-092 Protective glass for stereohead

ST-201 Iris aperture diaphragm for darkfield

ST-203 Glass immersion cell

ST-204 Translating cell holder

ST-205 Vacuum pick-up (with electric vacuum pump)

ST-207 Iris aperture diaphragm for brightfield

Antibacterial surface treatment, only for newly purchased microscope AB-030





How to connect the cameras to our microscopes.

Please refer to the Adapter reference list on Digital section.



Inspection & Industrial

v.7.5-OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# MET SERIES



**Metallurgical Microscopes** 

## **Metallurgical** Microscopy

Metallography is the study of the physical structure and components of metals, by using microscopy. Many different microscopy techniques are used in metallographic analysis.

Prepared specimens should be examined with the unaided eye after etching to detect any visible areas that have responded to the etchant differently from the norm as a guide to where microscopical examination should be employed. Light optical microscopy (LOM) examination should always be performed prior to any electron metallographic (EM) technique, as these are more time-consuming to perform and the instruments are much more expensive.

Further, certain features can be best observed with the LOM, e.g., the natural color of a constituent can be seen with the LOM but not with EM systems. Also, image contrast of microstructures at relatively low magnifications, e.g., <500X, is far better with the LOM than with the scanning electron microscope (SEM), while transmission electron microscopes (TEM) generally cannot be utilized at magnifications below about 2000 to 3000X. LOM examination is fast and can cover a large area. Thus, the analysis can determine if the more expensive, more time-consuming examination techniques using the SEM or the TEM are required and where on the specimen the work should be concentrated.

### **Brightfield and darkfield microscopy**

Most LOM observations are conducted using bright-field (BF) illumination, where the image of any flat feature perpendicular to the incident light path is bright, or appears to be white. But, other illumination methods can be used and, in some cases, may provide superior images with greater detail. Dark-field microscopy (DF), is an alternative method of observation that provides high-contrast images and actually greater resolution than bright-field. In dark-field illumination, the light from features perpendicular to the optical axis is blocked and appears dark while the light from features inclined to the surface, which look dark in BF, appear bright, or "self-luminous" in DF. Grain boundaries, for example, are more vivid in DF than BF.

### Polarized light microscopy

Polarized light (PL) is very useful when studying the structure of metals with non-cubic crystal structures (mainly metals with hexagonal close-packed (hcp) crystal structures). If the specimen is prepared with minimal damage to the surface, the structure can be seen vividly in cross-polarized light (the optic axis of the polarizer and analyzer are 90 degrees to each other, i.e., crossed). In some cases, an hcp metal can be chemically etched and then examined more effectively with PL. Tint etched surfaces, where a thin film (such as a sulfide, molybdate, chromate or elemental selenium film) is grown epitaxially on the surface to a depth where interference effects are created when examined with BF producing color images, can be improved with PL. If it is difficult to get a good interference film with good coloration, the colors can be improved by examination in PL using a sensitive tint (ST) filter.

### **Differential interference contrast microscopy**

Another useful imaging mode is differential interference contrast (DIC), which is usually obtained with a system designed by the Polish physicist Georges Nomarski. This system gives the best detail. DIC converts minor height differences on the plane-of-polish, invisible in BF, into visible detail. The detail in some cases can be quite striking and very useful. If an ST filter is used along with a Wollaston prism, color is introduced. The colors are controlled by the adjustment of the Wollaston prism, and have no specific physical meaning, per se. But, visibility may be better.



### **Oblique illumination**

DIC has largely replaced the older oblique illumination (OI) technique, which was available on reflected light microscopes prior to about 1975. In OI, the vertical illuminator is offset from perpendicular, producing shading effects that reveal height differences. This procedure reduces resolution and yields uneven illumination across the field of view. Nevertheless, OI was useful when people needed to know if a second phase particle was standing above or was recessed below the plane-of-polish, and is still available on a few microscopes. OI can be created on any microscope by placing a piece of paper under one corner of the mount so that the plane-of-polish is no longer perpendicular to the optical axis.

## **B-383MET** - Metallurgical Microscope

Brightfield upright microscope with IOS W-PLAN MET objectives and metallurgical attachment combining the exclusive **X-LED³** lighting source both for incident and transmitted illumination. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.
























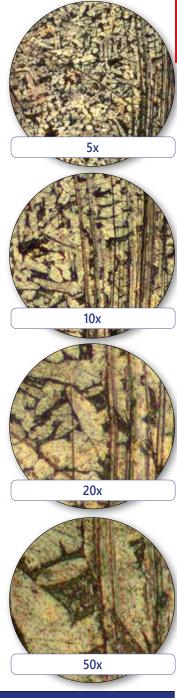

Part	Description
<b>Observation mode:</b>	Brightfield, incident polarized light.
Epi-illumination and polarizing filters:	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Field and aperture diaphragms, polarizer & analyzer filters.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 48 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/20 mm, high eye-point and secured by screw.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.

Part	Description				
Objectives (strain-free):	IOS W-PLAN MET 5x/0.12 IOS W-PLAN MET 10x/0.25 IOS W-PLAN MET 20x/0.40 IOS W-PLAN MET 50x/0.75 All with anti-fungus treatment.				
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range. With tempered glass plate.				
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.				
Condenser:	Abbe N.A. 1.25, with objective-coded iris diaphragm, focusable and centerable.				
Transmitted illumination (Fixed Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.				

## **B-510MET** - Metallurgical Microscope

Advanced routine laboratory microscope with IOS W-PLAN MET objectives and metallurgical attachment with the exclusive **X-LED³** lighting source for incident illumination only. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.























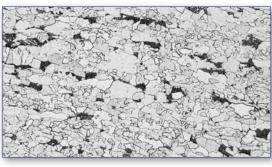



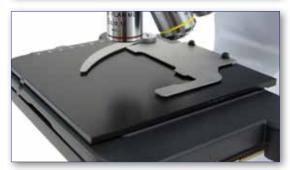



Part	Description
Observation mode:	Brightfield, simple polarized light, oblique illumination on incident light.
Epi-illumination and polarizing filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With aperture and field diaphragms, and oblique illumination system. With polarizer and analyzer. Multi-plug 100-240Vac/6Vdc external power supply.
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.
Interpupillary distance:	Adjustable between 50 and 75 mm.
Dioptric adjustment:	On the left eyepiece tube.
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.


Part	Description
Objectives (strain-free):	IOS W-PLAN MET 5x/0.12 IOS W-PLAN MET 10x/0.25 IOS W-PLAN MET 20x/0.40 IOS W-PLAN MET 50x/0.75 All with anti-fungus treatment.
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.


## **B-510METR** - Metallurgical Microscope

Advanced routine laboratory microscope with IOS W-PLAN MET objectives and metallurgical attachment with the exclusive **X-LED³** lighting source for both transmitted and incident illumination. The NCG (no cover glass) objectives are especially designed for microscopy use without a cover slip ideal for metallographic samples and other opaque specimens.





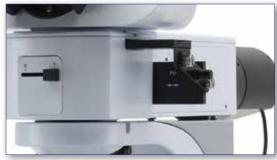

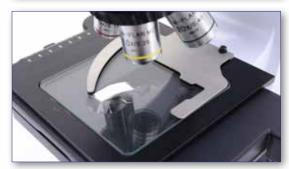



Part	Description	
Observation mode:	Brightfield on transmitted light. Brightfield, simple polarized light, oblique illumination on incident light.	
Epi-illumination and polarizing filters:	X-LED <sup>8</sup> with white 8 W LED (6.300 K) with brightness control. With aperture and field diaphragms, and oblique illumination system. With polarizer and analyzer.	
Head:	Trinocular (fixed 50/50), 30° inclined, 360° rotating.	
Interpupillary distance:	: Adjustable between 50 and 75 mm.	
Dioptric adjustment:	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	

Part	Description	
Objectives (strain-free):	IOS W-PLAN MET 5x/0.12 IOS W-PLAN MET 10x/0.25 IOS W-PLAN MET 20x/0.40 IOS W-PLAN MET 50x/0.75 All with anti-fungus treatment.	
Specimen stage:	Double layer rackless mechanical stage, 233x147 mm, 78x54 mm X-Y range. With tempered glass plate.	
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism with limit stop to prevent the contact between objective and specimen.	
Condenser:	Swing-out N.A. 0.2/0.9, with iris diaphragm, focusable and centerable.	
Transmitted illumination (Full Koehler type):	X-LED <sup>3</sup> with white 3.6 W LED (6.300 K) with brightness control. Multi-plug 100-240Vac/6Vdc external power supply.	

## **B-1000METBF** - Brightfield Metallurgical Microscope

The modular OPTIKA B-1000METBF offers superior quality **brightfield incident light**, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive X-LED<sup>8</sup> (8 W) transmitted illumination (Koehler system). The incident light relies on an incredibly bright 18 W LED illumination, designed by OPTIKA. B-1000 gives multiple options as manual or motorized configuration.



## **B-1000METBF** - Configuration Chart

Build the microscope that suites your needs by choosing among the components



M-1021M

Main body with focus system and X-LED8 illumination, for metallurgical model



M-1021M+M-1156\*+M-1149

Main body with motorized focus system and X-LED<sup>8</sup> illumination, for metallurgical model



M-1022M

Main body with focus system, for metallurgical model with incident light only



### M-1022M+M-1156\*+M-1149

Main body with motorized focus system, for metallurgical model with incident light only



M-1042

Sextuple reversed nosepiece, for RMS objectives with DIC slot



### M-1043+M-1156\*


Sextuple motorized reversed nosepiece, for RMS objectives with DIC slot





### IOS LWD U-PLAN MET

### IOS LWD U-PLAN F MET

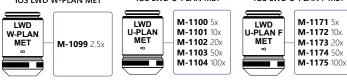

### M-1148

Mechanical stage with glass, for metallurgical model



### M-1147+M-1156\*

Motorized mechanical stage



### M-1154

0.70 N.A. swing-out MET condenser



<sup>\*</sup> Code M-1156 must be added only once for any motorized configuration

## **B-1000METDK** - Darkfield Metallurgical Microscope

The modular OPTIKA B-1000METDK offers superior quality **brightfield and darkfield incident light**, helping you working in a comfortable way during extended periods of use and performing reliable, accurate and rapid diagnosis benefiting from modularity, which gives the chance to create customized configurations tailored on customer needs. Versatile, robust, durable and sturdy, B-1000 offers premium quality optics, the state-of-the-art, exclusive **X-LED**<sup>8</sup> (8 W) transmitted illumination (Koehler system). The incident light relies on an incredibly bright 18 W LED illumination, designed by OPTIKA. B-1000 gives multiple options as manual or motorized configuration.



## **B-1000METDK** - Configuration Chart

Build the microscope that suites your needs by choosing among the components





M-1154 0.70 N.A. swing-out MET condenser



## IM-300METLD- Metallurgical Microscope

LED routine inverted microscope with IOS LWD U-PLAN MET objectives for material science and metallographic applications, combining a sturdy yet compact structure with dedicated components required in this field, like the NCG (no cover glass) objectives working without cover slide ideal for metallographic samples and other opaque specimens. A particularly simple and ingenious optical design allows stable alignments and smooth and accurate movements. This model is equipped with an 18W LED lighting system.



# IM-300METLD - Specifications



Part	Description	
Observation mode:	Brightfield, simple polarized light.	
<b>Epi-illumination and polarizing filters:</b> LED 18 W with brightness control.  With centrable aperture and field diaphragms.  With polarizer and 360° analyzer.  Supplied with blue (LBD) filter.		
Head:	Trinocular (2-position 100/0, 50/50), 45° inclined.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Diopter adjustment:	On the left eyepiece tube.	
Eyepieces:	WF10x/22 mm, high eye-point and with rubber cups.	
Nosepiece:	Quintuple revolving nosepiece, rotation on ball bearings.	

Part	Description
Objectives:	IOS LWD U-PLAN MET 5x/0.15 IOS LWD U-PLAN MET 10x/0.30 IOS LWD U-PLAN MET 20x/0.45 IOS LWD U-PLAN MET 50x/0.55 All with anti-fungus treatment.
Specimen stage:	Fixed stage, 250x160 mm, with round metal stage insert.
Focusing:	Coaxial coarse (adjustable tension) and fine focusing mechanism.





### 3

## **IM-5MET** - Metallurgical Microscope

Industrial and materials science inverted microscope especially designed for opaque specimens (including metals microstructure investigation and studies such as grain size, grain boundaries, phases, transformation, inclusions, and non-metals, as well as sample preparation and treatment) in metallography labs. Freely configurable lenses according to customer's preferences, FN 24 high eyepoint, infinity corrected optical system, coaxial focusing, mechanical stage, and epi-illumination attachment powered by halogen 12 V/100 W with brightness control. Sturdy and incredibly reliable, it is equipped with all the main controls in ergonomic position and with long lasting, efficient LED illumination to provide over 20 years of use.



## **IM-5MET** - Specifications



Part	Description	
Head:	Trinocular (split ratio: 100/0, 50/50), 45° inclined.	
Interpupillary distance:	Adjustable between 50 and 75 mm.	
Dioptric adjustment:	Both eyepieces.	
Eyepieces:	WF10x/24 mm, high eyepoint, secured by screw and with retractable rubber cups.	
Epi-illumination & filters:	Halogen 12 V/100 W with brightness control. With field and aperture diaphragms, polarizer and analyzer filters.	
Nosepiece:	Quintuple ball bearings revolving nosepiece, reversed.	
Objectives:	Selectable according to customer's preferences. All with anti-fungus treatment.	
Specimen stage:	Mechanical stage, 240x250 mm.	
Focusing:	Coaxial coarse and fine focusing mechanism with limit stop to prevent the contact between objective and specimen. Adjustable tension of coarse focusing knob.	

### IM-5MET is freely configurable in terms of objectives, by choosing among:

Included ■ Optional □

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1100	IOS LWD U-PLAN MET objective 5x/0.15	
M-1101	IOS LWD U-PLAN MET objective 10x/0.30	
M-1102	IOS LWD U-PLAN MET objective 20x/0.45	
M-1103	IOS LWD U-PLAN MET objective 50x/0.55	
M-1104	IOS LWD U-PLAN MET objective 100x/0.80 (dry)	

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, field flatness up to F.N. 25:		
M-1171	IOS LWD U-PLAN F MET objective 5x/0.15	
M-1172	IOS LWD U-PLAN F MET objective 10x/0.30	
M-1173	IOS LWD U-PLAN F MET objective 20x/0.50	
M-1174	IOS LWD U-PLAN F MET objective 50x/0.80	
M-1175	IOS LWD U-PLAN F MET objective 100x/0.90 (dry)	

MET Infinity-corrected Plan-Achromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:		
M-1094	IOS LWD U-PLAN MET BD objective 5x/0.15	
M-1095	IOS LWD U-PLAN MET BD objective 10x/0.30	
M-1096	IOS LWD U-PLAN MET BD objective 20x/0.45	
M-1097	IOS LWD U-PLAN MET BD objective 50x/0.55	
M-1098	IOS LWD U-PLAN MET BD objective 100x/0.80 (dry)	

MET Infinity-corrected Semi-Apochromatic, Long Working Distance objectives, for brightfield and darkfield, field flatness up to F.N. 25:		
M-1180	IOS LWD U-PLAN F MET BD objective 5x/0.15	
M-1181	IOS LWD U-PLAN F MET BD objective 10x/0.30	
M-1182	IOS LWD U-PLAN F MET BD objective 20x/0.50	
M-1183	IOS LWD U-PLAN F MET BD objective 50x/0.80	
M-1184	IOS LWD U-PLAN F MET BD objective 100x/0.90 (dry)	



v.7.5-OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

### **Headquarters and Manufacturing Facilities**

**OPTIKA® S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com


**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



# IS SERIES

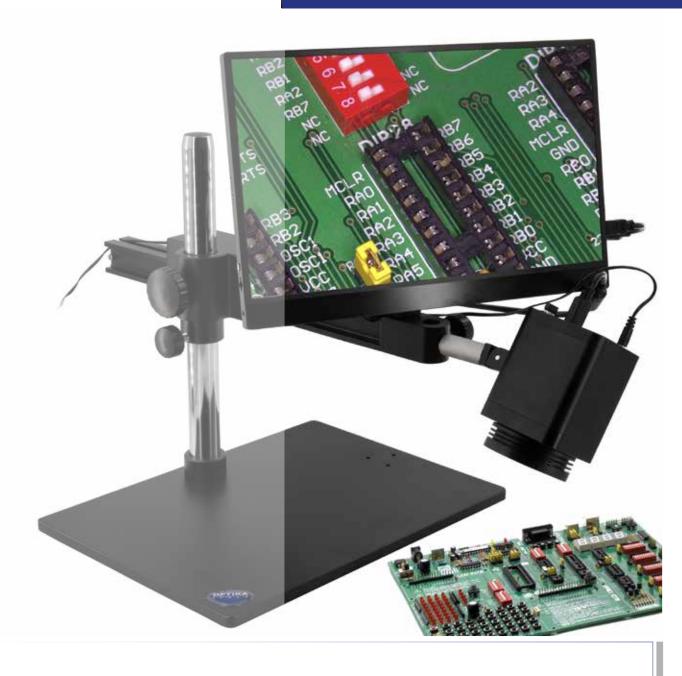

**Inspection Video Microscopes** 

## **IS-01** - Overview

OPTIKA **IS-01** is a compact, all-in-one video microscope with **integrated HDMI camera and built-in software**. Ideal for achieving up to 50X magnifications and performing measurements.

It's recommended when repeatability is required thanks to the click-stop mechanism at different magnifications steps.

With **60 frames** per second at full resolution and **1080 HD** image quality, you can see, on a plugged HDMI monitor, the same as what you see through the microscope eyepieces.



## **IS-01SMD** - Overview

OPTIKA **IS-01SMD** is a compact, all-in-one video microscope with **integrated HDMI camera and built-in software**. Ideal for achieving high magnifications and performing measurements.


It's recommended when repeatability is required thanks to the click-stop mechanismat different magnifications steps.

With **60 frames** per second at full resolution and **1080 HD** image quality, you can see, on a plugged HDMI monitor, the same as what you see through the microscope eyepieces. **The 3D rotational angle view attachment and the 7x to 50x zoom of IS-01SMD provide an unparalleled real-time 3D view ideal for PCB and Medical Devices inspection that is simple to use and remarkable in image quality.** 

Long working distance allows soldering rework and assembly by watching LCD screen. The Built-in camera is fitten with USB2.0 slots for storing image and video without any computer using a USB flash disk. Thanks to the measurement software, it is possible to perform dimensional measurements on screen.




## IS-01/IS-01SMD - Built-it Software Overview

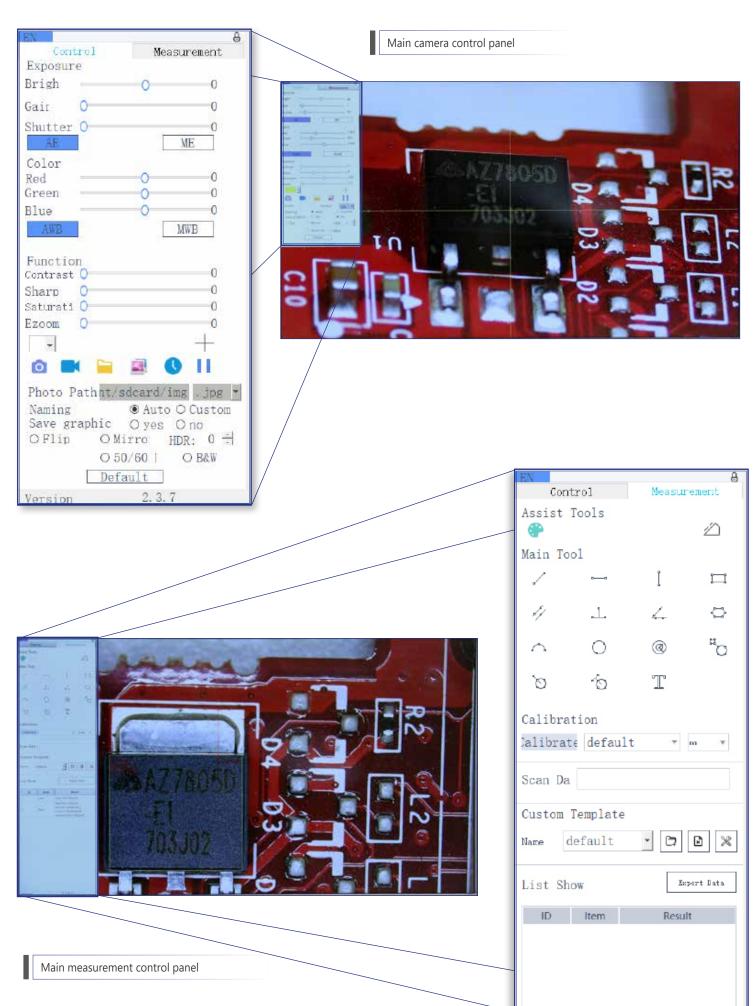


# **IS-01SMD** - Straight/3D view quick switching





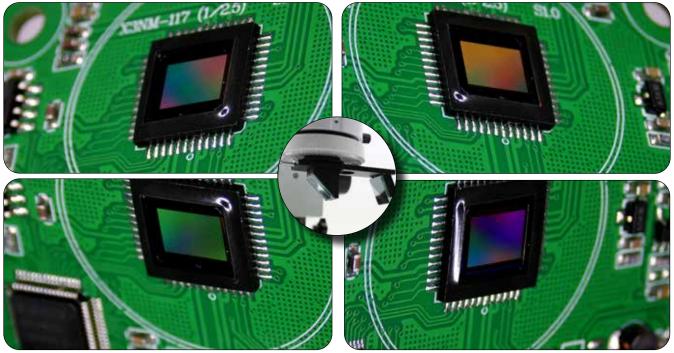


Captures of samples obtained by rotating the 3D attachment of the IS-01SMD model (optional on the IS-01 model)

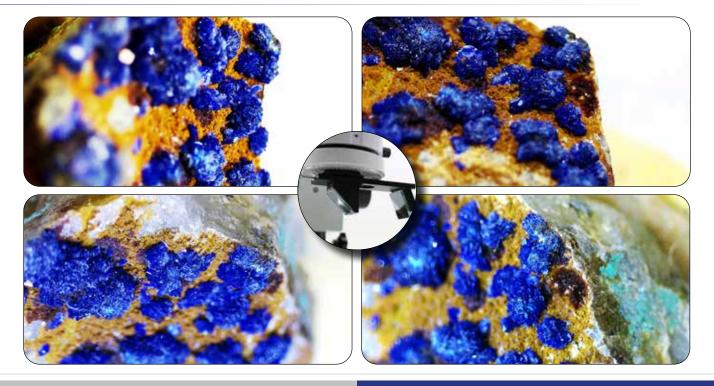


## **IS-01** - Straight View Zoom Factor


Straight View 0.7X

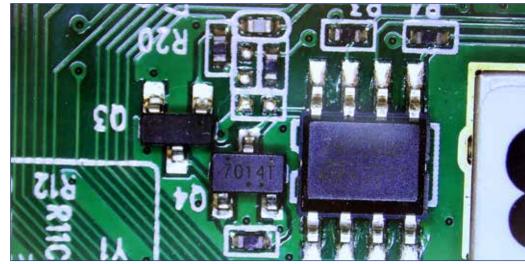

Straight View 2X

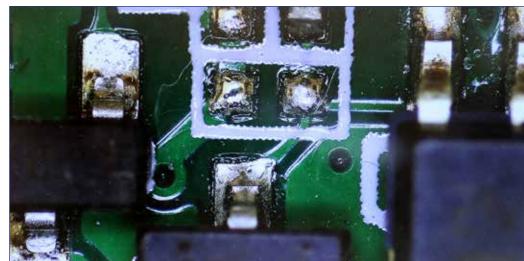

Straight View 5X

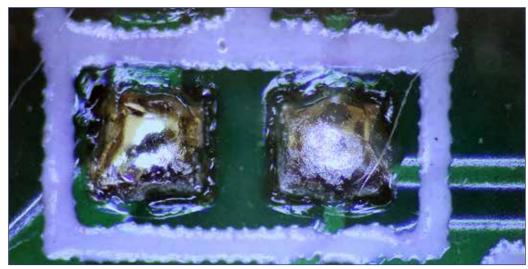



# IS-01/IS-01SMD - Technical Specifications

IS-01/IS-01SMD Main Body	
<b>Observation Method</b>	
Brightfield	Yes
Head	
Construction material	Plastic mold / Aluminum gears
Objective	
Optical system	160 mm
Anti-fungus treatment	Yes
Objective Type	Achromatic
Zoom type	Parfocal achromatic
Working distance (mm)	110
Standard magnifications	0.7x-5x
Zoom ratio	7.14
Zoom click stops	At 0.7x, 1x, 1.5x, 2x, 2.5x, 3x, 4x, 5x
Focusing System	
Туре	Coarse
Coarse total travel (mm)	50
Adjustable tension	Yes
Head holder internal diameter (mm)	76
Head	
Туре	Pillar
Pillar diameter (mm)	32

Incident Illumination	
Туре	LED
Light source power (W)	1
Illuminance (lux)	60000
<b>Brightness control</b>	Manual
Lifetime (hours)	> 65000
Temperature (K)	6300
Max. required power (W)	3
<b>Power Supply for Illumination</b>	
Туре	External
Microscope connector	Jack, 2.1 mm
Power plug type	Multi-plug (EU, UK, US)
Input voltage	110/240 Vac, 50/60 Hz
Output voltage	12 Vdc
<b>Additional Information</b>	
Maximum sample height (mm):	150
<b>Product Dimensions</b>	
Total height (mm)	510
Total width (mm)	320
Total depth (mm)	260
Weight	4.3 kg

### IS-01/IS-01SMD Digital Head

Resolution (MP)	2
Sensor type	CMOS
Sensor size	1/2.8"
Sensor aspect ratio	16/9
Full image size	1920x1080
Pixel size (µm)	2.9x2.9
Frame rate	60 fps (1920x1080)

S/N ratio (dB)	>50
Dynamic range (dB)	>50
Sensitivity (V/lux*s @550nm)	0.4
Exposure time (s)	0.02 - 1/10000
USB type	2.0
Software	built-in

IS-01SMD 3D Attachmen	t
Swing-out	Yes
360° Rotatable	Yes



## IS-4K2 - Overview

Advanced inspection system suitable for intensive use, ideal for video inspection in quality control of electronics, mechanics, and other industrial applications.

**IS-4K2** includes the **real-time full HD auto-focus camera** with **relevant zooming capalibilities** (optical zoom is 1x...18x).

Crystal clear 4K live view is shown trough a large HD monitor, with incredibly fast connection (30 fps). Its angle of view is **fully adjustable** and it enables **instant focus** in less than 1 sec. with no need to constantly adjust lens position. In addition, the working distance goes to infinity, making it the **recommended solution to inspect multi-layered objects**.

All functions are controlled directly from the screen via a wireless mouse. Images and videos can be saved on the SD card and easily transferred on any device. The system is completed by the SZ-STLX **boom stand** with a **special joint to enable any rotation** (transversal and longitudinal): the longitudinal movement is through the sliding of the horizontal rail (lockable); whilst also the overall height is both adjustable and lockable.



## **IS-4K2** - Technical Specifications



HDMI MONITOR	
Туре	LCD screen 360° rotating, freely inclinable for ergonomic vision
Size	15.6"
Power supply	5V / 4A
CAMERA	
<b>HDMI</b> camera resolution (MP)	8
Camera resolution (n° of pixels: W x H)	3840x2160
HDMI signal output	Yes
WiFi signal output	Yes
Sensor size	1/2.8"
Sensor technology	CMOS
Sensor type	SONY STARVIS
Optical zoom	1x-18x
Rolling shutter	Yes
Autofocus	Yes
Image format	16/9
Pixel size (mm)	1.45x1.45
Frame rate full resolution (fps)	30@3840x2160
Camera power	12V / 3A
Dimensions (mm)	80 x 80 x 116
Weight (Kg)	0,68

STAND	
Туре	Coarse simple overhanging
Coarse total travel (mm)	50
Adjustable tension	Yes
Head holder internal diameter (mm)	76
Pillar diameter (mm)	32
Maximum sample height (mm)	270 (from bench)
Total height (mm)	430
Base height (mm)	58
Base width (mm)	210
Base depth (mm)	255
Horizontal arm (mm)	790
Product weight (kg)	16.3
ACCESSORIES	
Accessories Included:	LED ringlight, HDMI cable (0,3m), Wireless Mouse, SD card and WiFi Adapter.

On-board controls and cross grid



### Variable working distance from 160 to 250mm.

By using a standard W.D. of 160-255 mm the camera can be calibrated and used for linear measurement regardless the zoom factor used. The system automatically "knows" the current zoom factor and performs the correct

Sliding horizontal rail for longitudinal positioning; adjustable longitudinal limit; adjustable and lockable height

High-performance boom stand with special joint for customized camera positioning

Wireless mouse and SD card are included

## IS-4K3 - Overview

Advanced inspection system suitable for intensive use, ideal for video inspection in quality control of electronics, mechanics, and other industrial applications.

IS-4K3 includes the real-time full HD auto-focus camera with relevant zooming capalibilities (optical zoom is 1x...18x).

Crystal clear 4K live view is shown trough a large HD monitor, with **incredibly fast connection** (30 fps). Its angle of view is **fully** adjustable and it enables instant focus in less than 1 sec. with no need to constantly adjust lens position. In addition, the working distance goes to infinity, making it the recommended solution to inspect multi-layered objects.


All functions are controlled directly from the screen via a wireless mouse.

Images and videos can be saved on the SD card and easily transferred on any device.

The system is completed by the **table clamping** SZ-STL5 **highly versatile flexible arm stand**, 360° rotating, ideal for high and large samples and on every bench thanks to its compact footprint, saving valuable space on the bench.





## **IS-4K3** - Technical Specifications

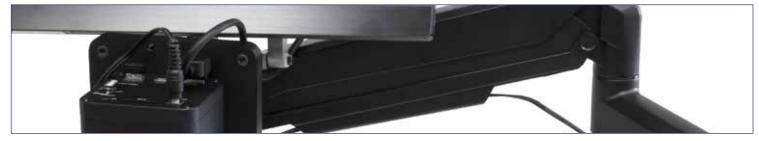

HDMI MONITOR	
Туре	LCD screen 360° rotating, freely inclinable for ergonomic vision
Size	15.6"
Power supply	5V / 4A
CAMERA	
HDMI camera resolution (MP)	8
Camera resolution (n° of pixels: W x H)	3840x2160
HDMI signal output	Yes
WiFi signal output	Yes
Sensor size	1/2.8"
Sensor technology	CMOS
Sensor type	SONY STARVIS
Optical zoom	1x-18x
Rolling shutter	Yes
Autofocus	Yes
Image format	16/9
Pixel size (mm)	1.45x1.45
Frame rate full resolution (fps)	30@3840x2160
Camera power	12V / 3A
Dimensions (mm)	80 x 80 x 116
Weight (Kg)	0,68
·	·

STAND	
Туре	Coarse, Pantograph, with table
Coarse total travel (mm)	50
Adjustable tension	Yes
Head holder internal diameter (mm)	76
Pillar diameter (mm)	35
Total height (mm)	327
Total width (mm)	300
Total depth (mm)	800
Horizontal arm (mm)	737
Product Weight (kg)	5.1
ACCESSORIES	
Accessories Included:	LED ringlight, HDMI cable (0,3m), Wireless Mouse, SD card and WiFi Adapter.

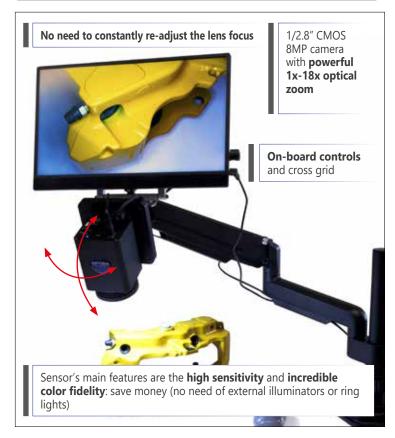

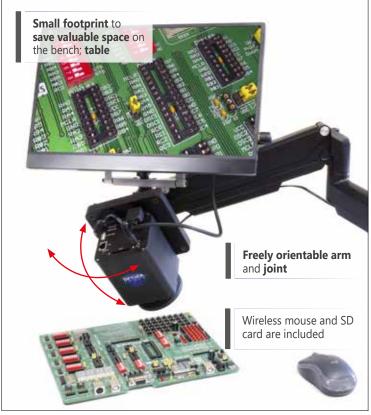

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA\* S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com OPTIKA<sup>®</sup> NSAth America OPTIKA<sup>®</sup> Central America OPTIKA<sup>®</sup> Africa

naa@optikasoopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com





## **CAMERAS & DIGITAL Solutions**

Cameras for Microscopy & Digital Scanner

## **CAMERAS & DIGITAL** Solutions

The most flexible way to create a digital microscope is to combine a digital camera to a standard microscope. You can then move the camera on different microscopes and different cameras on any microscope for an incredible versatility. Moreover, since microscope camera technology is improving so rapidly, it is easy to replace a camera as new ones are released.

OPTIKA offers an impressive amount of microscope camera, all of which are easy-to-use, affordable and professionally selected for excellent color resolution and rapid transmission. Downloadable, free of charge software is always available to enable the latest updates.

Applications include standard image capture and documentation for education, schools, home as well as professional image analysis for laboratories and industrial inspections, including very advanced solutions for critical applications like low light fluorescence imaging and material science applications.

### **Entry-Level Microscope Cameras**

page 357

### Smart & Affordable Solutions - Not Only For Schools & Education...

A wide range of affordable solutions with diversified live resolution for clear and crisp images to be combined via USB, HDMI or Wi-Fi.

### **Professional Microscope Cameras**

page 381

### Comprehensive Range, Remarkable Performance


An impressive offering compatible with any microscope brand to ensure the highest performance available for routine applications in professional environments, with USB, HDMI, 4K, Wi-Fi connections.

### **High Microscope Cameras**

page 407

### High Sensitivity Fluorescence Cameras

In order to detect the often low levels of fluorescence emitted by specimess, cameras used in fluorescence microscopy must have particular features, including high sensitivity and low noise, in order to capture as many photons as possible. These cameras are typically equipped with CCD sensor, although nowadays there is always a larger selection of scientific-grade CMOS.

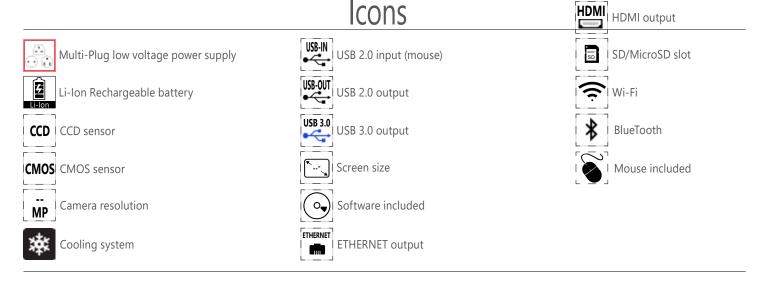

# Entry-Level Microscope Cameras



## Smart & Affordable Solutions Not Only For Schools & Education...

A wide range of affordable solutions with diversified live resolution for clear and crisp images to be combined via USB, HDMI or Wi-Fi ensuring smooth and productive teaching experiences on PC, Digital Imaging Units, projectors and other devices.

These models can also being used for basic requirements in laboratories adn industries.

Compatible with any microscope brand, thanks to the projection lens and rings (included in most of the cases).

## **Entry-Level Microscope Cameras**



### **E like Essential**


Simply the most essential eyepiece camera (C-E2). That's it!

**C-E2** E2 eyepiece camera, 2 MP CMOS, USB2.0


### **B like Basic - C-B Series**

The cost-effective OPTIKA C-B cameras are generally recommended for basic/general applications in education and home use. With C-B+ models, faster transmission speed is achieved through USB3.0 connection, being ideal especially on moving specimens.

C-B1B1 camera, 1.3 MP CMOS, USB2.0C-B3AB3 camera, 3.1 MP CMOS, USB2.0C-B5B5 camera, 5.1 MP CMOS, USB2.0C-B10+B10+ camera, 10 MP CMOS, USB3.0B18+ camera, 18 MP CMOS, USB3.0

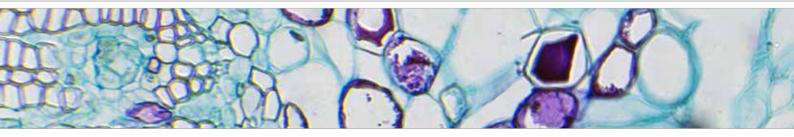

## **Entry-Level Microscope Cameras**

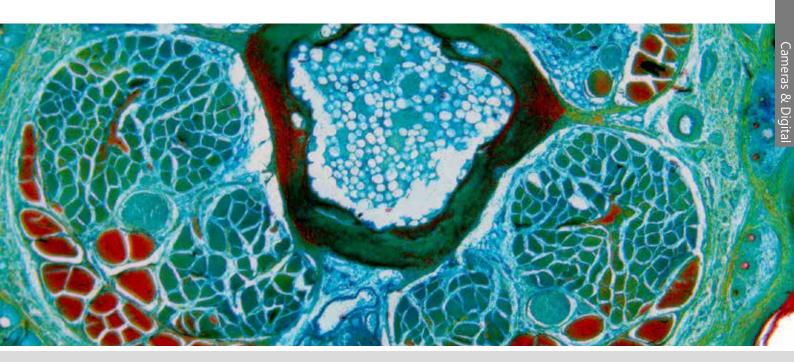

### **TB like Digital Imaging Unit - TB Series**

The combination of OPTIKA C-B cameras with Digital Imaging Unit for a completely new, revolutionary experience. Not a simple Digital Imaging Unit but a real PC with large touch screen for smooth and responsive control, representing an extremely comfortable solution for open discussions.

**TB-3WA** Digital Imaging Unit with B3A camera, 3.1 MP CMOS, USB2.0, EU **TB-5W** Digital Imaging Unit with B5 camera, 5.1 MP CMOS, USB2.0, EU

### **HE like HDMI Essential - C-HB & C-HBSC**

Recommended for its easy operation, no software installation is required with image captured on SD card. SC version includes a 11.5" Full HD screen for an all-in-one solution, with tilting features and saving space on the bench.


C-HB HB Camera, 1080p, 2 MP CMOS, HDMI, multi-plug

**C-HBSC** HB Camera, 1080p, 2 MP CMOS, HDMI, with screen, multi-plug

### WIFI Cameras - WiFi Series

The most flexible and versatile camera on the market! Allowing the camera to be moved from one microscope to another, whilst transferring the live view on any device (using Windows, Android or IOS).

**C-WF2** WF camera, 1 MP CMOS, Wi-Fi, multi-plug

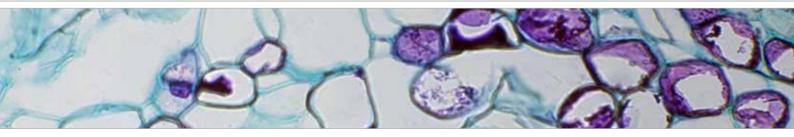

## **Essential**

## C-E2



Simply the most essential, user-friendly and handy eyepiece camera for low budget, with 2 MP resolution, CMOS sensor and USB2.0 connection.

Ready to use on any microscope with direct eyepiece tube connection, thanks to the rings included. Connectable also on trinocular tube.

Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.

### Model:

C-E2: E2 eyepiece camera, 2 MP CMOS, USB2.0

	C-E2
Digital camera resolution	2 MP
Analog camera resolution	NO
Signal output	USB 2.0
Audio Signal	NO
Sensor Size	1\3.2"
Sensor technology	CMOS
Image format	4\3
Full Image size	1600 x 1200
Frame rate full resolution	5 fps (1600x1200) / 7,5 fps (1280x1024) / 20fps (800x600) / 22fps (640x480)
Max Exposure time	Auto
ON board Memory	NO
<b>External Memory Card</b>	NO
External camera power	PC USB
White Balance	Auto
Gain Control	Auto
Back light control	Auto
Exposure control	Auto
C-Mount connection	NO
CS-Mount connection	NO
Arm length	-
8mm objective	NO

### Accessories included:

C-E2: 30 mm / 30.5 mm ring adapters, 1.8 m USB cable.

### **Essential**




- » Simple operation, driver-free
- » Universal connection to any microscope brand
- » Direct eyepiece & trino port connection
- » Crisp 2 MP images
- » High frame rate
- » Reliable color fidelity
- » Rings included
- » USB cable included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux



# Cameras & Digital

### C-B Series

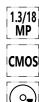



adapter.



Cost-effective, user-friendly cameras with several resolutions (from 1.3 up to 18 MP), CMOS sensor and USB2.0 or USB3.0 connection, recommended for basic/general applications in education and home use especially on moving specimens. Ready to use on any microscope with direct eyepiece tube connection, thanks to the C-mount projection lens and rings included. Connectable also to all the trinocular tube of different brands using the C-mount projection lens included or additional focusable C-Mount

Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.



- » Simple operation, driver-free
- » Universal connection to any microscope brand
- » Direct eyepiece & trino port connection
- » Crisp 1.3 MP crisp images
- » Reliable color fidelity
- » C-Mount projection lens and rings included
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux





### C-B Series

#### **C-B Models:**

**C-B1:** B1 camera, 1.3 MP CMOS, USB2.0 **C-B3:** B3 camera, 3.1 MP CMOS, USB2.0 **C-B5:** B5 camera, 5.1 MP CMOS, USB2.0

#### C-B+ Models:

**C-B10+:** B10+ camera, 10 MP CMOS, USB3.0 **C-B18+:** B18+ camera, 18 MP CMOS, USB3.0





### **C-B Series** - Specifications

	C-B1	C-B3A	C-B5
Digital camera recolution	1.3 MP (1280 x 1024)	3.1 MP (2048 x 1536)	
Digital camera resolution			5.1 MP (2592 x 1944)
Signal output	USB 2.0	USB 2.0	USB 2.0
Sensor Size	1/3"	1/2.8"	1/2.5″
Sensor technology	CMOS	CMOS	CMOS
Sensor type	Aptina CMOS	Aptina CMOS	Aptina CMOS
lmage format	5/4	4/3	4/3
Pixel size	3.6 x 3.6 µm	2.5 x 2.5 μm	2.2 x 2.2 μm
Frame rate full resolution	15 fps (1280 x 1024)	10.5 fps (2048x1536)	7 fps (2592 x 1944)
Frame rate other resolutions	50 fps (320 x 256)	15 fps (1920x1080)	27 fps (1280 x 960); 90fps (640 x 480)
Sensitivity	1 V/lux-second	600mV at 1/30sec	0.53 V/lux-second
Signal / noise ratio	44 dB	40.5 dB	40.5 dB
Dynamic range	74 dB	66.5 dB	66.5 dB
ADC conversion	8 Bit	8 Bit	8 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit ; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.14 msec - 2 sec	0.244 msec - 2 sec	0.294 msec - 2 sec
Binning	1x1; 2x2; 4x4	1x1; 2x2; 3x3	1x1; 2x2; 4x4
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter
Camera power	PC USB	PC USB	PC USB
C-mount	YES	YES	YES

### **C-B Series**



# C-B Contents: USB camera USB 2.0 cable 30 & 30.5 mm rings C-mount projection lens Micrometric slide





C-B10+	C-B18+
10 MP (3584 x 2748)	18 MP (4912 x 3684)
USB 3.0	USB 3.0
1/2.3"	1/2.3"
CMOS	CMOS
Aptina CMOS	Aptina CMOS
4/3	4/3
1.67 x 1.67 µm	1.25 x 1.25 μm
7.2 fps (3584 x 2746)	5.6 fps (4912 x 3684)
24.5 fps (1792 x 1372);	18.1 fps (2456 x 1842); 32.2 fps (1228 x 922)
0.31 V/lux-second	0.62 V/lux-second
34 dB	36.3 dB
65.2 dB	65.8 dB
8 Bit - 12 Bit	8 Bit - 12 Bit
1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
0.4 msec - 2 sec	0.1 msec - 2 sec
1x1; 2x2; 4x4	1x1; 2x2; 4x4
380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
PC USB	PC USB
YES	YES

Cameras & Digital

### **TB Series**



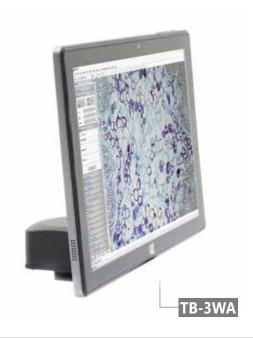

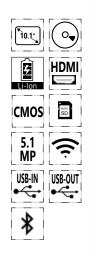




Exclusive, powerful Digital Imaging Unit combined to a CMOS sensor cameras with USB2.0 connection, recommended for discussion groups and educational purposes thanks to the easy operation, space-saving features and unparalleled comfort.

The unique holding solution for open discussion is 360° rotating and tilting for any adjustement, whilst the large touch screen provides fast, responsive and smooth control.

At any time, the Digital Imaging Unit can be easily detached to be used as a portable device.

Non-stop operation is granted by the simultaneous camera and power connection for long-term use and class/lesson alternation. Connectable also to all the trinocular tube of different brands using the projection lens included or additional focusable C-Mount adapter. Downloadable, free of charge software is always available to enable the latest updates.



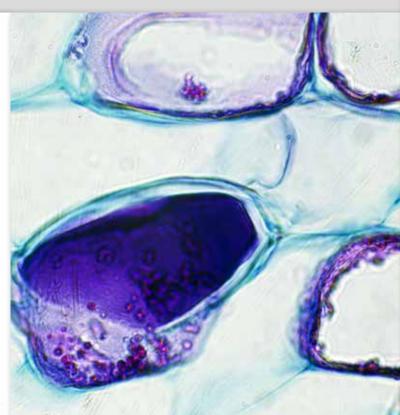

### **TB Series**

- » External digital camera connected to Digital Imaging Unit
- » Large touch screen with fast, responsive and smooth control
- » Easily detachable, can be used as a laptop
- » A 2-in-1 solution that you can use like a PC, being Windows-based
- » Simultaneous camera and power connection for long-term operation
- » Powerful Intel processor ensuring top performance and speed
- » High-resolution, vivid color graphic display
- » Crisp 3.1 or 5.1 MP images
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (direct or via dedicated adapter, on any microscope)
- » C-mount projection lens and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView for Windows






### **TB Series**

#### **Digital Imaging Unit Models:**

TB-3WA: Digital Imaging Unit with B3 camera, 3.1 MP CMOS, USB2.0, EU TB-5W: Digital Imaging Unit with B5 camera, 5.1 MP CMOS, USB2.0, EU

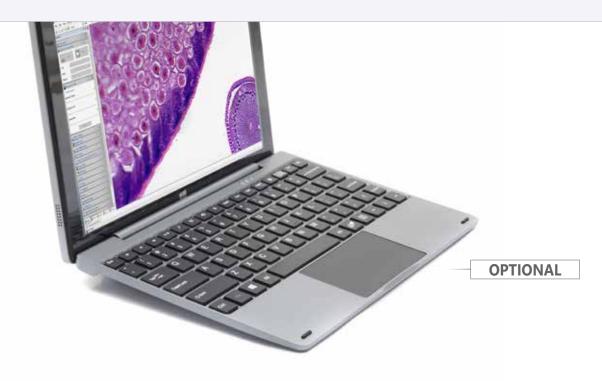

### **TB Series** - Specifications

AMERA TECHNICAL SPECIFICATIONS	TB-3WA	TB-5W
Digital camera resolution	3.1 MP (2048 x 1536)	5.1 MP (2592 x 1944)
ignal output	USB 2.0	USB 2.0
ensor Size	1/2"	1/2.5″
ensor technology	CMOS	CMOS
ensor type	Aptina CMOS	Aptina CMOS
mage format	4/3	4/3
ixel size	3.2 x 3.2 µm	2.2 x 2.2 μm
rame rate full resolution	12 fps (2048 x 1536)	7 fps (2592 x 1944)
rame rate other resolutions	32 fps (1024 x 768); 45 fps (680 x 510)	27 fps (1280 x 960); 90fps (640 x 480)
ensitivity	1 V/lux-second	0.53 V/lux-second
gnal / noise ratio	43 dB	40.5 dB
ynamic range	61 dB	66.5 dB
DC conversion	8 Bit	8 Bit
olor Depth	1 Bit ; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
xposure Time	0.244 msec - 2 sec	0.294 msec - 2 sec
inning	1x1; 2x2; 3x3	1x1; 2x2; 4x4
R filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
amera power	PC USB	PC USB
-mount	YES	YES

### **TB Series**



#### **TB Contents:**

USB camera
Digital Imaging Unit
C-mount projection lens
Micrometric slide
OTG cable (micro USB-C to USB-A)
USB cable USB-B to USB-A (0.5m)
External power supply

### **TB Series** - Digital Imaging Unit specifications

DIGITAL IMAGING UNIT TECH	HNICAL SPECIFICATIONS	
Operating system	Windows 11	
CPU	Intel® Celeron	
CPU speed	1.10 GHz	
Graphic card	Intel® HD Graphics 500	
RAM	Ram 6 GB DDR3	
Display size	LED 10.1" IPS Multi Touch Screen	
Display resolution	1920x1200	
Storage	Hdd 128 GB (EMMC)	
Network	802.11b/g/n/ac,2.4GHz	
Input ports	USB 3.0, Type C, Micro SD card reader, Microphone, DC in	
Output ports	Headphone - Micro HDMI	
Battery Type	Lithium-ion	
Battery capacity	3500 mAh	
Power consumption	24.05W	
Power supply	12V 2A EU / US	
Dimensions (mm)	247 x 165 X 9	
Weight (Kg)	0.57	
Language	Multilanguage	



### <sup>4</sup> C-HB & C-HBSC





Cost-effective, entry-level HD camera with 1080p, 2 MP resolution, CMOS sensor and HDMI connection, recommended for its easy operation, no software installation is required with image capturing on SD card.


Connection into the eyepiece tube or to the trinocular tube of any microscope brand via dedicated adapter (to be purchased separately). Wireless mouse, SD card and built-in software included.



### **C-HB & C-HBSC**

#### C-HB

- » Simple operation, built-in software
- » Crisp 2 MP images
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (direct or via dedicated adapter, on any microscope)
- » Mouse and SD card included

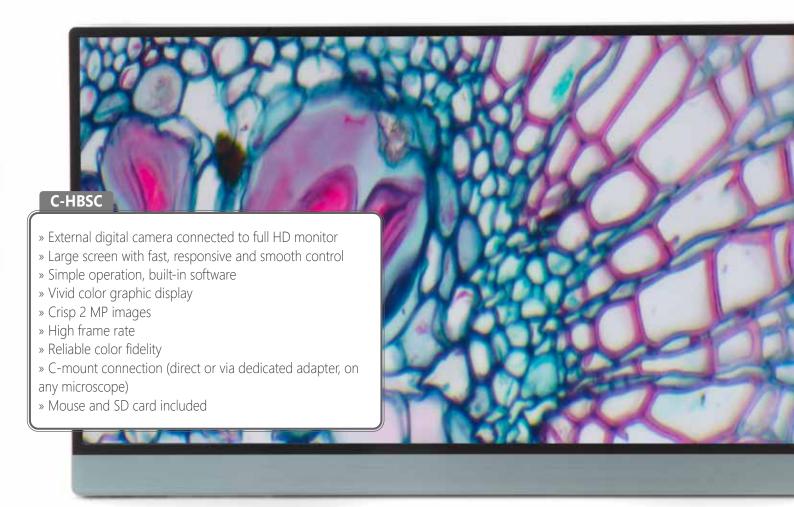

### <sup>4</sup> C-HB & C-HBSC

#### **C-HB Models:**

C-HB: HB camera, 1080p, 2 MP CMOS, HDMI, multi-plug

C-HBSC: HB camera, 1080p, 2 MP CMOS, HDMI, with screen, multi-plug



### C-HB & C-HBSC - Specifications

	С-НВ	С-НВЅС
Video resolution (USB output)	-	-
Video resolution (HDMI output)	HD 1080p	HD 1080p
Digital camera resolution	2 MP (1280 x 720)	2 MP (1280 x 720)
Signal output	HDMI	HDMI
Sensor Size	1/2.8"	1/2.8"
Sensor technology	CMOS	CMOS
Sensor type	sony starvis	SONY STARVIS
Image format	16/9	16/9
Pixel size	2.9 x 2.9 μm	2.9 x 2.9 μm
Frame rate (HDMI)	60@1920X1080	60@1920X1080
Sensitivity	1300 mV at 1/30sec	1300 mV at 1/30sec
Dark Signal	0.15mV at 1/30sec	0.15mV at 1/30sec
Exposure Time	0.01 msec - 1 sec	0.01 msec - 1 sec
Binning	1x1	1x1
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Camera power	DC 12V 1A	DC 12V 1A
C-mount	YES	YES
White balance	Auto/Manual//ROI	Auto/Manual//ROI
Live HDMI measurement	Yes	Yes
Exposure control	Auto / Manual	Auto / Manual

#### **C-HB Contents:**

HDMI camera

HDMI cable

SD card

Wireless mouse

Multi-plug external power supply

Calibration slide



#### **C-HBSC Contents:**

HDMI camera **HDMI** monitor HDMI cable

SD card

Wireless mouse

Multi-plug external power supply (2 pcs.)

calibration slide




### C-HB & C-HBSC - Monitor Specifications

MONITOR TECHNICAL SPECIFICATIONS		
Size	11,5"	
Power supply	12V / 2,5A	
HDMI cable	150 cm	





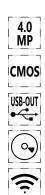




Cost-effective and user-friendly Wi-Fi camera with good resolution, CMOS sensor and Wi-Fi connection, recommended for basic/general applications in education and home use.

Direct Wi-Fi as no router is required thanks to the direct remote application for simplified use.

Ready to use on any microscope with direct eyepiece tube connection, thanks to the C-mount projection lens and rings included. Connectable also to all the trinocular tube of different brands using the C-mount projection lens included or additional focusable C-Mount adapter.

Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.

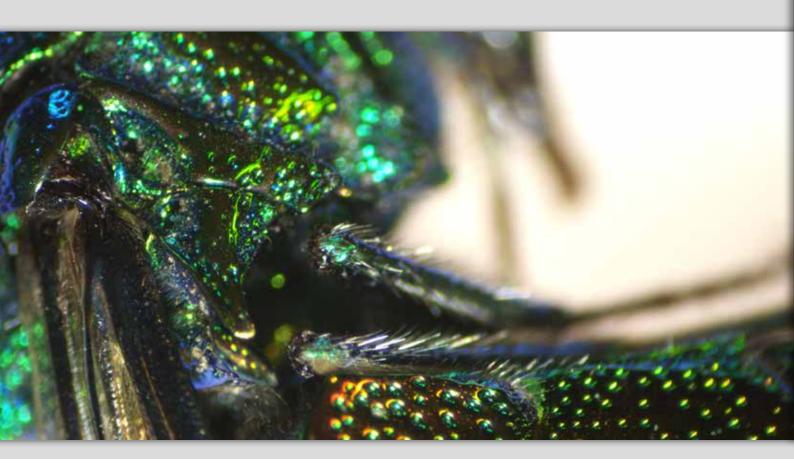



- » Simple, intuitive operation with powerful software
- » Universal connection to any microscope brand
- » Direct eyepiece & trino port connection
- » Wi-Fi interface
- » Direct Wi-Fi connection (no router is required)
- » Image and video capturing function when used in Wi-Fi mode
- » Reliable colour fidelity
- » C-mount projection lens and rings included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux





#### WiFi Series model:

C-WF2: WF camera, 4 MP CMOS, Wi-Fi, multi-plug



### WiFi Series - Specifications

CAMERA TECHNICAL SPECIFICATIONS	C-WF2
Camera sensor resolution	2688x1512
Signal output	USB 2.0
Sensor Size	1/1.8″
Sensor technology	CMOS
Sensor type	SONY EXMOR
Image format	16/9
Pixel size	2.9x2.9 mm
Frame rate full resolution	30fps@2688*1512
Frame rate other resolutions	30@2688*1512
Sensitivity	1 V/lux-second
Signal / noise ratio	43 dB
Dynamic range	61 dB
ADC conversion	8 Bit - 12Bit
Color Depth	8 Bit
Exposure Time	Auto
Binning	1x1
IR filter	380-650nm (IR CUT)
Camera power	Multiplug 100-240Vac/6Vdc external power supply
C-mount	YES
-	

#### **C-WF Contents:**

Wi-Fi camera
30 & 30.5 mm rings
C-mount projection lens
Micrometric slide
Multi-plug external power supply





### Recommended Camera Adapters

		Upright			
		Monocular Binocular (Ø 23 mm)	Trinocular (Ø 23 mm)	Binocular (Ø 30 mm)	Trinocular
Camera model	Sensor size	Ecovision / B-60 / B-150 B-190-290 / B-380 (ALC)	B-190 / B-290 B-380 (with M-699)	B-510 / B-810 / B-1000	B-380 / B-510 B-810 / B-1000
C-B1	1/3″	Included with the camera	Included with the camera	Included with the camera	M-620
C-B3A	1/2.8"	Included with the camera	Included with the camera	Included with the camera	M-620
C-B5	1/2.5"	Included with the camera	Included with the camera	Included with the camera	M-620.1
C-B10+	1/2.3"	Included with the camera	Included with the camera	Included with the camera	M-620.1
C-B18+	1/2.3"	Included with the camera	Included with the camera	Included with the camera	M-620.1
TB-3WA	1/2.8"	-	Included with the camera	M-115 + M-113.1	M-620
TB-5W	1/2.5"	-	Included with the camera	M-114 + M-113.1	M-620.1
С-НВ	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620
C-HBSC	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620
C-WF2	1/2.5"	Included with the camera	Included with the camera	M-114 + M-113.1	M-620.1









### **Recommended Camera Adapters**

v7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities


OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com

### Professional Microscope Cameras

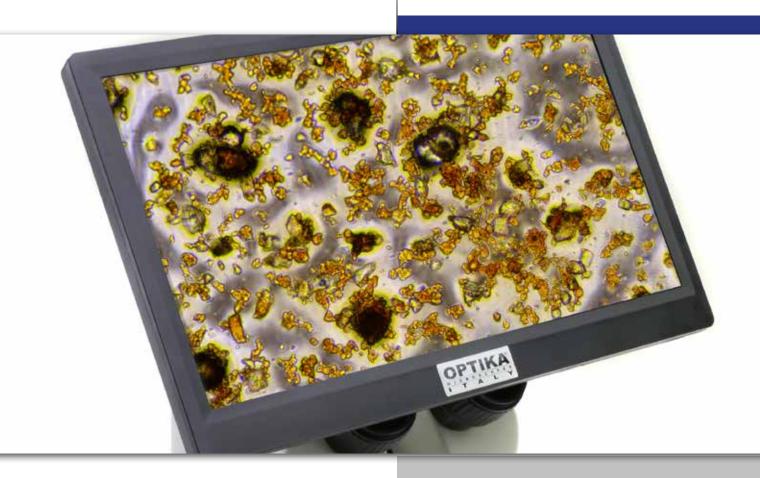


#### Professional Microscope Cameras Comprehensive Range, Remarkable Performance

An impressive offering compatible with any microscope brand to ensure the highest performance available for routine applications in professional environments, with USB, HDMI, 4K, Wi-Fi connections. Images and videos will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

The compact and elegantly designed housing conceals the very latest in camera technology.

### **Professional Microscope Cameras**

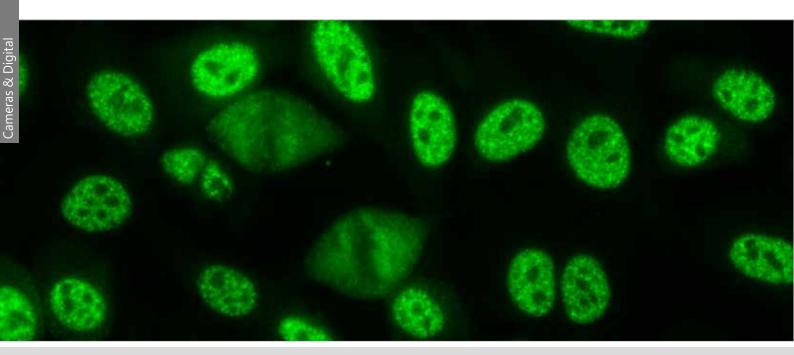

#### P as Professional

Recommended for professional use in laboratory and industrial field, the valuable yet affordable OPTIKA C-P cameras are equipped with topclass SONY EXMOR sensors and USB3.0 connection for premium features and faster transmission speed. Perfect for most of the brightfield, darkfield, phase contrast and metallographic applications when requiring PC/Laptop operation.

**C-P3** P3 Pro camera, 3.1 MP CMOS, USB3.0 P6 Pro camera, 6.3 MP CMOS, USB3.0

**C-P6AR** P6 Pro camera, 6.3 MP CMOS, USB3.0, with AR glass for IR applications

**C-P8** P8 Pro camera, 8.3 MP CMOS, USB3.0 P20 Pro camera, 20 MP CMOS, USB3.0

#### **GS like Global Shutter**

Global shutter camera are designed when particularly high-motion captures are needed, being perfect for enabling 'freeze frame' of fast changing events, exposing each and every pixel simultaneously on PC/Laptop.

Superb with moving specimens in brightfield, darkfield, phase contrast and even in polarized light observations thanks to the generous dynamic range which gives a great response to light and dark at the same time.


**C-P5GS** P5GS Pro global shutter camera, 5 MP CMOS, USB3.0

#### WH like Wi-fi & HDMI

The most versatile cameras with endless possibilities perfect for most of the brightfield, darkfield, phase contrast and material science applications. Take benefit from the on-board imaging software to display live view directly on monitor and projector, with data storage on SD card; or download the professional imaging software for PC.

**C-WH5** WH5 camera, 1080p, 5 MP CMOS, Wi-Fi/HDMI, multi-plug

**C-WH5SC** WH5 camera, 1080p, 5 MP CMOS, Wi-Fi/HDMI, with screen, multi-plug

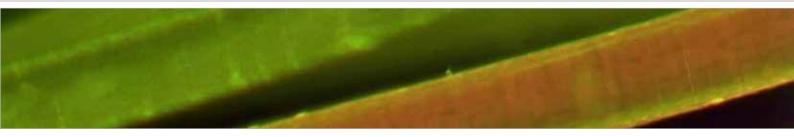

### **Professional Microscope Cameras**



#### **HA like HDMI Autofocus - C-HA**

Providing precise and ultra-fast automatic focus adjustment in any condition and in real time, perfect to compensate the lack of parfocality of the microscope or poor sample preparation, without any user effort.

Ideal for every use in brightfield, darkfield, phase contrast, polarized light and material science applications for easy connection to any monitor and projector.

**C-HA** HA autofocus camera, 2 MP CMOS, HDMI, multi-plug

#### 4K like Ultra-HD - C-HP4 / C-HUB4K

This is much more detailed than anything you're likely to have seen before: an immersive experience is created, improving detail and general sharpness. Ideal for every use in brightfield, darkfield, phase contrast, polarized light and material science applications for easy connection to Ultra-HD monitors, providing more pixels, more details and colors.

**C-HP4** HP4 camera, 4K/8MP CMOS, HDMI/USB

**C-HP4SC** HP4 camera, 4K/8MP CMOS, HDMI/USB, with 4K monitors

**C-HUB4K** HUB4K Camera,8Mp CMOS, 4K/USB/ETHERNET/WIFI multi video output

**(4**)

### **C-P Series**










Professional yet very easy to use cameras with several resolutions (from 3.1 up to 20 MP), large SONY EXMOR CMOS sensor and USB3.0 connection, recommended for general scientific or industrial purposes requiring rapid speed transmission.

The compact and elegantly designed housing conceals the very latest in camera technology.

Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.

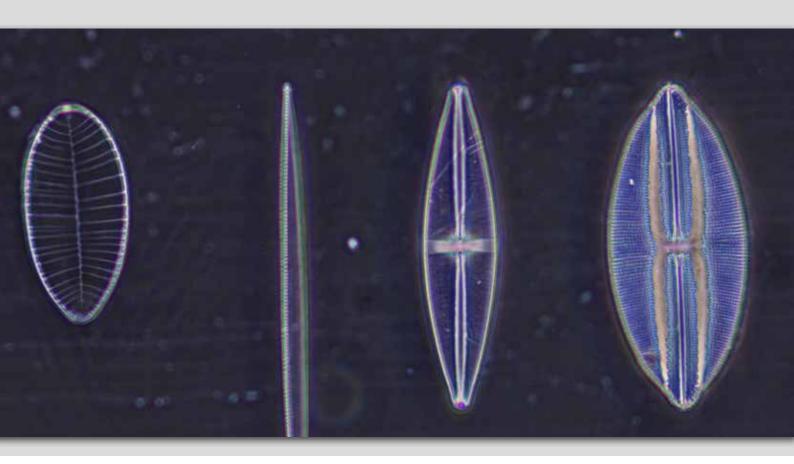

- » Simple operation, driver-free
- » Top-class SONY EXMOR sensor
- » Universal connection to any microscope brand
- » Crisp 3.1 to 20 MP images
- » USB3.0 for impressive high frame rate
- » Incredibly accurate colors
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux"




### C-P Series

#### **C-P Models:**

C-P3: P3 Pro camera, 3.1 MP CMOS, USB3.0 C-P6: P6 Pro camera, 6.3 MP CMOS, USB3.0

C-P6AR: P6 Pro camera, 6.3 MP CMOS, USB3.0 (AR glass)

C-P8: P8 Pro camera, 8.3 MP CMOS, USB3.0 C-P20: P20 Pro camera, 20 MP CMOS, USB3.0

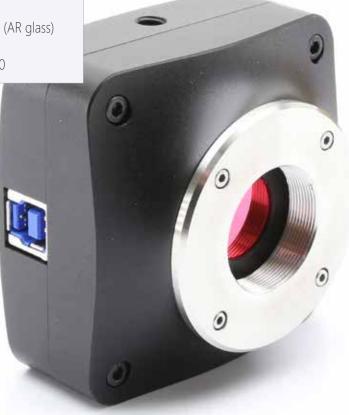

### **C-P Series** - Specifications

	C-P3	C-P6	C-P6AR	C-P8
Digital camera resolution	3.1 MP (2048 x 1536)	6.3 MP (3072 x 2048)	6.3 MP (3072 x 2048)	8.3 MP (3840 x 2160)
Signal output	USB 3.0	USB 3.0	USB 3.0	USB 3.0
Sensor Size	1/2.8″	1/1.8"	1/1.8"	1/1.8"
Sensor technology	CMOS	CMOS	CMOS	CMOS
Sensor type	SONY EXMOR	SONY EXMOR	SONY EXMOR	SONY EXMOR
Image format	4/3	3/2	3/2	16/9
Pixel size	2.5 x 2.5 μm	2.4 x 2.4 µm	2.4 x 2.4 µm	1.62 x 1.62 µm
Frame rate full resolution	50 fps (2048 x 1536)	30 fps (3072 x 2048)	30 fps (3072 x 2048)	32 fps (3840 x 2160)
Frame rate other resolutions	50 fps (1920 x 1080)	38 fps (1536 x 1024)	38 fps (1536 x 1024)	65 fps (1920 x 1080)
G Sensitivity	600mV at 1/30s	425mV at 1/30s	425mV at 1/30s	236mV at 1/30s
Dark Signal	0.15mV at 1/30s	0.15mV at 1/30s	0.15mV at 1/30s	0.1mV at 1/30s
Dynamic range	66 dB	66.8 dB	66.8 dB	65 dB
ADC conversion	8 Bit - 12Bit	8 Bit - 12Bit	8 Bit - 12Bit	8 Bit - 12Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit	1 Bit ; 4 Bit; 8 Bit; 24 Bit	1 Bit ; 4 Bit; 8 Bit; 24 Bit	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.1 msec - 15 sec	0.1 msec - 15 sec	0.1 msec - 15 sec	0.244 msec - 15 sec
Binning	1x1	1x1; 2x2	1x1; 2x2	1x1; 2x2
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)	380-1050 nm (AR glass)	380-650 nm (IR-cut filter)
Camera power	PC USB	PC USB	PC USB	PC USB
C-mount	YES	YES	YES	YES

### **C-P Series**



#### **C-P Contents:**

USB camera USB 3.0 cable



#### C-P20

20 MP (5440 x 3648)

USB 3.0

1"

CMOS

SONY EXMOR

3/2

2.4 x 2.4 µm

15 fps (5440 x 3648)

50 fps (2736 x 1824); 60 fps (1824 x 1216)

462mV at 1/30s

0.21mV at 1/30s

66.3 dB

8 Bit - 12Bit

1 Bit; 4 Bit; 8 Bit; 24 Bit

0.1 msec - 15 sec

1x1; 2x2; 3x3

380-650 nm (IR-cut filter)

PC USB

YES



### **C-PGS Model**











Professional yet very easy to use cameras with several resolutions, large SONY EXMOR CMOS sensor, USB3.0 connection and Global Shutter mode. Global Shutter mode can be easily thought of as a 'Snapshot' exposure mode, perfect for capturing images of moving objects and enabling 'freeze frame' capture of fast changing events, exposing each and every pixel simultaneously.

In addition, the generous dynamic range gives a great response to light and dark simultaneously, being recommened for polarizing light applications. The compact and elegantly designed housing conceals the very latest in camera technology.

Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.



- » Simple operation, driver-free
- » Top-class SONY EXMOR sensor
- » Universal connection to any microscope brand
- » Crisp 5MP images
- » Global Shutter technology for ""freeze frame"" capture
- » Generous dynamic range, recommened for polarizing light
- » USB3.0 for impressive high frame rate
- » Crystal-clear images, even for very short exposures
- » Outstanding noise characteristics, even in low lighting conditions
- » Incredibly accurate colors
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux"



### C-PGS Model

#### **C-PGS Model:**

C-P5GS: P5GS Pro global shutter camera, 5 MP CMOS, USB3.0



### **C-PGS Model** - Specifications


	C-P5GS
Digital camera resolution	5 MP (2448 x 2048)
Signal output	USB 3.0
Sensor Size	2/3"
Sensor technology	CMOS
Sensor type	SONY EXMOR
Image format	5/4
Pixel size	3.45 x 3.45 μm
Frame rate full resolution	35 fps (2448 x 2048)
Frame rate other resolutions	50 fps (1224 x 1024)
G Sensitivity	1146mV at 1/30s
Dark Signal	0.15mV at 1/30s
Dynamic range	70.60 dB
ADC conversion	8 Bit - 12 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.1 msec - 15 sec
Binning	1x1;
IR filter	380-650 nm (IR-cut filter)
Camera power	PC USB
C-mount	YES

Cameras & Digital

### **C-PGS** Contents:

USB camera



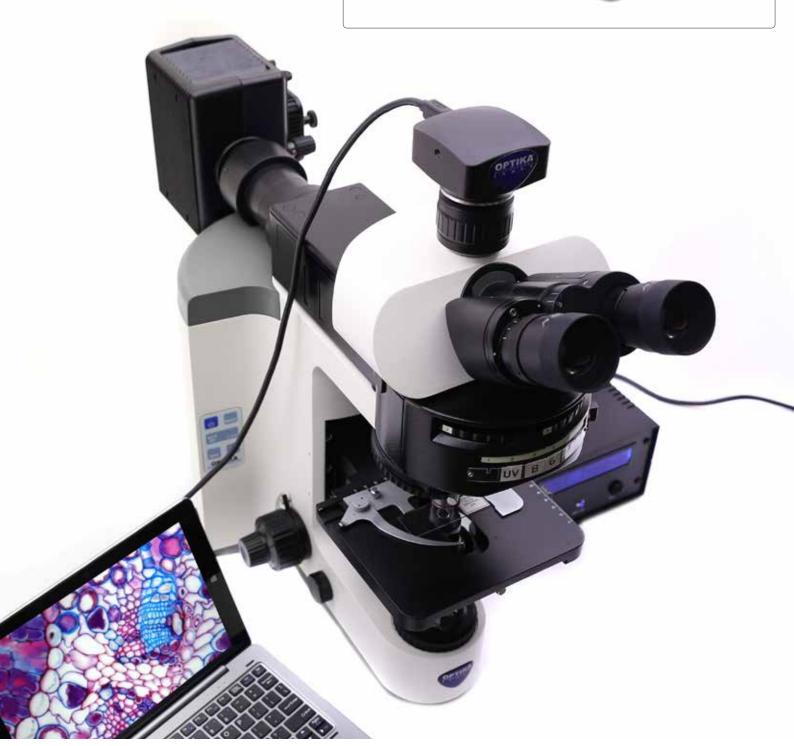

# Cameras & Digital

### C-WH5 & C-WH5SC








Smart and user-friendly dual output (HDMI & Wi-Fi) camera with good resolution (up to 5 MP), high-grade SONY CMOS sensor and HDMI/Wi-Fi connection, recommended for routine operations and whenever measurements are required. No software installation is required with image and video capturing on SD card when in HDMI mode. No router or external applications are required thanks to the quick and simple camera connection. At any time, it can be connected to PC and used via the downloadable, free of charge software (Windows), which is always available to enable the latest updates. The compact and elegantly designed housing conceals the very latest in camera technology. Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately).

Connection into the eyepiece tube or to the trinocular tube of any microscope brand via dedicated adapter (to be purchased separately). Available also as an all-in-one, space saving package including a 11.5" full HD monitor with compact footprint, enabling screen adjustment to ensure correct posture and eliminate fatigue during observation connectable to trinocular tube only via dedicated adapter (to be purchased separately). Wireless mouse, SD card and built-in software included.

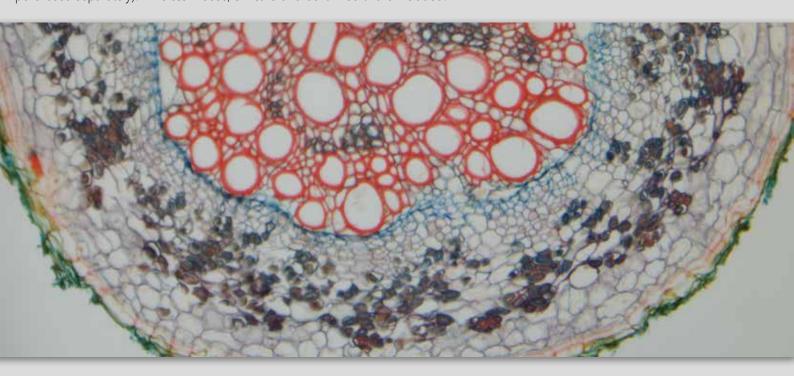

### C-WH5 & C-WH5SC

- » Simple operation, built-in software
- » Universal connection to any microscope brand
- » Crisp 1080p images and videos
- » Live measurements function
- » Dual output mode (HDMI/WiFi)
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (via dedicated adapter, on any microscope)
- » HDMI cable, Wireless adapter, mouse and SD card included
- » Downloadable, free of charge software
- » OPTIKA ProView for Windows



### © C-WH5 & C-WH5SC

#### **Models:**

C-WH5: WH5 camera, 1080p, 5 MP CMOS, Wi-Fi/HDMI, multi-plug

C-WH5SC: WH5 camera, 1080p, 5 MP CMOS, Wi-Fi/HDMI, with screen, multi-plug



### C-WH5 & C-WH5SC - Specifications

	C-WH5 & C-WH5SC
PC Camera resolution (MP)	5 MP
HDMI Camera resolution (MP)	2 MP
Digital camera resolution	1920 x 1080
HDMI Signal output	Yes
Sensor Size	1/1.8"
Sensor technology	CMOS
Sensor type	SONY
Image format	16/9
Pixel size	2.4x2.4 μm
Frame rate (HDMI)	60 fps (1920 x 1080 HDMI); 25fps (1920x1080) (WiFi)
Dynamic range	66 dB
Sensitivity	1120mV at 1/30s
Dark Signal	0.15mV at 1/30s
ADC conversion	8 Bit - 12Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.03 msec - 918 msec
Binning	1x1
IR filter	380-650 nm (IR-cut filter)
Camera power	DC 12V 1A
C-mount	YES

### C-WH5 & C-WH5SC





#### **C-WH5 Contents:**

HDMI camera HDMI cable WiFi adapter SD card Wireless mouse Multi-plug external power supply




#### **C-WH5SC Contents:**

HDMI camera
HDMI monitor
HDMI cable
WiFi adapter
SD card
Wireless mouse
Multi-plug external power
supply (2 pcs.)



### **C-WH5SC** - Monitor Specifications


## MONITOR TECHNICAL SPECIFICATIONS (C-WH5SC) Size 11,5" Power supply 12V / 2,5A HDMI cable 50 cm

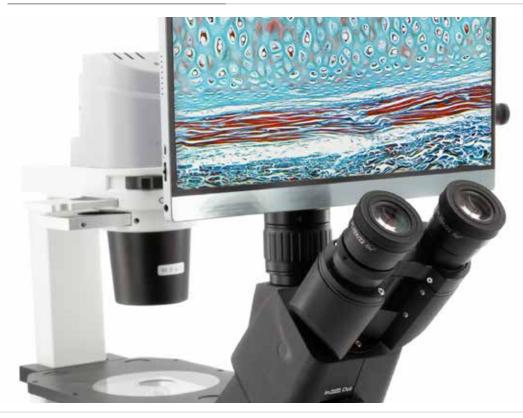



Cameras & Digital

### C-HA

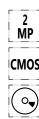






Impressive autofocusing FULL HD camera with 1080p, 2 MP resolution, CMOS sensor and HDMI connection, ensuring precise and ultra-fast automatic focus adjustment in any condition and in real time. Recommended for routine operations and perfect to compensate the lack of parfocality of the microscope without any user effort.

No software installation is required with image and video capturing on SD card.

Camera control panel shows exposure, white balance, color adjustment and sharpness when using mouse control.

Connection to the trinocular tube of any microscope brand via dedicated adapter (to be purchased separately).

Wireless mouse, SD card and built-in software included.

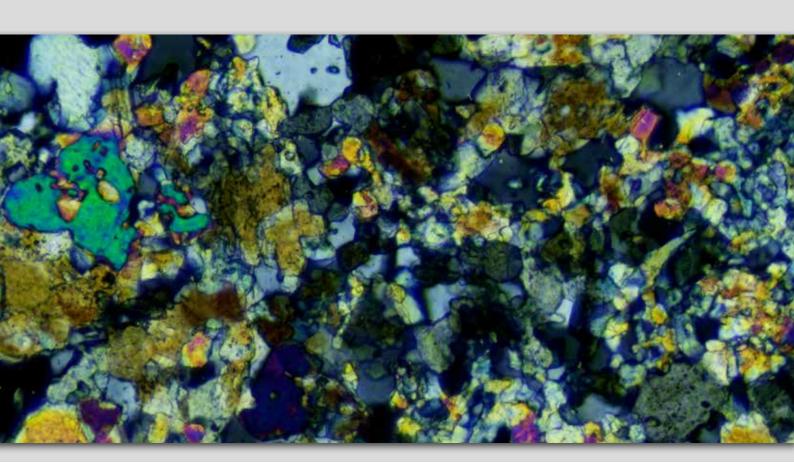

# С-НА

- » Simple operation, built-in software
- » Universal connection to any microscope brand
- » Rapid autofocusing system
- » Crisp 1080p images and videos
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (via dedicated adapter, on any microscope)
- » HDMI cable, mouse and SD card included



### <sup>4</sup> C-HA

#### Model:

C-HA: HA autofocus camera, 2 MP CMOS, HDMI, multi-plug



### **C-HA Series** - Specifications

	с-на
HDMI camera resolution	2 MP (1920 x 1080)
HDMI Signal output	Yes
Sensor Size	1/2.8"
Sensor technology	CMOS
Sensor type	Aptina CMOS
Image format	16/9
Pixel size	2.9 x 2.9 μm
Frame rate full resolution	50 fps (1920 x 1080)
G Sensitivity	510mV at 1/30sec
Dark Signal	0.15mV at 1/30sec
ADC conversion	8 Bit
Color Depth	8 Bit
Exposure Time	Auto
Binning	1x1
Cooling Temperature	None
Camera Power	5V 2A
IR filter	380-650 nm (IR-cut filter)
CS-mount	YES
C-mount	YES

### C-HA



#### **C-PGS Contents:**

USB camera USB 2.0 cable 30 & 30.5 mm rings C-mount projection lens Micrometric slide Wireless mouse SD card



OPTIKA



### C-HP4 & C-HP4SC







Professional ultra-high definition yet very easy to use 4K camera with 2160p, 8 MP resolution, large SONY CMOS sensor and HDMI connection to clearly reveal the sample's finest details on-screen and perform measurements.

No software installation is required with image and video capturing on SD card when in HDMI mode.

At any time, it can be connected to PC and used via the downloadable, free of charge software (Windows), which is always available to enable the latest updates. The compact and elegantly designed housing conceals the very latest in camera technology. Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Wireless mouse, SD card and built-in software included.

- » Simple operation, built-in software
- » Universal connection to any microscope brand
- » Ultra HD 4K 2160p images and videos
- » Dual output mode (HDMI/USB)
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (via dedicated adapter, on any microscope)
- » HDMI cable, mouse and SD card included
- » Downloadable, free of charge software
- » OPTIKA ProView for Windows



### **C-HP4** - Specifications

	C-HP4
PC Camera resolution (MP)	8 MP
HDMI Camera resolution (MP)	8 MP
Camera resolution (n° of pixels: W x H)	3840x2160
Color / Monochrome	Color
Sensor Size	1/1.8"
Sensor technology	CMOS
Sensor type	SONY
Image format	16/9
Pixel size	2.0 x 2.0 μm
Frame rate full resolution	30@ 3840x2160
G Sensitivity	505mV at 1/30s
Dark Signal	0.13mV at 1/30s
ADC conversion	8 Bit - 12 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.04 msec- 2 sec
Binning	1x1
IR filter	380-650 nm (IR-cut filter)
Camera power	12V 1A
C-mount	YES

### C-HP4 & C-HP4SC



#### **C-HP4 Contents:**

4k Camera USB Cable 180cm HDMI Cable 150cm Wireless Mouse 16gb SD Card Calibration slide





# C-HP4SC Contents: 4k Camera USB Cable 180cm HDMI Cable 150cm Wireless Mouse 16gb SD Card Calibration slide

#### Model:

C-HP4: HP4 camera, 4K/8 MP CMOS, USB/HDMI, multi-plug

C-HP4SC: HP4 camera, 4K/8MP CMOS, HDMI/USB, with 4K monitor

### **C-HP4SC** - Monitor Specifications

## MONITOR TECHNICAL SPECIFICATIONS (C-HP4SC) Size 11,5" Power supply 12V / 2,5A HDMI cable 50 cm



### C-HUB4K





Professional ultra-high definition yet very easy to use 4K camera with 2160p, 8 MP resolution, large SONY CMOS sensor and HDMI connection to clearly reveal the sample's finest details on-screen and perform measurements.


No software installation is required with image and video capturing on SD card when in HDMI mode.


At any time, it can be connected to PC and used via the downloadable, free of charge software (Windows), which is always available to enable the latest updates.


The compact and elegantly designed housing conceals the very latest in camera technology.

Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.


Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Wireless mouse, SD card and built-in software included.

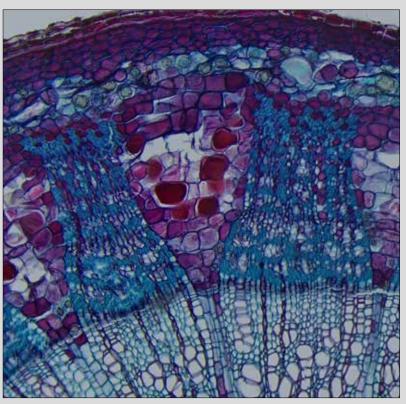

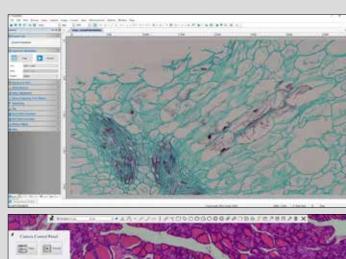

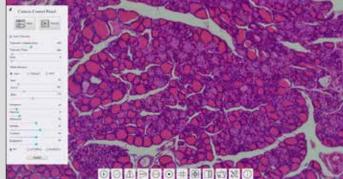

### C-HUB4K

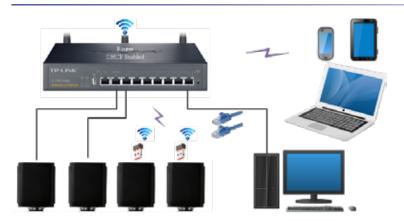

# C-HUB4K Contents: 4k Camera USB Cable 180cm HDMI Cable 150cm Wireless Mouse 16gb SD Card Calibration slide

- » Simple operation, built-in software
- » Universal connection to any microscope brand
- » Ultra HD 4K 2160p images and videos
- » Dual output mode (HDMI/Ethernet)
- » High frame rate
- » Reliable color fidelity
- » C-mount connection (via dedicated adapter, on any microscope)
- » HDMI cable, mouse and SD card included
- » Downloadable, free of charge software
- » OPTIKA ProView for Windows



#### Model:

C-HUB4K: HUB4K Camera,8Mp CMOS, 4K/USB/ETHERNET/WIFI multi video output

### **C-HUB4K** - Specifications

	C-HUB4K
PC Camera resolution (MP)	8 MP
HDMI Camera resolution (MP)	8 MP
Camera resolution (n° of pixels: W x H)	3840x2160
Color / Monochrome	Color
Sensor Size	1/1.8"
Sensor technology	CMOS
Sensor type	SONY
Image format	16/9
Pixel size	2.0 x 2.0 µm
Frame rate full resolution	30@ 3840x2160
G Sensitivity	505mV at 1/30s
Dark Signal	0.13mV at 1/30s
ADC conversion	8 Bit - 12 Bit
Color Depth	1 Bit; 4 Bit; 8 Bit; 24 Bit
Exposure Time	0.04 msec- 2 sec
Binning	1x1
IR filter	380-650 nm (IR-cut filter)
Camera power	12V 1A
C-mount	YES

### Camera Adapters Charts

			Upri	ight	
		Monocular Binocular (Ø 23 mm)	Trinocular (Ø 23 mm)	Binocular (Ø 30 mm)	Trinocular
Camera model	Sensor size	Ecovision / B-60 / B-150 B-190-290 / B-380 (ALC)	B-190 / B-290 B-380 (with M-699)	B-510 / B-810 / B-1000	B-380 / B-510 B-810 / B-1000
C-P3	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620
C-P6	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1
C-P8	1/1.8"	M-115	M-115	M-115 + M-113.1	M-620.1
C-P20	1"	-	-	-	M-620.3
C-P5GS	2/3"	M-118	M-118	M-118 + M-113.1	M-620.2
С-НР	1/1.9"	M-114	M-114	M-114 + M-113.1	M-620.1
C-HPSC	1/1.9"	M-114	M-114	M-114 + M-113.1	M-620.1
C-WH5	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1
C-WH5SC	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1
C-HA	1/2.8"	M-115	M-115	M-115 + M-113.1	M-620
C-HP4	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1
C-HUB4K	1/1.8"	M-114	M-114	M-114 + M-113.1	M-620.1















### **Camera Adapters Charts**

Inverted		Stereo	
Trinocular	Binocular (Ø 30.5 mm)	Binocular (Ø 30 mm)	Trinocular
IM-3 / IM-5	SFX	SLX / SZ / SZP	SLX / SZ / SZP
M-620	M-115 + M-113.2	M-115 + M-113.1	M-620
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620.1
M-620.1	M-115 + M-113.2	M-115 + M-113.1	M-620.1
M-620.3	-	-	M-620.3
M-620.2	M-118 + M-113.2	M-118 + M-113.1	M-620.2
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620.1
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620.1
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620.1
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620.1
M-620	M-115 + M-113.2	M-115 + M-113.1	M-620.1
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620
M-620.1	M-114 + M-113.2	M-114 + M-113.1	M-620



v7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities


OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

Optika Sales branches

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com

# High-End Microscope Cameras for Fluorescence Microscopy

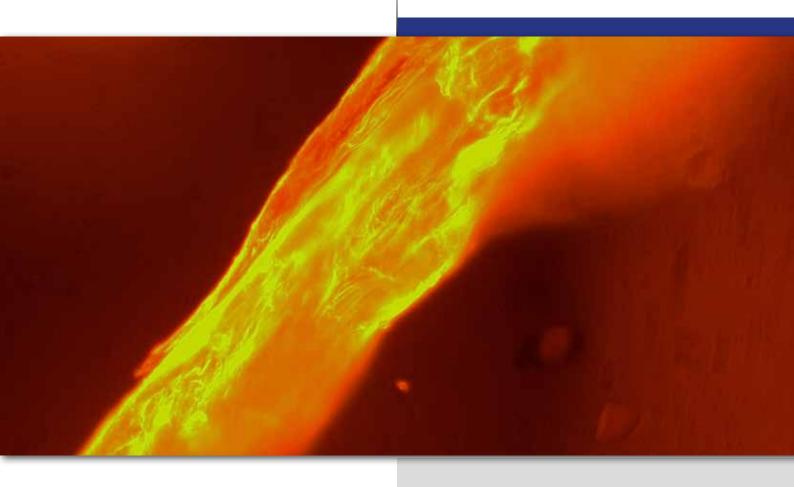

#### Large Selection of Best-in-class High Sensitivity Fluorescence Cameras


In order to detect the often low levels of fluorescence emitted by specimess, cameras used in fluorescence microscopy must have particular features, including high sensitivity and low noise, in order to capture as many photons as possible.

These cameras are typically equipped with CCD sensor, although nowadays there is always a larger selection of scientific-grade CMOS. Monochrome cameras are usually more suited to fluorescence imaging as they do not have a colour filter array, and enable more photons to reach the sensor, increasing their sensitivity very significantly compared to the color sensors.

The most indicated camera depends case by case, as it is of fundamental importance to consider the sample being imaged, the fluorochromes used, the required frame-rate, field of view, resolution and sensitivity.


All these elements drive to the selection of the right camera for a specific use.

### Migh-Level Microscope Cameras

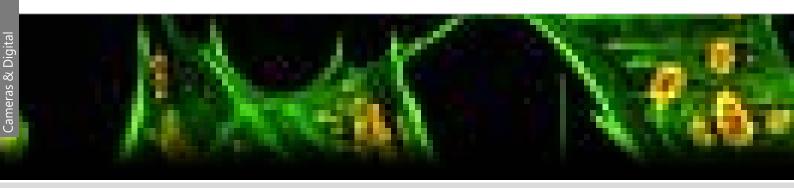

Cooled Scientific-grade CMOS Sensor, Rolling Shutter - C-P20CC & C-P20CM Cooled Scientific-grade CMOS Sensor, Global Shutter - C-P1CCGS & C-P1CMGS CCD Sensor, Rolling Shutter - C-P6FL

Model	C-P20CC	C-P20CM	C-P1CCGS	C-P1CMGS	C-P6FL
Sensor technology	Scientific-grade CMOS	Scientific-grade CMOS	Scientific-grade CMOS	Scientific-grade CMOS	CCD
Color / Monochrome	Color	Monochrome	Color	Monochrome	Color
Global / Rolling shutter	Rolling shutter	Rolling shutter	Global shutter	Global shutter	Rolling shutter
Resolution	20MP (5440 x 3648)	20MP (5440 x 3648)	1.7MP (1600 x 1100)	1.7MP (1600 x 1100)	6MP (2748 x 2200)
Frame rate	5 fps/10 fps/15 fps 30 fps	17.8 fps/41 fps/51 fps/64 fps	33 fps	94 fps	7.5 fps/14 fps
Sensitivity	426 mV at 1/30s	388 mV at 1/30s	4910 mV at 1/30s	8100 mV at 1/30s	1000 mV at 1/30s
Cooling system	Yes	Yes	Yes	Yes	No

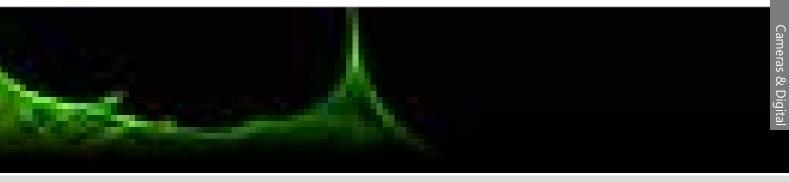

### C-P1CCGS & C-P1CMGS













Superb, stunning global shutter SONY EXMOR CMOS cameras with low resolution (1.7 MP), cooled large sensor and USB3.0 connection, recommended for specific scientific applications especially connected to fluorescence microscopy. Main key-ponits are the global shutter featuring its sensitivity, which makes this series impressive both with fast-moving specimens and in low-light fluorescence, especially thanks to the cooling function. Choose the monochrome version for superb sensitivity, being perfect for fluorescence imaging. The compact and elegantly designed housing conceals the very latest in camera technology. Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them. Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately).

Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.



### **C-P1CCGS & C-P1CMGS**

#### C-P1CCGS

- » Recommended for special applications, including fluorescence
- » Simple operation, driver-free
- » Top-class, large SONY EXMOR sensor
- » Cooling system for enhanced sensitivity
- » Global shutter for impressive performance on moving samples
- » Universal connection to any microscope brand
- » Crisp 1.7 MP images
- » USB3.0 for impressive high frame rate
- » Incredibly accurate colors
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux



### © C-P1CCGS & C-P1CMGS

#### **Models:**

C-P1CCGS: High Performance USB 3.0 C-mount Cooled Color Microscope Camera

C-P1CMGS: High Performance USB 3.0 C-mount Cooled Monochrome Microscope Camera



### **C-P1CCGS & C-P1CMGS** - Specifications

CAMERA TECHNICAL SPECIFICATIONS	C-P1CCGS	C-P1CMGS
Digital camera resolution	1.7 MP (1600 x 1100)	1.7 MP (1600 x 1100)
Signal output	USB 3.0	USB 3.0
Color / Monochrome	Color	Monochrome
Sensor Size	1.1"	1.1"
Sensor technology	CMOS	CMOS
Sensor type	SONY EXMOR	SONY EXMOR
Image format	3/2	3/2
Pixel size	9.0 x 9.0 μm	9.0 x 9.0 μm
Frame rate full resolution	33 fps (1600 x 1100)	94 fps (1600 x 1100)
G Sensitivity	4910mV at 1/30s	8100mV at 1/30s
Dark Signal	0.3mV at 1/30s	0.3mV at 1/30s
ADC conversion	8 Bit - 12 Bit	14 Bit
Color Depth	1 Bit ; 4 Bit; 8 Bit; 24 Bit	-
Exposure Time	0.1 msec - 3600 sec	0.1 msec - 3600 sec
Binning	1x1	1x1
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Cooling Temperature	-45°C	-45°C
Cooling power	12V 3A	12V 3A
Camera power	PC USB	PC USB
C-mount	YES	YES

### C-P1CCGS & C-P1CMGS



#### **C-P1CCGS & C-P1CMGS Contents:**

USB camera USB 3.0 cable Calibration slide External power supply



### C-P6FL






Top-class, easy to operate SONY EXVIEW CCD camera with high resolution (6 MP), large sensor and USB3.0 connection, recommended for specific scientific applications especially connected to fluorescence microscopy. Its particular sensitivity is relevantly important in low-light fluorescence, and if small changes in fluorescence need to be detected, determining how clear the image produced is. The compact and elegantly designed housing conceals the very latest in camera technology.

Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.

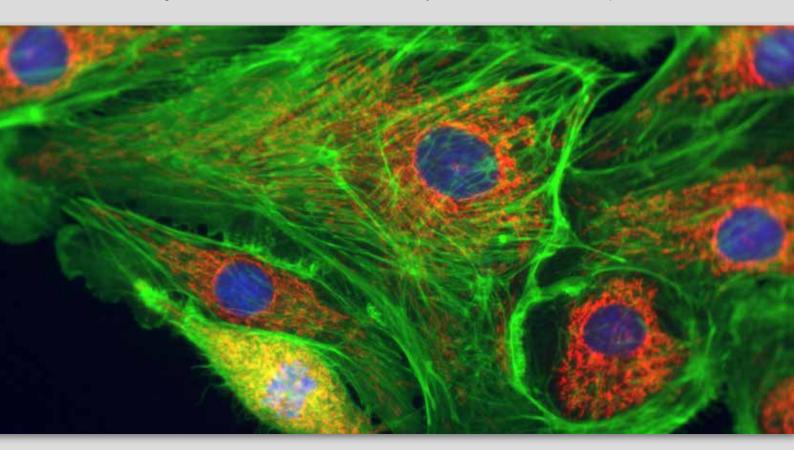

### C-P6FL

- » Recommended for special applications, including fluorescence
- » Simple operation, driver-free
- » Top-class, large SONY EXVIEW CCD sensor
- » Significant sensitivity for a non-cooled camera
- » Universal connection to any microscope brand
- » Crisp 6 MP images
- » USB3.0 for impressive high frame rate
- » Incredibly accurate colors
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux
- » C-mount connection (via dedicated adapter, on any microscope)



### <sup>4</sup> C-P6FL

#### Model:

**C-P6FL:** High Performance USB 3.0 C-mount Microscope Camera



### **C-P6FL** - Specifications


CAMERA TECHNICAL SPECIFICATIONSC-P6FLDigital camera resolution6 MP (2748 x 2200)Signal outputUSB 3.0Color / MonochromeColorSensor Size1"Sensor technologyCCDSensor typeSONY EXVIEWImage format5/4Pixel size4.54 x 4.54 μmFrame rate full resolution7.5 fps (2748 x 2200)Frame rate other resolution14 fps (2748 x 1092)Dynamic range (DB)62G Sensitivity1000mV at 1/30sDark Signal8mV at 1/30sADC conversion8 Bit - 12 BitColor Depth1 Bit; 4 Bit; 8 Bit; 24 BitExposure Time0.06 msec - 1000 secBinning1x1IR filter380-650 nm (IR-cut filter)Camera powerPC USBC-mountYES		
Signal outputUSB 3.0Color / MonochromeColorSensor Size1"Sensor technologyCCDSensor typeSONY EXVIEWImage format5/4Pixel size4.54 x 4.54 μmFrame rate full resolution7.5 fps (2748 x 2200)Frame rate other resolution14 fps (2748 x 1092)Dynamic range (DB)62G Sensitivity1000mV at 1/30sDark Signal8mV at 1/30sADC conversion8 Bit - 12 BitColor Depth1 Bit; 4 Bit; 8 Bit; 24 BitExposure Time0.06 msec - 1000 secBinning1x1IR filter380-650 nm (IR-cut filter)Camera powerPC USB	CAMERA TECHNICAL SPECIFICATIONS	C-P6FL
Color / Monochrome Sensor Size 1' Sensor technology CCD Sensor type SONY EXVIEW Image format 5/4 Pixel size 4.54 x 4.54 µm Frame rate full resolution 7.5 fps (2748 x 2200) Frame rate other resolution 14 fps (2748 x 1092) Dynamic range (DB) 62 G Sensitivity 1000mV at 1/30s Dark Signal 8mV at 1/30s ADC conversion 8 Bit - 12 Bit Color Depth 1 Bit; 4 Bit; 8 Bit; 24 Bit Exposure Time 0.06 msec - 1000 sec Binning 1x1 IR filter 380-650 nm (IR-cut filter) Camera power	Digital camera resolution	6 MP (2748 x 2200)
Sensor Size         1"           Sensor technology         CCD           Sensor type         SONY EXVIEW           Image format         5/4           Pixel size         4.54 x 4.54 μm           Frame rate full resolution         7.5 fps (2748 x 2200)           Frame rate other resolution         14 fps (2748 x 1092)           Dynamic range (DB)         62           G Sensitivity         1000mV at 1/30s           Dark Signal         8mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Signal output	USB 3.0
Sensor technology         CCD           Sensor type         SONY EXVIEW           Image format         5/4           Pixel size         4.54 x 4.54 μm           Frame rate full resolution         7.5 fps (2748 x 2200)           Frame rate other resolution         14 fps (2748 x 1092)           Dynamic range (DB)         62           G Sensitivity         1000mV at 1/30s           Dark Signal         8mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Color / Monochrome	Color
Sensor typeSONY EXVIEWImage format5/4Pixel size4.54 x 4.54 μmFrame rate full resolution7.5 fps (2748 x 2200)Frame rate other resolution14 fps (2748 x 1092)Dynamic range (DB)62G Sensitivity1000mV at 1/30sDark Signal8mV at 1/30sADC conversion8 Bit - 12 BitColor Depth1 Bit ; 4 Bit; 8 Bit; 24 BitExposure Time0.06 msec - 1000 secBinning1x1IR filter380-650 nm (IR-cut filter)Camera powerPC USB	Sensor Size	1"
Image format         5/4           Pixel size         4.54 x 4.54 μm           Frame rate full resolution         7.5 fps (2748 x 2200)           Frame rate other resolution         14 fps (2748 x 1092)           Dynamic range (DB)         62           G Sensitivity         1000mV at 1/30s           Dark Signal         8 mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Sensor technology	CCD
Pixel size       4.54 x 4.54 μm         Frame rate full resolution       7.5 fps (2748 x 2200)         Frame rate other resolution       14 fps (2748 x 1092)         Dynamic range (DB)       62         G Sensitivity       1000mV at 1/30s         Dark Signal       8mV at 1/30s         ADC conversion       8 Bit - 12 Bit         Color Depth       1 Bit; 4 Bit; 8 Bit; 24 Bit         Exposure Time       0.06 msec - 1000 sec         Binning       1x1         IR filter       380-650 nm (IR-cut filter)         Camera power       PC USB	Sensor type	SONY EXVIEW
Frame rate full resolution 7.5 fps (2748 x 2200)  Frame rate other resolution 14 fps (2748 x 1092)  Dynamic range (DB) 62  G Sensitivity 1000mV at 1/30s  Dark Signal 8mV at 1/30s  ADC conversion 8 Bit - 12 Bit  Color Depth 1 Bit ; 4 Bit; 8 Bit; 24 Bit  Exposure Time 0.06 msec - 1000 sec  Binning 1x1  IR filter 380-650 nm (IR-cut filter)  Camera power	Image format	5/4
Frame rate other resolution  14 fps (2748 x 1092)  Dynamic range (DB)  62  G Sensitivity  1000mV at 1/30s  Dark Signal  8mV at 1/30s  ADC conversion  8 Bit - 12 Bit  Color Depth  1 Bit ; 4 Bit; 8 Bit; 24 Bit  Exposure Time  0.06 msec - 1000 sec  Binning  1x1  IR filter  380-650 nm (IR-cut filter)  Camera power	Pixel size	4.54 x 4.54 μm
Dynamic range (DB)         62           G Sensitivity         1000mV at 1/30s           Dark Signal         8mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Frame rate full resolution	7.5 fps (2748 x 2200)
G Sensitivity         1000mV at 1/30s           Dark Signal         8mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Frame rate other resolution	14 fps (2748 x 1092)
Dark Signal         8mV at 1/30s           ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Dynamic range (DB)	62
ADC conversion         8 Bit - 12 Bit           Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	G Sensitivity	1000mV at 1/30s
Color Depth         1 Bit ; 4 Bit; 8 Bit; 24 Bit           Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	Dark Signal	8mV at 1/30s
Exposure Time         0.06 msec - 1000 sec           Binning         1x1           IR filter         380-650 nm (IR-cut filter)           Camera power         PC USB	ADC conversion	8 Bit - 12 Bit
Binning 1x1  IR filter 380-650 nm (IR-cut filter)  Camera power PC USB	Color Depth	1 Bit ; 4 Bit; 8 Bit; 24 Bit
IR filter 380-650 nm (IR-cut filter)  Camera power PC USB	Exposure Time	0.06 msec - 1000 sec
Camera power PC USB	Binning	1x1
The stage of the s	IR filter	380-650 nm (IR-cut filter)
<b>C-mount</b> YES	Camera power	PC USB
	C-mount	YES

4

# Cameras & Digital





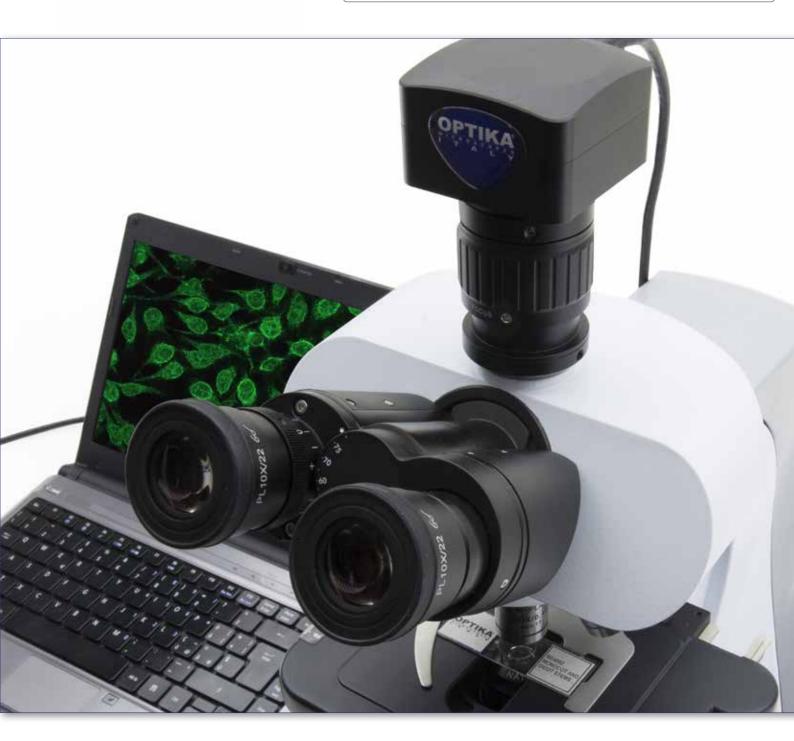

### C-P20CC & C-P20CM
















Ultra-professional yet intuitive cameras with incredibly high resolution (20 MP), cooled large SONY EXMOR CMOS sensor and USB3.0 connection, recommended for specific scientific applications especially connected to fluorescence microscopy. The cooling system affects sensitivity and therefore image clarity. Choose the monochrome version for superb sensitivity, being perfect for fluorescence imaging. The high resolution makes these models interesting also for morphological imaging, patch clamping and network studies.

The compact and elegantly designed housing conceals the very latest in camera technology.

Images will be of the highest quality and rich in contrast and detail with the top-class SONY sensors, worldwide recognized, ensuring beautiful true-to-life color and delivering incredibly accurate colors just as you see them.

Ideal to be connected to all the trinocular tube of different brands using the focusable C-Mount adapter (to be purchased separately). Downloadable, free of charge software (Windows, Mac OS or Linux) is always available to enable the latest updates.

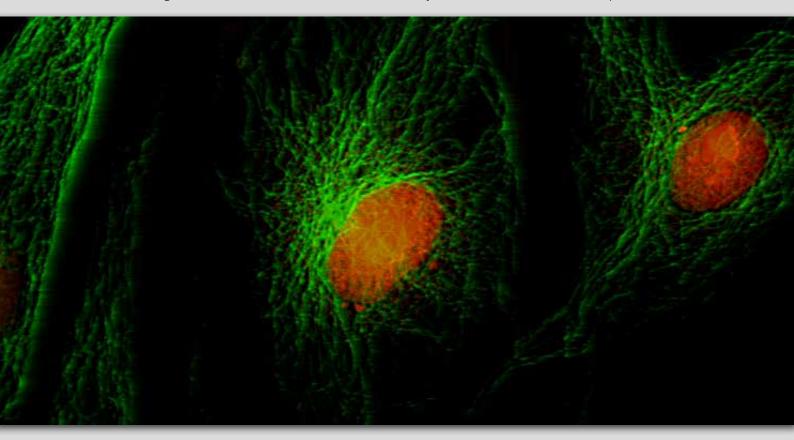

### **C-P20CC & C-P20CM**

#### C-P20CC

- » Recommended for special applications, including fluorescence
- » Simple operation, driver-free
- » Top-class, large SONY EXMOR sensor
- » Cooling system for enhanced sensitivity
- » Universal connection to any microscope brand
- » Crisp 20 MP images
- » USB3.0 for impressive high frame rate
- » Incredibly accurate colors
- » USB cable and calibration slide included
- » Downloadable, free of charge software
- » OPTIKA ProView & LiteView for Windows
- » OPTIKA LiteView for Mac OS or Linux

### 

#### **Models:**

C-P20CC: Pro Cooled Color camera, 20 MP CMOS, USB3.0

C-P20CM: Pro Cooled Monochromatic camera, 20 MP CMOS, USB3.0



### C-P20CC & C-P20CM - Specifications

CAMERA TECHNICAL SPECIFICATIONS	C-P20CC	C-P20CM
Digital camera resolution	20 MP (5440 x 3648)	20 MP (5440 x 3648)
Signal output	USB 3.0	USB 3.0
Color / Monochrome	Color	Monochrome
Sensor Size	1"	1"
Sensor technology	CMOS	CMOS
Sensor type	SONY EXMOR	SONY EXMOR
Image format	3/2	3/2
Pixel size	2.4 x 2.4 µm	2.4 x 2.4 μm
Frame rate full resolution	5 fps (5440 x 3648)	17.8 fps (5440 x 3648)
Frame rate other resolution	10 fps (4096x2160); 15 fps (2736x1824); 30fps (1824x1216)	41 fps (4096 x 2160); 51 fps (2736x1824); 64 fps (1824x1216)
G Sensitivity	426mV at 1/30s	388mV at 1/30s
Dark Signal	0.21mV at 1/30s	0.21mV at 1/30s
ADC conversion	8 Bit - 12 Bit	14 Bit
Color Depth	1 Bit ; 4 Bit; 8 Bit; 24 Bit	-
Exposure Time	0.1 msec - 3600 sec	0.1 msec - 3600 sec
Binning	1x1; 2x2; 3x3	1x1; 2x2; 3x3
IR filter	380-650 nm (IR-cut filter)	380-650 nm (IR-cut filter)
Cooling Temperature	-45°C	-45°C
Cooling power	12V 3A	12V 3A
Camera power	PC USB	PC USB
C-mount	YES	YES

### **C-P20CC & C-P20CM**



#### C-P20CC & C-P20CM Contents:

USB camera USB 3.0 cable Calibration slide External power supply





### Camera Adapters Charts

			<b>Upright</b>					
		Monocular Binocular (O 23 mm)	Trinocular (O 23 mm)	Binocular (O 30 mm)	Trinocular			
Camera model	Sensor size	Ecovision / B-60 / B-150 B-190-290 / B-380 (ALC)	B-190 / B-290 B-380 (with M-699)	B-510 / B-810 / B-1000	B-380 / B-510 B-810 / B-1000			
C-P20CC	1"	-	-	-	M-620.3			
C-P20CM	1"	-	-	-	M-620.3			
C-P1CCGS	1,1"	-	-	-	M-620.3			
C-P1CMGS	1,1"	-	-	-	M-620.3			
C-P6FL	1"	-	-	-	M-620.3			



### **Camera Adapters Charts**

Inverted		Stereo	
Trinocular	Binocular Binocular (O 30.5 mm) (O 30 mm)		Trinocular
IM-3 / IM-5	SFX	SLX / SZM / SZ / SZP	SLX / SZM / SZ / SZP
M-620.3	-	-	M-620.3
M-620.3	-	-	M-620.3
M-620.3	-	-	M-620.3
M-620.3	-	-	M-620.3
M-620.3	-	-	M-620.3



v7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA' S.r.I. Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com


Optika Sales branches

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



### OPTIKA SOFTWARE

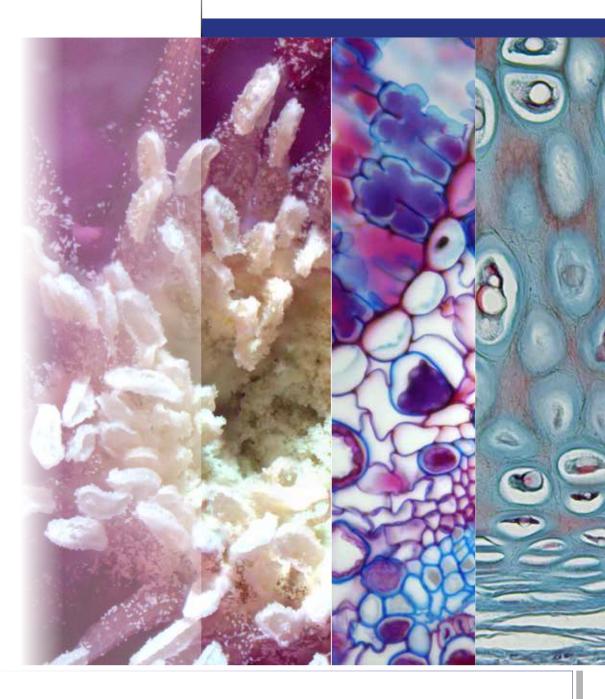

**OPTIKA SOFTWARE SUITES** 

### **OPTIKA SOFTWARE - Comparison chart**

#### **Software**

• Before proceeding with the SW installation, please check the table below "Software Function list" to identify the most suitable software.







#### **SOFTWARE FUNCTION LIST**

	FUNCTION		OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
	Simultaneous management of several cameras  GUI (Graphical User Interface)		X	X	X
			X		
	Report generator		X		X
	Archiving		X	Х	X
		Catalan	X	X	
		Chinese (simpl.)	X	X	
		Chinese (trad.)	X	X	
		Korean	X	X	
		English	X	X	X
GENERA		French	X	X	X
Z		German	X	X	X
5	Language	Indonesian	X	Х	
	Language	Italian	X	X	X
		Japanese	X	Х	
		Polish	X	Х	X
		Russian	X	Х	
		Spanish	X	X	X
		Swedish			Х
		Thai	Х	X	
		Turkish	X	Χ	

	FUNCTION		OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
	Measurements on "live"		X		
	Measurements on "captured"		X		X
	2D Measurements	Line	X		Х
S		Angle	X		
Ę		Parallel lines	X		
		Rectangle	X		
MEASUREMENT		Ellipse	X		
		Circle	X		
		Annulus	X		
		Arc	X		
		Curve	X		
		Polygon	Х		
	Particle count		X		
	Export to Excel		X		Х

# Cameras & Digital

### **OPTIKA SOFTWARE -** Comparison chart

#### **SOFTWARE FUNCTION LIST**

<b>FUNCTION</b> Simultaneous management of several cameras		OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
	tiff	X	Х	Х
	jpg	X	Х	X
	bmp	X	Х	X
Image formats	png	X	Х	
	рсх	X	Х	
	jp2	X	Х	
	dcm	X	Х	
IMAGE acquisition		X	X	X
	avi	X	Х	X
	wmv	X	Х	X
	mp4	X	Х	X
VIDEO formats	asf	X	X	X
VIDEO IOITIAIS	3gp	X	Х	X
	mov	X	Х	X
	h264	X	Х	X
	h265	X	Х	X
Continuous automatic exposure		X	X	X
Manual Exposure		X	X	X
Mobile spot for exposure		X	X	X
Resizable spot for exposure		X	X	X
Colour acquisition		X	X	X
Grey-scale acquisition		X	X	Х
Manual Time-Lapse		X		X
Automatic Time-Lapse		X		
Fast Image Acquisition		X	X	X
Focus Indicator		X		
White Balance		X	X	X
Black balance		X		
Background correction		X		
Dark Field Correction		X	X	
Image Enhancement		X	X	X
Live Histogram		X	X	X
Flip	Horizontal	X	X	Х
ΠΡ	Vertical	Х	X	Х
Rotate		X		

O
Z
10
G
S
П
C
0
02

FUNCTION	OPTIKA PRO VIEW	OPTIKA LITE VIEW	OPTIKA VISION LITE
Several function of image processing (filters)	X		
Multiple image combining	X		
EDF (Extended Depth of Focus)	X		
Colour Combine (Multi-Fluorescence Imaging)	X		
Shift Correction	X		
HDR (High Dynamic Range)	X		
Layer Management	X		X
Text Overlay	X		
Ruler Overlay	X		
Measurement Overlay	X		
Grids	X		X

4083.Wifi, 4083.4 and 4083.EC2 work with Vision Lite only. Cameras with HDMI connection only, do not require any software.

### OPTIKA Vision Lite - Extremely Intuitive Software


Optika Vision Lite has been designed and developed to be incredibly intuitive, simple and easy to use for customers needing a convenient solution to be combined with OPTIKAM cameras.

- » Friendly interface, multilanguage
- » Capture still images & stream live videos
- » Perform linear measurements
- » Export comprehensive reports

#### Friendly interface, multilanguage

Engineered for easy user interaction and optimized image acquisition, the main purpose of OPTIKA Vision Lite is ensure clear communication.

- •An efficient means to efficiently completing your jobs
- •Pleasant, easy-to-navigate menus
- Eight languages pre-installed, others upgreadable

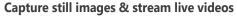

Use the live preview to accurately focus your image and change parameters to obtain the perfect final result you are looking for. Images can be saved in different formats and even as test reports, including personal comments.

Additional features:

- Image stack acquisition
- Grid addition for rapid considerations
- Image flipping option available


#### **Perform linear measurements**


Perform linear measurements in an extremely way just by drawing a line after creating your preferred calibration based on the magnification.


- Accurate measurements through simple calibration
- Comprehensive data export (notes & measures included)
- Indicate particular objects in the image to add persona comments

#### **Export comprehensive reports**

Detailed test reports can be generated, printed and saved. Reports can be also customized with company logos.

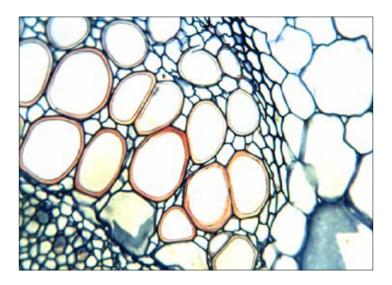

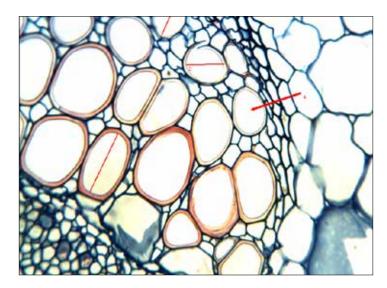

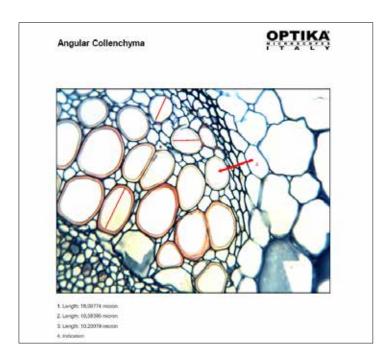

### **OPTIKA LITEView - Life is Easier**

**OPTIKA LITEView** is a basic image acquisition software. The user who simply wants acquire a still image or a video, with no no need to perform measurements, has, with this powerful and intuitive software, the perfect solution.

- -) Simple management of «live» image
- -) Acquisition of still images or video
- -) Basic imaging functions
- -) Background correction

#### Simple management of «live» image

Image preview is freely customizable by the user. A simple White Balance function with a mobile spot allows to perform the balance even on very small areas, once the specimen has been framed and focused.


#### Basic functions:

- Automatic or manual acquisition
- Possibility to have «live» and «capture» at different resolutions
- White Balance with mobile spot
- Background correction for the acquisition of perfectly illuminated images.


#### Capturing still images or video

Just select the option and the software performs: acquiring still images or videos is simply and intuitive.



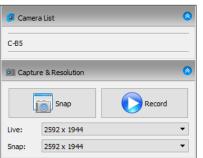



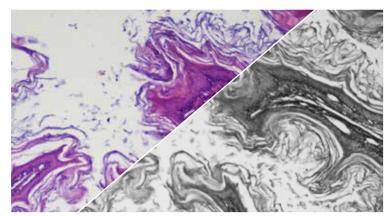

Color / Grey scales


#### **Basic imaging functions**

Image parameters can be modified according user's needs. Color, Contrast and Gamma can be chaned in real time. More, it is possible to use a color camera in «SGrey Scales» modo in order to increase the camera sensitivity.

#### **Background Correction**

Any inhomogeneity of illumination of the microscope can be corrected by using the background correction function.

This allows to obtain a faithful reproduction of the image without annoying inhomogeneity due to a not perfect illumination.

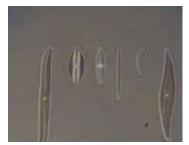

No Background correction



With Background correction

### OPTIKA PROView - Professional Image Analysis

**OPTIKA PROVIEW** is a professional image analysis software. The user who needs to acquire an image or video and to perform a series of processings or measurements, can easily achieve incredible results thanks to this software. PROView incorporates all the functions of the LITEView package, but in addition allows:

- White Balance and Black Balance
- Simultaneous management of several cameras
- Graphical User Interface fully customizable
- Imaging of Multichannel Fluorescence Images with «pixel shift» function
- Multilanguage Software

#### **Beginners? Experts?**


An «On-line» manual will help any user (no matter on how expert he can be) to get the best from the software


#### **Images always perfect**

The management of the acquisition parameters allows to get always the best from your camera. White balance, black balance, background correction, «live» management of Colors, Contrast, Gamma, Gain and Exposure Time ensure to obtain a faithful image. A numerical focus indicator will ensure an optimal focusing, also on specimaens with different focal planes.



It is possible to obtain the balance either on the whole frame or on a small ROI (Region Of Interest) of the image simply resizing and moving the spot in one part of the specimen

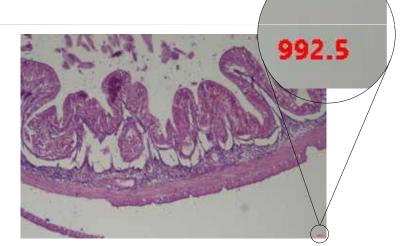

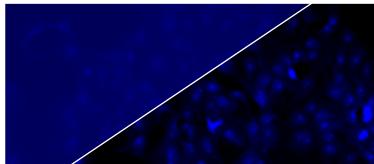


No black correction / Black correction

#### **Multichannel Fluorescence Image processing**


Acquire fluorescence images with a specific filtercube, use a false color for the used fluorochrome, get a single multichannel image is simply and intuitive.




DAPI (UV)

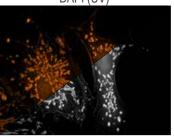

TEXAS RED (G)

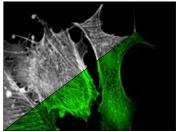

FITC (B)

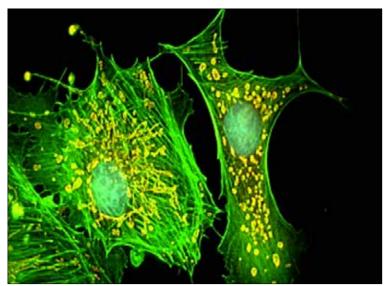

Combined multichannel image


### **OPTIKA PROView - Professional Image Analysis**

#### «Pixel Shift» function


Fluorescence ilter cubes, sometimes, are not perfectly aligned.

During acquisition of multichannel luorescence images, this can cause a non perfect overlapping of the different signals, making the colocalization calculation almost impossible.

«Pixel Shift» function allows to correct these small misalignments:

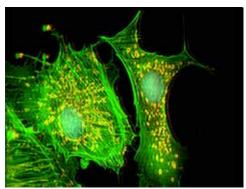

Orignal image

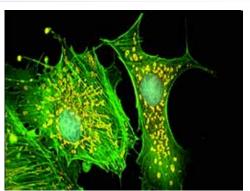


Corrected image

#### HDR (High Dynamic Range) acquisition

Acquisition of different images with different exposure times allows this function to create a final image where bright and dark zones of the specimen are perfectly displayed.



Standard Dynamic Range

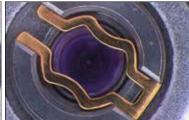


High Dynamic Range

#### **Extended Depth of Focus (EDF)**

Acquire images with different focal planes, specially on specimens observed under a stereomicroscope, and to obtain a focused final image with a theoretical infinite focus. EDF function (also known as «Z-stack») allows a very refined image processing.





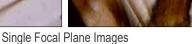





#### Stitching & Tiling


Get an image with high resolution but, at the same time, have a wide view of the specimen under observation. Impossible? No. The multiple image alignment function allows to get a singe image starting from adjacent images of the specimen.



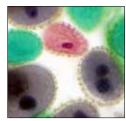

Separate Images

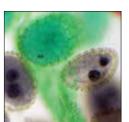

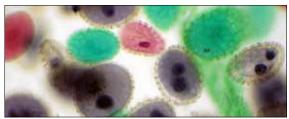


Stitched image

### **OPTIKA PROView - Professional Image Analysis**

#### Measurements

User can perform measurements on the «live» image (no need to capture an image) and on captured images.

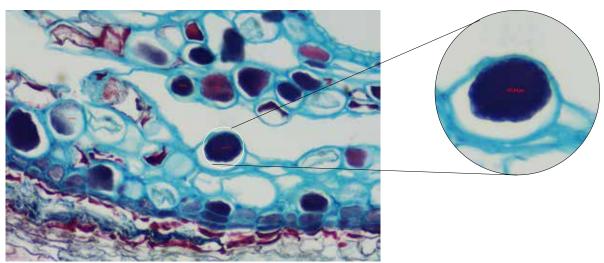

#### **From Beginners To Experts**


Measurements available:

- · linear measurements
- angles
- circles
- annuli
- poligons
- touch count

#### **Report Generator**

At the end of the analysis it is possible to export images and measurement results either on a Excel sheet and on a Report Generator in MS Word format.

The template is freely configurable and can be modified according to laboratory standards.

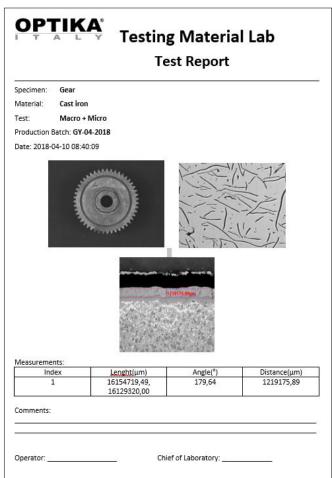

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

#### **Optika Sales branches**

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com





Cameras & Digital

 $v\,7.5 - OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**


**OPTIKA® S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com



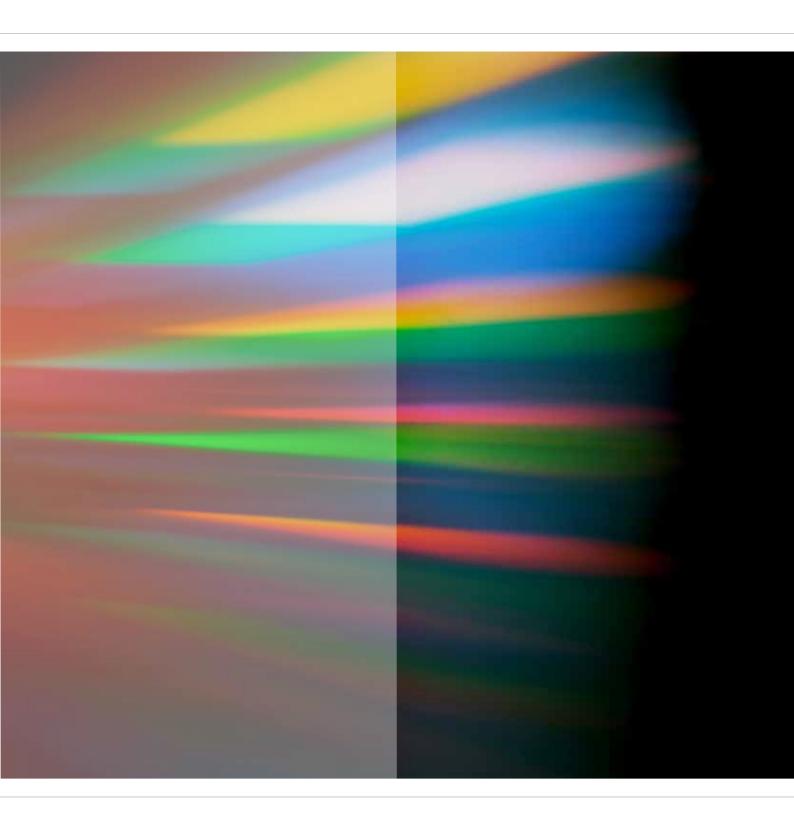

### **POLARIMETRY AND REFRACTOMETRY**

# Extremely Realiable Polarimeter for Practical Experience

# THE STATE-OF-THE-ART POLARIMETER Combining wide measuring range with simplicity

- » Entirely Made in Italy
- » Outstanding quality/price ratio
- » New design for a compact and space-saving solution

# MONOCHROMATIC LED LIGHT SOURCE Instead of traditional sodium bulb

- » Unparalleled illumination
- » 30x longer lifespan than conventional polarimeters
- » Plug & play, no heating time required



### **Polarimeter - POL-X**

The POL-X polarimeter from Optika precisely determines the optical rotation of substances, quickly and easily. POL-X measures the optical rotation over the entire measuring range at the same high accuracy. It's suitable for many applications, from routine measurements to demanding experiments for innovative projects. With this new polarimeter, OPTIKA meets the requirements of the pharmaceutical, cosmetics, chemical and medical industries and creates an ideal tool for R&D applications.



### **Polarimeter - POL-X**

### **Technical Specifications**

Measuring range of optical rotation: ± 180°

**Resolution:** 1° **Accuracy:** 0.05°

Magnification factor of the glass: 4x

**Light source:** Monochromatic LED, 1.2 W, Id = 590 nm (equivalent to sodium lamp)

Length of test tubes: 100 mm and 200 mm

Weight: 1.7 Kg

Overall dimensions: 450x180x320 mm

**Accessories**:

POL-1.2 - 100 mm polarimeter tube POL-1.3 - 200 mm polarimeter tube



# **Target Industries**



### Food

Concentration and purity of compounds (carbohydrates, lactose, raffinose, various starches, fructose, levulose, sucrose, naturalmonosaccharides, glucose, maltrose, xylose) in sugar based foods, cereals and syrups.

### Chemical

Identification and characterization of polymers (bio-, natural or synthetic), foods, cereals and syrups.





### Pharmaceutical

Determination of product purity by measuring specific and optical rotation of amino acids, antibiotics, dextrose, steroids, amino sugars, cocaine, diuretics, tranquilizers, analgesics, codeine, serums, vitamins, etc.

# Fragrance, Flavor & Essential Oils

Inspection of incoming raw materials, such as camphors, gums, orange oil, citric acid, lavender oil, spearmint oil, glygericacid, lemon oil, etc.



# Refractometers - HR & HRD Series




OPTIKA offers a wide range of hand refractometers to ensure incredibly rapid and convenient measurement of concentration in liquid and semi-solid samples, combining accurate performance with excellent repeatability.

The refractive index of a substance is related to its specific density; a refractometer is used to measure the purity or the concentration of a sample when mixed.

Refractometers are suitable for a wide range of applications, which includes the control of blend ratios in light industrial applications (such as glycols, battery acid, heat exchange fluids, coolants, quenchants and hydraulic oils) and represents the ideal solution for users working in the food industry (fruit, beverages, confectionery, jam, honey and other sugar based products).

- » Easy operation
- » Built-in LED illumination
- » Brix or scale specific (including ATC)
- » Sturdy construction with rubber handgrip

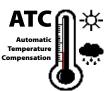

**ATC** - Automatic Temperature Compensation
No need to worry about temperature change during your measurement.



## **HR** Series - Hand Refractometers

#### HR-110N - Hand Refractometer, 0-20% Brix, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy	
0-20% Brix	0.1 Brix	±0.10 Brix	





#### HR-130N - Hand Refractometer, 0-32% Brix, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy
0-32% Brix	0.2 Brix	±0.2 Brix





#### HR-140N - Hand Refractometer, 0-32% Brix, ATC, Triple Scale, Built-in LED Illuminato

Measure Range	Min. Div	Accuracy
0-32% Brix	0.2 Brix	±0.2 Brix
0-140 °Oe	1°Oe	±1°Oe
0-27 KMW(Babo)	0.2 KMW(Babo)	±0.2KMW(Babo)





#### HR-150N - Hand Refractometer, 0-80% Brix, Single Scale, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy	
0-80% Brix	0.5 Brix	±0.5 Brix	





#### HR-160N - Hand Refractometer for Urine and Protein, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy
1.000-1.050 sg	0.002 sg	± 0.002 sg
1.3300 RI - 1.3600 RI	0.0005 RI	± 0.0005 RI
0-12 g/dl	0.2 g/dl	± 0.2 g/dl





#### HR-170N - Hand Refractometer for Battery and Antifreeze, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy
E: -60° F - 32° F	10° F	± 10°F
P: -50° F - 32° F	10° F	± 10°F
B: 1.15-1.30 sg	0.01 sg	± 0.1sg





#### HR-180N - Hand Refractometer for Wine and Grape, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy	
0-80% v/v	1,0% v/v	±1,0 v/v	





#### HR-190N - Hand Refractometer, 0-28% Salinity, ATC, Built-in LED Illuminator

Measure Range	Min. Div	Accuracy	
0-28% Salinity	0.2%	±0.2%	





# **HRD** Series - Digital Refractometers

#### HRD-300N - Digital Refractometer, Brix and Refractive Index

Measure Range	Min. Div	Accuracy
Brix 0-50%	0.10%	±0.2 %
Refractive index: 1.3330-1.4200nD	0.0001 nd	± 0.0003 nD



#### HRD-400N - Digital Refractometer, Brix, Salinity and Refractive Index

Measure Range	Min. Div	Accuracy
Brix 0-50%	± 0.1 %	± 0.2 %
Salinity 0.0-28.0%	± 0.1%	± 0.2 %
Refractive index: 1.3330-1.4200 nD	± 0.0001 nD	± 0.0003 nD



#### HRD-500N - Digital Refractometer for Urine, Serum Protein and Refractive Index

Measure Range	Min. Div	Accuracy
Urine SP.G 1.000-1050	0.001	±0.002
SERUM P. 0-12 g/dl	0.1	±0.2
Refractive index:	±0.0001 nD	±0.0003 nD
1.3330-1.3900 nD		







# **Refractometers - 2WAJ ABBE Bench Refractometer**



Abbe bench-top refractometer measures refractive index ND and average color dispersion NF-NC of transparent and semi-transparent liquid or solid samples. It is a must-have equipment in factories, teaching institutes and science research centers related to oil, grease, pharmacy, painting, food, sugar-refining and geological industries, among the others. It is unique thanks to its ability in measuring also solid samples, such as film, glass, and other transparent materials.

Main prism: horizontal

**Secondary prism:** fitted on a hinge **Refraction index scale:** Nd 1.300 – 1.700

**Precision:** Nd ± 0.0003 **Division:** Nd 0.0005

**Sugars scale:** 0-95% from Nd 1.300 – 1.530 **Precision:** 0-50% = 0.2%; 51-95% = 0.1%

Division: 0.25%

Thermometer scale: 0°C - 70°C, div. 0.1°C

Weight: 4 Kg

**Size:** 140x100x235 mm



 $v\,7.5-OPTIKA\ reserves\ the\ right\ to\ make\ corrections,\ modifications,\ enhancements,\ improvements\ and\ other\ changes\ to\ its\ products\ at\ any\ time\ without\ notice.$ 

#### **Headquarters and Manufacturing Facilities**

**OPTIKA**° **S.r.I.** Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com









**OPTIKA Microscopes** is committed to pursuing the full realization and continuous evolution of the Quality Management System in order to consolidate and improve the Company image, also through the commitment and professionalism of all Company staff.

v 7.5 - OPTIKA reserves the right to make corrections, modifications, enhancements, improvements and other changes to its products at any time without notice.

Headquarters and Manufacturing Facilities

OPTIKA S.r.I.

Via Rigla, 30 - 24010 Ponteranica (BG) - ITALY - Tel.: +39 035.571.392 - info@optikamicroscopes.com

**Optika Sales branches** 

OPTIKA° Spain OPTIKA° China OPTIKA° India spain@optikamicroscopes.com china@optikamicroscopes.com india@optikamicroscopes.com **OPTIKA**° North America **OPTIKA**° Central America **OPTIKA**° Africa

namerica@optikamicroscopes.com camerica@optikamicroscopes.com africa@optikamicroscopes.com